

TS 144 – HAEMOPHILUS GROWTH SUPPLEMENT

INTENDED USE

For cultivation of *Haemophilus influenzae*.

COMPOSITION

Ingredients	Concentration
NAD	7.500mg
Haematin	7.500mg

(per vial sufficient for 500 ml medium)

Applied in Medium

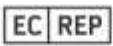
TM 1352 - HAEMOPHILUS TEST AGAR BASE

INSTRUCTION FOR USE

Rehydrate the contents of one vial aseptically with 5.0 ml sterile distilled water. Mix well and aseptically add to 500 ml sterile, molten, cooled (45-50°C) TM 1352 - HAEMOPHILUS TEST AGAR BASE. Mix well and aseptically pour into sterile petri plates.

STORAGE

Vials should be stored in sealed container at 2-8°C.

 GMP Good Manufacturing Practices Certified	 IVD For In Vitro Diagnostic Use	 QTY. Quantity	 LOT/ B. NO. Lot / Batch Number	 REF Catalogue Number	 Manufacturer
 Temperature Unit	 EC REP Authorized Representative Mallin+Gottlieb Buckhorn 10, 40142 Düsseldorf, Germany	 European Conformity	 QR Code	 Consults Instructions for Use	 Best Before

NOTE: Please consult the Material Safety Data Sheet for information regarding hazards and safe handling Practices.

***For Lab Use Only**
Revision: 2 March., 2022