



Chromogenic Culture Media

For Rapid Differentiation & Accurate Identification

Traditional Culture Media supports microbial growth. **Chromogenic Culture Media** go a step further — enabling rapid differentiation and presumptive identification of microorganisms through distinct colony coloration.


Chromogenic Culture Media combines high-performance nutrients with specialized chromogenic substrates such as X-Gal, X-Glucuronides, Salmon-GAL that react with microbial enzymes to produce characteristic coloured colonies — simplifying routine microbiological analysis.

Key Benefits

- Faster identification within 18–24 hours
- Reduces biochemical confirmation testing
- Saves laboratory time & operational cost
- Supports detection of multiple organisms on one medium
- Simplifies interpretation and reporting
- High selectivity, sensitivity & consistency
- Consistent batch-to-batch performance
- Quality-compliant manufacturing standards

Culture Media & Supplements


Chromogenic Culture Media (Most Acceptable Products)

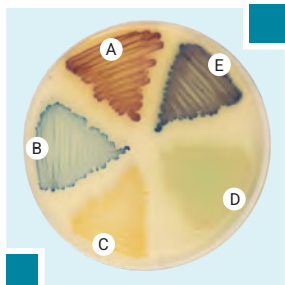

TM 1199 | CHROMOGENIC UTI AGAR

For identification and confirmation of microorganisms causing urinary tract infections


Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Escherichia coli</i>	25922	Pink- Purple
B	<i>Enterococcus faecalis</i>	29212	Blue
C	<i>Staphylococcus aureus</i>	25923	Cream yellow
D	<i>Pseudomonas aeruginosa</i>	27853	Colourless (Greenish pigment may be observed)
E	<i>Klebsiella pneumoniae</i>	13883	Bluish Purple

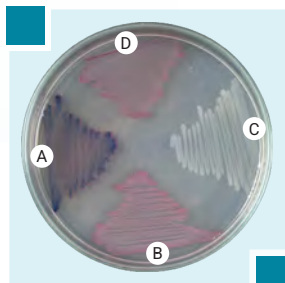

TM 1338 | CHROMOGENIC COLIFORM AGAR W/SLS

For simultaneous detection of total coliforms & *Escherichia coli* in water and foods

Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Escherichia coli</i>	25922	Dark blue to violet (Red color appears around the colony)
B	<i>Klebsiella pneumoniae</i>	13076	Colorless
C	<i>Salmonella enteritidis</i>	13883	Light pink
D	<i>Enterobacter aerogenes</i>	13048	Salmon to Red

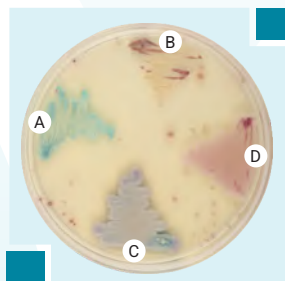

TM 1197 | CHROMOGENIC CANDIDA AGAR

For fast isolation and identification of *Candida* species from mixed flora

Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Candida albicans</i>	10231	Light green
B	<i>Candida glabrata</i>	2001	Cream to white
C	<i>Candida tropicalis</i>	1369	Blue to purple
D	<i>Candida krusei</i>	34135	Purple, fuzzy



TM 1339 | CHROMOGENIC E. COLI AGAR

For easy enumeration of *E. coli* without membrane filtration, or pre-incubation on Mineral Modified Glutamate Medium Base

Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Escherichia coli</i>	25922	Bluish green
B	<i>Salmonella enteritidis</i>	13076	Colourless



TM 1858 | CHROMOGENIC COLIFORM AGAR (CCA)

For determination of coliform bacteria particularly *Enterobacteriaceae* on the basis of their ability to ferment lactose

Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Escherichia coli</i>	25922	Dark blue to violet colonies
B	<i>Citrobacter freundii</i>	8090	Pink to red colonies
C	<i>Salmonella enteritidis</i>	13076	Colourless colonies
D	<i>Klebsiella aerogenes</i>	13048	Pink to red colonies

Culture Media & Supplements


Chromogenic Culture Media



TM 1639 | CHROMOGENIC UTI AGAR, MODIFIED

For enumeration and differentiation of enteric pathogens in urinary tract infections


Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Escherichia coli</i>	25922	Pink/Purple Colonies
B	<i>Enterococcus faecalis</i>	25923	Golden Yellow
C	<i>Staphylococcus aureus</i>	29212	Blue-Green (Small)
D	<i>Klebsiella pneumoniae</i>	13883	Blue to Purple Mucoid
E	<i>Pseudomonas aeruginosa</i>	27853	Colorless (Greenish pigment may be observed)
F	<i>Proteus mirabilis</i>	12453	Light Brown

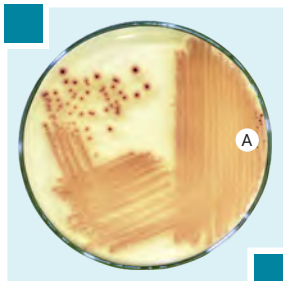

TM 1640 | CHROMOGENIC VIBRIO AGAR

For selective isolation and differentiation of *Vibrio* species

Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Vibrio cholerae</i>	15748	Salmon to Purple

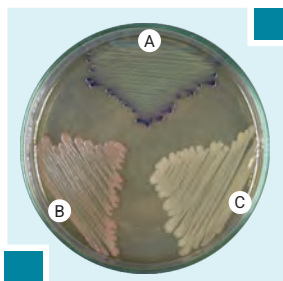

TM 1340 | CHROMOGENIC ECC AGAR

For presumptive identification of *Escherichia coli* and other coliforms

Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Escherichia coli</i>	25922	Blue/Purple
B	<i>Klebsiella pneumoniae</i>	13883	Rose/Pink
C	<i>S. typhimurium</i>	14028	Colourless



TM 1634 | CHROMOGENIC LISTERIA AGAR BASE (MODIFIED)

For selective identification and differentiation of *Listeria monocytogenes*

Symbol	Culture Name	ATCC	Appearance of Colony
A	<i>Listeria monocytogenes</i>	19111	Bluish Green

Culture Media & Supplements

Chromogenic Dehydrated Culture Media

CODE	PRODUCT NAME		F.Wt. Gm/Ltr	PACK SIZE	
TM 1841	CHROMOGENIC A. RAMBACH AGAR * for detection and isolation of <i>Salmonella</i> species in clinical sample	(Part I) (Part II)	30.70 10 ml	500 gm -	
TM 1523	CHROMOGENIC BACILLUS AGAR * for isolation & differentiation between various species of <i>Bacillus</i> using chromogenic substrates		49.20	100 gm 500 gm	
TS 058	POLYMYXIN B SELECTIVE SUPPLEMENT *		#11 vl	5 vl 25 vl	
TM 2110	CHROMOGENIC BACILLUS AGAR BASE for isolation and differentiation between various species of <i>Bacillus</i> by chromogenic method		49.22	100 gm 500 gm	
TS 281	BACILLUS SELECTIVE SUPPLEMENT		#11 vl	5 vl	
TM 1885	CHROMOGENIC BACILLUS CEREUS AGAR * for isolation and differentiation of <i>Bacillus cereus</i> from food		40.98	100 gm 500 gm	
TS 255	BACILLUS CEREUS SELECTIVE SUPPLEMENT *		#25 vl	5 vl	
TM 2111	CHROMOGENIC BIFIDOBACTERIUM AGAR * for the differentiation of <i>Bifidobacterium</i> and <i>Lactobacillus</i> species		59.48	100 gm 500 gm	
TM 2112	CHROMOGENIC CAMPYLOBACTER AGAR BASE * for selective isolation and presumptive identification of <i>Campylobacter</i> species		59.53	100 gm 500 gm	
TS 145	CAMPYLOBACTER SELECTIVE SUPPLEMENT (KARMALI) *		#01 vl	5 vl 25 vl	
TM 1197	CHROMOGENIC CANDIDA AGAR (CHROMOGENIC CANDIDA DIFFERENTIAL AGAR) * for fast isolation and identification of <i>Candida</i> species from mixed flora		42.72	100 gm 500 gm	
TM 1588	CHROMOGENIC CANDIDA DIFFERENTIAL AGAR MODIFIED * for fast isolation and identification of <i>Candida</i> species from clinical and non-clinical specimens		42.05	100 gm 500 gm	
TS 213	CHROMOGENIC CANDIDA SELECTIVE SUPPLEMENT *		#24 vl	5 vl 25 vl	
TM 1833	CHROMOGENIC CLED AGAR BASE * for isolation and differentiation of UTI pathogens		43.00	500 gm	
TM 1858	CHROMOGENIC COLIFORM AGAR (CCA) (ISO 9308-1:2014, ISO 11133:2014) * for determination of coliforms and <i>Escherichia coli</i> in water samples		26.45	100 gm 500 gm	
TM 1338	CHROMOGENIC COLIFORM AGAR (W/ SLS) * for simultaneous detection of total coliforms and <i>Escherichia coli</i> in water and foods		27.00	100 gm 500 gm	
TM 2113	CHROMOGENIC CANDIDA DIFFERENTIAL AGAR BASE * selective and differential medium for rapid isolation and identification of <i>Candida</i> species from mixed cultures		31.20	100 gm 500 gm	
TS 282	CHROMOGENIC CANDIDA DIFFERENTIAL SELECTIVE SUPPLEMENT *		#32 vl	5 vl	
TM 2114	CHROMOGENIC CLOSTRIDIAL AGAR BASE * for selective isolation and presumptive identification of <i>Clostridium</i> species		47.81	100 gm 500 gm	
TM 2115	CHROMOGENIC COLIFORM AGAR ((CCA) W/ 1% AGAR (ISO 9308-1:2014)) * for recommended for detection of total coliforms and <i>Escherichia coli</i> in water samples		25.92	100 gm 500 gm	
TM 2116	CHROMOGENIC COLIFORM AGAR MODIFIED * recommended for the simultaneous detection of <i>Escherichia coli</i> and thermotolerant coliforms in water, milk, dairy products and other food samples		24.00	100 gm 500 gm	
TM 2117	CHROMOGENIC COLISTIN RESISTANT AGAR BASE * recommended for isolation and differentiation of gram negative colistin resistant microorganisms		38.10	500 gm	
TS 283	CHROMOGENIC COLISTIN RESISTANT SELECTIVE SUPPLEMENT *		#13 vl	5 vl	

Culture Media & Supplements

Chromogenic Dehydrated Culture Media

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 2118	CHROMOGENIC CRONOBACTER ISOLATION AGAR (CCI AGAR) (ISO /TS 22964: 2017) * recommended for the isolation and identification of <i>Cronobacter sakazakii</i> from food product	32.40	500 gm	
TM 1854	CHROMOGENIC EC 0157:H7 AGAR (ISO 16654:2001) * for isolation and differentiation of <i>Escherichia coli</i> 0157:H7 from food and environmental samples	28.85	100 gm 500 gm	
TS 248	CHROMOGENIC EC 0157:H7 SELECTIVE SUPPLEMENT *	#18 vl	5 vl 25 vl	
TM 2119	CHROMOGENIC E. COLI AGAR * for the detection and enumeration of <i>Escherichia coli</i> in foods without further confirmation on membrane filter or by indole reagent	36.57	100 gm 500 gm	
TM 2120	CHROMOGENIC EC BROTH (W/ RUG) * recommended for detection of <i>Escherichia coli</i> in water and food samples by a chromogenic and fluorogenic method	10.56	100 gm	
TM 2121	CHROMOGENIC EC 0157:H7 AGAR, MODIFIED * recommended for isolation and differentiation of <i>Escherichia coli</i> 0157:H7 from food, environmental and clinical samples	28.85	100 gm 500 gm	
TS 005	POTASSIUM TELLURITE 1% (1 ml/vl) *	#5 vl	5 vl 25 vl	
TM 2441	CHROMOGENIC ECO157:H7 SELECTIVE AGAR BASE, MODIFIED * for presumptive enumeration of <i>Escherichia coli</i> 0157:H7 by membrane filtration technique from food samples	66.82	100 gm 500 gm	
TS 313	CHROMOGENIC EC 0157:H7 SELECTIVE SUPPLEMENT, MODIFIED *	#8 vl	5 vl	
TM 1340	CHROMOGENIC ECC AGAR * for presumptive identification of <i>Escherichia coli</i> and other coliforms in food and environmental samples	55.83	100 gm 500 gm	
TM 1341	CHROMOGENIC ECC SELECTIVE AGAR BASE * for detection of <i>Escherichia coli</i> and coliforms in water and food samples	26.48	100 gm 500 gm	
TS 187	CHROMOGENIC ECC SELECTIVE SUPPLEMENT *	#38 vl	5 vl	
TM 2122	CHROMOGENIC ECC SELECTIVE AGAR BASE, MODIFIED * for detection of <i>Escherichia coli</i> and coliforms in water and food samples	39.30	500 gm	
TM 1343	CHROMOGENIC ECD W/ MUG * for detection of <i>Escherichia coli</i> by chromogenic and fluorogenic substrates	53.17	100 gm 500 gm	
TM 2123	CHROMOGENIC ENRICHMENT BROTH BASE FOR EC 0157:H7* for isolation and selective differentiation of <i>Escherichia coli</i> 0157:H7 from food, environmental and clinical samples by chromogenic method	22.80	100 gm 500 gm	
TS 248	CHROMOGENIC EC 0157:H7 SELECTIVE SUPPLEMENT*	#44 vl	5 vl 25 vl	
TM 1631	CHROMOGENIC ENTEROBACTER SAKAZAKII AGAR * for isolation and identification of <i>Enterobacter sakazakii</i> from dairy and food products	51.67	100 gm 500 gm	
TM 1853	CHROMOGENIC ENTEROBACTER SAKAZAKII AGAR, MODIFIED (ISO 22964:2017, ISO 22964:2006) * for isolation and identification of <i>Cronobacter sakazakii</i> from milk and milk products	30.75	100 gm 500 gm	
TM 1344	CHROMOGENIC ENTEROCOCCI BROTH * for identification and differentiation of enterococci from water	18.59	100 gm 500 gm	
TM 1632	CHROMOGENIC ENTEROCOCCUS FAECIUM AGAR BASE for identification and differentiation of <i>Enterococcus faecium</i> from faeces, sewage and water supplies	29.00	100 gm 500 gm	
TS 215	CHROMOGENIC ENTEROCOCCUS FAECIUM SELECTIVE SUPPLEMENT*	#19 vl	5 vl	

Culture Media & Supplements

Chromogenic Dehydrated Culture Media

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE
TM 1633	CHROMOGENIC KLEBSIELLA SELECTIVE AGAR BASE * for selective isolation of <i>Klebsiella</i> species from water can be used by membrane filter technique	40.80	100 gm 500 gm
TS 204	CHROMOGENIC KLEBSIELLA SELECTIVE SUPPLEMENT *	#25 vl	5 vl
TM 1339	CHROMOGENIC TRYPTONE BILE GLUCURONIDE AGAR (CHROMOGENIC <i>E. COLI</i> AGAR) (TBX AGAR)(ISO 16649-1:2001 / 16649-2:2001 / 16649-3:2015) * for enumeration of <i>Escherichia coli</i> from food samples, animal feed and water samples	36.57	100 gm 500 gm
TM 2124	CHROMOGENIC L. MONO DIFFERENTIAL AGAR BASE * for the selective and differential isolation, enumeration and identification of <i>Listeria monocytogenes</i> and <i>Listeria</i> species based on PCPLC activity	67.20	500 gm
TS 285	LECITHIN SOLUTION*	#16 vl	5 vl
TS 286	MODIFIED L.MONO SELECTIVE SUPPLEMENT*		5 vl
TM 1845	CHROMOGENIC L-MONO LISTERIA DIFFERENTIAL AGAR * for selective identification and differentiation of <i>Listeria monocytogenes</i>	70.30	100 gm
TS 205	CHROMOGENIC LISTERIA SELECTIVE SUPPLEMENT *	#15 vl	5 vl
TS 031	CHROMOGENIC ENRICHMENT SUPPLEMENT *		5 vl
TM 1634	CHROMOGENIC LISTERIA AGAR BASE (Modified) * for selective identification and differentiation of <i>Listeria</i> species	67.25	100 gm 500 gm
TS 205	CHROMOGENIC LISTERIA SELECTIVE SUPPLEMENT *	#15 vl	5 vl
TM 004	CHROMOGENIC LISTERIA AGAR BASE (ISO 11290-1 & 2:2004, ISO 11133:2014) * for selective identification and differentiation of <i>Listeria</i> species	70.05	100 gm 500 gm
TS 205	CHROMOGENIC LISTERIA SELECTIVE SUPPLEMENT *	#15 vl	5 vl
TS 031	CHROMOGENIC ENRICHMENT SUPPLEMENT *		5 vl
TM 1635	CHROMOGENIC MeReSA AGAR BASE * for isolation and identification of Methicillin resistant <i>Staphylococcus aureus</i> from clinical samples	83.30	100 gm 500 gm
TS 206	CHROMOGENIC MeReSa SELECTIVE SUPPLEMENT *	#12 vl	5 vl
TS 219	CEFOXITIN SUPPLEMENT*	#12 vl	5 vl
TM 2127	CHROMOGENIC M-COLICONFIRM BROTH * recommended for detection of <i>E.coli</i> & other total coliforms in water samples by membrane filtration	17.43	500 gm
TS 287	ECC SELECTIVE SUPPLEMENT MODIFIED *	#29 vl	5 vl
TS 042	TTC SOLUTION 1% (10 ml/vl) *	#21 vl	5 vl 25 vl
TM 1636	CHROMOGENIC MM AGAR * for identification and differentiation of <i>Salmonella</i> and non- <i>Salmonella</i> from water and clinical samples	49.13	100 gm 500 gm
TM 1638	CHROMOGENIC M-LAURYL SULPHATE AGAR * for enumeration and differentiation of <i>E. coli</i> and other coliforms by membrane filter technique	88.00	100 gm 500 gm
TM 1198	CHROMOGENIC OGYE AGAR BASE * for isolation and enumeration of Yeasts and Moulds from milk and milk products using chromogenic substrates	37.10	100 gm 500 gm
TS 071	OXYTETRA SELECTIVE SUPPLEMENT *	#27 vl	5 vl 25 vl
TM 2125	CHROMOGENIC RAJHANS MEDIUM (SALMONELLA AGAR *) for identification and differentiation of <i>Salmonella</i> species from among the members of <i>Enterobacteriaceae</i> , especially <i>Proteus</i> species	46.82	100 gm 500 gm

Culture Media & Supplements

Chromogenic Dehydrated Culture Media

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 1426	CHROMOGENIC SALMONELLA AGAR * for isolation and differentiation of <i>Salmonella</i> species from coliforms	27.90	100 gm 500 gm	
TM 1839	CHROMOGENIC SALMONELLA AGAR * for identification of <i>Salmonella</i> species from other organisms in the family <i>Enterobacteriaceae</i>	55.20	100 gm 500 gm	
TS 252	Salmonella SELECTIVE SUPPLEMENT	#10 vl	5 vl	
TM 1824	CHROMOGENIC SALMONELLA AGAR, MODIFIED * for identification and differentiation of <i>Salmonella</i> species from among the members of <i>Enterobacteriaceae</i> , especially <i>Proteus</i> species	42.34	100 gm 500 gm	
TM 1337	CHROMOGENIC STAPHYLOCOCCUS AUREUS AGAR BASE * for isolation and identification of staphylococci	63.00	100 gm 500 gm	
TS 001	EGG YOLK TELLURITE EMULSION (100 ml/vl) *	#4 vl	5 vl	
TM 1199	CHROMOGENIC UTI AGAR * for identification and confirmation of microorganisms causing urinary tract infections	32.45	100 gm 500 gm	
TM 1909	CHROMOGENIC UTI AGAR* for presumptive identification of microorganisms mainly causing urinary tract infections	56.80	100 gm 500 gm	
TM 2126	CHROMOGENIC UNIVERSAL DIFFERENTIAL MEDIUM * it is a differential medium recommended for presumptive identification of microorganisms from clinical and non-clinical specimens	35.00	100 gm 500 gm	
TM 1639	CHROMOGENIC UTI AGAR, MODIFIED * for enumeration and differentiation of enteric pathogens in urinary tract infections	55.44	100 gm 500 gm	
TS 207	DMACA REAGENT (10 ml/vl)	-	1 vl	
TS 208	TDA REAGENT (10 ml/vl)	-	1 vl	
TM 1825	CHROMOGENIC UTI SELECTIVE AGAR * for identification, differentiation and confirmation of enteric bacteria from specimens such as urine	56.94	100 gm 500 gm	
TM 1640	CHROMOGENIC VIBRIO AGAR * for selective isolation and differentiation of <i>Vibrio</i> species	67.50	100 gm 500 gm	
TM 1907	CHROMOGENIC VRE AGAR BASE * for identification of Vancomycin Resistant Enterococci from clinical specimens	50.95	100 gm 500 gm	
TS 253	CHROMOGENIC VRE AGAR SUPPLEMENT *	#20 vl	5 vl	

Culture Media & Supplements

Dehydrated Culture Media (as per BIS)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 649	ALKALINE PEPTONE WATER (IS:5887 (Part V) 1976, reaffirmed 2005) for enrichment of <i>Vibrio</i> species	15.00	100 gm 500 gm	
TM 936	ANDRADE PEPTONE WATER (IS:5887 (Part I and IV) 1976, reaffirmed 2005) a basal medium to study fermentation reactions by adding carbohydrates	15.10	100 gm 500 gm	
TM 635	BAIRD PARKER AGAR BASE (IS:5887 (Part II) 1976, reaffirmed 2005) for isolation and enumeration of coagulase positive staphylococci from food and other products	65.00	100 gm 500 gm	
TS 001	EGG YOLK TELLURITE EMULSION (100 ml/vl) *	#4 vl	5 vl	
TS 002	EGG YOLK EMULSION (100 ml/vl) *		5 vl	
TS 003	POTASSIUM TELLURITE 3.5% (10 ml/vl) *	#3 vl	5 vl	
TM 674	BILE SALT AGAR (IS:5887 (Part V) 1976, reaffirmed 2005) for isolation and enumeration of bile tolerant enteric bacilli	40.00	100 gm 500 gm	
TM 952	BRILLIANT GREEN AGAR BASE W/ PHOSPHATES (IS:5887 (Part III) 1999, reaffirmed 2005) for selective isolation of <i>Salmonella</i> by inhibiting <i>E. coli</i> , <i>Proteus</i> and <i>Pseudomonas</i> species	52.00	100 gm 500 gm	
TS 013	SULPHA SUPPLEMENT *	#20 vl	5 vl	
TM 678	BRILLIANT GREEN BILE BROTH 2% (BIS 5401-2 : 2012) for detection and confirmation of coliform bacteria in water and foods	40. 00	100 gm 500 gm	
TM 686	BUFFERED PEPTONE WATER (IS:5887 (Part III) 1999, reaffirmed 2005) for pre-enrichment of injured <i>Salmonella</i> species prior to selective enrichment and isolation	20.07	100 gm 500 gm	
TM 703	COOKED MEAT MEDIUM (R.C. MEDIUM) (as per BIS) for cultivation of aerobes and anaerobes especially pathogenic <i>Clostridia</i> and also for maintenance of stock cultures. It is recommended by BIS committee under the specifications IS:5887(Part II)-1976	11.54	100 gm 500 gm	
TM 086	DECARBOXYLASE TEST MEDIUM BASE (FALKOW) (IS:5887 (Part V) 1976, reaffirmed 2005) for cultivation and differentiation of bacteria based on their decarboxylase activity	09.00	100 gm 500 gm	
TM 1870	DEOXYCHOLATE CITRATE AGAR MEDIUM (IS:5887 (Part VII) 1999, reaffirmed 2005) for isolation of <i>Shigella</i> species from food samples	55.45	100 gm 500 gm	
TM 1171	EMB AGAR, LEVINE (IS:5887 (Part I) 1976, reaffirmed 2005) for isolation, enumeration and differentiation of <i>Enterobacteriaceae</i>	37.50	100 gm 500 gm	
TM 1179	ETHYL VIOLET AZIDE BROTH (E.V.A. BROTH) (IS:5887 (Part II) 1976, reaffirmed 2005) for selective and confirmatory detection of enterococci as an indicator of faecal pollution in water	35.80	100 gm 500 gm	
TM 119	GLUCOSE SALT TEEPOL BROTH (DOUBLE PACK) (IS:5887 (Part V) 1976, reaffirmed 2005) for enrichment of <i>Vibrio parahaemolyticus</i> in drinking water	(Part I) 48.00 (Part II) 04.00	100 gm 500 gm	
TM 742	HUGH LEIFSON MEDIUM (IS:5887 (Part V) 1976, reaffirmed 2005) for detecting aerobic and anaerobic breakdown of glucose	20.33	100 gm 500 gm	

Culture Media & Supplements

Dehydrated Culture Media (as per BIS)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 2143	KAUFFMAN MUELLER'S TETRATHIONATE BROTH BASE (IS:5887 (Part I)-1999) recommended as selective enrichment medium for isolation of <i>Shigella</i> species from food samples	132.50	500 gm	
TM 1036	MacCONKEY AGAR (IS:5887 (Part I & II) 1976, reaffirmed 2005) for isolation and differentiation of lactose fermenting and non-fermenting lactose enteric bacteria	55.07 55.07	100 gm 500 gm	
TM 1238	MacCONKEY BROTH W/ NEUTRAL RED (IS:5887 (Part I & II) 1976, reaffirmed 2005) for selective enrichment and enumeration of coliforms	40.07 80.14	100 gm 500 gm	
TM 1239	MacCONKEY BROTH W/ NEUTRAL RED (DOUBLE STRENGTH) (IS:5887 (Part I) 1976, reaffirmed 2005) for primary isolation of coliforms from large samples like water	80.14	100 gm 500 gm	
TM 2230	MACCONKEY BROTH PURPLE (IS:5401-1:2012) for presumptive identification of coliforms from variety of specimens such as water, milk and food etc	40.02	100 gm 500 gm	
TM 787	MILK AGAR W/CETRIMIDE (DOUBLE PACK) (IS:13428:1998, reaffirmed 2005) for detection and enumeration of <i>Pseudomonas aeruginosa</i> in water	(Part I) 100.00 (Part II) 19.80	120 gm 500 gm	
TM 228	MOTILITY TEST MEDIUM (EDWARDS AND EWING) (IS:5887 (Part I & V) 1976, reaffirmed 2005) for testing motility of enteric bacteria	22.00	500 gm	
TM 685	MR-VP MEDIUM (BUFFERED GLUCOSE BROTH) (IS:5887 (Part I, III & IV) 1976, reaffirmed 2005) for differentiation of coli-aerogenes group by MR-VP test	15.00	100 gm 500 gm	
TM 197	MYP AGAR BASE (PHENOL RED EGG YOLK POLYMYXIN AGAR BASE) (IS:5887 (Part VI) 1976, reaffirmed 2005) for isolation and identification of pathogenic staphylococci and <i>Bacillus</i> species	46.00	100 gm 500 gm	
TS 058	POLYMYXIN B SELECTIVE SUPPLEMENT*	#22 vl	5 vl 25 vl	
TS 002	EGG YOLK EMULSION (100 ml/vl)*	#11 vl	5 vl	
TM 797	NITRATE BROTH (IS:5887 (Part IV)-1999) for detection of nitrate reduction and enumeration of <i>Bacillus cereus</i>	39.00	100 gm 500 gm	
TM 798	NUTRIENT AGAR NO. 2 (IS:5887 (Part I, II & V) 1976, reaffirmed 2005) general purpose culture medium	40.00	100 gm 500 gm	
TM 1270	NUTRIENT BROTH (W/ 1% PEPTONE) (IS:5887 (Part I & II) 1976, reaffirmed 2005, IS:5887 (Part IV) 1999) general purpose culture medium	25.00	100 gm 500 gm	
TM 1059	NUTRIENT GELATIN (IS:5887 (Part VII) 1999, reaffirmed 2005) for detection of gelatin liquefaction by proteolytic microorganisms	158.00	500 gm	
TM 805	PEPTONE WATER (IS:5887 (Part I) 1976, reaffirmed 2005) general purpose growth medium & as the base of carbohydrate fermentation media	25.00	100 gm 500 gm	
TM 828	PLATE COUNT AGAR (IS:5402-1969 First Reprint 1983) for determination of microbial counts in milk and other dairy products by pour plate method	30.00	100 gm 500 gm	
TM 389	SELENITE BROTH (SELENITE F BROTH) (DOUBLE PACK) (IS:5887 (Part I, III & IV) 1976, (reaffirmed 2005) for isolation and enrichment of <i>Salmonella</i> from faeces, urine or other pathological material	(Part I) 19.00 (Part II) 04.00	100 gm 500 gm	
TM 1291	SIMMONS CITRATE AGAR (IS:5887 (Part I) 1976, reaffirmed 2005) for differentiation between faecal coliforms and members of the aerogenes group on the basis of citrate utilization	24.28	100 gm 500 gm	

Approximate number of vials required for 500 gm of dehydrated medium. | * Store Between 2-8°C

Culture Media & Supplements

Dehydrated Culture Media (as per BIS)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE
TM 1299	TCBS AGAR (IS : 5887 (Part V) 1976, reaffirmed 2005) for selective isolation of <i>Vibrio cholerae</i> and other enteropathogenic <i>Vibrio</i>	89.00	100 gm 500 gm
TM 440	TERGITOL - 7 AGAR BASE (IS : 5887 (Part I) 1976, reaffirmed 2005) for selective isolation and identification of coliform bacteria from water	33.13	100 gm 500 gm
TS 042	T.T.C. SOLUTION 1% (10 ml/vl) *	5 vl	5 vl 25 vl
TM 1874	TRIPLE SUGAR IRON AGAR (IS : 5887 (Part I, III and V) 1976, reaffirmed 2005) for confirmation of gram negative enteric bacilli on basis of dextrose, lactose and sucrose fermentation and H ₂ S production	64.32	100 gm 500 gm
TM 891	TRYPTONE WATER W/O NaCl (IS : 5887 (Part V) 1976, reaffirmed 2005) for detection of <i>Vibrio cholerae</i> and <i>Vibrio parahaemolyticus</i>	10.00	100 gm 500 gm
TM 2410	UREA AGAR BASE (CHRISTENSEN) (IS : 5887 (Part-I)-1976) for detection of urease production, particularly by <i>Proteus vulgaris</i> , micrococci & <i>Paracolon</i> organisms	24.01	100 gm 500 gm
TS 030	UREA 40% (5 ml/vl) *	#268 vl	5 VL 25 vl
TM 1314	VIOLET-RED BILE AGAR (IS : 5401 (Part 1) : 2002) for selective isolation, detection and enumeration of coli-aerogenes in water, milk products	41.53	100 gm 500 gm
TM 1319	WILLIS AND HOBBS' MEDIUM BASE (IS : 5887 (Part 4) : 1999) for isolation and identification of <i>Clostridium</i> from foods	47.03	500 gm
TS 002	EGG YOLK EMULSION (100 ml/vl)*	#4vl	5 vl
TS 098	WILLIS AND HOBBS' SUPPLEMENT*	#22 vl	5 vl
TM 498	YEAST GLUCOSE CHLORAMPHENICOL AGAR (IS:5403:1999 reaffirmed 2005) for selective enumeration of Yeasts and Moulds in milk and milk products	40.00	100 gm 500 gm

Dehydrated Culture Media (as per ISO)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE
TM 1943	ACETAMIDE BROTH (DOUBLE PACK) (ISO 16266-2:2018) for confirmation of non-fermentative Gram-negative bacteria, particularly <i>Pseudomonas aeruginosa</i>	(Part I) 21.40 (Part II) 01.40	100 gm 500 gm
TM 008	ALKALINE PEPTONE WATER (pH 8.6) for detection and enrichment of <i>Vibrio</i> species	50.00	500 gm
TM 1878	ALKALINE SALINE PEPTONE WATER (ASPW) (ISO 8261, ISO 7218, ISO 6887, ISO 21872-1&2:2007) for enrichment of <i>Vibrio</i> species from food and water samples	40.00	100 gm 500 gm
TM 1505	L-ARGININE DIHYDROLASE MEDIUM, MODIFIED (ISO 22964:2006) for confirmation of <i>Enterococcus sakazakii</i> from milk and milk products	09.01	500 gm
TM 1881	BACILLUS CEREUS SELECTIVE AGAR BASE (MYP) (ISO 7932:2004) selective isolation and enumeration of <i>Bacillus cereus</i>	46.03	500 gm
TS 058	POLYMYXIN B SELECTIVE SUPPLEMENT *	#22 vl	5 vl 25 vl
TS 002	EGG YOLK EMULSION (100 ml/vl) *	#11 vl	5 vl
TM 1579	BAIRD PARKER AGAR BASE (RPF) (ISO 6888-1 & 2:1999) for isolation & enumeration of coagulase positive staphylococci from food & pharma products	56.00	500 gm
TS 001	EGG YOLK TELLURITE EMULSION (100 ml/vl) *	#10 vl	5 vl
TS 176	RPF SUPPLEMENT (FIBRINOGEN PLASMA TRYPSIN INHIBITOR SUPPLEMENT) *	#80 vl	5 vl 25 vl

Culture Media & Supplements

Dehydrated Culture Media (as per ISO)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 036	BILE ESCULIN AGAR (ISO 10273:1994) for isolation & identification of <i>Yersinia enterocolitica</i> from food and animal feeding products	64.50	100 gm 500 gm	
TM 1883	BILE ESCULIN AZIDE AGAR (ISO 7899-2:2000) for isolation and presumptive identification of faecal streptococci	56.65	500 gm	
TM 1842	BLOOD AGAR BASE NO. 2 (ISO 11290-1:2017) for isolation, cultivation and detection of haemolytic activity of streptococci, pneumococci and other fastidious microorganisms	42.50	100 gm 500 gm	
TS 006	BRUCELLA SELECTIVE SUPPLEMENT *	#26 vl	5 vl 25 vl	
TS 007	CAMPYLOBACTER SUPPLEMENT-I (BLASER-WANG)*	#26 vl	5 vl	
TS 008	CAMPYLOBACTER SUPPLEMENT-II (BUTZLER) *	#26 vl	5 vl 25vl	
TS 009	CAMPYLOBACTER SUPPLEMENT-III (SKIRROW) *	#26 vl	5 vl	
TS 010	CAMPYLOBACTER GROWTH SUPPLEMENT *	#26 vl	5 vl	
TS 011	STREPTO SUPPLEMENT *	#26 vl	5 vl 25 vl	
TM 1146	BLOOD FREE CAMPYLOBACTER SELECTIVITY AGAR BASE (ISO 10272-1 & 2:2017) for selective isolation and differentiation of <i>Campylobacter</i> species	45.50	500 gm	
TS 102	CAMPYLOBACTER SUPPLEMENT V (BFCSA) *	#22 vl	5 vl	
TS 103	CAT SELECTIVE SUPPLEMENT *	#22 vl	5 vl	
TM 1510	BOLTON BROTH BASE (ISO 10272-1:2017) for selective enrichment of <i>Campylobacter</i> species from foods	27.60	500 gm	
TS 179	BOLTON SELECTIVE SUPPLEMENT*	#37 vl	5 vl	
TM 1512	BRIGHT GREEN AGAR W/ PHOSPHATES (ISO 6785:2001, ISO 6579:1981, ISO 3565:1975) for selective isolation of <i>Salmonella</i>	54.69	500 gm	
TS 013	SULPHA SUPPLEMENT *	#20 vl	5 vl	
TM 365	BRIGHT GREEN BILE BROTH 2% (BRIGHT GREEN LACTOSE BILE BROTH 2%) (ISO 4831:2006, ISO 4832:2006) for detection and confirmation of coliform bacteria in water and foods	40.00	100 gm 500 gm	
TM 1894	BUFFERED CHARCOAL YEAST EXTRACT AGAR MEDIUM (ISO 11731-2:2017) for selective isolation and cultivation of <i>Legionella</i> species from cooling towers, water samples, clinical and other materials	35.00	100 gm 500 gm	
TS 018	LEGIONELLA SELECTIVE SUPPLEMENT IV (MWY)*	#5 vl	5 vl	
TS 019	LEGIONELLA SUPPLEMENT *	#5 vl	5 vl	
TS 115	LEGIONELLA (GVPC)SELECTIVE SUPPLEMENT*	#5 vl	5 vl 25 vl	
TS 254	PCP SUPPLEMENT*	#5 vl	5 vl	
TM 307	BUFFERED PEPTONE WATER (ISO 6579-1:2017, 11133:2014, 11290-2:2017, 21528:2017, 6887-1/2/3/4:2017) for pre enrichment of injured <i>Salmonella</i> species prior to selective enrichment and isolation	20.00	100 gm 500 gm	
TM 1522	CETRIMIDE BROTH BASE (ISO 8360-2:1988) for cultivation of <i>Pseudomonas aeruginosa</i> from water samples using membrane filter technique	31.90	500 gm	
TM 1858	CHROMOGENIC COLIFORM AGAR (CCA) (ISO 9308-1:2014, ISO 11133:2014) for determination of coliforms and <i>Escherichia coli</i> in water samples	26.45	100 gm 500 gm	
TM 2115	CHROMOGENIC COLIFORM AGAR (CCA) W/1% AGAR(ISO 9308-1:2014) for recommended for detection of total coliforms and <i>Escherichia coli</i> in water samples	25.92	100 gm 500 gm	

Approximate number of vials required for 500 gm of dehydrated medium. | * Store Between 2-8°C



Culture Media & Supplements

Dehydrated Culture Media (as per ISO)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 2118	CHROMOGENIC CRONOBACTER ISOLATION AGAR (CCI AGAR)(ISO /TS 22964: 2017) recommended for the isolation and identification of <i>Cronobacter sakazakii</i> from food products	32.40	500 gm	
TM 1854	CHROMOGENIC EC 0157:H7 AGAR (ISO 16654:2001) for isolation and differentiation of <i>Escherichia coli</i> 0157:H7 from food and environmental samples	28.85	100 gm	
TS 248	CHROMOGENIC EC 0157:H7 SELECTIVE SUPPLEMENT *	#18 vl	5 vl 25 vl	
TM 1853	CHROMOGENIC ENTEROBACTER SAKAZAKII AGAR, MODIFIED (ISO 22964:2017, ISO 22964:2006) for isolation and identification of <i>Cronobacter sakazakii</i> from milk and milk products	30.75	100 gm 500 gm	
TM 2045	COLUMBIA BLOOD AGAR BASE (ISO 10272-2:2017) recommended for selective detection and enumeration of <i>Campylobacter</i> species from food chain	44.00	500 gm	
TM 004	CHROMOGENIC LISTERIA AGAR BASE (ISO 11290-1 & 2:2004, ISO 11133:2014) * for selective identification and differentiation of <i>Listeria</i> species	70.05	100 gm 500 gm	
TS 205	CHROMOGENIC LISTERIA SELECTIVE SUPPLEMENT *	#15 vl	5 vl	
TS 031	CHROMOGENIC ENRICHMENT SUPPLEMENT *	#15 vl	5 vl	
TM 1339	CHROMOGENIC TRYPTONE BILE GLUCURONIDE AGAR (CHROMOGENIC E.COLI AGAR) (TBX AGAR)(ISO 16649-1 & 2 :2001, ISO 16649-3:2015, ISO 11133:2014) * for enumeration of <i>Escherichia coli</i> from food samples, animal feed and water samples	36.57	100 gm 500 gm	
TM 1852	DICHLORAN ROSE BENGAL CHLORAMPHENICOL AGAR (DRBC AGAR) (ISO 21527-1:2008, 11133:2014) for selective isolation of Yeasts and Moulds of significance in food spoilage	31.60	100 gm 500 gm	
TM 417	EC BROTH (ISO 7251:2005, 11866-1:1997) for selective enumeration of faecal and non faecal coliforms in water	37.00	100 gm 500 gm	
TM 1804	EE BROTH, MOSSEL (ISO 21528-1:2004, ISO 7402:1993) for selective enrichment of <i>Enterobacteriaceae</i> in bacteriological examination of foods	43.47	500 gm	
TM 371	EMB AGAR, LEVINE (ISO 21150:2015) for isolation, enumeration and differentiation of members of <i>Enterobacteriaceae</i> from pharma, dairy & food products	37.50	100 gm 500 gm	
TM 2073	EUGONIC LT 100 BROTH BASE W/O TWEEN 80 (ISO 21149:2017) for the enrichment and detection of mesophilic aerobic bacteria present in cosmetic products	32.40	500 gm	
TM 1472	FRASER BROTH BASE (ISO 11290-1:2017) for isolation and enumeration of <i>Listeria monocytogenes</i> from food and animal feeds	57.35	100 gm 500 gm	
TS 034	FRASER SUPPLEMENT *	#19 vl	5 vl	
TS 035	FRASER SELECTIVE SUPPLEMENT *	#19 vl	5 vl	
TM 118A	GIOLITTI-CANTONI BROTH BASE (ISO 6888-3:2003) for selective enrichment of <i>Staphylococcus aureus</i> from foods	55.20	500 gm	
TS 003	POTASSIUM TELLURITE 3.5% (10 ml/vl) *	#18 vl	5 vl	
TS 005	POTASSIUM TELLURITE 1% (1 ml/vl) *	#18 vl	5 vl 25 vl	

Culture Media & Supplements

Dehydrated Culture Media (as per ISO)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 1937	HALF FRASER BROTH BASE (ISO 11290-1:2017) for the selective enrichment of <i>Listeria</i> species from foods	55.00	500 gm	
TS 034	FRASER SUPPLEMENT *	#19 vl	5 vl	
TS 035	FRASER SELECTIVE SUPPLEMENT *	#10 vl	5 vl	
TM 121	HEKTOEN ENTERIC AGAR (ISO 21567:2004) for differential & selective isolation of <i>Salmonella</i> and <i>Shigella</i> from enteric pathological specimens	76.67	100 gm 500 gm	
TM 2138	IRON SULPHITE AGAR MODIFIED (ISO 15213:2003) for the enumeration of sulphite-reducing bacteria growing under anaerobic conditions	42.00	500 gm	
TM 1553	KLIGLER IRON AGAR (ISO 13720-1995, ISO 10273:2003) for identification of <i>Pseudomonas</i> species and for differential identification of gram-negative enteric bacilli based on dextrose and lactose fermentation and H ₂ S production	58.02	500 gm	
TM 2173	L-LYSINE DECARBOXYLASE SALINE BROTH (ISO 21872-1-2017) recommended for biochemical confirmation of <i>Vibrio parahaemolyticus</i> and <i>Vibrio cholerae</i> from food and animal feeding stuffs on the basis of lysine decarboxylation	19.01	500 gm	
TM 2265	L-ORNITHINE DECARBOXYLASE BROTH (ISO / TS 22964: 2017) for detection of the ability of microorganisms to decarboxylate ornithine	14.01	100 gm	
TM 1251	LACTOSE TTC AGAR (WITH SODIUM HEPTADECYL SULPHATE) (ISO 9308-1:2000) for detection and enumeration of <i>E. coli</i> and other coliforms in water by membrane filtration technique	57.15	100 gm 500 gm	
TS 042	TTC SOLUTION 1% (10 ml/vl) *	#3 vl	5 vl 25 vl	
TM 1556	LACTOBACILLUS MRS AGAR * (ISO 1995, DRAFT ISO/DIS 13720:2010) for the isolation and enumeration of Lactic Acid Bacteria from meat and meat products	65.30	500 gm	
TM 150	LAURYL SULPHATE BROTH (LAURYL TRYPTOSE BROTH) (ISO 11866-2, ISO 4831:2006) for detection and enumeration of coliform bacteria in water, waste water, dairy products and other food sample	35.60	100 gm 500 gm	
TM 1829	LAURYL SULPHATE TRYPTOSE BROTH, MODIFIED (ISO 22964:2006) for the selective pre-enrichment medium for <i>Cronobacter sakazakii</i> from milk and milk products	64.60	500 gm	
TM 1226	LISTERIA IDENTIFICATION AGAR BASE (PALCAM) (ISO 11290-2) for detection and enumeration of <i>Listeria monocytogens</i> from food and animal feeds	69.00	100 gm 500 gm	
TS 119	LISTERIA SELECTIVE SUPPLEMENT (PALCAM) *	#15 vl	5 vl 25 vl	
TM 1229	LISTERIA OXFORD MEDIUM BASE (ISO 11290-1) for isolation of <i>Listeria</i> species from pathological samples	55.50	100 gm 500 gm	
TS 120	OXFORD LISTERIA SUPPLEMENT *	#19 vl	5 vl	
TS 121	LISTERIA MOXALACTAM SUPPLEMENT *		5 vl	
TM 1830	LYSINE DECARBOXYLASE BROTH W/O PEPTONE (ISO 6579:2002, ISO 22964:2006) for differentiating <i>Salmonella</i> serotype arizonae from the <i>Bethesda Ballerup</i> group of <i>Enterobacteriaceae</i>	09.00	100 gm 500 gm	
TM 775	M-LAURYL SULPHATE BROTH (ISO 9308-1:2000) for enumeration of <i>Escherichia coli</i> in water by membrane filter technique	76.20	500 gm	

Culture Media & Supplements

Dehydrated Culture Media (as per ISO)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 1572	MUG LAURYL SULPHATE BROTH, MODIFIED (ISO 11866-2:1997) for selective enrichment of presumptive <i>Escherichia coli</i> and other coliforms from milk & milk products; also, recommended for enumeration purpose using of Most Probable Number (MPN) method	36.70	500 gm	
TM 1875	MacCONKEY BROTH PURPLE W/ BCP (ISO 9308-2:1990) for presumptive identification of coliforms from water	40.01	100 gm 500 gm	
TM 529	MacCONKEY BROTH W/ BCP & NaCl (ISO 9308-2:1990) for presumptive identification of coliforms from water, milk and foods etc	40.00	100 gm 500 gm	
TM 1915	MacCONKEY SORBITOL AGAR BASE (ISO 16654:2001) a selective medium for isolation and detection of <i>Escherichia coli</i> 0157:H7	50.03	500 gm	
TS 087	TELLURITE-CEFEXIME SUPPLEMENT *	#20 vl	5 vl	
TM 1843	MANNITOL EGG YOLK POLYMYXIN AGAR (ISO 7932:2004, ISO 21871:2006) for the enumeration of <i>B. cereus</i> in foodstuffs and other samples	46.03	500 gm	
TS 058	POLYMYXIN B SELECTIVE SUPPLEMENT *	#22 vl	5 vl 25 vl	
TS 002	EGG YOLK EMULSION (100 ml/vl) *	#11 vl	5 vl	
TM 1857	MINERAL MODIFIED GLUTAMATE BROTH BASE (ISO 16649-3:2005) for enumeration of coliform bacteria in water and wastewater samples	17.70	500 gm	
TM 2215	MRS AGAR, MODIFIED * (ISO 15214:1998) recommended for the selective cultivation of Lactic acid bacteria from food	69.21	500 gm	
TM 2217	MRS SELECTIVE AGAR BASE (ISO 20128:2006) recommended for the selective cultivation of Lactic acid bacteria from food	65.15	100 gm 500 gm	
TS 298	CIPROFLOXACIN CLINDAMYCIN SELECTIVE SUPPLEMENT*	#8 vl	5 vl	
TM 2246	MOTILITY NITRATE MEDIUM, BUFFERED (ISO 7937:2004) for isolation and detection of <i>Clostridium perfringens</i> on the basis of motility and nitrate test	19.50	500 gm	
TM 1399	MUELLER KAUFFMAN TETRATHIONATE NOVOBIOCIN BROTH BASE (ISO 6579-1:2017) for enrichment and isolation of <i>Salmonella</i> by suppressing <i>Proteus</i> species	89.50	500 gm	
TS 152	MKTT NOVOBIOCIN SUPPLEMENT *	#6 vl	5 vl	
TM 1038A	NUTRIENT AGAR (pH 6.8) (ISO 10273:2003, ISO 21528:2004, ISO 19250:2010, ISO 6579:2015) a general purpose nutrient medium for testing of water, sewage, faeces & other materials	23.00	100 gm 500 gm	
TM 2256	NUTRIENT AGAR 1.5% (ISO 1995, ISO/DIS 13720:2010) recommended general purpose nutrient medium which can be used for cultivation of fastidious microorganisms after appropriate enrichment	28.00	500 gm	
TM 2274	PARK & SANDER ENRICHMENT BROTH BASE (ISO/DIS 10272:1995) for selective enumeration of thermo-tolerant <i>Campylobacter</i> species from food	29.35	500 gm	
TS 073	PARK & SANDER SELECTIVE SUPPLEMENT A *	#15 vl	5 vl	
TS 074	PARK & SANDER SELECTIVE SUPPLEMENT B *	#15 vl	5 vl	
TM 1844	PEPTONE SORBITOL BROTH (W/Bile Salt) (ISO 10273:2017) for selective enrichment of <i>Yersinia enterocolitica</i> from dairy products	(Part I) (Part II) 21.00 10.00	500 gm -	
TM 2281	PHOSPHATE BUFFERED SALINE (FOR LISTERIA) (ISO 11290-2:2017) recommended for the preparation of dilutions for <i>Listeria</i> species for further testing from food sample	18.37	500 gm	
TM 2270	PSB BROTH MODIFIED (ISO 1994, ISO/DIS 10273) for primary enrichment and enumeration of <i>Yersinia enterocolitica</i> from food	30.77	500 gm	

Culture Media & Supplements

Dehydrated Culture Media (as per ISO)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TM 2300	RPF AGAR BASE (ISO 6888-2:1999) recommended for the enumeration of coagulase positive staphylococci from food and animal feeding stuff	58.00	500 gm	
TS 176	RPF SUPPLEMENT (FIBRINOGEN PLASMA TRYPSIN INHIBITOR SUPPLEMENT) *	#87 vl	5 vl 25 vl	
TM 2324	SALINE NUTRIENT AGAR FOR VIBRIO (ISO 21872-1-2017) recommended for enrichment of <i>Vibrio parahaemolyticus</i> and <i>Vibrio cholerae</i> species	43.00	500 gm	
TM 1847	SALINE MEAT YEAST AGAR (as per ISO 16649) for identification of <i>Vibrio parahaemolyticus</i> from food products and animal feeding products	58.30	500 gm	
TM 2325	SALINE PEPTONE WATER W/ 6% NaCl (ISO 21872-1-2017) recommended as a biochemical test to distinguish between <i>Vibrio</i> species based on salt tolerance	43.00	500 gm	
TM 2326	SALINE PEPTONE WATER W/ 10% NaCl (ISO 21872-1-2017) recommended for detection of <i>Vibrio parahaemolyticus</i> and <i>Vibrio cholerae</i> from food	110.00	500 gm	
TM 2327	SALINE TRYPTONE / TRYPTOPHAN MEDIUM (ISO 8914:1990) recommended for detection of indole production by <i>Vibrio parahaemolyticus</i>	41.00	500 gm	
TM 1848	SALT POLYMYXIN BROTH BASE (ISO 8914:1990) for detection and enumeration of <i>Vibrio</i> species	29.70	500 gm	
TS 058	POLYMYXIN B SELECTIVE SUPPLEMENT *	#34 vl	5 vl 25 vl	
TM 294	SELENITE CYSTINE BROTH (FLUID SELENITE CYSTINE BROTH) (DOUBLE PACK) (ISO 6579:1993) an enrichment medium for isolation of <i>Salmonella</i> species from food, dairy and clinical samples	(Part I) 19.01 (Part II) 04.00	100 gm 500 gm	
TM 2336	SLANETZ AND BARTLEY MEDIUM (ISO/DIS 7899-2: 2000) for detection and enumeration of faecal streptococci from water samples by membrane filtration technique	46.50	500 gm	
TM 1611	SOYABEAN BILE BROTH W/ NOVOBIOCIN (ISO 16654:2001) * for the detection and enrichment of <i>Escherichia coli</i> O157:H7 from foods	33.02	500 gm	
TM 544	STANDARD METHODS AGAR (PLATE COUNT AGAR) (ISO 4833-1 & 2:2013) for determination of plate counts of microorganisms in foods, water, waste water and also from clinical samples	23.50	100 gm 500 gm	
TM 1097	SYNTHETIC SEA SALT (ISO 9308-3:1998) for preparation of special diluents	19.19	500 gm	
TM 2374	TOLUIDINE BLUE DNA AGAR (ISO 8870:2006(E) 83:2006(E)) for detection of thermostable deoxyribonuclease activity to establish speciation of <i>S. aureus</i> in contaminated foods	26.56	100 gm	
TM 1819	TRIPLE SUGAR IRON AGAR (ISO 6785:2001 / ISO 6579:2017) for identification of gram negative enteric bacteria on the basis of sugar fermentation & H ₂ S production	65.00	100 gm 500 gm	

Approximate number of vials required for 500 gm of dehydrated medium. | * Store Between 2-8°C

Culture Media & Supplements

Dehydrated Culture Media (as per ISO)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE
TM 2385	TRYPTONE SALT BROTH (ISO:1999 ISO/DIS 6887-1:2017) for preparation of specimens, stock suspensions and decimal dilutions for microbiological tests	09.50	100 gm 500 gm
TM 549	TRYPTONE SOYA YEAST EXTRACT BROTH (ISO 11290-2:1998) for isolation and cultivation of <i>Listeria</i> from Henry's light	36.00	500 gm
TM 1849	TRYPTONE WATER (as per ISO 11866) for detection of indole production by coliforms	25.00	100 gm 500 gm
TM 2400	TRYPTONE YEAST SODIUM SULPHITE AGAR BASE (ISO 14189:2013) recommended for the enumeration of <i>Clostridium perfringens</i> from water	42.00	500 gm
TS 076	PERFRINGEN'S T.S.C. SUPPLEMENT (TSC SUPPLEMENT) *	#1 vI	5 vI 25 vI
TS 250	CLOSTRIDIUM PERFRINGENS SUPPLEMENT *	#22 vI	5 vI
TM 1826	TSC AGAR BASE (PERFRINGENS AGAR BASE) (ISO 7937: 2004, ISO 14189:2013) for enumeration of <i>Clostridium perfringens</i> from food	42.00	500 gm
TS 076	PERFRINGENS T.S.C. SUPPLEMENT *	#24 vI	5 vI 25 vI
TS 250	CLOSTRIDIUM PERFRINGENS SUPPLEMENT *	#24 vI	5 vI
TM 2406	TWEEN ESTERASE TEST AGAR BASE (ISO 10273:2017) recommended for confirmation of <i>Yersinia enterocolitica</i>	30.01	100 gm 500 gm
TM 2411	UREA INDOLE BROTH, MODIFIED (ISO 10273:2003) for confirmation of <i>Yersinia enterocolitica</i> by urease and indole test	30.03	100 gm 500 gm
TM 483	VIOLET RED BILE GLUCOSE AGAR W/O LACTOSE (ISO 21528-1 & 2:2017) for detection and enumeration of <i>Enterobacteriaceae</i> in raw foods	38.53	500 gm
TM 1860	VIOLET RED BILE LACTOSE AGAR (ISO 4832:2006) for detection and enumeration of coliform bacteria in food, water and dairy products	38.53	100 gm 500 gm
TM 2412	VP MEDIUM FOR LISTERIA (ISO 11290-2-2017) recommended for the distinction of <i>Listeria</i> species from other species based on Voges-Proskauer test from food samples	17.00	500 gm
TM 1621	XLD AGAR MODIFIED (ISO 6579-2002) for isolation and enumeration of <i>Salmonella typhi</i> and other <i>Salmonella</i> species	55.43	100 gm 500 gm
TM 1827	YEAST EXTRACT AGAR (ISO 6222:1999) for enumeration of microorganisms from water	24.00	100 gm 500 gm
TMH 119	BAIRD PARKER AGAR BASE (as per EP/IP/BP) for the isolation and enumeration of coagulase positive staphylococci from food, pharmaceutical and other materials	63.00	100 gm 500 gm
TS 001	EGG YOLK TELLURITE EMULSION (100 ml/vI) *	#4 vI	5 vI
TS 002	EGG YOLK EMULSION (100 ml/vI)	#4 vI	5 vI
TS 003	POTASSIUM TELLURITE 3.5% (10 ml/vI)	#4 vI	5 vI
TS 005	POTASSIUM TELLURITE 1% (1ml /vI) *	#3 vI	5 vI 25 vI

Culture Media & Supplements

Harmonized Culture Media - Pharma Microbiological Testing Media (USP/JP/EP/BP/IP)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TMH 101	BUFFERED NaCl-PEPTONE SOLUTION (as per USP/EP/JP/BP) dilution fluid for samples in case of microbiological contamination	16.10	100 gm 500 gm	
TMHV 101	BUFFERED NaCl-PEPTONE SOLUTION (as per USP/EP/JP/BP) (VEG.) dilution fluid for samples in case of microbiological contamination	16.10	100 gm 500 gm	
TMH 113	CETRIMIDE AGAR (as per USP/EP/JP/BP) for selective isolation of <i>Pseudomonas aeruginosa</i>	45.30	100 gm 500 gm	
TMH 116	COLUMBIA AGAR (as per USP/EP/JP/BP) for detection of <i>Clostridium sporogenes</i> from pharmaceutical products	44.00	100 gm 500 gm	
TS 217	GENTAMYCIN SUPPLEMENT *	#12 vI	5 vI	
TMH 107	EE BROTH MOSSEL (ENTEROBACTERIA ENRICHMENT BROTH, MOSSEL) (as per USP/BP/JP/EP) for selective enrichment of <i>Enterobacteriaceae</i>	45.01	100 gm 500 gm	
TMH 110	MacCONKEY AGAR (as per USP/EP/BP/JP) for isolation, enumeration and enrichment of <i>Enterobacteriaceae</i>	50.03	100 gm 500 gm	
TMH 118	MacCONKEY AGAR (as per USP/EP/BP/JP) for selective isolation and differentiation of <i>E. coli</i> and other enteric bacteria from pharmaceutical products in accordance with the microbial limit testing by harmonized method	50.03	100 gm 500 gm	
TMH 378	MacCONKEY AGAR W/O CV, NaCl, w/ 0.5% SODIUM TAUROCHOLATE for the selection and recovery of the <i>Enterobacteriaceae</i> and related enteric gram-negative bacilli from clinical, food and water samples	55.00	100 gm 500 gm	
TMH 109	MacCONKEY BROTH (as per USP/EP/JP/BP) for detection and enumeration of coliform bacteria	35.01	100 gm 500 gm	
TMH 114	MANNITOL SALT AGAR BASE (as per USP/BP/EP/JP) for selective isolation and enumeration of <i>Staphylococcus</i> species	111.02	100 gm 500 gm	
TS 002	EGG YOLK EMULSION (100 ml/vI) *		5 vI	
TMHV 114	MANNITOL SALT AGAR BASE (as per USP/BP/EP/JP) (VEG.) for selective isolation and enumeration of <i>Staphylococcus</i> species	111.02	100 gm 500 gm	
TS 002	EGG YOLK EMULSION (100 ml/vI) *		5 vI	
TMH 105	POTATO DEXTROSE AGAR (as per USP/EP/JP/BP) for the cultivation of Yeasts and Moulds	39.00	100 gm 500 gm	
TMH 111	RAPPAPORT VASSILIADIS SALMONELLA ENRICHMENT BROTH (as per USP/EP/JP/BP) for selective enrichment & isolation of <i>Salmonella</i> species from pharmaceutical products	42.50	100 gm 500 gm	
TMH 115	REINFORCED CLOSTRIDIAL BROTH (as per USP/BP/JP/EP) for isolation, cultivation and enumeration of <i>Clostridia</i> species, highly nutritive for <i>Clostridium sporogenes</i> and other anaerobes	38.00	100 gm 500 gm	
TMH 104	SABOURAUD DEXTROSE AGAR (as per USP/EP/JP/BP) for selection, isolation and cultivation of Yeasts and Fungi	65.00	100 gm 500 gm	
TMH 106	SABOURAUD DEXTROSE BROTH (as per USP/EP/JP/BP) for selection, isolation and cultivation of Yeasts, Moulds and aciduric microorganisms	30.00	100 gm 500 gm	

Approximate number of vials required for 500 gm of dehydrated medium. | * Store Between 2-8°C



Culture Media & Supplements

Dehydrated Culture Media - Pharma Microbiological Testing Media (USP/JP/EP/BP/IP)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TMH 103	SOYA CASEIN DIGEST AGAR (TRYPTONE SOYA AGAR) (CASO AGAR) (as per USP/EP/JP/BP/IP) for the cultivation of various microorganisms from pharmaceuticals products in accordance with harmonized method	40.00	100 gm 500 gm	
TMH 102	SOYBEAN CASEIN DIGEST MEDIUM (TRYPTONE SOYA BROTH) (CASO BROTH) (as per USP/EP/JP/BP) for cultivation of wide variety of microorganisms recommended for sterility testing of Moulds and lower bacteria	30.00	100 gm 500 gm	
TMHV 102	SOYBEAN CASEIN DIGEST MEDIUM (TRYPTONE SOYA BROTH) (CASO BROTH) (as per USP/EP/JP/BP) (VEG.) for cultivation of wide variety of microorganisms recommended for sterility testing of Moulds and lower bacteria	30.00	100 gm 500 gm	
TMH 108	VIOLET RED BILE DEXTROSE AGAR (as per USP/BP/EP/JP) for detection and enumeration of coliform bacteria	41.53	100 gm 500 gm	
TMH 117	VIOLET RED BILE GLUCOSE AGAR (USP/EP/JP/BP/IP) for detection and enumeration of <i>Enterobacteriaceae</i> especially subculturing of bile tolerant gram negative bacteria from pharmaceutical products in accordance with microbial limit test	41.53	100 gm 500 gm	
TMH 112	XLD AGAR (XYLOSE LYSINE DEOXYCHOLATE AGAR) (as per USP/BP/JP/EP) for selective differentiation and enrichment medium of <i>Salmonella</i> and <i>Shigella</i> species	55.18	100 gm 500 gm	
TSM 301	ALTERNATIVE THIOGLYCOLLATE MEDIUM (NIH THIOGLYCOLLATE MEDIUM) (as per USP) for sterility testing of biological products	29.00	500 gm	
TSMV 301	ALTERNATIVE THIOGLYCOLLATE MEDIUM (NIH THIOGLYCOLLATE MEDIUM) (as per USP) (VEG.) for sterility testing of biological products	29.00	500 gm	
TSMV 301D	ALTERNATIVE THIOGLYCOLLATE MEDIUM (as per USP) (VEG.), (NIH THIOGLYCOLLATE MEDIUM) (Double wrapped) for evaluation of sterility in manufacturing process	29.00	500 gm	
TSM 2129	FILL TEST MEDIUM recommended for the evaluation of sterility in manufacturing process for easy detection of contamination by Media Fill Test	30.10	500 gm	
TSMV 2129	FILL TEST MEDIUM (VEG.) recommended for the evaluation of sterility in manufacturing process for easy detection of contamination by Media Fill Test	30.10	500 gm	

Culture Media & Supplements

Sterile Dehydrated Culture Media (γ -Irradiated)

CODE	PRODUCT NAME	F.Wt. Gm/Ltr	PACK SIZE	
TSMV 2129D	FILL TEST MEDIUM (VEG.) (Double wrapped) for the evaluation of sterility in manufacturing process for easy detection of contamination by Media Fill Test	30.10	500 gm	
TSM 318	FLUID THIOGLYCOLLATE MEDIUM (as per USP) for sterility testing of biologicals and for cultivation of aerobic, anaerobic and microaerophilic organisms	29.75	500 gm	
TSM 377	LURIA BROTH for cultivation and maintenance of recombinant strains of <i>Escherichia coli</i>	20.00	500 gm	
TSMH 2228	MacCONKEY BROTH (as per USP/EP/JP/BP) for detection enumeration of coliform bacteria	35.01	500 gm	
TSM 350	NUTRIENT BROTH for general cultivation of less fastidious microorganisms, can be enriched with blood or other biological fluids	13.00	500 gm	
TSM 332	SOYA CASEIN DIGEST MEDIUM (TRYPTONE SOYA BROTH) (CASO BROTH) for the evaluation of sterility in manufacturing processes	30.00	500 gm	
TSMV 332	SOYA CASEIN DIGEST MEDIUM (TRYPTONE SOYA BROTH) (CASO BROTH) (VEG.) for the evaluation of sterility in manufacturing processes	30.00	500 gm	
TSM 332K	SOYA CASEIN DIGEST MEDIUM WITH BCP routine use media for multiplication of large number of microorganisms with indicator with addition of sugar can be used	30.01	500 gm	
TSMH 332	SOYBEAN CASEIN DIGEST MEDIUM (as per USP/EP/JP/BP) for the evaluation of sterility in manufacturing processes	30.00	500 gm	
TSMVH 332	SOYBEAN CASEIN DIGEST MEDIUM (as per USP/EP/JP/BP) (VEG.) for the evaluation of sterility in manufacturing processes	30.00	500 gm	
TSMH 345	SOYA CASEIN DIGEST AGAR (as per USP/IP) general purpose medium used with or without blood for enrichment and isolation of various fastidious microorganisms	40.00	500 gm	
TSM 1938	SOYABEAN CASEIN DIGEST MEDIUM W/ MANNITOL, STERILE POWDER is γ -irradiated sterile powder recommended for the evaluation of sterility in manufacturing process. It can also be used for cultivation of a wide variety of microorganisms	30.02	100 gm	
TSMH 2454	STERILE TRYPTONE SOYA BROTH W/COLOUR INDICATOR (AS PER USP/EP/JP) for cultivation of a wide variety of microorganisms and recommended for sterility testing of Moulds and lower bacteria	27.80	500 gm	

Approximate number of vials required for 500 gm of dehydrated medium. | * Store Between 2-8°C





Media Supplement

Supporting Accurate Growth & Differentiation

Culture Media forms the base for microbial cultivation. **Media Supplements** enhance its capability by improving growth promotion, selectivity, enrichment, and accurate identification of microorganisms.

Specially formulated with selective agents, enrichment factors, chromogenic components, blood, serum, and growth-promoting substances, Media Supplements optimize microbiological performance for targeted isolation and dependable test results.

Key Benefits

- Enhances growth of fastidious microorganisms
- Improves selectivity and microbial recovery
- Supports accurate identification and differentiation
- Compatible with a wide range of Culture Media
- Supports efficient routine microbiological analysis

β-Lactamase Supplement

Key Benefits

- Broad-spectrum β-lactam inactivation
- Ensures effective microbial recovery
- High enzymatic activity for consistent performance
- Supplied as lyophilized vials for enhanced stability and shelf life

Applications

- Sterility Testing
- Environmental Monitoring
- Antibiotic Residue Detection

Culture Media & Supplements

Media Supplements

CODE	PRODUCT NAME	PACK SIZE	
TS 267	ACRIFLAVIN-CEFSULODIN-VANCOMYCIN SUPPLEMENT (ACV SUPPLEMENT) for the isolation of <i>Enterohemorrhagic E. coli</i> (EHEC) Recommended With : TM 2020	5 vl	
TS 099	AEROMONAS SELECTIVE SUPPLEMENT * for the selective isolation of <i>Aeromonas</i> species Recommended With : TM 1136	5 vl	
TS 178	ALBUMIN GLUCOSE SUPPLEMENT * for rapid cultivation of <i>Mycobacterium tuberculosis</i> Recommended With : TM 099	5 vl	
TS 097	AMPICILLIN SUPPLEMENT for isolation of <i>Aeromonas hydrophila</i> from foods Recommended With : TM 1074 / TM 1137	5 vl	
TS 319	AMPICILLIN DEXTRIN SELECTIVE SUPPLEMENT for selective isolation and differentiation of <i>Aeromonas</i> species from water samples using membrane filter technique Recommended With : TM 1958, TM 1959	5 vl	
TS 129	ANTHRACIS SELECTIVE SUPPLEMENT for selective isolation of <i>Bacillus anthracis</i> Recommended With : TM 1271 / TM 2269	5 vl	
TS 255	BACILLUS CEREUS SELECTIVE SUPPLEMENT for selective isolation of <i>Bacillus species</i> Recommended With : TM 1885	5 vl	
TS 281	BACILLUS SELECTIVE SUPPLEMENT recommended for the selective isolation of <i>Bacillus</i> species Recommended With : TM 2110	5 vl	

* Store Between 2-8°C



CODE	PRODUCT NAME	PACK SIZE	
TS 100	BACTEROIDES SELECTIVE SUPPLEMENT * for selective isolation of <i>Bacteriodes</i> species Recommended With : TM 1144	5 vl	
TS 211	BASIC FUCHSIN (6 g/vl) for standard test for lactose fermenting organisms Recommended With : TM 373 / TMV 373	1 vl	
TS 322	BCSA SELECTIVE SUPPLEMENT for isolation of <i>Burkholderia cepacia</i> from the respiratory secretions of patients with cystic fibrosis and other non-clinical specimens Recommended With : TM 2022	5 vl 25 vl	
TS 261	β LACTAMASE I SUPPLEMENT (W/ 3300 UNITS/VIAL OF PENICILLINASE) for inactivation of <i>Penicillin</i> group of antibiotics	1 vl	
TS 262	β LACTAMASE II NON STERILE POWDER SUPPLEMENT for inactivation of <i>Penicillins</i> and <i>Cephalosporins</i> of first, second, third, fourth & fifth generation and <i>Carbapenem</i>	1 vl	
TS 263	β LACTAMASE MIXTURE SUPPLEMENT (W/ (>500 β I UNITS/VIAL) & (>50 β II UNITS/VIAL)) for inactivation of <i>Penicillins</i> and <i>Cephalosporins</i> of first, second, third, fourth generation and <i>Carbapenem</i>	1 vl	
TS 263G	β LACTAMASE MIXTURE SUPPLEMENT (STERILE) (W/ (>500 β I UNITS/VIAL) & (>50 β II UNITS/VIAL)) for inactivation of <i>Penicillins</i> and <i>Cephalosporins</i> of first, second, third, fourth generation and <i>Carbapenem</i>	1 vl	
TS 264	β LACTAMASE MIXTURE SUPPLEMENT (W/ 2 β II UNITS/MG & 20 β I UNITS/MG) for inactivation of <i>Penicillins</i> and <i>Cephalosporins</i> of first, second, third, fourth generation and <i>Carbapenem</i>	1 vl	
TS 326	β LACTAMASE II SUPPLEMENT (W/ >2 β II UNITS/MG) for inactivation of <i>Penicillins</i> and <i>Cephalosporins</i> of first, second, third, fourth generation and <i>Carbapenem</i>	1 vl	
TS 179	BOLTON SELECTIVE SUPPLEMENT for selective isolation of <i>Campylobacter</i> species from foods Recommended With : TM 1510	5 vl	
TS 012	BORDETELLA SELECTIVE SUPPLEMENT * for selective isolation of <i>Bordetella pertussis</i> Recommended With : TM 044 / TM 1147 / TM 062 / TM 063	5 vl	
TS 004	B.P. SULPHA SUPPLEMENT for suppression of growth and swarming of <i>Proteus</i> species Recommended With : TM 358 / TMV 358	5 vl	
TS 045	BROMO CRESOL PURPLE (15 mg/vl) for selective isolation and enumeration of faecal streptococci Recommended With : TM 136	5 vl	
TS 006	BRUCELLA SELECTIVE SUPPLEMENT for selective isolation of <i>Brucella</i> species from milk Recommended With : TM 041 / TM 514 / TM 680 / TM 050 / TM 682 / TM 1513 / TM 071 / TM 1842 / TM 2011	5 vl 25 vl	
TS 269	BROMO THYMOL BLUE SUPPLEMENT recommended for the selective cultivation of pathogenic fungi Recommended With : TM 090	5 vl	
TS 323	BURKHOLDERIA SELECTIVE SUPPLEMENT for isolation of <i>Burkholderia cepacia</i> from the respiratory secretions of patients with cystic fibrosis and other nonclinical specimens Recommended With : TM 1934	5 vl	
TS 010	CAMPYLOBACTER GROWTH SUPPLEMENT for the enhancement of growth & aerotolerance of <i>Campylobacter fetus</i> Recommended With : TM 041 / TM 071 / TMV 071 / TM 514 / TM 680 / TM 050 / TM 2044 / TM 2105 / TM 1842 / TM 2011	5 vl	

Culture Media & Supplements

Media Supplements

CODE	PRODUCT NAME	PACK SIZE	
TS 278	CAMPYLOBACTER SELECTIVE SUPPLEMENT, ABEYTA recommended by FDA BAM for the selective isolation of <i>Campylobacter</i> species Recommended With : TM 2105	5 vl	
TS 007	CAMPYLOBACTER SUPPLEMENT-I (Blaser Wang) for selective isolation of <i>Campylobacter</i> species Recommended With : TM 041 / TM 056 / TM 071 / TMV 071 / TM 1842 / TM 2011 / TM 514 / TM 1157 / TM 2044	5 vl	
TS 008	CAMPYLOBACTER SUPPLEMENT-II (Butzler) for selective isolation of <i>Campylobacter</i> species Recommended With : TM 041 / TM 514 / TM 071 / TM 2011 / / TMV 071 / TM 2044 / TM 1842	5 vl 25 vl	
TS 009	CAMPYLOBACTER SUPPLEMENT-III (Skirrow) for selective isolation of <i>Campylobacter</i> species Recommended With : TM 041 / TM 514 / TM 680 / TM 056 / TM 071 / TM 2044 / TM 1842 / TM 2011 / TMV 071	5 vl	
TS 067	CAMPYLOBACTER SUPPLEMENT-IV (Preston) for isolation of <i>Campylobacter</i> species Recommended With : TM 603 / TM 956 / TM 1590 / TMV 1590	5 vl	
TS 102	CAMPYLOBACTER SUPPLEMENT-V (BFCSA) for selective isolation of <i>Campylobacter</i> species Recommended With : TM 1146	5 vl	
TS 026	CAMPYLOBACTER SELECTIVE SUPPLEMENT for the selective isolation of <i>Helicobacter pylori</i> (<i>Campylobacter pylori</i>) Recommended With : TM 071 / TM 2044 / TMV 071	5 vl	
TS 081	CAMPYLOBACTER SELECTIVE SUPPLEMENT (Wesley) for the selective isolation of <i>Campylobacter</i> species Recommended With : TM 928	5 vl	
TS 145	CAMPYLOBACTER SELECTIVE SUPPLEMENT (Karmali) for isolating thermotolerant <i>Campylobacter</i> species Recommended With : TM 1361 / TM 2112	5 vl 25 vl	
TS 146	CAMPYLOBACTER SELECTIVE SUPPLEMENT (With Hemin) (Karmali) for isolating thermotolerant <i>Campylobacter</i> species Recommended With : TM 1361	5 vl 25 vl	
TS 027	CAMPYLOBACTER SUPPLEMENT-VI (Butzler) for selective isolation of thermotolerant <i>Campylobacter</i> species Recommended With : TM 071 / TMV 071 / / TM 2044	5 vl 25 vl	
TS 103	CAT SELECTIVE SUPPLEMENT for selective isolation of thermophilic <i>Campylobacter</i> species Recommended With : TM 1146	5 vl	
TS 273	CATC SUPPLEMENT recommended for the selective isolation and detection of <i>Enterococcus faecalis</i> by means of TTC reduction Recommended With : TM 2041	5 vl	
TS 104	C.B.I. SUPPLEMENT for the selective isolation of <i>Clostridium botulinum</i> Recommended With : TM 1154	5 vl	
TS 101	CCDA SELECTIVE SUPPLEMENT for the selective cultivation of <i>Campylobacter</i> species Recommended With : TM 1145	5 vl	
TS 275	CC SUPPLEMENT * recommended for the selective cultivation of pathogenic fungi Recommended With : TM 090	5 vl	

CODE	PRODUCT NAME	PACK SIZE	
TS 301	CEFIXIME SUPPLEMENT recommended for selective isolation of <i>E.coli</i> O157:H7 Recommended With : TM 2231	5 vI	
TS 219	CEFOXITIN SUPPLEMENT for selective isolation of Methicillin Resistant <i>Staphylococcus aureus</i> from clinical specimens Recommended With : TM 1635 / TM 1561 / TM 1635	5 vI	
TS 075	CETRINIX SUPPLEMENT for the selective isolation of <i>Pseudomonas</i> species Recommended With : TM 626	5 vI 25 vI	
TS 077	CFC SUPPLEMENT for the selective isolation of <i>Pseudomonas</i> species Recommended With : TM 626	5 vI 25 vI	
TS 053	CHLORAMPHENICAL SELECTIVE SUPPLEMENT for the selective isolation of Yeasts and Moulds Recommended With : TM 563 / TM 667 / TM 982 / TMV 982 / TM 276	5 vI 25 vI	
TS 138	CHLORTETRACYCLINE SELECTIVE SUPPLEMENT for selective isolation of <i>Penicillium viridicatum</i> and related species from foods Recommended With : TM 1330 / TM 1584 / TM 2306 / TM 2319 / TM 2320 / TM 1951 / TM 2438	5 vI 25 vI	
TS 213	CHROMOGENIC CANDIDA SELECTIVE SUPPLEMENT for isolation and identification of <i>Candida</i> species from mixed cultures Recommended With : TM 1588	5 vI 25 vI	
TS 282	CHROMOGENIC CANDIDA DIFFERENTIAL SELECTIVE SUPPLEMENT recommended for rapid and direct isolation and identification of <i>Candida</i> species from mixed cultures Recommended With : TM 2113	5 vI	
TS 283	CHROMOGENIC COLISTIN RESISTANT SELECTIVE SUPPLEMENT for isolation & differentiation of colistin-resistant pathogens Recommended With : TM 2117	5 vI	
TS 248	CHROMOGENIC EC O157:H7 SELECTIVE SUPPLEMENT for selective isolation and detection of <i>Escherichia coli</i> O157:H7 from food samples Recommended With : TM 1854 / TM 2123 / TM 2017	5 vI 25 vI	
TS 313	CHROMOGENIC EC O157:H7 SELECTIVE SUPPLEMENT, MODIFIED for selective isolation and detection of <i>Escherichia coli</i> O157:H7 from food samples Recommended With : TM 2441	5 vI	
TS 187	CHROMOGENIC ECC SELECTIVE SUPPLEMENT for the selective isolation of <i>Escherichia coli</i> and coliforms from water and food samples Recommended With : TM 1341	5 vI	
TS 215	CHROMOGENIC ENTEROCOCCUS FAECIUM SELECTIVE SUPPLEMENT for selective isolation and detection of <i>Enterococcus faecium</i> Recommended With : TM 1632	5 vI	
TS 031	CHROMOGENIC ENRICHMENT SUPPLEMENT for the isolation, detection and enumeration of <i>Listeria monocytogenes</i> Recommended With : TM 1845 / TM 004	5 vI	
TS 204	CHROMOGENIC KLEBSIELLA SELECTIVE SUPPLEMENT for the selective isolation of <i>Klebsiella</i> species from water and other sources Recommended With : TM 1633	5 vI	
TS 205	CHROMOGENIC LISTERIA SELECTIVE SUPPLEMENT for the isolation, detection and enumeration of <i>Listeria monocytogenes</i> Recommended With : TM 1845 / TM 004 / TM 1634	5 vI	

Culture Media & Supplements

Media Supplements

CODE	PRODUCT NAME	PACK SIZE	
TS 206	CHROMOGENIC MeReSa SELECTIVE SUPPLEMENT for the selective isolation of Methicillin Resistant <i>Staphylococcus aureus</i> from clinical specimens Recommended With : TM 1635	5 vl	
TS 302	CHLORTETRACYCLINE SELECTIVE SUPPLEMENT, MODIFIED recommended for selective isolation of <i>Penicillium viridicatum</i> and related species from foods Recommended With : TM 2232	5 vl	
TS 259	CHLORTETRACYCLINE SELECTIVE SUPPLEMENT, MODIFIED for selective isolation of <i>Penicillium viridicatum</i> and related species isolated from food, in accordance with FDA BAM, 2001, Chapter 18 Yeasts, Moulds and Mycotoxins Recommended With : TM 116	5 vl	
TS 295	CHROMOGENIC SUPPLEMENT for the enumeration of faecal coliform by membrane filter technique Recommended With : TM 2186	5 vl	
TS 253	CHROMOGENIC VRE AGAR SUPPLEMENT for selective isolation of Vancomycin Resistant enterococci (VRE) Recommended With : TM 1907	5 vl	
TS 298	CIPROFLOXACIN CLINDAMYCIN SELECTIVE SUPPLEMENT recommended for selective cultivation of Lactic acid bacteria from food Recommended With : TM 2217	5 vl	
TS 024	CLOSTRIDIUM DIFFICILE SUPPLEMENT for selective isolation of <i>Clostridium difficile</i> Recommended With : TM 068 / TM 2042	5 vl 25 vl	
TS 274	CLOSTRIDIUM DIFFICILE SELECTIVE SUPPLEMENT recommended for the selective isolation of <i>Clostridium difficile</i> Recommended With : TM 2043	5 vl	
TS 250	CLOSTRIDIUM PERFRINGENS SUPPLEMENT for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 1826 / TM 2400	5 vl	
TS 271	COLISTIN SELECTIVE SUPPLEMENT for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 2027	5 vl	
TS 141	CPC SUPPLEMENT for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 1422	5 vl	
TS 315	CYCLOSERINE-CEFOXITIN SUPPLEMENT for the isolation and presumptive identification of <i>Clostridium difficile</i> , a recognized cause of pseudomembranous colitis (antimicrobial agent-associated) Recommended With : TM 2446	5 vl 25 vl	
TS 207	DMACA REAGENT (10 ml/vl) for enumeration and differentiation of enteric pathogens in urinary tract infections Recommended With : TM 1639	1 vl	
TS 028	DERMATO SUPPLEMENT for the determine the tryptophan-deaminase activity Recommended With : TM 083	5 vl	
TS 059	DIPHThERIA VIRULENCE SUPPLEMENT (Part A & B) (DOUBLE PACK) for isolation and presumptive identification of <i>Corynebacterium diphtheriae</i> Recommended With : TM 450	1 Kit	
TS 135	DOYLE'S ANTIBIOTIC SUPPLEMENT for selective isolation of <i>Campylobacter</i> species Recommended With : TM 987	5 vl 25 vl	

CODE	PRODUCT NAME	PACK SIZE	
TS 287	ECC SELECTIVE SUPPLEMENT MODIFIED recommended for detection of <i>Escherichia coli</i> and other coliform Recommended With : TM 2127	5 vl	
TS 001	EGG YOLK TELLURITE EMULSION (100 ml/vl) for identification of staphylococci Recommended With : TM 1579 / TM 635 / TM 636 / TM 943 / TMV 358 / TM 1337 / TM 358 / TMV 358 / TMH 119	5 vl	
TS 290	EGG YOLK EMULSION, 50% for use in various Culture Media in accordance with FDA BAM 1998 Recommended With : TM 2151, TM 2311	5 vl	
TS 002	EGG YOLK EMULSION (100 ml/vl) for use in various Culture Media Recommended With : TM 1120 / TM 358 / TM 1881 / TM 635 / TM 943 / TMV 358 / TM 1154 / TM 1886 / TM 724 / TM 1843 / TMV 615 / TM 197 / TM 206 / TMV 206 / TM 1206 / TM 1125 / TM 1207 / TM 1221 / TMV 197 / TM 197 / TM 210 / TMHV 114 / TMH 114	5 vl	
TS 168	ENTEROCOCCUS FAECIUM SELECTIVE SUPPLEMENT for selective isolation and detection of <i>Enterococcus faecium</i> Recommended With : TM 1470	5 vl	
TS 157	ENTEROCUCUS SELECTIVE SUPPLEMENT for enumeration of intestinal enterococci from surface and waste water Recommended With : TM 1388 / TM 2220	5 vl	
TS 090	ESCULIN SUPPLEMENT (0.5 g/vl) for detection of group D streptococci by means of esculin hydrolysis Recommended With : TM 037	5 vl	
TS 033	FRASER ENRICHMENT SUPPLEMENT for isolation and enrichment of <i>Listeria monocytogenes</i> from food and environmental samples Recommended With : TM 113	5 vl	
TS 035	FRASER SELECTIVE SUPPLEMENT for selective isolation and enrichment of <i>Listeria monocytogenes</i> from food and environmental samples Recommended With : TM 113 / TM 1472 / TM 1937 / TM 2092	5 vl	
TS 034	FRASER SUPPLEMENT for selective isolation and enumeration of <i>Listeria monocytogenes</i> from food, animal feeds etc Recommended With : TM 1472 / TM 1937 / TM 2092 / TM 2093	5 vl	
TS 277	FRIIS SUPPLEMENT a selective supplement for detection of non-avian <i>Mycoplasmas</i> in pharmaceutical products in accordance with European Pharmacopeia Recommended With : TM 2095	5 vl	
TS 184	GBS SUPPLEMENT for selective isolation and presumptive identification of group B streptococci Recommended With : TM 1545	5 vl	
TS 036	GC SUPPLEMENT W/ANTIBIOTICS for selective isolation and cultivation of pathogenic <i>Neisseria</i> Recommended With : TM 116 / TMV 116 / TM 933	5 vl 25 vl	
TS 056	G.N. SPORE ANAEROBIC SUPPLEMENT for the selective isolation of gram-negative anaerobes Recommended With : TM 915 / TM 916	5 vl 25 vl	
TS 299	GROWTH SUPPLEMENT-I FOR MSM a growth supplement for the enrichment of <i>Salmonella</i> species Recommended With : TM 2218	5 vl	
TS 300	GROWTH SUPPLEMENT-II FOR MSM a growth supplement for the enrichment of <i>Salmonella</i> species Recommended With : TM 2218	5 vl	

Culture Media & Supplements

Media Supplements

CODE	PRODUCT NAME	PACK SIZE	
TS 109	GTC SUPPLEMENT for selective isolation of enterococci Recommended With : TM 1190	5 vl	
TS 025	G. VAGINALIS SELECTIVE SUPPLEMENT for selective isolation of <i>Gardnerella vaginalis</i> Recommended With : TM 071 / TMV 071 / TM 2044	5 vl	
TS 217	GENTAMYCIN SUPPLEMENT for detection of <i>Clostridium perfringens</i> from pharma products Recommended With : TM 701 / TMH 116	5 vl	
TS 128	GENTA-OXY SELECTIVE SUPPLEMENT for selective isolation of Vancomycin Resistant enterococci Recommended With : TM 1061	5 vl 25 vl	
TS 111	GEORGE KIMMIG SELECTIVE SUPPLEMENT for selective isolation and cultivation of fungi Recommended With : TM 1211	5 vl	
TS 047	GRUFT MYCOBACTERIAL SUPPLEMENT for selective cultivation of <i>Mycobacteria</i> Recommended With : TM 375 / TM 526 / TM 2167	5 vl	
TS 021	HEAMOGLOBIN POWDER for isolation of <i>Neisseria spp.</i> Recommended With : TM 064 / TM 116 / TMV 116 / TM 879 / TM 933 / TM 1105 / TM 2264	100 gm	
TS 144	HAEMOPHILUS GROWTH SUPPLEMENT for cultivation of <i>Haemophilus influenzae</i> Recommended With : TM 1352	5 vl	
TS 014	HORSE SERUM (Store at 2-20°C) for isolation and cultivation of <i>Mycoplasma</i> or <i>Trichomonas</i> or <i>Streptococcus</i> species Recommended With : TM 514 / TM / 680 / TM 050 / TM 1513 / TM 1544 / TM 1545 / TM 1366 / TM 768 / TM 1265 / TM 233 / TM 235 / TM 794 / TM 1266 / TM 1267 / TM 233 / TM 1308 / TM 1307	100 ml	
TS 068	KANAMYCIN SULPHATE SELECTIVE SUPPLEMENT for selective isolation of enterococci Recommended With : TM 748 / TM 749	5 vl	
TS 314	KANAMYCIN AND VANCOMYCIN SUPPLEMENT for the isolation and partial identification of obligately anaerobic gram-negative bacilli Recommended With : TM 2445	5 vl	
TS 110	KIMMIG SELECTIVE SUPPLEMENT for selective isolation and cultivation of fungi Recommended With : TM 1211	5 vl	
TS 134	KL VIRULENCE ENRICHMENT (20 ml/vl) for cultivation and in-vitro toxicity testing of <i>Corynebacterium diphtheriae</i> Recommended With : TM 984	5 vl 25 vl	
TS 257	KLEBSIELLA SELECTIVE SUPPLEMENT for selective isolation and easy detection of <i>Klebsiella</i> species from water and other sources Recommended With : TM 582 / TM 2192	5 vl	
TS 227	L. MONO SELECTIVE SUPPLEMENT-I (Earlier TS 147) for isolation of <i>Listeria</i> species Recommended With : TM 1443 / TM 1806 / TM 2148 / TM 1375	5 vl	

* Store Between 2-8°C



CODE	PRODUCT NAME	PACK SIZE	
TS 228	L. MONO SELECTIVE SUPPLEMENT-II (Earlier TS 148) for isolation of <i>Listeria</i> species Recommended With : TM 1443 / TM 1806 / TM 1375	5 vI	
TS 321	L. MONO ENRICHMENT SUPPLEMENT-I for selective differentiation of <i>Listeria monocytogenes</i> from other <i>Listeria</i> species Recommended With : TM 1443 / TM 1806	5 vI	
TS 229	L. MONO ENRICHMENT SUPPLEMENT-II (Earlier TS 149) * for selective differentiation of <i>Listeria monocytogenes</i> from other <i>Listeria</i> species Recommended With : TM 1443 / TM 1806 / TM 1375	5 vI	
TS 161	LACHICA SUPPLEMENT for selective isolation of <i>Aeromonas hydrophilla</i> Recommended With : TM 1369	1 vI	
TS 049	10% LACTIC ACID SOLUTION (10 ml/vI) for adjustment of acidic pH Recommended With : TM 080 / TM 340 / TM 326 / TM 397 / TM 204 / TMV 204 / TM 205 / TMV 205 / TM 114 / TM 115 / TM 398 / TM 397	5 vI 25 vI	
TS 069	LACTIC SUPPLEMENT (10 ml/vI) for selective isolation of Lactic acid bacteria in beer and brewing processes Recommended With : TM 752 / TM 753	5 vI	
TS 308	LACTOBACILLI SUPPLEMENT recommended for the selective isolation of lactobacilli from wine Recommended With : TM 2379	5 vI 25 vI	
TS 154	LANDER VPT SUPPLEMENT-A for selective enrichment and transport of <i>Campylobacter</i> spp.	2 vI	
TS 160	LANDER AF SUPPLEMENT-B for selective enrichment and transport of <i>Campylobacter</i> spp.	1 vI	
TS 294	LCN SUPPLEMENT used for the selective isolation of <i>Mycobacterium</i> from specimens containing mixed flora Recommended With : TM 2168	5 vI 25 vI	
TS 285	LECITHIN SOLUTION for selective differentiation of <i>Listeria monocytogenes</i> and <i>Listeria</i> species based on PCPLC activity Recommended With : TM 2124	5 vI	
TS 044	LEGIONELLA GROWTH SUPPLEMENT (Twin Pack) (Part A & B) for enhancing growth of <i>Legionella</i> species Recommended With : TM 760 / TM 759	5 vI	
TS 195	LEGIONELLA GROWTH SUPPLEMEN (W/O L-CYSTEINE) for enhancing growth of <i>Legionella</i> species Recommended With : TM 1218 / TM 2156	5 vI 25 vI	
TS 015	LEGIONELLA SELECTIVE SUPPLEMENT for selective isolation of <i>Legionella</i> species Recommended With : TM 684 / TM 759	5 vI 25 vI	
TS 016	LEGIONELLA SELECTIVE SUPPLEMENT-II for selective isolation of <i>Legionella</i> species from mixed cultures Recommended With : TM 684	5 vI	
TS 078	LEPTOSPIRA ENRICHMENT SUPPLEMENT for cultivation of <i>Leptospira</i> species	5 vI	

Culture Media & Supplements

Media Supplements

CODE	PRODUCT NAME	PACK SIZE	
TS 017	LEGIONELLA SELECTIVE SUPPLEMENT-III for selective isolation of <i>Legionella</i> species from mixed cultures Recommended With : TM 684	5 vl	
TS 018	LEGIONELLA SELECTIVE SUPPLEMENT-IV (MWY) for selective isolation of <i>Legionella</i> species Recommended With : TM 684 / TM 1894	5 vl	
TS 019	LEGIONELLA SUPPLEMENT for enhanced growth of <i>Legionella</i> species Recommended With : TM 684 / TM 1894 / TM 759	5 vl	
TS 114	LEGIONELLA GROWTH SUPPLEMENT (BCYE) for enhancing growth of <i>Legionella</i> species Recommended With : TM 1218 / TM 2156	5 vl	
TS 114A	LEGIONELLA SELECTIVE SUPPLEMENT (GVPN) for selective isolation of <i>Legionella</i> species Recommended With : TM 1218 / TM 2156	5 vl	
TS 115	LEGIONELLA (GVPC) SELECTIVE SUPPLEMENT for isolation of <i>Legionella</i> species Recommended With : TM 1218 / TM 1894 / TM 2156	5 vl 25 vl	
TS 116	LEGIONELLA BMPS SELECTIVE SUPPLEMENT for selective isolation of <i>Legionella</i> species Recommended With : TM 1218 / TM 2156	5 vl 25 vl	
TS 041	LINCO T SUPPLEMENT (LINCOMYCIN-COLISTIN-AMPHOTERICIN-B TRIMETHOPRIM) * contains antibiotics, which are targeted to selective isolation and cultivation of <i>Neisseria</i> species in GC AGAR BASE Recommended With : TM 116 / TMV 116	5 vl 25 vl	
TS 121	LISTERIA MOXALACTAM SUPPLEMENT for selective isolation and cultivation of <i>Listeria monocytogenes</i> from food and environmental specimens Recommended With : TM 1229	5 vl	
TS 292	LISTERIA MOXALACTAM SUPPLEMENT MODIFIED recommended for the selective isolation of <i>Listeria</i> species from food samples in accordance with FDA BAM,1998 Recommended With : TM 2162, TM 2163	5 vl	
TS 092	LISTERIA SELECTIVE SUPPLEMENT-II for selective isolation and cultivation of <i>Listeria monocytogenes</i> Recommended With : TM 1333 / TM 1444 / TM 761 / TMV 761	5 vl 25 vl	
TS 119	LISTERIA SELECTIVE SUPPLEMENT (PALCAM) for the selective isolation and identification of <i>Listeria monocytogenes</i> Recommended With : TM 1226 / TM 1227	5 vl 25 vl	
TS 117	LISTERIA UVM SUPPLEMENT-I for selective isolation and cultivation of <i>Listeria monocytogenes</i> Recommended With : TM 1225	5 vl	
TS 118	LISTERIA UVM SUPPLEMENT-II for selective isolation and cultivation of <i>Listeria monocytogenes</i> from clinical specimens Recommended With : TM 1225	5 vl	
TS 289	LM SELECTIVE SUPPLEMENT recommended for presumptive enumeration of <i>Listeria</i> species Recommended With : TM 2149	5 vl	
TS 231	MeReSa SELECTIVE SUPPLEMENT for selective isolation of Methicillin Resistant <i>Staphylococcus aureus</i> from clinical specimens Recommended With : TM 1561	5 vl	

* Store Between 2-8°C



CODE	PRODUCT NAME	PACK SIZE	
TS 318	M-CP SELECTIVE SUPPLEMENT-I recommended for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 2130	5 vl	
TS 320	M-CP SELECTIVE SUPPLEMENT-II recommended for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 2130	5 vl	
TS 317	M-CP SELECTIVE SUPPLEMENT, MODIFIED recommended for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 2130	5 vl	
TS 297	M-PA SELECTIVE SUPPLEMENT recommended for selective isolation of <i>Pseudomonas aeruginosa</i> Recommended With : TM 2197	5 vl	
TS 096	MCBRIDE LISTERIA SUPPLEMENT for selective isolation of <i>Listeria</i> species Recommended With : TM 905	5 vl	
TS 216	MEROPENEM SUPPLEMENT for enrichment of enterococci Recommended With : TM 1643 / TM 1642	5 vl	
TS 050	MIDDLEBROOK ADC GROWTH SUPPLEMENT for cultivation of <i>Mycobacteria</i> Recommended With : TM 213	5 vl	
TS 060	MIDDLEBROOK OADC GROWTH SUPPLEMENT for isolation and cultivation of <i>Mycobacteria</i> Recommended With : TM 215 / TM 786 / TM 618 / TM 214 / TM 590	5 vl	
TS 152	MKTT NOVOBIOCIN SUPPLEMENT for selective enrichment and isolation of <i>Salmonella</i> Recommended With : TM 1399	5 vl	
TS 268	MODIFIED CFC SELECTIVE SUPPLEMENT for selective isolation of <i>Pseudomonas</i> species Recommended With : TM 2023 / TM 2024	5 vl	
TS 270	MODIFIED CPC SUPPLEMENT for selective isolation of <i>Pseudomonas</i> species Recommended With : TM 2027	5 vl	
TS 286	MODIFIED L. MONO SELECTIVE SUPPLEMENT for selective isolation of <i>L. monocytogenes</i> and <i>Listeria</i> species based on PCPLC activity Recommended With : TM 2124	5 vl	
TS 113	MOXALACTAM SUPPLEMENT for isolation and cultivation of <i>Listeria monocytogenes</i> Recommended With : TM 1213	5 vl	
TS 123	MUELLER TELLURITE SERUM (25 ml/vl) for isolation, cultivation and differentiation of <i>Corynebacterium diphtheriae</i> Recommended With : TM 1050	5 vl	
TS 296	MUG SUPPLEMENT (50 MG PER VIAL) for measuring β - glucuronidase activity, for rapid and sensitive identification of <i>Escherichia coli</i> Recommended With : TM 2186	5 vl	
TS 124	MYCOPLASMA CULTIVATION SUPPLEMENT for isolation and cultivation of <i>Mycoplasma</i> Recommended With : TM 1265	5 vl	
TS 052	MYCOPLASMA ENRICHMENT SUPPLEMENT for isolation of <i>Mycoplasma</i> Recommended With : TM 233 / TM 794 / TM 235	5 vl 25 vl	

Culture Media & Supplements

Media Supplements

CODE	PRODUCT NAME	PACK SIZE	
TS 310	MYCOPLASMA SELECTIVE SUPPLEMENT supplement for the isolation of <i>Mycoplasma</i> Recommended With : TM 2405	5 vl	
TS 249	MYCOPLASMA UROGENITAL SELECTIVE SUPPLEMENT for selective cultivation of <i>Mycoplasma</i> associated with urogenital infections Recommended With : TM 1267	5 vl 25 vl	
TS 020	NALIDIXIC SELECTIVE SUPPLEMENT for selective isolation of <i>Pseudomonas aeruginosa</i> from clinical specimens Recommended With : TM 060 / TMV 060 / TM 2035	5 vl	
TS 293	NAMC LISTERIA SELECTIVE SUPPLEMENT recommended for the selective enrichment of <i>Listeria</i> species from food Recommended With : TM 2165	5 vl	
TS 303	NEO ENRICHMENT SELECTIVE SUPPLEMENT recommended for the selective isolation of <i>Listeria</i> species from food samples in 24 hours Recommended With : TM 2253	5 vl	
TS 095	NEOMYCIN SUPPLEMENT for selective isolation of <i>Streptococcus</i> species Recommended With : TM 927 / TM 1274	5 vl	
TS 055	NON SPORE ANAEROBIC SUPPLEMENT for selective isolation or cultivation of nonsporing anaerobic bacteria Recommended With : TM 915 / TM 916	5 vl	
TS 051	NOVOBIOCIN SUPPLEMENT for presumptive identification of <i>Aeromonas hydrophila</i> Recommended With : TM 1400 / TM 1127 / TM 546	5 vl 25 vl	
TS 070	NOVOBIOCIN SELECTIVE SUPPLEMENT for the presumptive detection of <i>Salmonella</i> species in foods and feed materials Recommended With : TM 770	5 vl	
TS 304	OFPBL SELECTIVE SUPPLEMENT recommended for selective isolation of <i>Burkholderia cepacia</i> from clinical specimens as well as non-clinical samples Recommended With : TM 2262	5 vl	
TS 062	OLEIC ALBUMIN SUPPLEMENT for cultivation of <i>Mycobacterium tuberculosis</i> Recommended With : TM 100 / TM 720	5 vl 25 vl	
TS 153	OXACILLIN RESISTANCE SELECTIVE SUPPLEMENT for screening Oxacillin Resistant microorganisms Recommended With : TM 1402	5 vl	
TS 120	OXFORD LISTERIA SUPPLEMENT for selective isolation of <i>Listeria</i> species Recommended With : TM 1229 / TM 2163	5 vl	
TS 071	OXYTETRA SELECTIVE SUPPLEMENT for selective isolation and cultivation of Yeasts and Moulds Recommended With : TM 1198 / TM 799 / TM 1061 / TM 800	5 vl 25 vl	
TS 073	PARK AND SANDERS SELECTIVE SUPPLEMENT-A for selective isolation and enumeration of thermotolerant <i>Campylobacter</i> species Recommended With : TM 1335 / TM 1403 / TM 2274	5 vl	
TS 074	PARK AND SANDERS SELECTIVE SUPPLEMENT-B for selective isolation and enumeration of thermotolerant <i>Campylobacter</i> species Recommended With : TM 1158 / TM 1335 / TM 1403 / TM 2274	5 vl	
TS 254	PCP SUPPLEMENT for selective isolation and cultivation of <i>Legionella</i> species from cooling towers, clinical and other material Recommended With : TM 1894	5 vl	

CODE	PRODUCT NAME	PACK SIZE	
TS 203	PEMBA SUPPLEMENT for selective cultivation of <i>Bacillus cereus</i> Recommended With : TM 1587	5 vL 25 vL	
TS 054	PERFRINGEN'S S.F.P. SUPPLEMENT (SFP SUPPLEMENT) for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 615 / TMV 615 / TM 846	5 vL 25 vL	
TS 158	PERFRINGENS SUPPLEMENT A * for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 808	5 vL 25 vL	
TS 159	PERFRINGENS SUPPLEMENT B * for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 808	5 vL 25 vL	
TS 076	PERFRINGEN'S T.S.C. SUPPLEMENT (TSC SUPPLEMENT) for the selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 615 / TMV 615 / TM 1826 / TM 809 / TM 902 / TM 2311 / TM 2400	5 vL 25 vL	
TS 311	PIG SERUM for culturing of <i>mycoplasmas</i> and viruses Recommended With : TM 2405	100 ml	
TS 058	POLYMYXIN B SELECTIVE SUPPLEMENT for selective isolation of various microorganisms Recommended With : TM 1120 / TM 1881 / TM 1523 / TM 1207 / TM 197 / TMV 197 / TM 1843 / TM 281 / TM 851 / TM 1848 / TM 1931 / TM 2221	5 vL 25 vL	
TS 132	POTASSIUM CHLORATE SUPPLEMENT for selective enrichment of <i>Yersinia enterocolitica</i> Recommended With : TM 1357	5 vL	
TS 232	POTASSIUM LACTATE 50% (10 ml/vL) for isolation and enumeration of wild yeast and pitching yeasts Recommended With : TM 172	5 vL	
TS 137	POTASSIUM SORBATE 10% (10 ml/vL) for the detection and isolation of acid resistant yeasts, <i>Zygosaccharomyces bailii</i> and <i>Zygosaccharomyces rouxii</i> in salads, sauces and dressings Recommended With : TM 1329	5 vL	
TS 003	POTASSIUM TELLURITE 3.5% (10 ml/vL) for selective isolation of staphylococci and <i>Corynebacteria</i> Recommended With : TM 358 / TMV 358 / TM 635 / TM 943 / TM 118 / TM 118A / TM 522	5 vL	
TS 005	POTASSIUM TELLURITE 1% (1 ml/vL) for the selective isolation of staphylococci and <i>Corynebacteria</i> Recommended With : TM 636 / TM 670 / TM 411 / TM 707 / TM 1165 / TM 984 / TMV 221 / TM 1048 / TM 221 / TM 235 / TM 438 / TM 879 / TM 439 / TM 1109 / TM 395 / TMV 395 / TM 1722 / TM 1820 / TM 2121 / TM 2420 / TMH 119	5 vL 25 vL	
TS 139	PROPIONIBACTERIA GROWTH SUPPLEMENT for selective cultivation of <i>Propionibacteria</i> Recommended With : TM 1331	5 vL	
TS 170	PYR REAGENT (10 ml/vL) for detection of Pyrrolidonyl arylamidase enzyme activity of β -hemolytic streptococci	1 vL	
TS 288	RAPID LISTERIA SELECTIVE SUPPLEMENT for selective enrichment of <i>Listeria</i> species from food samples Recommended With : TM 2128	5 vL	

Culture Media & Supplements

Media Supplements

CODE	PRODUCT NAME	PACK SIZE	
TS 176	RPF SUPPLEMENT (FIBRINOGEN PLASMA TRYPSIN INHIBITOR SUPPLEMENT) for isolation and enumeration of coagulase positive staphylococci from food and clinical samples Recommended With : TM 358 / TM 1579 / TMV 358 / TM 2300	5 vl 25 vl	
TS 171	RIPPEY CABELLI SUPPLEMENT for cultivation of <i>Aeromonas hydrophila</i> Recommended With : TM 1480	5 vl	
TS 048	ROSOLIC ACID (0.1 g/vl) for isolation of coliform bacteria Recommended With : TM 178 / TMV 178 / TM 179 / TM 180	5 vl	
TS 305	S.T.A. SELECTIVE SUPPLEMENT recommended for the selective isolation of <i>Brochothrix thermosphacta</i> from meat products Recommended With : TM 2315	5 vl 25 vl	
TS 212	SALENRICH SELECTIVE SUPPLEMENT for two-step enrichment of sub-lethally injured <i>Salmonella</i> from food and feeds Recommended With : TM 1606	5 vl 25 vl	
TS 306	SALMONELLA SELECTIVE ENRICHMENT SUPPLEMENT recommended for selective isolation and differentiation of <i>Salmonella</i> species from coliforms by chromogenic method Recommended With : TM 2328	5 vl	
TS 252	SALMONELLA SELECTIVE SUPPLEMENT for isolation and identification of <i>Salmonella</i> spp. Recommended With : TM 1839	5 vl	
TS 307	SELECTIVE SUPPLEMENT FOR MRSA recommended for the selective detection of Methicillin Resistance <i>Staphylococcus aureus</i> from clinical specimens Recommended With : TM 2332	5 vl 25 vl	
TS 167	SM SELECTIVE SUPPLEMENT for selective isolation and cultivation of <i>Pseudomonas solanacearum</i> Recommended With : TM 1076	5 vl	
TS 054	S.F. P. SUPPLEMENT (Perfringens S.F.P. Supplement) for selective isolation of <i>Clostridium perfringens</i> Recommended With : TM 615 / TMV 615 / TM 846	5 vl 25 vl	
TS 082	SHIGELLA SELECTIVE SUPPLEMENT for selective isolation of <i>Shigella</i> species Recommended With : TM 857	5 vl	
TS 108	SODIUM BICARBONATE SOLUTION (20 ml/vl) for selective isolation of <i>Enterococcus</i> species Recommended With : TM 1190	5 vl	
TS 089	STAPH STREPTO SUPPLEMENT for selective isolation of staphylococci and streptococci Recommended With : TM 071 / TMV 071 / TM 2044	5 vl	
TS 265	STERILE BETA LACTAMASE-I SUPPLEMENT (READY-TO-USE) for inactivation of Penicillin	6 vl	
TS 291	STERILE CHARCOAL SUPPLEMENT FOR LEGIONELLA AGAR recommended for growth of <i>Legionella</i> species Recommended With : TM 2156	5 vl 25 vl	
TS 256	STREPTOCOCCUS SELECTIVE SUPPLEMENT recommended for the selective isolation and cultivation of streptococci of medical and veterinary importance Recommended With : TM 071 / TMV 071 / TM 2044	5 vl 25 vl	
TS 011	STREPTO SUPPLEMENT for selective cultivation of <i>Streptococcus</i> species Recommended With : TM 041 / TM 071 / TM 1842 / TM 2011 / TM 2044	5 vl 25 vl	

CODE	PRODUCT NAME	PACK SIZE	
TS 013	SULPHA SUPPLEMENT for selective cultivation of <i>Salmonella</i> species Recommended With : TM 364 / TM 047 / TM 046 / TM 952 / TM 1512	5 vI	
TS 208	TDA REAGENT (10 ml/vI) for determination of tryptophan-deaminase activity Recommended With : TM 1639	1 vI	
TS 087	TELLURITE CEFIXIME SUPPLEMENT for isolation of <i>Escherichia coli</i> 0157:H7 Recommended With : TM 1393 / TM 1915	5 vI	
TS 131	TICARCILLIN SUPPLEMENT for selective enrichment of <i>Yersinia enterocolitica</i> Recommended With : TM 1357	5 vI	
TS 063	TOLUIDINE BLUE (0.1 g/vI) for detection of deoxyribonuclease activity Recommended With : TM 081	5 vI 25 vI	
TS 080	TRIBUTYRIN (10 ml/vI) for detection of lipolytic microorganisms Recommended With : TM 885	5 vI	
TS 251	TRICHODERMA HARZIANUM SELECTIVE SUPPLEMENT for selective isolation of <i>Trichoderma harzianum</i> Recommended With : TM 1873	5 vI	
TS 112	TRICHOMONAS SELECTIVE SUPPLEMENT-I for selective isolation of <i>Trichoderma</i> species from clinical specimens Recommended With : TM 1212	5 vI	
TS 133	TRICHOMONAS SELECTIVE SUPPLEMENT-II for selective isolation of <i>Trichoderma</i> species Recommended With : TM 1307	5 vI	
TS 042	TTC SOLUTION 1% (10 ml/vI) for detection of microbial growth by means of TTC reduction Recommended With : TM 1151 / TM 136 / TM 747 / TM 138 / TM 751 / TM 139 / TMV 751 / TM 1251 / TM 388 / TM 226 / TM 1274 / TM 440 / TMV 440 / TM 441 / TM 1311	5 vI 25 vI	
TS 309	TYCSB SUPPLEMENT recommended for the selective isolation of <i>S.mutans</i> Recommended With : TM 2398	5 vI	
TS 030	UREA 40% (5 ml/vI) for detection of urease activity Recommended With : TM 477 / TM 1196 / TM 1554 / TM 781 / TM 394 / TM 1851	5 vI 25 vI	
TS 126	UREA SOLUTION 5% (5 ml/vI) for selective cultivation of <i>Mycoplasma hominis</i> and <i>Ureaplasma urealyticum</i> Recommended With : TM 1267	5 vI 25 vI	
TS 043	VANCLO-T SUPPLEMENT (Vancomycin-Colistin-Amphotericin-B-Trimethoprim) for selective isolation of <i>Neisseria</i> species Recommended With : TM 116 / TMV 116	5 vI 25 vI	
TS 218	VANCOMYCIN SUPPLEMENT for isolation of <i>Cronobacter sakazakii</i> Recommended With : TM 1642 / TM 2047	5 vI 25 vI	
TS 280	VANCOMYCIN SUPPLEMENT for isolation of <i>Cronobacter sakazakii</i> Recommended With : TM 1829	5 vI	

Culture Media & Supplements

Media Supplements

CODE	PRODUCT NAME	PACK SIZE	
TS 038	V.C.N. SUPPLEMENT for selective isolation of <i>Neisseria gonorrhoeae</i> and <i>Neisseria meningitidis</i> Recommended With : TM 116 / TMV 116 / TM 933 / TM 1105	5 vl 25 vl	
TS 039	V.C.N.T. SUPPLEMENT for selective isolation of <i>Neisseria gonorrhoeae</i> and <i>Neisseria meningitidis</i> Recommended With : TM 116 / TMV 116 / TM 933 / TM 1105	5 vl 25 vl	
TS 022	VITAMINS AMINO GROWTH SUPPLEMENT (VITAMINS & AMINO ACIDS MIXTURE) (Double pack) for cultivation of wide variety of microorganisms Recommended With : TM 1267 / TM 879 / TM 933 / TM 1105 / TM 064 / TM 116 / TMV 116	5 vl 25 vl	
TS 272	VITAMINO GROWTH SUPPLEMENT, MODIFIED for the isolation and cultivation of fastidious microorganisms, especially <i>Neisseria</i> and <i>Haemophilus</i> species, from a variety of clinical specimens Recommended With : TM 2264	5 vl	
TS 266	VITAMINS K1 SUPPLEMENT a vitamin growth supplement used for the isolation of anaerobic organisms Recommended With : TM 1960	5 vl	
TS 276	V.P. SUPPLEMENT for selective isolation & differentiation of <i>Legionella</i> species Recommended With : TM 983	5 vl	
TS 098	WILLIS AND HOBBS'S SUPPLEMENT for isolation of <i>Clostridium perfringens</i> and <i>Clostridium botulinum</i> species from food Recommended With : TM 1319	5 vl	
TS 072	XLT4 SUPPLEMENT for isolation of non-typhi <i>Salmonella</i> Recommended With : TM 493	1 vl	
TS 023	YEAST AUTOLYSATE SUPPLEMENT for the cultivation of <i>Neisseria</i> species Recommended With : TM 064 / TM 116 / TMV 116	5 vl 25 vl	
TS 057	YERSINIA SELECTIVE SUPPLEMENT for selective isolation of <i>Yersinia enterocolitica</i> Recommended With : TM 501 / TM 2439 / TM 2440	5 vl 25 vl	
TS 321	YERSINIA SELECTIVE SUPPLEMENT-2 for use, in selective isolation and enumeration of <i>Yersinia enterocolitica</i> from food samples, in accordance with FDA BAM, 1998 Recommended With : TM 2439	5 vl	



Biological Media Bases

Reliable Support for Microbial Growth

It's provide the essential components needed to prepare high-quality **Culture Media**, supporting reliable microbial growth and consistent laboratory results across applications.

What sets Us Apart



Pharmaceutical-grade
raw materials



Stringent quality
control



Batch-to-batch
consistency



Optimized
formulations



Superior solubility,
clarity, and
growth performance



Customizable
blends



Wide product
range

Culture Media & Supplements

Biological Media Bases

CODE	PRODUCT NAME	PACK SIZE	
242 M	AGAR AGAR TYPE-I For general purpose, Bacteriological Grade	100 gm 500 gm 5 Kg	
1201	AGAR AGAR POWDER Bacteriological Grade	100 gm 500 gm	
1202	AGAR AGAR POWDER Highly Purified	100 gm 500 gm	
1203	AGAR AGAR POWDER (as per IP) Bacteriological Grade	100 gm 500 gm	
1243	AGAR AGAR POWDER HIGH GEL At low temp. like NMT 38°C	100 gm 500 gm	
1228	AGAR SPECIAL Equivalent to Agar Noble	100 gm 500 gm	
3540	AGAR AGAR For Molecular Biology	100 gm 500 gm	
1525	B. MEAT EXTRACT PASTE (BOVINE) For general purpose, Bacteriological Grade	500 gm 5 Kg	
1204	B. MEAT EXTRACT POWDER (TYPE-I) (BOVINE) For general purpose, Bacteriological Grade	500 gm 5 Kg	
1205 V	VEG. BE EXTRACT POWDER Culture Media Ingredient, Equivalent to Beef Extract Powder	500 gm	
1205	B. MEAT EXTRACT POWDER (BOVINE) For Microbiology	500 gm	
1206	B. MEAT EXTRACT (STD) TBL POWDER (Bovine) Culture Media Ingredient	500 gm	
1207	BILE SALT POWDER Culture Media Ingredient, Bacteriological Grade	100 gm 500 gm	
1234	BILE SALT (STD) TBL POWDER Culture Media Ingredient	500 gm	
1208	BILE SALT MIXTURE Equivalent to Bile Salt No.3	100 gm	
1235	BILE SALT MIXTURE (STD) TBL POWDER Culture Media Ingredient	100 gm	
1236 V	VEG. BIOPEPTONE Culture Media Ingredient, Equivalent to Biopeptone	500 gm	
1236	BIOPEPTONE (MIXTURE OF CASEIN & MEAT PEPTONE) Used as additional enrichment of microorganisms	500 gm	
1238	BIO PEPTONE (STD) TBL POWDER Culture Media Ingredient	500 gm	
1240 V	VEG. BH INFUSION POWDER Culture Media Ingredient, Equivalent to Brain Heart Infusion Powder	500 gm	
1240	BRAIN HEART INFUSION POWDER Culture Media Ingredient used for cultivation of fastidious microorganisms	500 gm	
1262 V	VEG. CA HYDROLYSATE Culture Media Ingredient, Equivalent to Casein Acid Hydrolysate	500 gm	
1262	CASEIN ACID HYDROLYSATE (TECHNICAL) Culture Media Ingredient	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
1230	CASEIN ACID HYDROLYSATE (STD) TBL POWDER Culture Media Ingredient	500 gm	
1231	CASEIN ACID HYDROLYSATE Less than 3% NaCl	500 gm	
1232	CASEIN ACID HYDROLYSATE Special for Pertussis Vaccine production	500 gm	
1233	CASEIN ACID HYDROLYSATE (Vitamin Free) Culture Media Ingredient	500 gm	
1501	CASITONE Peptic Digest Casein With Dipeptides and Tripeptides	500 gm	
1512 V	VEG. CEH ENZYMATIC HYDROLYSATE (TYPE-I) Culture Media Ingredient, Equivalent to Casein Enzymatic Hydrolysate	500 gm	
1512	CASEIN ENZYMATIC HYDROLYSATE (TYPE-I) (TRYPTONE TYPE-I) Tryptic Digest Casein With Dipeptides	500 gm	
1513 V	VEG. CEH ENZYMATIC HYDROLYSATE (STD) TBL POWDER Culture Media Ingredient	100 gm 500 gm	
1513	CASEIN ENZYMATIC HYDROLYSATE (STD) TBL POWDER Culture Media Ingredient	100 gm 500 gm	
1580	COTTON SEED PEPTONE Bacteriological Grade	500 gm	
1216	CORN MEAL EXTRACT POWDER Culture Media Ingredient	500 gm	
3531	CORN STEEP LIQUOR	250 Kg	
3533	DRIED YEAST POWDER (ACTIVE)	25 Kg	
4501	FISH PEPTONE Bacteriological Grade	500 gm	
1210	GELATIN CRYSTAL, Bloom Type B Bacteriological Grade	500 gm	
1502 V	VEG. GEL PEPTONE Culture Media Ingredient, Equivalent to Gelatin Peptone	500 gm	
1502	GELATONE (GELATIN PEPTONE) Culture Media Ingredient	500 gm	
3534	GELATONE (STD.) TBL POWDER Culture Media Ingredient	500 gm	
4502	GLUTEN RICE HYDROLYSATE Culture Media Ingredient	500 gm	
4503	GLUTEN MAIZE HYDROLYSATE Culture Media Ingredient	500 gm	
1211	HAEMOGLOBIN POWDER Culture Media Supplement	100 gm	
1251 V	VEG. H INFUSION POWDER Culture Media Ingredient, Equivalent to Heart Infusion Powder	500 gm	
1251	HEART INFUSION POWDER A media base for cultivation of fastidious Microorganisms	500 gm	

* Store Between 2-8°C



Culture Media & Supplements

Biological Media Bases

CODE	PRODUCT NAME	PACK SIZE	
1242 V	VEG. LB HYDROLYSATE Culture Media Ingredient, Equivalent to Lactalbumin Hydrolysate	500 gm	
1242	LACTALBUMIN HYDROLYSATE Culture Media Ingredient with Sugar Amino Acids	500 gm	
1302	INACTIVE DRIED YEAST Inactive dry yeast is a natural dough relaxation agent that can improve dough quality of many different baking products	500 gm 25 Kg	
4504	LACTALBUMIN HYDROLYSATE TBL Rich in all essential amino acids and suitable for vaccine productions	500 gm	
3535	LACTOSE (MONO) (γ-IRRADIATED) Culture Media Ingredient	500 gm 5 Kg 50 Kg	
4509	LACTOSE (MONO), STERILE (γ-IRRADIATED) (TRIPLE PACK) Culture Media Ingredient	500 gm 5 Kg 50 Kg	
1297	LIVER EXTRACT PASTE Bacteriological Grade	500 gm	
1212	LIVER EXTRACT POWDER Bacteriological Grade	500 gm	
1212 V	VEG. Liv EXTRACT POWDER Culture Media Ingredient, Equivalent to Liver Extract Powder	500 gm	
1241	LIVER EXTRACT POWDER (PROTOLYSED) For cultivation of fastidious anaerobes and bulk production of vaccines, steroids and enzymes	500 gm	
1239 V	VEG. Liv HYDROLYSATE POWDER Enzymatic Veg. Hydrolysate, equivalent to Liver Hydrolysate Powder used for cultivation of anaerobes	500 gm	
1239	LIVER HYDROLYSATE POWDER Hydrolysis of liver by enzyme, ingredients of culture media used for cultivation of anaerobes	500 gm	
1213 V	VEG. Liv INFUSION POWDER Equivalent to Liver Infusion Powder Suitable for Vaccine Production	500 gm	
1213	LIVER INFUSION POWDER Suitable for Vaccine Production	500 gm	
1298	MALT EXTRACT PASTE Bacteriological Grade	500 gm	
1214	MALT EXTRACT POWDER Bacteriological Grade	500 gm	
1215	MALT EXTRACT (STD) TBL POWDER Culture Media Ingredient	500 gm	
3536	MALTOSE (MONOHYDRATE) (γ-IRRADIATED) Culture Media Ingredient	500 gm 5 Kg	
4505	MALTOSE (Monohydrate) NRC GRADE For Vaccine Production	500 gm	
3537	D-MANNITOL (γ-IRRADIATED) Culture Media Ingredient	500 gm 5 Kg	
1299	MEAT EXTRACT PASTE Bacteriological Grade	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
4508	MEAT EXTRACT TYPE-I Bacteriological Grade	500 gm	
1217 V	VEG. M EXTRACT POWDER Culture Media Ingredient, Equivalent to Meat Extract Powder	500 gm	
1217	MEAT EXTRACT POWDER Culture Media Ingredient	500 gm	
1218	MEAT EXTRACT (STD) TBL POWDER Culture Media Ingredient	500 gm	
1219 V	VEG. M INFUSION POWDER Culture Media Ingredient, Equivalent to Meat Infusion Powder	500 gm	
1219	MEAT INFUSION POWDER For Vaccine Production	500 gm	
1504 V	VEG. MP PEPTONE Culture Media Ingredient, Equivalent to Meat Peptone	500 gm	
1504	MEAT PEPTONE Culture Media Ingredient	500 gm	
1259	MEAT PEPTONE-P Peptic Digest Animal Meat Tissue	500 gm	
1260	MEAT PEPTONE-T Tryptic Digest Animal Meat Tissue	500 gm	
1507 V	VEG. MYCO PEPTONE Culture Media Ingredient, Equivalent to Mycological Peptone	500 gm	
1507	MYCOLOGICAL PEPTONE Suitable for Cultivation of Yeasts & Moulds	500 gm	
1518	MYCOLOGICAL PEPTONE (STD) TBL POWDER Culture Media Ingredient	500 gm	
1524	OAT MEAL POWDER Culture Media Ingredient	500 gm	
1220	OX BILE POWDER (DRIED) Bacteriological Grade, General Purpose, Culture Media Ingredient	100 gm 500 gm	
4506	PEA PROTEIN HYDROLYSATE (PEA PEPTONE) Bacteriological Grade	500 gm	
3542	PEPTONE-F Bacteriological Grade	25 Kg	
1527	PEPTONE PASTE Bacteriological Grade	500 gm 5 Kg	
1506	PEPTONE-R Regular Bacteriological Grade	500 gm 5 Kg	
1578	PEPTONE-RG (GRANULAR) Bacteriological Grade	500 gm 5 Kg	
1505 V	VEG. PEPTONE-TBL Culture Media Ingredient, Equivalent to Peptone TBL	500 gm 5 Kg	
1505	PEPTONE-TBL For cultivation of fastidious microorganisms	500 gm 5 Kg	
1519	PEPTONE (STD) TBL POWDER Culture Media Ingredient	500 gm	
1508	PEPTONE SPECIAL Equivalent to Neopeptone	500 gm	

Culture Media & Supplements

Biological Media Bases

CODE	PRODUCT NAME	PACK SIZE	
1581	PEPTONE TYPE-III Bacteriological Grade	500 gm	
1221 V	VEG. PM HYDROLYSATE Culture Media Ingredient, Equivalent Peptonized Milk	500 gm	
1221	PEPTONIZED MILK Suitable for Lactobacilli, Yeast & Moulds	500 gm	
1509 V	VEG. PP PEPTONE Culture Media Ingredient, Equivalent to Proteose Peptone	500 gm 5 Kg	
1509	PROTEOSE PEPTONE (Tryptic Digest Animal Tissue) Rich in Proteose's Peptides	500 gm 5 Kg	
1319	POTATO PEPTONE For microbial cultures in industrial fermentation processes & cell culture supplementation	500 gm 25 Kg	
1523	PROTEOSE PEPTONE, (STD) TBL POWDER Culture Media Ingredient	500 gm	
221 M	SODIUM CHOLATE (CHOLIC ACID SODIUM) A selective inhibitory agent used in Culture Media	25 gm 100 gm	
220 M	SODIUM DEOXYCHOLATE A selective inhibitory agent used in Culture Media	25 gm 100 gm	
212 M	SODIUM TAUROCHOLATE A selective inhibitory agent used in Culture Media	100 gm 500 gm	
1261	SODIUM TAUROGLYCOCHOLATE A selective inhibitory agent used in Culture Media	100 gm 500 gm	
3543	SOLUBLE STARCH Culture Media Ingredient	500 gm	
1528	SOYATONE-R (SOYA PEPTONE-R) Bacteriological Grade	500 gm 5 Kg	
3546	SOYATONE TYPE-I (SOYA PEPTONE TYPE-I) Bacteriological Grade	500 gm 5 Kg	
1510	SOYATONE (SOYA PEPTONE) Peptic Digest of Soyabean Meal	500 gm 5 Kg	
1511	SOYATONE (STD), TBL POWDER (Soya Peptone) Culture Media Ingredient	500 gm	
3548	SOYATONE HN	25 Kg	
3544	SOYATONE PASTE	25 Kg	
1522	TRYPTONE-D Tryptic Digest Casein for Vaccine Production	500 gm	
1520	TRYPTONE-T For Tetanus Toxin	500 gm	
1529	TRYPTONE-R General purpose for Bacteriology	500 gm 5 Kg	

CODE	PRODUCT NAME	PACK SIZE	
1512 V	VEG. TRYPTONE (TYPE-I) Veg Enzymatic Digest with Dipeptides	500 gm	
1512	TRYPTONE (TYPE-I) Tryptic enzymatic digest of casein with dipeptides	500 gm	
1513 V	VEG. TRYPTONE (STD) TBL POWDER Culture Media Ingredient	100 gm 500 gm	
1513	TRYPTONE (STD) TBL POWDER Culture Media Ingredient	100 gm 500 gm	
1521 V	VEG. TRYP HYDROLYSATE Enzymatic Hydrolysate of protein that Replace Meat Infusion	500 gm 5 Kg	
1521	TRYPTOSE Enzymatic hydrolysate of Protein that replace Meat Infusion	500 gm 5 Kg	
1265	TRYPTOSE (STD) TBL POWDER To promote excellent growth of highly fastidious microorganisms	500 gm	
1629	WHEAT PEPTONE Culture Media Ingredient	500 gm	
1631	WHEAT EXTRACT TYPE-I Culture Media Ingredient	500 gm	
1637	WHEAT EXTRACT TBL POWDER Culture Media Ingredient	500 gm	
1300	YEAST AUTOLYSATE POWDER For Bacteriology	500 gm	
1531	YEAST EXTRACT PASTE For Bacteriology	500 gm 5 Kg	
3550	YEAST EXTRACT POWDER-F	25 Kg	
1532	YEAST EXTRACT POWDER-R Bacteriological Grade	500 gm 5 Kg	
1224	YEAST EXTRACT POWDER TYPE-I Bacteriological Grade	500 gm 5 Kg	
1225	YEAST EXTRACT TBL POWDER Low salt content refined rich in B-group of Vitamins	500 gm 5 Kg	
1264	YEAST EXTRACT (STD) TBL POWDER Culture Media Ingredient	500 gm	
4141	YEAST EXTRACT (ULTRA PURIFIED) Culture Media Ingredient	500 gm 5 Kg	
1628	YEAST PEPTONE Culture Media Ingredient	500 gm	

* Store Between 2-8°C





Ready-to-Use Culture Media

Open, Inoculate, and Incubate

- Meticulous packaging
- Favourably convenient
- Approved Quality
- Reliable growth formulations
- High reproducibility
- Triple-layer packaging
- Gamma-irradiated packs
- Proven Sterility
- Complies with Harmonized pharmacopeia (USP/IP/JP/EP/BP) & ISO regulations
- Industry-tailored range
- Large variety of formats in various sizes
- Sterility ensured
- Easy-open sterile packs

Sub-Categories


Plates | Bottles | Tubes | Kits | Bags | Slants



Ready-to-Use Culture Media Plates

Simple, High-end, and Reliable

- **Standard Sizes:** Pick from 90 mm, 55 mm, and 100 mm plates to fit your testing needs.
- **Antibiotic Susceptibility Testing:** Get precise results with 100 mm plates made for this specific analysis.
- **Partitioned Plates:** Boost efficiency with 90 mm plates offering Bi-partitioning

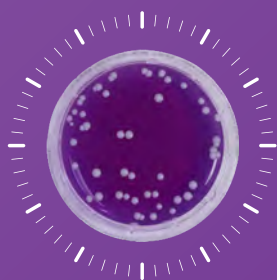


55 mm Plates



90 mm Plates



100 mm Plates

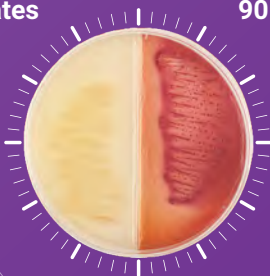

Bi-Plates



Chromogenic Plates

CODE	PRODUCT NAME		PACK SIZE
TMP 033	BAIRD PARKER AGAR PLATE for the isolation and enumeration of coagulase positive staphylococci from food and clinical sample	 	50 Plts
TMP 034	BLOOD AGAR PLATE W/ GENTAMICIN (5MCG / PLATE) * for selective isolation of <i>Streptococcus pneumoniae</i>		20 Plts 50 Plts
TMP 001	BRAIN HEART INFUSION AGAR PLATE for cultivation of fastidious pathogenic Bacteria, Yeasts and Moulds		50 Plts
TMP 035	BRAIN HEART INFUSION AGAR PLATE W/ BLOOD for cultivation of fastidious pathogenic bacteria requiring blood for growth		50 Plts
TMP 066	BURKHOLDERIA CEPACIA SELECTIVE AGAR (BCSA) (AS PER USP) for use in qualitative procedures for the selective and differential isolation of <i>Burkholderia cepacia</i> complex from respiratory secretions of patients with cystic fibrosis	 	50 Plts
TMP 002	CETRIMIDE AGAR PLATE for selective isolation of <i>Pseudomonas aeruginosa</i> from pharmaceutical product and clinical specimens	 	50 Plts
TMP 002GT	CETRIMIDE AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for selective isolation of <i>Pseudomonas aeruginosa</i> from pharmaceutical product (as per USP/EP/BP/JP/IP) and clinical specimens	 	50 Plts
TMP 020	CHOCOLATE AGAR PLATE * for isolation and cultivation of fastidious microorganisms like <i>Neisseria gonorrhoeae</i>		20 Plts 50 Plts
TMP 1858	CHROMOGENIC COLIFORM AGAR PLATE* for determination of coliform bacteria particularly <i>Enterobacteriaceae</i> on the basis of their ability to ferment lactose		20 Plts 50 Plts
TMP 073T	CHROMOGENIC E.COLI AGAR PLATE (TRIPLE PACK) for the detection and enumeration of <i>Escherichia coli</i> in foods without further confirmation on membrane filter or by indole reagent		20 Plts 50 Plts
TMP 1199	CHROMOGENIC UTI AGAR PLATE * for identification and confirmation of microorganisms causing urinary tract infections		20 Plts 50 Plts
TMP 063	CLED AGAR W/ ANDRADE INDICATOR PLATE for isolation and differentiation of microorganisms based on lactose fermentation		50 Plts
TMP 014	CLED AGAR W/ BROMOTHYMOL BLUE PLATE for isolation and differentiation of urinary pathogens on the basis of lactose fermentation		50 Plts
TMP 032GT	COLUMBIA AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for the selection and subculture of <i>Clostridium</i> species in accordance with harmonized method of USP/EP/BP/JP/IP		50 Plts
TMP 031	COLUMBIA 5% BLOOD AGAR PLATE for isolation and cultivation of fastidious organisms		20 Plts 50 Plts
TMP 031A	COLUMBIA BLOOD AGAR W/ 10% SHEEP BLOOD for isolation and cultivation of fastidious organisms		50 Plts
TMP 074	DEOXYCHOLATE CITRATE AGAR (DCA) PLATE for selective isolation of enteric pathogens		50 Plts
TMP 003	DEY-ENGLEY NEUTRALIZING AGAR PLATE for disinfectant testing, where neutralization agent is important for determining its bactericidal activity		50 Plts
TMP 003GT	DEY-ENGLEY NEUTRALIZING AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for disinfectant testing, where neutralization agent is important for determining its bactericidal activity		50 Plts
TMP 025	EMB AGAR, LEVINE PLATE for isolation, enumeration and differentiation of <i>Enterobacteriaceae</i>	 	50 Plts
TMP 037	HEKTOEN ENTERIC AGAR PLATE for differential & selective isolation of <i>Salmonella</i> and <i>Shigella</i> from enteric pathological specimens	 	50 Plts
TMP 004	LURIA BERTANI AGAR PLATE for the cultivation and maintenance of recombinant strains of <i>E.coli</i> for genetic and molecular studies, may be used for routine isolation, cultivation of not particularly fastidious microorganisms		50 Plts

Ready-to-Use Culture Media

90 mm Plates

CODE	PRODUCT NAME		PACK SIZE
TMP 004A	LURIA BERTANI AGAR PLATE (W/ 100ug/ml AMPICILLIN) * for the cultivation and maintenance of recombinant strains of <i>E.coli</i> for genetic and molecular biology studies		50 Plts
TMP 004 KI	LURIA BERTANI AGAR PLATE (W/ 100ug/ml KANAMYCIN & IPTG) * for the cultivation and maintenance of recombinant strains of <i>E.coli</i> for genetic and molecular biology studies		50 Plts
TMP 004K	LURIA BERTANI AGAR PLATE (W/ 100ug/ml KANAMYCIN) * for the cultivation and maintenance of recombinant strains of <i>E.coli</i> for genetic and molecular biology studies		50 Plts
TMP 004C	LURIA BERTANI AGAR PLATE (W/ 50ug/ml CHLORAMPHENICOL) for the cultivation and maintenance of recombinant strains of <i>E.coli</i> for genetic and molecular biology studies		50 Plts
TMP 1898	MacCONKEY AGAR "TBL" for cultivation and differentiation of enteric bacteria and gram-positive microorganisms from clinical samples especially <i>Enterococcus faecalis</i> within 18 hours		50 Plts
TMP 015	MacCONKEY AGAR PLATE (W/ 0.15% BILE SALTS, CV and NaCl) for selective isolation and differentiation of coliform organisms and other enteric pathogens	 	50 Plts
TMPH 015T	MacCONKEY AGAR PLATE (W/ 0.15% BILE SALTS, CV and NaCl) (TRIPLE PACK) for selective isolation and differentiation of coliform organisms and other enteric bacteria in accordance with the harmonized method of USP/EP/BP/JP/IP	  	50 Plts
TMP 016	MacCONKEY AGAR PLATE (W/O CV, NaCl W/ 0.5% SODIUM TAUROCHOLATE) for the selection and recovery of the <i>Enterobacteriaceae</i> and related enteric gram-negative bacilli from clinical, food and water samples	 	50 Plts
TMP 072T	MALT EXTRACT AGAR PLATE (TRIPLE PACK) for detection, isolation and enumeration of Yeasts and Moulds		50 Plts
TMP 021	MANNITOL SALT AGAR PLATE for selective isolation of pathogenic staphylococci	 	50 Plts
TMPH 021	MANNITOL SALT AGAR PLATE for selection and subculture of <i>Staphylococcus aureus</i> in accordance with harmonized method USP/EP/BP/JP/IP	 	50 Plts
TMPH 021T	MANNITOL SALT AGAR PLATE (TRIPLE PACK) for selection and subculture of <i>Staphylococcus aureus</i> in accordance with harmonized method USP/EP/BP/JP/IP	 	50 Plts
TMP 005	MUELLER HINTON AGAR PLATE for determination of susceptibility of microorganisms to antimicrobial agents		50 Plts
TMP 061	MUELLER HINTON AGAR PLATE W/ 5% SHEEP BLOOD for determination of susceptibility of <i>Streptococcus species</i> to antimicrobial agents		20 Plts 50 Plts
TMP 080	MUELLER HINTON AGAR NO. 2 + 5% SHEEP BLOOD for testing susceptibility of common and rapidly growing bacteria using antimicrobial discs by using Kirby bauer technique		20 Plts 50 Plts
TMP 006	NUTRIENT AGAR PLATE a general purpose medium used for cultivation of wide variety of microorganisms		50 Plts
TMP 006GT	NUTRIENT AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) a general purpose medium used for cultivation of wide variety of microorganisms		50 Plts
TMP 007	PLATE COUNT AGAR for determination of plate counts of microorganisms in food, water and waste water		50 Plts
TMP 007GT	PLATE COUNT AGAR (γ-IRRADIATED) (TRIPLE PACK) for determination of plate counts of microorganisms in food, water and waste water		50 Plts
TMP 008	POTATO DEXTROSE AGAR PLATE for the isolation and enumeration of Yeasts & Moulds from water, dairy, other food products and clinical samples		50 Plts

CODE	PRODUCT NAME		PACK SIZE	
TMPH 008	POTATO DEXTROSE AGAR PLATE for the subculture of fungi in accordance with harmonized method of USP/EP/BP/JP		50 Plts	
TMP 081GT	PSEUDOMONAS AGAR-F (FOR FLUORESCIN) (as per IP) (γ-IRRADIATED) (TRIPLE PACK) for detection of fluorescein production by <i>Pseudomonas</i> species		50 Plts	
TMP 082GT	PSEUDOMONAS AGAR-P (FOR PYOCYANIN) (as per IP) (γ-IRRADIATED) (TRIPLE PACK) for detection of Pyocyanin production by <i>Pseudomonas</i> species		50 Plts	
TMP 009	R2A AGAR PLATE for heterotrophic plate count of treated potable water, using longer incubation period		50 Plts	
TMP 009G	R2A AGAR PLATE (γ-IRRADIATED) for heterotrophic plate count of treated potable water, using longer incubation period		50 Plts	
TMP 009GT	R2A AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for heterotrophic plate count of treated potable water, using longer incubation period		50 Plts	
TMP 069GT	REINFORCED CLOSTRIDIAL AGAR (γ-IRRADIATED) (TRIPLE PACK) Used for the cultivation and enumeration of Clostridia and other anaerobes		50 Plts	
TMP 010	SABOURAUD CHLORAMPHENICOL AGAR PLATE for selective cultivation of Yeasts and Moulds	  	50 Plts	
TMP 010G	SABOURAUD CHLORAMPHENICOL AGAR PLATE (γ-IRRADIATED) for selective cultivation of Yeasts and Moulds	 	50 Plts	
TMP 010GT	SABOURAUD CHLORAMPHENICOL AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for selective cultivation of Yeasts and Moulds	 	50 Plts	
TMP 068	SABOURAUD CYCLOHEXIMIDE CHLORAMPHENICOL AGAR PLATE for selective isolation and cultivation of pathogenic Fungi		50 Plts	
TMP 011	SABOURAUD DEXTROSE AGAR PLATE for cultivation of Yeasts, Moulds and aciduric microorganisms		50 Plts	
TMP 011G	SABOURAUD DEXTROSE AGAR PLATE (γ-IRRADIATED) for cultivation of Yeasts, Moulds and aciduric microorganisms	  	50 Plts	
TMP 011GT	SABOURAUD DEXTROSE AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for the subculture of <i>Candida albicans</i> in accordance with harmonized method of USP/EP/BP/JP/IP	 	50 Plts	
TMP 086GT	SABOURAUD DEXTROSE AGAR PLATE W/ 0.3% SOYA LECITHIN AND 3% TWEEN 80 (γ-IRRADIATED) (TRIPLE PACK) for determining the efficiency of sanitization of containers, equipment surfaces, water miscible cosmetics, etc.	  	50 Plts	
TMP 040	SABOURAUD DEXTROSE AGAR PLATE W/ CYCLOHEXIMIDE for selective isolation of Fungi		50 Plts	
TMP 087GT	SABOURAUD DEXTROSE AGAR PLATE W/ LTHTh (γ-IRRADIATED) (TRIPLE PACK) for determining the efficiency of sanitization of containers, equipment surfaces, water miscible cosmetics, etc.	  	50 Plts	
TMP 058	SALMONELLA SHIGELLA AGAR PLATE (SS AGAR PLATE) for differentiation and selective isolation of <i>Salmonella</i> and <i>Shigella</i> species from pathological specimens, suspected foodstuff etc	 	50 Plts	
TMP 017	SHEEP BLOOD AGAR PLATE * for cultivation of fastidious organisms and studying haemolytic reactions		20 Plts 50 Plts	
TMP 018	SHEEP BLOOD AGAR BASE, MODIFIED PLATE * for cultivation and studying haemolytic reaction of <i>Bacillus cereus</i>	 	20 Plts 50 Plts	
TMP 019	SHEEP BLOOD SOYA BEAN CASEIN DIGEST AGAR PLATE * for cultivation of wide variety of microorganisms and studying haemolytic reactions	 	20 Plts 50 Plts	
TMP 012	SOYA BEAN CASEIN DIGEST AGAR PLATE for the subculture of aerobic organisms	 	50 Plts	
TMP 012G	SOYA BEAN CASEIN DIGEST AGAR PLATE (γ-IRRADIATED) for the subculture of aerobic organisms in accordance with the harmonized method of USP/EP/BP/JP/IP		50 Plts	
TMP 012GT	SOYA BEAN CASEIN DIGEST AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for the subculture of aerobic organisms in accordance with the harmonized method of USP/EP/BP/JP/IP		50 Plts	

Ready-to-Use Culture Media

90 mm Plates

CODE	PRODUCT NAME		PACK SIZE	
TMP 065GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.5% GLYCEROL, 0.5% POLYSORBATE 80, 0.07% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc		50 Plts	
TMP 064GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.5% GLYCEROL & 1% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc		50 Plts	
TMP 042GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCERINE & 1% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		50 Plts	
TMP 067GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.5% POLYSORBATE 80 & 0.1% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK) for determining of sanitization of containers, equipment's, surfaces, water miscible cosmetics etc		50 Plts	
TMP 043GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.1% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of microorganisms		50 Plts	
TMP 044GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.05% TWEEN 80 (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of microorganisms		50 Plts	
TMP 026GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of microorganisms		50 Plts	
TMP 045GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.5% LECITHIN & 4% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		50 Plts	
TMP 027GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.5% LECITHIN & 4% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.	 	50 Plts	
TMP 028GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 2% GLYCEROL (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of microorganisms		50 Plts	
TMP 030GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.5% LECITHIN & 2% Polysorbate 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.	 	50 Plts	
TMP 049GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ β-LACTAMASE (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of organisms and for inactivation of β -lactam antibiotics		50 Plts	
TMP 050GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ β-LACTAMASE II (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of organisms and for inactivation of cephalosporins of first, second, third and fourth generation		50 Plts	
TMP 051GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ β-LACTAMASE MIXTURE (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of organisms and for inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems		50 Plts	
TMP 052GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ β-LACTAMASE MIXTURE (25 IU/PLATE) (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of organisms and for inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems		50 Plts	
TMP 053GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN AND TWEEN 80 W/ β-LACTAMASE MIXTURE (γ-IRRADIATED) (TRIPLE PACK) recommended for determining efficiency of containers, equipment surfaces, water miscible cosmetics and for inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems		50 Plts	
TMP 054GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN AND TWEEN 80 W/ β-LACTAMASE II (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of containers, equipment surfaces, water miscible cosmetics and for inactivation of cephalosporins of first, second, third and fourth generation and penems		50 Plts	
TMP 055GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN AND TWEEN 80 W/ β-LACTAMASE (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of containers, equipment surfaces, water miscible cosmetics and inactivation of β -lactam antibiotics		50 Plts	
TMP 056GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.5% POLYSORBATE 80, 0.07% SOYA LECITHIN AND 5 IU/PLATE β-LACTAMASE MIXTURE (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of aerobes and fungi and for inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems		50 Plts	
TMP 029GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL & 4% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc		50 Plts	

CODE	PRODUCT NAME		PACK SIZE	
TMP 013	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN AND POLYSORBATE 80 (TYPTONE SOYA AGAR PLATE W/LECITHIN & POLYSORBATE 80) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		50 Plts	
TMP 013GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN AND POLYSORBATE 80 (TYPTONE SOYA AGAR PLATE W/LECITHIN & POLYSORBATE 80) (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		50 Plts	
TMPV 013GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN AND POLYSORBATE 80 (VEG.) [TYPTONE SOYA AGAR PLATE W/LECITHIN & POLYSORBATE 80 (VEG.)] (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		50 Plts	
TMP 057GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LTHTh (γ-IRRADIATED) (TRIPLE PACK) for determining the efficiency of sanitization of containers, equipment surfaces etc and for enumeration of organisms from water insoluble & fatty products containing antimicrobial or preservatives		50 Plts	
TMP 062GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.5% POLYSORBATE 80, 0.07% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK) for determining the efficiency of sanitization of containers, equipment, surfaces, water-miscible cosmetics, etc.		50 Plts	
TMP 070GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 4% GLYCEROL, 0.5% POLYSORBATE 80, 0.07% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc.		50 Plts	
TMP 084GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.5% POLYSORBATE 80 and 0.07% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of aerobes and fungi and for inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems	 	50 Plts	
TMP 085GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.3% SOYA LECITHIN AND 3% TWEEN 80 (γ-IRRADIATED) (TRIPLE PACK) for determining of sanitization of containers, equipment's, surfaces, water miscible cosmetics etc.	 	50 Plts	
TMP 088GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.05% POLYSORBATE 80, 0.07% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc.	 	50 Plts	
TMP 022	TCBS AGAR PLATE for selective isolation of <i>Vibrio cholerae</i> and enteropathogenic <i>Vibrios</i>	 	50 Plts	
TMP 083GT	VOGEL JOHNSON AGAR MEDIUM (as per IP) (γ-IRRADIATED) (TRIPLE PACK) for selective isolation of coagulase positive, mannitol fermenting <i>S. aureus</i> from foods & clinical samples		50 Plts	
TMP 023	VIOLET RED BILE AGAR PLATE for selective isolation, detection and enumeration of coli-aerogenes bacteria in water, milk and other dairy, food products		50 Plts	
TMPH 023	VIOLET RED BILE GLUCOSE AGAR PLATE for selection and subculture of bile tolerant organisms in accordance with the harmonized method of USP/EP/BP/JP/IP		50 Plts	
TMP 024	XYLOSE LYSINE DEOXYCHOLATE AGAR (XLD AGAR) PLATE for selective isolation and enumeration of <i>Salmonella typhi</i> and other <i>Salmonella</i> species	 	50 Plts	
TMPH 024	XYLOSE LYSINE DEOXYCHOLATE AGAR (XLD AGAR) PLATE for selective isolation and enumeration of <i>Salmonella typhi</i> and other <i>Salmonella</i> species with the harmonized method of USP/EP/BP/JP/IP	 	50 Plts	
TMPH 024T	XYLOSE LYSINE DEOXYCHOLATE AGAR (XLD AGAR) PLATE (TRIPLE PACK) for selective isolation and enumeration of <i>Salmonella typhi</i> and other <i>Salmonella</i> species with the harmonized method	 	50 Plts	

55 mm Contact Plate Equivalent to RODAC Plate

CODE	PRODUCT NAME		PACK SIZE	
TSP 001GT	CETRIMIDE AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for the isolation and subculture of <i>Pseudomonas aeruginosa</i> in accordance with harmonized method of USP/EP/BP/JP/IP	 	100 Plts	
TSP 1858	CHROMOGENIC COLIFORM AGAR PLATE for determination of coliform bacteria particularly Enterobacteriaceae on the basis of their ability to ferment lactose		100 Plts	
TSP 002	DEY-ENGLEY NEUTRALIZING AGAR PLATE for disinfectant testing, where neutralization agent is important for determining its bactericidal activity		100 Plts	

Ready-to-Use Culture Media

55 mm Contact Plate Equivalent to RODAC Plate

CODE	PRODUCT NAME		PACK SIZE	
TSP 003GT	DEY-ENGLEY NEUTRALIZING AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for disinfectant testing, where neutralization agent is important for determining its bactericidal activity		100 Plts	
TSP 004G	POTATO DEXTROSE AGAR (γ-IRRADIATED) for the subculture of Fungi in accordance with harmonized method of USP/EP/BP/JP		100 Plts	
TSP 005G	SABOURAUD CHLORAMPHENICOL AGAR PLATE (γ-IRRADIATED) for selective cultivation of Yeasts and Moulds		100 Plts	
TSP 060GT	SABOURAUD DEXTROSE AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for cultivation of Yeasts, Moulds and aciduric microorganisms		100 Plts	
TSP 082GT	SABOURAUD DEXTROSE AGAR PLATE W/ LTHTh (γ-IRRADIATED) (TRIPLE PACK) for determining the efficiency of sanitization of containers, equipment surfaces, water miscible cosmetics, etc.		100 Plts	
TSP 016GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.5% GLYCEROL, 0.5% POLYSORBATE 80, 0.07% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc		100 Plts	
TSP 006GT	SOYA BEAN CASEIN DIGEST AGAR PLATE (γ-IRRADIATED) (TRIPLE PACK) for the subculture of aerobic organisms in accordance with the harmonized method of USP/EP/BP/JP/IP		100 Plts	
TSP 007GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of microorganisms in accordance with harmonized method of USP/EP/BP/JP		100 Plts	
TSP 008GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN AND POLYSORBATE 80 (TYPTONE SOYA AGAR PLATE W/ LECITHIN & POLYSORBATE 80) (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		100 Plts	
TSPV 008GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN AND POLYSORBATE 80 (Veg) [TYPTONE SOYA AGAR PLATE W/LECITHIN & POLYSORBATE 80 (Veg) (γ-IRRADIATED) (TRIPLE PACK)] for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		100 Plts	
TSP 079GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN & TWEEN 80, W/ β LACTAMASE, W/ 1% GLYCEROL (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of containers, equipment surfaces, water miscible cosmetics		100 Plts	
TSP 009GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.5% LECITHIN & 4% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc.		100 Plts	
TSP 010GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.5% LECITHIN & 2% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc.		100 Plts	
TSP 011GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL & 4% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		100 Plts	
TSP 063GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL & 0.5% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for cultivation of wide variety of microorganisms		100 Plts	
TSP 012GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.5% LECITHIN & 4% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		100 Plts	
TSP 064GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LTHTh (γ-IRRADIATED) (TRIPLE PACK) for determining the efficiency of sanitization of containers, equipment surfaces etc. and for enumeration of organisms from water insoluble & fatty products containing antimicrobial or preservatives		100 Plts	
TSP 065GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCERINE AND 1% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		100 Plts	
TSP 066GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.05% TWEEN 80 (γ-IRRADIATED) (TRIPLE PACK) sterility test medium for cultivation of wide variety of microorganisms		100 Plts	
TSP 080GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.07% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.		100 Plts	

Ready-to-Use Culture Media

55 mm Contact Plate Equivalent to RODAC Plate

CODE	PRODUCT NAME	PACK SIZE
TSP 067GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCERINE & 1% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK)  for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc.	100 Plts
TSP 068GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% LECITHIN & 4% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK)  for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc.	100 Plts
TSP 069GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.1% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK)  for cultivation of wide variety of microorganisms	100 Plts
TSP 070G	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ β-LACTAMASE (γ-IRRADIATED)  for cultivation of wide variety of organisms and for inactivation of β-lactam antibiotics	100 Plts
TSP 070GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ β-LACTAMASE (γ-IRRADIATED) (TRIPLE PACK)  for cultivation of wide variety of organisms and for inactivation of β-lactam antibiotics	100 Plts
TSP 071GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ β-LACTAMASE MIXTURE (γ-IRRADIATED) (TRIPLE PACK)  for cultivation of wide variety of organisms and for inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems	100 Plts
TSP 072GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ β-LACTAMASE II (γ-IRRADIATED) (TRIPLE PACK)  for cultivation of wide variety of organisms and for inactivation of cephalosporins of first, second, third and fourth generation	100 Plts
TSP 073GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ TWEEN 80, SOYA LECITHIN & β-LACTAMASE MIXTURE (γ-IRRADIATED) (TRIPLE PACK)  recommended for determining efficiency of containers, equipment, surfaces, water miscible cosmetics and for inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems	100 Plts
TSP 015GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.5% POLYSORBATE 80, 0.07% SOYA LECITHIN (γ-IRRADIATED) (TRIPLE PACK)  for determining the efficiency of sanitization of containers, equipment, surfaces, water-miscible cosmetics, etc.	100 Plts
TSP 074GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN, TWEEN 80 & β-LACTAMASE-II (γ-IRRADIATED) (TRIPLE PACK)  for determining efficiency of containers, equipment, surfaces, water miscible cosmetics and for inactivation of cephalosporins of first, second, third and fourth generation and penems	100 Plts
TSP 074AGT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN, TWEEN 80, 1% GLYCEROL & β-LACTAMASE-II (γ-IRRADIATED) (TRIPLE PACK)  for determining efficiency of containers, equipment surfaces, water miscible cosmetics and for inactivation of cephalosporins of first, second, third and fourth generation and penems	100 Plts
TSP 075GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN, TWEEN 80 & β-LACTAMASE (γ-IRRADIATED) (TRIPLE PACK)  for determining efficiency of containers, equipment, surfaces, water miscible cosmetics and inactivation of β-lactam antibiotics	100 Plts
TSP 075AGT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN, TWEEN 80, 1% GLYCEROL & β-LACTAMASE (γ-IRRADIATED) (TRIPLE PACK)  for determining efficiency of containers, equipment surfaces, water miscible cosmetics	100 Plts
TSP 076GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 1% GLYCEROL, 0.5% POLYSORBATE 80, 0.07% SOYA LECITHIN & 5 IU/PLATE β-LACTAMASE MIXTURE (γ-IRRADIATED) (TRIPLE PACK)  for determining efficiency of containers, equipment, surfaces, water miscible cosmetics and inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems	100 Plts
TSP 077GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ POLYSORBATE 80, GLYCEROL W/ 5 IU OF β-LACTAMASE II & 50 IU OF β-LACTAMASE 1/100 ML (γ-IRRADIATED) (TRIPLE PACK)  for determining efficiency of containers, equipment, surfaces, water miscible cosmetics and inactivation of penicillins, cephalosporins of first, second, third and fourth generation and penems	100 Plts
TSP 081GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ 0.5% GLYCEROL & 1% POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK)   for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc	100 Plts
TSP 078GT	TOTAL PLATE COUNT AGAR PLATE W/ β-LACTAMASE (γ-IRRADIATED) (TRIPLE PACK)  for determining of plate counts of microorganisms in foods, water and wastewater and inactivation of β-lactam antibiotics	100 Plts

100 mm Plates

CODE	PRODUCT NAME	PACK SIZE
TMPC 001	DICHLORAN ROSE BENGAL CHLORAMPHENICOL (DRBC) AGAR PLATE for selective isolation of Fungi-Yeasts and Moulds of significance in food spoilage. The composition and performance criteria are in accordance with ISO 21527-1:2008	50 Plts

Ready-to-Use Culture Media

100 mm Plates

CODE	PRODUCT NAME	PACK SIZE	
TMPC 004	MUELLER HINTON AGAR NO. 2 for determination of susceptibility of microorganisms to antimicrobial agents	50 Plts	
TMPC 002	MYP AGAR PLATE for isolation and identification of pathogenic staphylococci and <i>Bacillus</i> species	50 Plts	
TMPC 013GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN & POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) to determine the efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc	50 Plts	
TMPC 003GT	SOYA BEAN CASEIN DIGEST AGAR PLATE W/ LECITHIN & POLYSORBATE 80 (TRYPTONE SOYA AGAR PLATE W/ LECITHIN & POLYSORBATE 80) (γ-IRRADIATED) (TRIPLE PACK) for determining efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics etc	50 Plts	

Bi-Plates

CODE	PRODUCT NAME	PACK SIZE	
TMPB 001	NUTRIENT-MacCONKEY AGAR PLATE for selective isolation and differentiation of coliform and other enteric pathogens Combination of Nutrient Agar (TMP 006) + MacConkey Agar w/ 0.15% Bile Salts, CV and NaCl (TMP 015)	50 Plts	
TMPB 002	CLED-MacCONKEY AGAR PLATE for isolation and differentiation of urinary pathogens on the basis of lactose fermentation Combination of CLED Agar W/ Bromothymol Blue (TMP 014) + MacConkey Agar w/ 0.15% Bile Salts, CV and NaCl (TMP 015)	50 Plts	
TMPB 003	XLD-MacCONKEY AGAR PLATE for selective isolation and enumeration of <i>Salmonella</i> species and differentiation of enteric pathogens Combination of XLD Agar (TMP 024)+ MacConkey Agar w/ 0.15% Bile Salts, CV and NaCl (TMP 015)	50 Plts	
TMPB 004	CETRIMIDE-MacCONKEY AGAR PLATE for selective isolation of <i>Pseudomonas</i> and differentiation of coliform and other enteric pathogens Combination of Cetrimide Agar (TMP 002)+ MacConkey Agar w/ 0.15% Bile Salts, CV and NaCl (TMP 015)	50 Plts	
TMPB 005	BLOOD-MacCONKEY AGAR PLATE for isolation and cultivation of fastidious organisms and differentiation of coliforms and other enteric pathogens Combination of Blood Agar (TMP 017) + MacConkey Agar w/ 0.15% Bile Salts, CV and NaCl (TMP 015)	50 Plts	
TMPB 006	MacCONKEY-MANNITOL SALT AGAR PLATE for cultivation and differentiation of enteric bacteria, Restricting the swarming of <i>Proteus</i> species along with potentially pathogenic Gram positive organisms especially pathogenic staphylococci Combination of MacConkey Agar w/ 0.15% Bile Salts, CV and NaCl (TMP 015) + Mannitol Salt Agar (TMP 021)	50 Plts	
TMPB 007	BLOOD-CHOCOLATE AGAR PLATE for isolation of <i>Neisseria</i> and other fastidious microorganisms Combination of Sheep Blood Agar (TMP 017) + Chocolate Agar (TMP 020)	50 Plts	
TMPB 008	BLOOD-MANNITOL SALT AGAR PLATE for isolation of <i>Neisseria</i> and other fastidious microorganisms along with potentially pathogenic Gram positive organisms especially pathogenic staphylococci Combination of Sheep Blood Agar (TMP 017) + Mannitol Salt Agar (TMP 021)	50 Plts	
TMPB 009	CHOCOLATE-MacCONKEY AGAR PLATE for the isolation and cultivation of fastidious organisms and differentiation of coliforms and other enteric pathogens Combination of Chocolate Agar (TMP 020) + MacConkey Agar W/O CV, NaCl W/ 0.5% Sodium Taurocholate (TMP 016)	50 Plts	

Glass Bottles

CODE	PRODUCT NAME	PACK SIZE	
TRM 360	BLOOD AGAR BASE (INFUSION AGAR) W/O BLOOD for selective isolation of fastidious pathogenic microorganisms after addition of blood	 100 ml*25	
TRM 361	BRAIN HEART INFUSION AGAR for cultivation of fastidious microorganisms like Bacteria, Yeasts and Moulds	 100 ml*25	
TRM 405	CLED AGAR (W/ BROMO THYMOL BLUE) (BROLACIN AGAR) for isolation and differentiation of urinary pathogen by lactose fermentation	 100 ml*25	
TRM 054	CLED AGAR W/ ANDRADE INDICATOR (CYSTINE LACTOSE ELECTROLYTE DEFICIENT AGAR) for isolation and differentiation of microorganisms based on lactose fermentation	 100 ml*25	
TRM 1199	CHROMOGENIC UTI AGAR * for presumptive identification of microorganisms mainly causing urinary tract infections	 100 ml*25	
TRM 1858	CHROMOGENIC COLIFORM AGAR * for determination of coliform bacteria particularly <i>Enterobacteriaceae</i> on the basis of their ability to ferment lactose	 100 ml*25	
TRM 349	MacCONKEY AGAR W/O CV, W/ 0.15% BILE SALTS for selective isolation and differentiation of lactose fermenting and lactose non-fermenting enteric pathogens	  100 ml*25	

Ready-to-Use Culture Media

Glass Bottles

CODE	PRODUCT NAME		PACK SIZE
TRMH 110	MacCONKEY AGAR W/ 0.15% BILE SALTS, CV AND NACL (USP/EP/BP/JP/IP) for the selection and subculture of <i>E.coli</i> in accordance with harmonized method		100 ml*25
TRM 378	MacCONKEY AGAR W/ SODIUM TAUROCHOLATE W/O CV & NaCl for the selection and recovery of the <i>Enterobacteriaceae</i> and related enteric gram-negative bacilli from clinical, food, water samples		100 ml*25
TRM 339	MUELLER HINTON AGAR for cultivation of <i>Neisseria</i> and for determination of susceptibility of microorganisms		100 ml*25
TRM 206	MANNITOL SALT AGAR for selective isolation of pathogenic staphylococci from clinical and non-clinical samples		100 ml*25
TRMH 114	MANNITOL SALT AGAR (USP/EP/BP/JP/IP) for selective isolation of pathogenic staphylococci from pharmaceutical products in accordance with microbial limit test		100 ml*25
TRM 341	NUTRIENT AGAR general purpose medium for cultivation of less fastidious microorganisms		100 ml*25
TRM 544	PLATE COUNT AGAR (STANDARD PLATE AGAR) for determination of plate count of microorganisms in food, waste water and clinical samples		100 ml*25
TRM 344	POTATO DEXTROSE AGAR for isolation and enumeration of Yeasts and Moulds from clinical, dairy and other food products		100 ml*25
TRMH 344	POTATO DEXTROSE AGAR (USP/EP/BP/JP/IP) for isolation and enumeration of Yeasts and Moulds from clinical, dairy and other food products		100 ml*25
TRM 269	R2A AGAR for heterotrophic plate count of treated potable water using longer incubation time		100 ml*25
TRM 622	SABOURAUD CHLORAMPHENICOL AGAR for selective isolation and cultivation of Yeasts and Moulds		100 ml*25
TRM 387	SABOURAUD DEXTROSE AGAR for cultivation of Yeasts, Moulds and aciduric bacteria from clinical and non-clinical samples		100 ml*25
TRMH 387	SABOURAUD DEXTROSE AGAR (USP/EP/JP/BP/IP) for cultivation of Yeasts, Moulds and aciduric bacteria from pharmaceutical products in accordance with microbial limit testing		100 ml*25 500 ml*8
TRM 002	SABOURAUD DEXTROSE AGAR W/ 0.5% SOYA LECITHIN & 4% TWEEN 80 used for cultivation of Yeasts, Moulds and aciduric bacteria		100 ml*25
TRMH 345	SOYA CASEIN DIGEST AGAR (USP/EP/JP/BP/IP) for cultivation of variety of microorganisms from pharmaceutical products accordance with harmonized method		100 ml*25 500 ml*8
TRM 001	SOYA CASEIN DIGEST AGAR W/ 0.5% SOYA LECITHIN & 4% TWEEN 80 for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc		100 ml*25
TRM 1927	VIOLET RED GLUCOSE AGAR for selective isolation, detection and enumeration of coli-aerogenes bacteria in water, milk other dairy products and clinical samples		100 ml*25
TRMH 117	VIOLET RED BILE GLUCOSE AGAR (USP/EP/JP/BP/IP) for detection and enumeration of <i>Enterobacteriaceae</i> especially subculturing of bile tolerant gram negative bacteria from pharmaceutical products in accordance with microbial limit test		100 ml*25
TRM 492	XLD AGAR for selective isolation and enumeration of <i>Salmonella typhi</i> and other <i>Salmonella</i> species		100 ml*25

Glass Bottles-Wash Solutions

CODE	PRODUCT NAME		PACK SIZE
TMKH 006	BUFFERED NaCl-PEPTONE SOLUTION, pH 7.0 (USP/EP/JP/BP/IP) dilution fluid for samples in case of microbiological contamination		100 ml*25 500 ml*8
TMKH 009	BUFFERED NaCl-PEPTONE SOLUTION, pH 7.0 W/ SOYA LECITHIN (USP/EP/JP/BP/IP) for the preparation of test suspension		100 ml*25
TMKH 004	DILUTING FLUID A (USP) used as diluent for sterility testing of pharma products		100 ml*25 500 ml*8
TMKH 011	DILUTING FLUID D (USP) used as diluent for sterility testing of pharma products		100 ml*25

Ready-to-Use Culture Media

Glass Bottles-Wash Solutions

CODE	PRODUCT NAME		PACK SIZE
TMKH 012	DILUTING FLUID K (USP) used as diluent for sterility testing of pharma products		100 ml*25
TMKH 013	PEPTONE WATER 0.1% used as diluent for microbial enumeration purposes	 	100 ml*25 500 ml*8
TMKH 013S	PEPTONE WATER 0.1% used as diluent for microbial enumeration purposes		100 ml*25
TMKH 027	PHOSPHATE BUFFER pH 7.0 used as a diluent		100 ml*25
TMKH 014	STERILE SALINE 0.85% used as a diluent		100 ml*25 500 ml*8

Glass Bottles-Microbial Limit Test

CODE	PRODUCT NAME		PACK SIZE
TMKH 005	EE BROTH MOSSEL (ENTEROBACTERIA ENRICHMENT BROTH, MOSSEL) (USP/EP/JP/BP/IP) for selective enrichment of <i>Enterobacteriaceae</i>		100 ml*25 500 ml*8
TMKH 008	MacCONKEY BROTH (USP/EP/JP/BP/IP) for the selective enrichment of <i>Escherichia coli</i> in accordance with harmonized method		100 ml*25 500 ml*8
TMKH 015	RAPPAPORT VASSILIADIS SALMONELLA ENRICHMENT BROTH (USP/EP/JP/BP/IP) for selective enrichment of <i>Salmonella</i> species from pharmaceutical and clinical samples in accordance with harmonized method		100 ml*25
TMKH 025	REINFORCED CLOSTRIDIAL BROTH (as per USP/BP/JP/EP) for isolation, cultivation and enumeration of <i>Clostridia</i> species, highly nutritive for <i>Clostridium sporogenes</i> and other anaerobes		100 ml*25
TMKH 007	SABOURAUD DEXTROSE BROTH (USP/EP/JP/BP/IP) for the enrichment of <i>Candida albicans</i> in accordance with harmonized method		100 ml*25 500 ml*8
TMKH 030S	T.A.T. BROTH WITH TWEEN 20 (as per USP) for sterility testing of highly viscous or gelatinous substances in cosmetic and pharma		100 ml*25

Glass Bottles-Miscellaneous Use

CODE	PRODUCT NAME		PACK SIZE
TMK 016	GN BROTH, HAJNA for the selective enrichment of gram negative organisms of the enteric group from clinical & non clinical samples		100 ml*25
TMK 020	GN BROTH MEDIUM 11 (as per IP 2018) for enrichment of <i>Shigella</i> from pharmaceutical products	 	100 ml*25
TMK 017	0.1% PEPTONE WATER W/ 0.1% TWEEN 80 used as a growth medium and for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc.	 	100 ml*25
TMK 018	0.1% PEPTONE WATER W/ 0.5% SOYA LECITHIN AND 4% TWEEN 80 used as a growth medium and for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc		100 ml*25
TMK 019	0.1% PEPTONE WATER W/ 0.5% SOYA LECITHIN AND 4% TWEEN 20 used as a growth medium and for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc.		100 ml*25

Glass Bottles-Sterility Test Media

CODE	PRODUCT NAME		PACK SIZE
TMKH 010	ALTERNATIVE THIOGLYCOLLATE MEDIUM (USP) for sterility testing of biological products		100 ml*25
TMKH 028PS	BUFFERED NaCl-PEPTONE SOLUTION, pH 7.8 (USP/EP/JP/BP/IP) dilution fluid for samples in case of microbiological contamination		200 ml*25

 200 ml Bottle with Screw Cap

 Recommended For Pharmaceuticals,  Clinical,  Food | * Store Between 2-8°C





Blood Culture Bottles

Diagnostic Tool for the Detection of Bacteremia & Fungemia

Significance

- Detect microorganisms in bloodstream infections
- Identify causative agents and infection source
- Support antimicrobial susceptibility testing for therapy

How to use Blood Culture Bottles?

- Label the bottle and remove the cap seal
- Disinfect the exposed rubber stopper
- Collect blood using a sterile needle and syringe
- Inject the sample into the bottle aseptically
- Incubate at 35-37°C for 24-48 hours (extend up to 7 days if required)

Ready-to-Use Culture Media

#Glass Bottles-Sterility Test Media

CODE	PRODUCT NAME		PACK SIZE	
TMKH 001	FLUID THIOGLYCOLLATE MEDIUM (USP/EP/JP/BP/IP) for sterility testing of biological and for cultivation of aerobes, anaerobes and microaerophilic organisms		200 ml*25 500 ml*8	
TMKH 001N	FLUID THIOGLYCOLLATE MEDIUM (USP/EP/JP/BP/IP) for sterility testing of biological and for cultivation of aerobes, anaerobes and microaerophilic organisms (narrow mouth bottles)		100 ml*25	
TMKH 001S	FLUID THIOGLYCOLLATE MEDIUM (USP/EP/JP/BP/IP) for sterility testing of biological and for cultivation of aerobes, anaerobes and microaerophilic organisms		100 ml*25	
TMKH 021	FLUID THIOGLYCOLLATE MEDIUM W/ 0.5% SOYA LECITHIN & 4% TWEEN 80 for sterility testing of biological and for cultivation of aerobes, anaerobes and microaerophilic organisms		100 ml*25	
TMKH 023	SCDM W/ 0.5% SOYA LECITHIN & 1% TWEEN 80 for cultivation of wide variety of microorganisms and for determining efficiency of sanitization of containers, equipment, surfaces, water miscible cosmetics etc		100 ml*25	
TMKH 003N	SOYA CASEIN DIGEST MEDIUM (USP/EP/JP/BP/IP) for the evaluation of sterility in manufacturing processes (narrow mouth bottles)		100 ml*25	
TMKH 003	SOYA CASEIN DIGEST MEDIUM (USP/EP/JP/BP/IP) for the evaluation of sterility in manufacturing processes		100 ml*25 200 ml*25 500 ml*8	
TMKH 029	SOYA CASEIN DIGEST MEDIUM (USP/EP/JP/BP/IP) (DOUBLE STRENGTH) for the evaluation of sterility in manufacturing processes		100 ml*25	
TMKH 003S	SOYA CASEIN DIGEST MEDIUM (USP/EP/JP/BP/IP) for the evaluation of sterility in manufacturing processes		90 ml*25 100 ml*25	
TMKH 022	SOYA CASEIN DIGEST MEDIUM W/ 0.5% SOYA LECITHIN & 4% TWEEN 80 for the evaluation of sterility in manufacturing processes		100 ml*25	
TMKH 026S	SOYA CASEIN DIGEST MEDIUM W/ 0.5% SOYA LECITHIN & 5% POLYSORBATE 80 for the evaluation of sterility in manufacturing processes		100 ml*25	

Blood Culturing System (Ready-to-Use Broth Media) (Aluminium-Capped Bottle)

CODE	PRODUCT NAME		PACK SIZE	
TMK 362	BRAIN HEART INFUSION BROTH for detection of fastidious microorganisms in blood		25 ml*50 50 ml*25	
TMK 362S	BHI SUPPLEMENTED W/ 0.05% SPS for detection of microorganisms associated with blood culture		25 ml*50 50 ml*25	
TMK 303	BILE BROTH BASE for cultivation of <i>Enterobacteriaceae</i> group		25 ml*50 50 ml*25	
TMK 319S	FLUID THIOGLYCOLLATE MEDIUM W/ 0.05% SPS for cultivation of aerobes, anaerobes and microaerophiles		25 ml*50 50 ml*25	
TMK 308S	GLUCOSE BROTH W/ 0.05% SPS (BLOOD CULTURE BOTTLE) for detection of microorganisms in blood		25 ml*50 50 ml*25	
TMK 374	HARTLEY DIGEST BROTH for isolation of various bacteria from blood especially streptococci and <i>Corynebacterium diphtheriae</i>		25 ml*50 50 ml*25	
TMK 374S	HARTLEY DIGEST BROTH W/ 0.05% SPS for isolation of various bacteria from blood especially streptococci and <i>Corynebacterium diphtheriae</i>		25 ml*50 50 ml*25	
TMK 350	NUTRIENT BROTH for general cultivation of microorganisms		25 ml*50 50 ml*25	
TMK 332	TRYPTONE SOYA BROTH (SOYA CASEIN DIGEST MEDIUM) for detection of microorganisms in blood		25 ml*50 50 ml*25	
TMK 332S	TRYPTONE SOYA BROTH W/ 0.05% SPS for detection of microorganisms in blood		25 ml*50 50 ml*25	

25 ml pack size called Paediatric & 50 ml pack size called Adult | # Once open, recommended to consume immediately

Recommended For Pharmaceuticals, Clinical, Food | * Store Between 2-8°C





Ready-to-Use Culture Media Bags

Reducing the Workflow!

Open ➡ Dispense ➡ Use

Experience the convenience and efficiency of **Ready-to-Use Culture Media Bags** for your microbial testing needs

Benefits

- Versatile
- Flexible
- Time-Saving
- Easy-to-Use
- Easy Disposal
- Cost-Effective
- Safe Handling
- Space-Saving
- Energy-Efficient
- Quality Guaranteed
- Comply as per ISO testing
- Resilient Hanging Design
- Contamination-Free System
- 12 months Longer Shelf Life



TMF 001: Buffered Peptone Water (BPW)

For Food & Beverage Industry
(*Salmonella* testing)



TMF 002: Half Fraser Broth Bag

Clinical Testing (Clinical specimens),
Food Industry (*Listeria* testing)



TMF 003: Tryptone Soya Broth/ Soya Casein Digest Medium

For Pharmaceutical Industry
(Media fill activity)

Ready-to-Use Culture Media

#Ready-to-Use Broth Media (Bags)

CODE	PRODUCT NAME		PACK SIZE
TMF 001	BUFFERED PEPTONE WATER for pre-enrichment of injured <i>Salmonella</i> species prior to selective enrichment and isolation	 	225 ml*20 5 Ltr*2
TMF 002	HALF FRASER BROTH for collection and shipment of clinical specimens	 	5 Ltr*2
TMF 003	TRYPTONE SOYA BROTH / SOYA CASEIN DIGEST MEDIUM for sterility testing and cultivation of fastidious and non-fastidious microorganisms	 	5 Ltr*2

Ready-to-Use Liquid Medium (Tubes)

CODE	PRODUCT NAME		PACK SIZE
TMT 033	BRAIN HEART INFUSION BROTH for the cultivation of fastidious and non-fastidious microorganisms, including aerobic and anaerobic bacteria, from clinical and non-clinical specimens	 	5 ml*25 5 ml*50
TMTH 001	BUFFERED SODIUM CHLORIDE PEPTONE SOLUTION pH 7.0 (USP/EP/BP/JP/IP) for the preparation of test suspension in accordance with harmonized methods		9 ml*25 9 ml*50 10 ml*25 10 ml*50
TMTH 002	BUFFERED SODIUM CHLORIDE PEPTONE SOLUTION pH 7.0 W/ 0.05% POLYSORBATE 80 (USP/EP/BP/JP/IP) for the preparation of test suspension		10 ml*25 10 ml*50
TMT 021	DEY-ENGLEY NEUTRALISING BROTH * for neutralising and testing antiseptics and disinfectants	 	10 ml*25 10 ml*50
TMT 015	DILUENT TSB WITH CAP 4 for total aerobic microbial, count of water soluble specimens to be tested in pharmaceuticals, cosmetic etc		10 ml*25 10 ml*50
TMTH 003	ENTEROBACTERIA ENRICHMENT BROTH, MOSSEL (USP/EP/BP/JP/IP) for the enrichment of bile tolerant organisms in accordance with the harmonized methods		10 ml*25 10 ml*50
TMTH 019	FLUID THIOGLYCOLLATE MEDIUM (USP/EP/BP/JP/IP) for sterility testing of biologicals and cultivation of aerobes and microphiles in accordance with harmonized method		5 ml*25 9 ml*25 9 ml*50 10 ml*25 10 ml*50
TMT 042	GLUCOSE BROTH for detection of microorganisms in blood		5 ml*50
TMTH 022	GN BROTH (IP 2018) for enrichment of <i>Shigella</i> from pharmaceutical products in accordance with Indian Pharmacopeia		10 ml*25 10 ml*50
TMT 023	LURIA BERTANI BROTH for cultivation and maintenance of recombinant strains of <i>E.coli</i> and may be used for routine cultivation of not particularly fastidious microorganisms		10 ml*25 10 ml*50
TMTH 024	MacCONKEY BROTH (USP/EP/BP/JP/IP) for the selective enrichment of <i>Escherichia coli</i> in accordance with the harmonized method		5 ml*50 10 ml*50
TMTH 043	MacCONKEY BROTH (USP/EP/BP/JP/IP) (DOUBLE STRENGTH) for selective enrichment of <i>Escherichia coli</i> in accordance with harmonized method		5 ml*50
TMT 040	MOLECULAR TRANSPORT MEDIUM (VIRAL LYSIS TRANSPORT MEDIUM) for the stabilization, transportation and inactivation of infectious viral agents Quantity: 1.5 ml medium in 10 to 15 ml self stand tube		1.5 ml*50

Once open, recommended to consume immediately

 Recommended For Pharmaceuticals,  Clinical,  Food | * Store Between 2-8°C



Ready-to-Use Culture Media

#Ready-to-Use Liquid Medium (Tubes)

CODE	PRODUCT NAME		PACK SIZE
TMT 004	MR-VP MEDIUM to perform Methyl Red & Voges-Proskauer tests for enteric gram negative bacilli		10 ml*25 10 ml*50
TMT 005	PEPTONE WATER general purpose growth medium and used as the base of carbohydrate fermentation media	 	5 ml*50 10 ml*25 10 ml*50
TMTH 006	PHOSPHATE BUFFER pH 7.0 (USP) used as a diluent		9 ml*25 9 ml*50
TMT 007	PHOSPHATE BUFFER pH 7.2 (FDA BAM) used as a diluent		9 ml*25 9 ml*50
TMT 025	R2A BROTH for cultivation and maintenance of heterotropic bacteria from potable water		10 ml*25 10 ml*50
TMT 008	RAPPAPORT VASSILIADIS SALMONELLA ENRICHMENT BROTH for selective enrichment of <i>Salmonella</i> species from clinical samples		10 ml*25 10 ml*50
TMTH 009	RAPPAPORT VASSILIADIS SALMONELLA ENRICHMENT BROTH (IP) for selective enrichment of <i>Salmonella</i> species as per Indian Pharmacopeia		10 ml*25 10 ml*50
TMTH 010	RAPPAPORT VASSILIADIS SALMONELLA ENRICHMENT BROTH (USP/EP/BP/JP/IP) for selective enrichment of <i>Salmonella</i> species with harmonized Pharmacopeia		10 ml*25 10 ml*50
TMT 026	REHYDRATION FLUID recommended as diluting fluid for performing growth promotion test use with longer stability		9 ml*25 9 ml*50
TMTH 027	SABOURAUD DEXTROSE BROTH (USP/EP/BP/JP/IP) for the enrichment of <i>Clostridia</i> from pharmaceutical products in accordance with the microbial limit testing by harmonized methods	 	10 ml*25 10 ml*50
TMTH 011	SOYA BEAN CASEIN DIGEST MEDIUM (USP/EP/BP/JP/IP) sterility test media prepared in accordance with harmonized method		9 ml*25 9 ml*50 10 ml*25 10 ml*50
TMTH 013	SOYA BEAN CASEIN DIGEST MEDIUM W/ LTHTh for determining the efficiency of sanitization of containers, equipments, surfaces, water miscible cosmetics		10 ml*25 10 ml*50
TMT 028	STERILE SALINE 0.85% used as a diluent	  	10 ml*50
TMT 014	STERILE SALINE 0.9% used as a diluent	  	9 ml*50 10 ml*25 10 ml*50
TMT 029	STERILE SALINE 0.9% W/ 0.05% TWEEN 80 used as a diluent		9 ml*50 10 ml*25
TMT 031	STERILE SALINE 0.9% W/ 0.05% SOYA LECITHIN & 0.5% POLYSORBATE 80 used as a diluent		10 ml*25
TMT 016	TETRATHIONATE BROTH for selective isolation of <i>Salmonella</i> from food and other pathological materials	 	10 ml*25 10 ml*50
TMT 032	TRYPTONE BROTH W/ 10% NaCl recommended for enumeration of <i>Staphylococcus aureus</i>		10 ml*25
TMT 017	TRYPTOPHAN MEDIUM for detection of indole production	 	10 ml*25 10 ml*50
TMT 030	TWEEN 80 recommended as a surfactant and emulsifier		10 ml*25
TMT 012	UNIVERSAL TRANSPORT MEDIUM for collection and transport of clinical specimen for recovery of viral agents Quantity: 3 ml medium in 10 to 15 ml self stand tube		3 ml*50

Once open, recommended to consume immediately |  200 ml Bottle with Screw Cap



 Recommended For Pharmaceuticals,  Clinical,  Food | * Store Between 2-8°C

Ready-to-Use Culture Media

#Ready-to-Use Liquid Medium (Tubes)

CODE	PRODUCT NAME		PACK SIZE	
TMT 034	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS for collection and transport of clinical specimen for recovery of viral agents Quantity: 3 ml medium in 10 to 15 ml self stand tube		3 ml*50	
TMT 044	UREA BROTH for identification of bacteria on the basis of urea utilization, especially for the differentiation of <i>Proteus</i> , <i>Salmonella</i> and <i>Shigella</i> species in accordance with IP		10 ml*50	
TMT 018	UREA INDOLE MEDIUM * for differentiation of microorganism especially <i>Enterobacteriaceae</i> on the basis of their ability to hydrolyze urea and indole production	 	10 ml*25 10 ml*50	
TMT 020	VIRAL TRANSPORT MEDIUM for collection and transport of clinical specimen for recovery of viral agents Quantity: 3 ml medium in 10 to 15 ml self stand tube		3 ml*50	
TMT 035	VIRAL TRANSPORT MEDIUM W/ GLASS BEADS for collection and transport of clinical specimen for recovery of viral agents Quantity: 3 ml medium in 10 to 15 ml self stand tube		3 ml*50	
TMT 036	UNIVERSAL TRANSPORT MEDIUM for collection and transport of clinical specimen for recovery of viral agents Quantity: 1.3 ml medium in 5 to 10 ml self stand tube		1.3 ml*50	
TMT 037	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS for collection and transport of clinical specimen for recovery of viral agents Quantity: 1.3 ml medium in 5 to 10 ml self stand tube		1.3 ml*50	
TMT 038	VIRAL TRANSPORT MEDIUM for collection and transport of clinical specimen for recovery of viral agents Quantity: 1.3 ml medium in 5 to 10 ml self stand tube		1.3 ml*50	
TMT 039	VIRAL TRANSPORT MEDIUM W/ GLASS BEADS for collection and transport of clinical specimen for recovery of viral agents Quantity: 1.3 ml medium in 5 to 10 ml self stand tube		1.3 ml*50	

#Ready-to-Use Liquid Medium (Tubes)

This Ready-to-Use Liquid Media Tube is pre-measured, presterilized media ideal to be used with Filtration Funnel and Membrane Filtration Funnel for detection and enumeration of microorganisms in water as well as product in pharmaceutical, food and beverages processing industries

CODE	PRODUCT NAME		PACK SIZE	
TMA 001	MF-ENDO BROTH for total coliforms		20 Tubes 50 Tubes	
TMA 002	M-FC BROTH for fecal organism		20 Tubes 50 Tubes	
TMA 003	M-FC WITH ROSOLIC ACID BROTH for fecal organism		20 Tubes 50 Tubes	
TMA 004	M-TGE BROTH for total bacteria		20 Tubes 50 Tubes	
TMA 005	M-TGE WITH TTC INDICATOR BROTH for total bacteria		20 Tubes 50 Tubes	
TMA 006	TRYPTICASE SOYA BROTH-USP for total bacteria		20 Tubes 50 Tubes	
TMA 007	KF-STREPTOCOCCAL BROTH for total bacteria		20 Tubes 50 Tubes	
TMA 008	PSEUDOMONAS BROTH for <i>Pseudomonas</i> sp		20 Tubes 50 Tubes	
TMA 009	M-GREEN YM BROTH for Yeasts and Moulds		20 Tubes 50 Tubes	
TMA 010	ORANGE SERUM BROTH for <i>Lactobacillus</i> , acid resistant bacteria		20 Tubes 50 Tubes	

Once open, recommended to consume immediately

 Recommended For Pharmaceuticals,  Clinical,  Food | * Store Between 2-8°C



Ready-to-Use Culture Media

#Ready-to-Use Liquid Medium (Tubes)

This Ready-to-Use Liquid Media Tube is pre-measured, presterilized media ideal to be used with Filtration Funnel and Membrane Filtration Funnel for detection and enumeration of microorganisms in water as well as product in pharmaceutical, food and beverages processing industries

CODE	PRODUCT NAME		PACK SIZE
TMA 011	HPC MEDIA WITH TTC INDICATOR for total bacteria		20 Tubes 50 Tubes
TMA 012	R2A BROTH for total bacteria		20 Tubes 50 Tubes

Ready-to-Use Slants (For Tuberculosis Tests)

CODE	PRODUCT NAME		PACK SIZE
TMS 12	BILE ESCULIN AGAR SLANTS (10 Slants) * for isolation & identification of <i>Yersinia enterocolitica</i> from food and animal feeding products		1 Kit
TMS 01	L.J. MEDIUM SLANT (10 Slants) * for cultivation of <i>Mycobacterium tuberculosis</i>		1 Kit
TMS 01L	L.J. MEDIUM SLANT (10 Long Slants) for cultivation of <i>Mycobacterium tuberculosis</i>		1 Kit
TMS 07	TUBERCULOSIS FIRST LINE KIT (7 Slants) * with five antitubular drugs (Isoniazide, Streptomycin, Ethambutol, Rifampicin, Pyrazinamide) + 2 controls		1 Kit
TMS 08	TUBERCULOSIS SECOND LINE KIT (10 Slants) * with seven antitubular drugs (Kanamycin, Amikacin, Ethionamide, D-Cycloserine, Clarithromycin, Ciprofloxacin, p-Amino salicylic acid) + 3 controls		1 Kit
TMS 02	KLIGLER IRON AGAR SLANT (10 Slants) * for differential identification of gram-negative enteric bacilli on the basis of the fermentation of dextrose, lactose and H ₂ S Production	 	1 Kit

Ready-to-Use Slants/Butts (Glass Tube)

CODE	PRODUCT NAME		PACK SIZE
TMS 03	LYSINE IRON AGAR SLANT (10 Slants) * for detection of enteric organism especially <i>Salmonella arizonae</i> , based on their ability to decarboxylate or deaminate lysine and to form hydrogen sulphide (H ₂ S)	 	1 Kit
TMS 10	MIU MEDIUM (10 Butts) * for detection of motility, urease and indole production		1 Kit
TMS 04	MOTILITY INDOLE LYSINE IRON AGAR BUTT SLANT (10 Slants) for identification of members of <i>Enterobacteriaceae</i> on the basis of motility, lysine decarboxylase, lysine deaminase and indole production		1 Kit
TMS 11	SOYA CASEIN DIGEST AGAR SLANT (10 Slants) * general purpose medium used for enrichment, isolation of fastidious microorganisms and also for sterility testing	 	1 Kit
TMS 05	SIMMON'S CITRATE AGAR SLANT (10 Slants) * to differentiate gram-negative bacteria on the basis of citrate utilization	 	1 Kit
TMS 06	TRIPLE SUGAR IRON AGAR SLANT (10 Slants) * to differentiate gram-negative enteric bacilli based on carbohydrate fermentation	 	1 Kit
TMS 09	UREA AGAR SLANT (10 Slants) * to differentiate <i>Enterobacteriaceae</i> on the basis of their ability to produce urease	 	1 Kit

Ready-to-Use Kits

CODE	PRODUCT NAME		PACK SIZE
TMK 021	CHROMOGENIC ECC BROTH KIT for rapid detection and confirmation of presence/absence <i>Escherichia coli</i> and coliforms from water sample		1 Kit
TMK 07	H₂S TEST POWDER KIT (10 TEST) (1 Kit Contains: 10 Bottles) for detection of <i>Salmonella</i> , <i>Citrobacter</i> and <i>E. coli</i> from water		1 Kit

Once open, recommended to consume immediately



 Recommended For Pharmaceuticals,  Clinical,  Food | * Store Between 2-8°C



Combo Sterility Kit

For Sterility Testing

Exceptionally made to meet the strict requirements of pharmacopeia like **USP/EP/JP/IP**, that ensure the highest standards of quality and accuracy in sterility testing.

Offer unmatched convenience and versatility, which cater to various testing methods employed in the pharmaceutical industry.

Variants



TMCK 001 (TMCKH 001/001N + TMCKH 003/003N)

- Designed for closed system methods
- Features an aluminium- rubber sealed glass bottle
- Ensures a controlled and sterile testing environment



TMCK 001S (TMCKH 001S + TMCKH 003S)

- Tailored for open system methods
- Comes with a screw cap for direct inoculations
- Simplifies the testing process while maintaining sterility

Applications



Raw Material
Validation



Finished Product
Quality Control



Testing Other
Articles

Ready-to-Use Culture Media

Ready-to-Use Kits

CODE	PRODUCT NAME		PACK SIZE	
TMK 01	WATER TESTING KIT (COLI-CHECK) for detection of coliforms in potable water		1 Kit 10 Kit	
TMK 02	PA E. COLI KIT (100 ml Capacity & Dehydrated Culture Medium-3X Concentration) for detection of presence and absence of coliform bacteria in water		10 Kit	
TMK 03	RAPID COLIFORM KIT * for rapid detection and confirmation of <i>Escherichia coli</i> and other coliforms from water on the basis of enzyme substrate reaction		1 Kit	
TMK 04	RAPID ENTEROCOCCI TEST KIT for rapid and easy identification & differentiation of enterococci from water		1 Kit	
TMK 05	SALT TESTING KIT iodine test solution		10 Kit	
TMK 06	H₂S TEST STRIP KIT (10 TEST) (1 Kit Contains: 10 Bottles) for detection of <i>Salmonella</i> , <i>Citrobacter</i> and <i>E. coli</i> from water		1 Kit	
TMK 08L	FOOD PATHOGEN DETECTION KIT for rapid detection of food pathogens such as <i>E. coli</i> , <i>E. coli</i> O157: H7, <i>Salmonella</i> , <i>Listeria</i> , <i>Staphylococcus</i> and <i>Clostridium</i> species		5 Kit 10 Kit	
TMK 08S	FOOD PATHOGEN DETECTION KIT for rapid detection of food pathogens such as <i>E. coli</i> , <i>E. coli</i> O157: H7, <i>Salmonella</i> , <i>Listeria</i> , <i>Staphylococcus</i> and <i>Clostridium</i> species		5 Kit 10 Kit	
TMK 08M	FOOD PATHOGEN DETECTION KIT for rapid detection of food pathogens such as <i>E. coli</i> , <i>E. coli</i> O157: H7, <i>Salmonella</i> , <i>Listeria</i> , <i>Staphylococcus</i> and <i>Clostridium</i> species		5 Kit 10 Kit	
TMK 09	H₂S TEST MEDIUM for the detection of <i>Salmonella</i> species and <i>Citrobacter</i> species from water samples		1 Kit	
TMK 10	SULPHATE REDUCING BACTERIA TEST KIT (1 Kit Contains: 10 Bottles) for the detection of Sulphate Reducing Bacteria for testing 20 ml water sample		1 Kit	
TMK 11	SULPHATE REDUCING BACTERIA TEST KIT (1 Kit Contains: 10 Bottles) for the detection of Sulphate Reducing Bacteria for testing 20 ml water sample		1 Kit	
TMK 12	WATER TEST KIT for rapid and simultaneous detection of <i>Salmonella</i> species, <i>E. coli</i> , <i>Citrobacter</i> species and <i>Vibrio</i> species		1 Kit	
TMK 14	WATER TEST KIT (1 Kit Contains: 10 Bottles) for primary detection of <i>Salmonella</i> , <i>Citrobacter</i> and <i>E. coli</i> based on H ₂ S production in plastic bottles		1 Kit	
TMK 15	SELECTIVE E. COLI TEST KIT for the rapid detection and confirmation of <i>E. coli</i> based on enzyme-substrate reaction from water samples		1 Kit	

Combo Sterility Kits

CODE	PRODUCT NAME		PACK SIZE	
TMCK 001	COMBO STERILITY KIT (USP/EP/JP/IP) each kit consists of 2 glass bottles (100 ml) of Fluid Thioglycollate Medium (TMKH 001/001N) and Soya Casein Digest Medium (TMKH 003/003N). Recommended for all purpose		5 Kit 20 Kit	
TMCK 001S	COMBO STERILITY KIT (USP/EP/JP/IP) each kit consists of 2 glass bottles (100 ml) of Fluid Thioglycollate Medium (TMKH 001S) and Soya Casein Digest Medium (TMKH 003s). Recommended for all purpose		5 Kit 20 Kit	

Once open, recommended to consume immediately



 Recommended For Pharmaceuticals,  Clinical,  Food | Store Between 2 -8°C *



Ready-to-Use Transport Swabs with Medium

Rest Easy, Samples are Protected

Capture & Preserve:

Allows the capture and preservation of microbial specimens with the utmost efficiency and reliability. Each swab is carefully designed to maximize sample collection and retention, ensuring accurate and robust downstream analysis.

Instant Activation:

Forget about the hassle of rehydration or additional steps. Ready-to-Use Transport Medium is instantly activated upon contact with the sample, providing immediate preservation and stabilization of the specimen. Save time and simplify your workflow!

Optimal Viability:

Formulated with proprietary ingredients that enhance the survival of a wide range of microorganisms during transportation, it guarantees reliable results.

Long Shelf-life:

With an extended shelf life, the medium remains stable, maintaining its integrity and reliability for longer periods.

Versatile:

Designed to accommodate diverse microbiological specimens, including bacteria, viruses, or fungi. It's the perfect solution for clinical, veterinary, environmental, and industrial applications.

Secure & Leak-proof:

Our innovative packaging ensures a secure and leak-proof transportation experience. The tubes are equipped with a tight-sealing cap, preventing leakage and cross-contamination.

Convenient:

Each tube is individually packaged, ensuring sterility and easy-to-use handling. The sealed packaging eliminates the risk of contamination and provides a sterile environment for the specimens.

Ready-to-Use Culture Media

Transport Swab with Medium (Disposable)

CODE	PRODUCT NAME		PACK SIZE
THTS 001	STERILE SWABS W/ 0.85% SALINE for propagation of pathogenic cocci and other fastidious associated with blood culture work and allied pathological investigations		10 Nos
THTS 001A	STERILE SWABS W/ 0.90% SALINE (COTTON SWAB) for propagation of pathogenic cocci and other fastidious associated with blood culture work and allied pathological investigations		10 Nos
THTS 011	STERILE SWABS W/ 0.90% SALINE (NYLON SWAB) for propagation of pathogenic cocci and other fastidious associated with blood culture work and allied pathological investigations		10 Nos
THTS 014GT	STERILE NYLON SWAB W/ NORMAL SALINE TUBE & STERILE SCDM, SOYA LECITHIN AND POLYSORBATE 80 TUBE (γ-IRRADIATED) (TRIPLE PACK) for isolating microbial specimens in presence of antiseptics and disinfectants		50 Nos
THTS 015	STERILE SWABS W/ 0.85% SALINE (NYLON SWAB) recommended for propagation of pathogenic cocci and other fastidious associated with blood culture work and allied pathological investigations		10 Nos
THTS 016	TRANSPORT SWABS W/ ALTERNATIVE THIOGLYCOLLATE MEDIUM (NIH THIOGLYCOLLATE MEDIUM) for transportation of aerobes, anaerobes and microaerophiles		10 Nos 50 Nos
THTS 415	TRANSPORT SWABS W/ CARY-BLAIR MEDIUM BASE (TRANSPORT MEDIUM W/O CHARCOAL) for holding clinical specimens during collection, like, <i>Enterobacter aerogenes</i> , <i>Klebsiella pneumoniae</i> , <i>Neisseria meningitidis</i> , <i>Salmonella typhimurium</i> , <i>Vibrio cholera</i> , etc.		10 Nos
THTS 715	TRANSPORT SWABS W/ DEY-ENGLEY NEUTRALIZING AGAR (D/E AGAR DISINFECTANT TESTING) for transportation of microorganisms from sanitized surface		10 Nos
THTS 341	TRANSPORT SWABS W/ NUTRIENT AGAR for enumeration of <i>Enterobacteriaceae</i>		10 Nos
THTS 420	TRANSPORT SWABS W/ STUART TRANSPORT MEDIUM for preservation and transport of <i>Neisseria</i> species and other fastidious organisms, like, <i>Haemophilus influenzae</i> , <i>Neisseria gonorrhoeae</i> , <i>Shigella flexneri</i> , <i>Streptococcus pneumoniae</i> , etc		10 Nos
THTS 456	TRANSPORT SWABS W/ TRANSPORT MEDIUM AMIES W/O CHARCOAL for transportation and preservation of clinical samples like, <i>Neisseria gonorrhoeae</i> , <i>Shigella flexneri</i> , <i>Streptococcus pneumoniae</i> , <i>Salmonella typhi</i> , etc		10 Nos
THTS 009	TRANSPORT SWABS W/ TRANSPORT MEDIUM AMIES W/ CHARCOAL for transportation and preservation of microbiological samples		10 Nos
THTS 010	TRANSPORT SWABS W/ BUFFERED PEPTONE SODIUM CHLORIDE 7.0 pH use for transport of specimen		10 Nos
THTS 002	TRANSPORT SWABS W/ AMIES MEDIUM (A) with 1.0 ml medium and one swab recommended for collection and transport of aerobic, anaerobic and fastidious organisms from wound, skin and throat		50 Nos
THTS 003	TRANSPORT SWABS W/ AMIES MEDIUM (B) with 1.0 ml medium and two swabs recommended for collection and transport of aerobic, anaerobic and fastidious organisms from wound, skin and throat and MRSA Screening		50 Nos
THTS 004	TRANSPORT SWABS W/ AMIES MEDIUM (C) with 1.0 ml medium and three swabs recommended for collection and transport of aerobic, anaerobic and fastidious organisms from wound, skin and throat and MRSA Screening and multiple body sides		50 Nos
THTS 005	TRANSPORT SWABS W/ AMIES MEDIUM (D) with 1.0 ml medium and one swab recommended for collection and transport of aerobic, anaerobic and fastidious organisms from nasopharyngeal, paediatric and urogenital		50 Nos

Ready-to-Use Culture Media

Transport Swab with Medium (Disposable)

CODE	PRODUCT NAME		PACK SIZE	
THTS 006	TRANSPORT SWAB W/ SOYA BEAN CASEIN DIGEST MEDIUM W/ 6.5% NaCl with 2.0 ml TSB medium and one swab recommended for collection & transport of aerobic, anaerobic and fastidious organisms from nose, throat, axilla, perneum, groin for MRSA Screening		50 Nos	
THTS 007	TRANSPORT SWAB W/ SOYA BEAN CASEIN DIGEST MEDIUM W/ 6.5% NaCl with 2.0 ml TSB medium and one swab recommended for collection & transport of aerobic, anaerobic and fastidious organisms from rectum		50 Nos	
THTS 008	TRANSPORT SWAB W/ SELENITE MEDIUM (A) With 2.0 ml medium recommended for enrichment of enteric organisms from fecal specimens		50 Nos	
THTS 012	TRANSPORT SWAB W/ CARY BLAIR MEDIUM for recovery of aerobic, anaerobic and fastidious bacteria from faecal specimen		50 Nos	
TMT 041GT	TRANSPORT SWABS W/ DEY-ENGLEY NEUTRALIZING BROTH (γ-IRRADIATED) (TRIPLE PACK) for transporting microbial specimens in presence of antiseptics and disinfectants		50 Nos	



Surface Sampling Swabs

Collect, Preserve, Test

Simplify the Precision Microbial Collection

Superior Sampling Precision

The Swabs are designed to collect the most accurate and representative samples, ensuring effective environmental monitoring and regulatory compliance.

Contamination Control

Safeguard your production environment against potential contaminants. The Swabs are manufactured under stringent quality standards to minimize the risk of cross-contamination.

User-friendly Design

Ergonomically designed Swab handles and easy-to-use instructions make sampling a breeze, reducing the chances of errors and ensuring consistent results.

Enhanced Detection

With these quality Swabs, you can detect even trace amounts of contaminants, allowing for swift corrective actions and preventing potential product recalls.

Versatility

Whether you need to sample flat surfaces, intricate equipment, or hard-to-reach corners, these Swabs are engineered for adaptability.

Global Compliance

Comply with international regulatory standards effortlessly. These Swabs are designed to meet the rigorous requirements of the pharmaceutical industry.

Ready-to-Use Culture Media

Surface Sampling Swabs

CODE	PRODUCT NAME	PACK SIZE	
TMST 001T	STERILE SWABS W/ 0.85% SALINE (TRIPLE PACK) for collection and transportation of microorganisms from surface	50 Nos	
TMST 001GT	STERILE SWABS W/ 0.85% SALINE (γ-IRRADIATED) (TRIPLE PACK) for collection and transportation of microorganisms from surface	50 Nos	
TMST 002T	STERILE SWABS W/ 0.9% SALINE (TRIPLE PACK) for collection and transportation of microorganisms from surface	50 Nos	
TMST 002GT	STERILE SWABS W/ 0.9% SALINE (γ-IRRADIATED) (TRIPLE PACK) for collection and transportation of microorganisms from surface	50 Nos	
TMST 003GT	STERILE SWABS W/ 0.9% SALINE & SECOND TUBE CONTAINING SCDM W/ BCP (γ-IRRADIATED) (TRIPLE PACK) for collection and transportation of microorganisms from surface	50 Nos	
TMST 004GT	PREWETTED SWAB IN NORMAL SALINE AND SECOND TUBE CONTAINING SCDM W/ SOYA LECITHIN, POLYSORBATE 80 (γ-IRRADIATED) (TRIPLE PACK) for collection and isolation of microbial specimens in presence of antiseptics and disinfectants	50 Nos	
TMST 005GT	STERILE SWABS W/ 0.9% SALINE, 0.5% POLYSORBATE 80, 0.05% SOYA LECITHIN & W/ STERILE 60 IU / TUBE β-LACTAMASE-II (γ-IRRADIATED) (TRIPLE PACK) for transportation of microbial specimens by neutralization and inactivation of first, second, third and fourth generation cephalosporins	50 Nos	
TMST 006GT	STERILE SWAB W/ TRYPTONE BROTH AND 10% NaCl (γ-IRRADIATED) (TRIPLE PACK) for the selective transportation of <i>Staphylococcus</i>	50 Nos	
TMST 007GT	TRANSPORT SWAB FOR ENVIRONMENTAL SAMPLING (γ-IRRADIATED) (TRIPLE PACK) for detection of microbial contamination on various surfaces	50 Nos	
TMST 008GT	TRANSPORT SWAB W/ SCDM W/ LECITHIN, POLYSORBATE 80 & NORMAL SALINE TUBE (γ-IRRADIATED) (TRIPLE PACK) for transporting microbial specimen	50 Nos	
TMST 009GT	TRANSPORT SWAB W/ SOYA BEAN CASEIN DIGEST MEDIUM (γ-IRRADIATED) (TRIPLE PACK) for transportation of microorganisms from surface	50 Nos	
TMST 010GT	TRANSPORT SWAB W/ DEY-ENGLEY NEUTRALIZING BROTH (γ-IRRADIATED) (TRIPLE PACK) for transportation of microorganisms from sanitized surface	50 Nos	
TMST 011GT	TRANSPORT SWABS W/ SOYA LECITHIN (0.07%W/V) & TWEEN 80 (0.5%V/V) (γ-IRRADIATED) (TRIPLE PACK) for transporting microbial specimen	50 Nos	
TMST 012	STERILE NEUTRALIZER SWABS W/ 25000 IU / LITRE β-LACTAMASE, 0.9% SALINE W/ 0.5% TWEEN 80 & 0.07% SOYA LECITHIN for isolating microbial specimens in presence of antiseptics and disinfectants	50 Nos	
TMST 013	STERILE VISCOSWABS W/ STERILE 0.9% SALINE W/ 0.5% POLYSORBATE 80 & 0.05% SOYA LECITHIN IN TUBES for used as rinse solution for surface, equipment sampling in clean rooms and isolators	50 Nos	
TMST 014GD	STERILE SWABS W/ RINSE SOLUTION (γ-IRRADIATED) (DOUBLE PACK) for propagation of pathogenic cocci and other fastidious associated with blood culture work and allied pathological investigations	50 Nos	
TMDS 001	TM DUO-SNAP SWAB (SCDM) for environmental monitoring and sterility testing in dry and hard to reach areas and isolators, as well as for evaluating sanitization efficiency using a single-tube swab system	 50 Nos	
TMDS 002	TM DUO-SNAP SWAB (SCDM W/ TWEEN 80 & LECITHIN) for monitoring cleanrooms and hard to reach areas where the neutralization of antiseptics and disinfectants is essential for assessing its bactericidal effectiveness	 50 Nos	
TMDS 003	TM DUO-SNAP SWAB (SCDM W/ 1% TWEEN 80 & 0.01% LECITHIN) for monitoring cleanrooms and hard to reach areas where the neutralization of antiseptics and disinfectants is essential for assessing its bactericidal effectiveness	 50 Nos	
TMDS 004	TM DUO-SNAP SWAB (W/ BUFFERED PEPTONE WATER) for microbiological sample collection from environmental surfaces, isolators, food processing equipment and other hard to reach areas using a single tube swab system	 50 Nos	
TMDS 005	TM DUO-SNAP SWAB W/ DEY ENGLEY NEUTRALIZING BROTH W/ BCP for monitoring cleanrooms and hard to reach areas where the neutralization of antiseptics and disinfectants is essential for assessing its bactericidal effectiveness	 50 Nos	
TMDS 006	TM DUO-SNAP SWAB W/ FLUID THIOGLYCOLLATE MEDIUM for environmental monitoring and sterility testing in dry and hard to reach areas and isolators, as well as for evaluating sanitization efficiency using a single-tube swab system	 50 Nos	

TM Duo-Snap Swabs


Simplifying Environmental Monitoring

TM Duo-Snap Swabs is an all-in-one system for efficient surface sampling—no pipettes, no broth prep!

Key Features

- Pre-filled, sterile and gamma-irradiated
- Snap-valve for contamination-free broth release
- Effective on wet and dry surfaces
- Consistent, reliable sampling results

How to Use

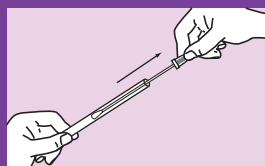


1. Open the tube and remove the swab




2. Swab the target area to collect the sample



3. Return the swab to the tube

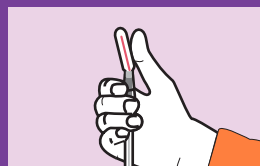

4. Bend the snap valve at a 45° angle until a snapping sound is heard

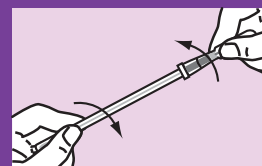

5. Gently squeeze the bulb to release the medium into the tube, then incubate the swab directly & observe for turbidity after the incubation period

Ready-to-Use Culture Media

🔍 Surface Sampling Swabs

CODE	PRODUCT NAME	PACK SIZE
TMDS 007	TM DUO-SNAP SWAB W/ LETHEEN BROTH for microbiological sample collection from food processing equipment and environmental surfaces, where neutralization of disinfectant residues with antimicrobial properties is required	NEW 50 Nos
TMDS 008	TM DUO-SNAP SWAB W/ SCDM W/ 2% GLYCEROL & 0.5% POLYSORBATE 80 for environmental monitoring and sterility testing in dry and hard to reach areas, as well as for evaluating sanitization efficiency using a single-tube swab system	NEW 50 Nos
TMDS 009	TM DUO-SNAP SWAB W/ SCDM W/ 2% GLYCEROL & 0.5% POLYSORBATE 80 W/ 50IU OF BL-I & BL-II for environmental monitoring and sterility testing of dry and hard to reach areas, and inactivation of first, second, third, and fourth generation penicillins, cephalosporins, and penems using a single tube swab system	NEW 50 Nos
TMDS 010	TM DUO-SNAP SWAB W/ SCDB W/ 1% TWEEN 80 & 0.1% LECITHIN for monitoring cleanrooms and hard to reach areas where the neutralization of antiseptics and disinfectants is essential for assessing its bactericidal effectiveness	NEW 50 Nos
TMDS 011	TM DUO-SNAP SWAB W/ SCDM W/ INDICATOR for environmental monitoring and sterility testing in dry and hard to reach areas and isolators, as well as for evaluating sanitization efficiency using a single tube swab system	NEW 50 Nos
TMDS 012	TM DUO-SNAP SWAB W/ SCDM W/ 0.5% TWEEN 80 & 0.05% SOYA LECITHIN for environmental monitoring and sterility testing in dry and hard to reach areas, as well as for evaluating sanitization efficiency using a single-tube swab system	NEW 50 Nos

Ready-to-Use Culture Media

Viral Transport Kits

CODE	PRODUCT NAME	PACK SIZE	
TMVT 001	VIRAL TRANSPORT KIT (SINGLE NYLON FLOCKED SWAB) w/ 3 ml Viral Transport Medium in 10-15 ml polypropylene tube and one sterile flocked nylon swab with breakpoint	50 Nos	
TMVT 002	VIRAL TRANSPORT KIT (WITH NASAL AND THROAT SWAB) w/ 3 ml Viral Transport Medium in 10-15 ml polypropylene tube, one Nasopharyngeal nylon flocked swab with breakpoint and Oropharyngeal viscose throat swab without breakpoint	50 Nos	
TMVT 003	VIRAL TRANSPORT KIT (DOUBLE NYLON FLOCKED SWAB) w/ 3 ml Viral Transport Medium in 10-15 ml polypropylene tube with one Nasopharyngeal and one Oropharyngeal nylon flocked swab with breakpoints	50 Nos	
TMVT 009	VIRAL TRANSPORT KIT (SINGLE VISCOSSE SWAB) w/ 3 ml Viral Transport Medium in 10-15 ml polypropylene tube, one sterile viscose swab without breakpoint	50 Nos	
TMVT 010	VIRAL TRANSPORT KIT (NYLON FLOCKED AND POLYESTER SWAB) w/ 3 ml Viral Transport Medium in 10-15 ml polypropylene tube, one Nasopharyngeal nylon flocked swab and one polyester swab with breakpoint	50 Nos	
TMVT 011	VIRAL TRANSPORT KIT (SINGLE POLYESTER SWAB) w/ 3 ml Viral Transport Medium in 10-15 ml polypropylene tube and one sterile polyester swab with breakpoint	50 Nos	
TMVT 006	VIRAL TRANSPORT KIT (NYLON FLOCKED SWAB AND VISCOSSE SWAB) w/ 1.3 ml Viral Transport Medium in 5-10 ml polypropylene tube with one sterile nylon flocked swab with breakpoint and one sterile viscose swab without breakpoint	50 Nos	
TMVT 007	VIRAL TRANSPORT KIT (DOUBLE NYLON FLOCKED SWAB) w/ 1.3 ml Viral Transport Medium in 5-10 ml polystyrene tube with one sterile Nasopharyngeal and one sterile Oropharyngeal nylon flocked swab with breakpoint	50 Nos	
TMVT 008	VIRAL TRANSPORT KIT (SINGLE NYLON FLOCKED SWAB) w/ 1.3 ml Viral Transport Medium in 5-10 ml polypropylene tube with one sterile nylon flocked swab with breakpoint	50 Nos	

Universal Transport Kits

CODE	PRODUCT NAME	PACK SIZE	
TMUT 001	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (SINGLE NYLON FLOCKED SWAB) w/ 3ml Universal Transport Medium with three glass beads in 10-15 ml polypropylene tube, one Nasopharyngeal nylon flocked swab with breakpoint and Oropharyngeal viscose throat swab without breakpoint	50 Nos	
TMUT 002	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (NASAL AND THROAT SWAB) w/ 3ml Universal Transport Medium with three glass beads in 10-15 ml polypropylene tube, one Nasopharyngeal nylon flocked swab with breakpoint and Oropharyngeal viscose throat swab without breakpoint	50 Nos	
TMUT 003	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (DOUBLE NYLON FLOCKED SWAB) w/ 3ml Universal Transport Medium with three glass beads in 10-15 ml polypropylene tube with one Nasopharyngeal and one Oropharyngeal nylon flocked swab with breakpoints	50 Nos	
TMUT 009	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (SINGLE VISCOSSE SWAB) w/ 3ml Universal Transport Medium with three glass beads in 10-15 ml polypropylene tube, one sterile viscose swab without breakpoint	50 Nos	
TMUT 010	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (NYLON FLOCKED AND POLYESTER SWAB) w/ 3ml Universal Transport Medium with three glass beads in 10-15 ml polypropylene tube, one Nasopharyngeal nylon flocked swab and one polyester swab with breakpoint	50 Nos	
TMUT 011	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (SINGLE POLYESTER SWAB) w/ 3ml Universal Transport Medium with three glass beads in 10-15 ml polypropylene tube and one sterile polyester swab with breakpoint	50 Nos	
TMUT 006	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (NYLON FLOCKED SWAB) w/ 1.3 ml Universal Transport Medium with three glass beads in 5-10 ml polypropylene tube with one sterile nylon flocked swab with breakpoint and one sterile viscose swab without breakpoint	50 Nos	
TMUT 007	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (DOUBLE NYLON FLOCKED SWAB) w/ 1.3 ml Universal Transport Medium with three glass beads in 5-10 ml polypropylene tube with one sterile Nasopharyngeal and one sterile oropharyngeal nylon flocked swab with breakpoint	50 Nos	
TMUT 008	UNIVERSAL TRANSPORT MEDIUM W/ GLASS BEADS (SINGLE NYLON FLOCKED SWAB) w/ 1.3 ml Universal Transport Medium with three glass beads in 5-10 ml polypropylene tube with one sterile nylon flocked swab	50 Nos	

CODE	PRODUCT NAME	PACK SIZE	
MTM 001	MOLECULAR TRANSPORT KIT (SINGLE NYLON FLOCKED SWAB) w/ 3 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polypropylene tube and one sterile nylon flocked swab with breakpoint	50 Nos	
MTM 002	MOLECULAR TRANSPORT KIT (NASAL AND THROAT SWAB) w/ 3 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polypropylene tube, one Nasopharyngeal nylon flocked swab with breakpoint & one Oropharyngeal viscose throat swab without breakpoint	50 Nos	
MTM 003	MOLECULAR TRANSPORT KIT (DOUBLE NYLON FLOCKED SWAB) w/ 3 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polystyrene tube with one Nasopharyngeal and one Oropharyngeal nylon flocked swab with breakpoints	50 Nos	
MTM 009	MOLECULAR TRANSPORT KIT (SINGLE VISCOSE SWAB) w/ 1.5 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polypropylene tube, one sterile viscose swab without breakpoint	50 Nos	
MTM 010	MOLECULAR TRANSPORT KIT (NYLON FLOCKED AND POLYESTER SWAB) w/ 1.5 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polypropylene tube, one Nasopharyngeal nylon flocked swab and one polyester swab with breakpoint	50 Nos	
MTM 011	MOLECULAR TRANSPORT KIT (SINGLE POLYESTER SWAB) w/ 1.5 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polypropylene tube and one sterile polyester swab with breakpoint	50 Nos	
MTM 006	MOLECULAR TRANSPORT KIT (NASAL AND THROAT SWAB) w/ 1.5 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polypropylene tube, one Nasopharyngeal nylon flocked swab with breakpoint & one Oropharyngeal viscose throat swab without breakpoint	50 Nos	
MTM 007	MOLECULAR TRANSPORT KIT (DOUBLE NYLON FLOCKED SWAB) w/ 1.5 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polypropylene tube with one Nasopharyngeal and one Oropharyngeal nylon flocked swab with breakpoints	50 Nos	
MTM 008	MOLECULAR TRANSPORT KIT (SINGLE NYLON FLOCKED SWAB) w/ 1.5 ml Molecular Transport Medium (Viral Lysis Transport Medium) in 10-15 ml polypropylene tube and one sterile nylon flocked swab with breakpoint	50 Nos	



Antibiotic Sensitivity Discs

CE Marking Certification

The assessment of antimicrobial susceptibility in bacterial and fungal isolates is a vital and widespread procedure for determining the most effective antimicrobial agent(s) against infectious organisms.

Antibiotic Sensitivity Discs are an indispensable tool for conducting Antimicrobial Susceptibility Testing. These discs are carefully crafted by impregnating specially selected filter paper with antimicrobial solutions that adhere to the rigorous standards set by WHO, FDA, and CLSI. The impregnation process ensures accurate and uniform application of the solution across the paper.

Applications

- Available in Blister Sets, Cartridges, and Vials
- Comply with International standards, including guidelines from the CLSI
- Enables precise determination of susceptibility patterns
- Comprehensive range of Discs

To know about "TM OCTAClick DISPENSER,"
Please refer to page number **249** for details.

Antibiotic Sensitivity Discs*

Concentration of Antibiotics as per CLSI (Formerly known as NCCLS) Standards

CODE	PRODUCT NAME	ABBR.	CON.	PACK SIZE	
TBD 001	AMIKACIN	AK	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 003	AMOXICLAV (AMOXICILLIN/ CLAVULANIC ACID)	AMC	30 mcg (20/10 mcg)	5 ct 1 B. set 1 vl 5 vl	
TBD 004	AMPICILLIN	AMP	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 005	AMPICILLIN/ SULBACTAM	A/S	20 mcg (10/10 mcg)	5 ct 1 B. set 1 vl 5 vl	
TBD 006	AZITHROMYCIN	AT	15 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 046	AZTREONAM	AZT	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 047	CARBENICILLIN	CB	100 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 007	CEFACTOR	CF	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 009	CEFAZOLIN	CZ	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 048	CEFDINIR	CDR	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 049	CEFEPIME	CPM	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 010	CEFIXIME	CFM	5 mcg	5 ct 1 B. set 1 vl 5 vl	

CODE	PRODUCT NAME	ABBR.	CON.	PACK SIZE	
TBD 011	CEFOPERAZONE	CS	75 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 012	CEFOTAXIME (CEPHOTAXIME)	CTX	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 075	CEFOXITIN	CXN	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 050	CEFPODOXIME	CPD	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 051	CEFPROZIL	CPR	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 014	CEFTAZIDIME	CAZ	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 052	CEFTAZIDIME/ CLAVULANIC ACID	CAC	30/10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 053	CEFTIZOXIME	CZX	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 013	CEFTRIAXONE	CTR	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 015	CEFUROXIME	CXM	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 016	CEPHALOTHIN	CEP	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 017	CHLORAMPHENICOL	C	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 018	CIPROFLOXACIN	CIP	5 mcg	5 ct 1 B. set 1 vl 5 vl	

Antibiotic Sensitivity Discs*

Concentration of Antibiotics as per CLSI (Formerly known as NCCLS) Standards

CODE	PRODUCT NAME	ABBR.	CON.	PACK SIZE	
TBD 019	CLARITHROMYCIN	CL	15 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 020	CLINDAMYCIN	CD	2 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 054	COLISTIN (METHANE SULPHONATE)	CLM	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 022	CO-TRIMOXAZOLE (SULPHAMETH- OXAZOLE/ TRIMETHOPRIM)	COT	25 mcg (23.75/ 1.25 mcg)	5 ct 1 B. set 1 vl 5 vl	
TBD 023	DOXYCYCLINE HCl	DO	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 055	ERTAPENEM	ETP	10 mcg	1 vl 5 vl	
TBD 024	ERYTHROMYCIN	E	15 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 056	FAROPENEM	FAR	5 mcg	1 vl 5 vl	
TBD 057	FLUCONAZOLE (ANTIFUNGAL)	FLC	25 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 076	FOSFOMYCIN 	FO	200 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 058	GATIFLOXACIN	GAT	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 059	GEMIFLOXACIN	GEM	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 025	GENTAMICIN	GEN	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 060	GENTAMICIN	HLG	120 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 061	IMIPENEM	IPM	10 mcg	1 vl 5 vl	

CODE	PRODUCT NAME	ABBR.	CON.	PACK SIZE	
TBD 026	KANAMYCIN	K	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 062	LEVOFLOXACIN	LE	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 063	LINEZOLID	LZ	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 027	LOMEFLOXACIN	LOM	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 028	MEROPENEM	MR	10 mcg	1 vl 5 vl	
TBD 064	MINOCYCLINE	MI	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 065	MOXIFLOXACIN	MO	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 029	NALIDIXIC ACID	NA	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 077	NEOMYCIN 	N	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 078	NEOMYCIN 	N	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 030	NETILLIN (NETIMICIN SULPHATE)	NET	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 031	NITROFURANTOIN	NIT	300 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 032	NORFLOXACIN	NX	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 079	NOVOBIOCIN 	NV	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 080	NOVOBIOCIN 	NV	30 mcg	5 ct 1 B. set 1 vl 5 vl	

ABBR.: Abbreviation | CON.: Concentration | Store between 2-8°C *
Pack Size: 5 ct = 5*50 Discs in Plastic Container | 1 Blister Set (B. set) = 5*50 Discs in Blister Pack
1 vl = 100 Discs in vial | 5 vl = 5*100 Discs in vial



Antibiotic Sensitivity Discs*

Concentration of Antibiotics as per CLSI (Formerly known as NCCLS) Standards

CODE	PRODUCT NAME	ABBR.	CON.	PACK SIZE	
TBD 081	NYSTATIN 	NS	100 units	5 ct 1 B. set 1 vl 5 vl	
TBD 033	OFLOXACIN	OF	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 066	OXACILLIN	OX	1 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 067	PEFLOXACIN	PF	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 034	PENICILLIN-G	P	10 units	5 ct 1 B. set 1 vl 5 vl	
TBD 035	PIPERACILLIN	PC	100 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 068	PIPERACILLIN / TAZOBACTAM	PIT	100/10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 036	POLYMYXIN-B	PB	300 units	5 ct 1 B. set 1 vl 5 vl	
TBD 069	PRULIFLOXACIN (ULIFLOXACIN)	PRU	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 037	RIFAMPICIN	R	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 038	SPARFLOXACIN	SPX	5 mcg	5 ct 1 B. set 1 vl 5 vl	

CODE	PRODUCT NAME	ABBR.	CON.	PACK SIZE	
TBD 070	SPECTINOMYCIN	SPT	100 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 039	STREPTOMYCIN	S	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 071	STREPTOMYCIN	HLS	300 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 072	TEICOPLANIN	TEI	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 040	TETRACYCLINE	TE	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 073	TICARCILLIN	TI	75 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 041	TICARCILLIN/ CLAVULANIC ACID	TCC	75/2.5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 042	TOBRAMYCIN	TOB	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 043	TRIMETHOPRIM	TR	5 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 044	VANCOMYCIN	VA	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 074	VORICONAZOLE (Antifungal)	VRC	1 mcg	5 ct 1 B. set 1 vl 5 vl	

Concentration of Antibiotics as per Non-CLSI (Formerly known as NCCLS) Standards

CODE	PRODUCT NAME	ABBR.	CON.	PACK SIZE	
TBD 002	AMOXICILLIN	AX	10 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 045	AMOXICILLIN	AX	30 mcg	5 ct 1 B. set 1 vl 5 vl	

CODE	PRODUCT NAME	ABBR.	CON.	PACK SIZE	
TBD 008	CEFADROXIL (CEPHADROXIL)	CFR	30 mcg	5 ct 1 B. set 1 vl 5 vl	
TBD 021	CLOXACILLIN	CX	5 mcg	5 ct 1 B. set 1 vl 5 vl	

Antibiotic Sensitivity Discs*

Concentration of Antibiotics as per CLSI (Formerly known as NCCLS) Standards

PRODUCT CODE	ANTIMICROBIAL AGENT	SYMBOL	DISC CONTENT	INTERPRETATIVE CRITERIA			QUALITY CONTROL LIMITS (mm)													
				SENSITIVE mm or more	INTERMEDIATE mm	RESISTANT mm or less	E.coli ATCC 25922	S.aureus ATCC 25923	Paeruginosa ATCC 27853	E.coli ATCC 35218	S.aureus ATCC 25913	E.faecalis ATCC 29212	H.influenzae ATCC 49247	H.influenzae ATCC 49766	K.pneumoniae ATCC 700603	N.gonorrhoeae ATCC 49266	S.pneumoniae ATCC 49619	C.jejuni ATCC 33560		
TBD 001	AMIKACIN <i>Enterobacteriaceae, Paeruginosa Acinetobacter & Staphylococcus</i>	AK	30 mcg	17	15-16	14	19-26	20-26	18-26	-	-	-	-	-	-	-	-	-		
TBD 003	AMOXICLAV (Amoxicillin/Clavulanic) <i>Enterobacteriaceae Staphylococcus spp. Haemophilus influenzae & haemophilus parainfluenzae</i>	AMC	30 mcg	18	14-17	13	18-24	28-36	-	17-22	-	-	-	-	-	-	-	-		
				20	-	19	-	-	-	-	-	15-23	-	-	-	-	-	-		
TBD 004	AMPICILLIN <i>Enterobacteriaceae Staphylococcus spp. Enterococcus spp. Haemophilus infleznaze & Haemophilus parainfluenzae Streptococcus spp. beta haemolytic group</i>	AMP	10 mcg	17	14-16	13	15-22	-	-	6	-	-	-	-	-	-	-	-		
				29	-	28	-	27-35	-	-	-	-	-	-	-	-	-	-		
				17	-	16	-	-	-	-	-	-	-	-	-	-	-	-	-	
				22	19-21	18	-	-	-	-	-	-	13-21	-	-	-	-	-	-	
TBD 005	AMPICILLIN/SULBACTAM <i>Enterobacteriaceae, Acinetobacter & Haemophilus infleznaze & Haemophilus parainfluenzae</i>	A/S	10/10 mcg	15	12-14	11	19-24	29-37	-	13-19	-	-	-	-	-	-	-	-		
				20	-	19	-	-	-	-	-	14-22	-	-	-	-	-	-	-	
TBD 006	AZITHROMYCIN <i>Staphylococcus, S.pneumoniae, Streptococcus spp., Viridians group & Streptococcus spp. Beta haemolytic group haemophilus infleznaze & Haemophilus parainfluenzae Neisseria meningitidis Salmonella Typhi</i>	AZM	15 mcg	18	14-17	13	-	21-26	-	-	-	-	-	-	-	-	19-25	-		
				12	-	-	-	-	-	-	-	13-21	-	-	-	-	-	-		
				20	-	-	-	-	-	-	-	-	-	-	-	30-38	-	-	-	
				13	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	
TBD 046	AZTREONAM <i>Enterobacteriaceae Paeruginosa haemophilus infleznaze & Haemophilus parainfluenzae</i>	AZT	30 mcg	21	18-20	17	28-36	-	-	31-38	-	-	-	-	-	10-16	-	-		
				22	16-21	15	-	-	23-29	-	-	-	-	-	-	-	-	-	-	-
TBD 047	CARBENICILLIN	CB	100 mcg	-	-	-	23-29	-	18-24	-	-	-	-	-	-	-	-	-		
				-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	
TBD 007	CEFACLOR <i>Enterobacteriaceae Staphylococcus spp. Haemophilus infleznaze & Haemophilus parainfluenzae</i>	CF	30 mcg	18	15-17	14	23-27	27-31	-	-	-	-	-	-	-	-	24-32	-		
				20	17-19	16	-	-	-	-	-	-	25-31	-	-	-	-	-	-	
TBD 009	CEFAZOLIN <i>Enterobacteriaceae Enterobacteriaceae (uncomplicated UTIs) Staphylococcus spp.</i>	CZ	30 mcg	23	20-22	19	21-27	-	-	-	-	-	-	-	-	-	-	-		
				15	-	14	-	-	-	-	-	-	-	-	-	-	-	-	-	
TBD 048	CEFDINIR <i>Enterobacteriaceae Staphylococcus spp. Haemophilus infleznaze & Haemophilus parainfluenzae</i>	CDR	5 mcg	20	17-19	16	24-28	25-32	-	-	-	-	-	-	-	-	40-49	26-31		
				20	-	-	-	-	-	-	-	24-31	-	-	-	-	-	-	-	
TBD 049	CEFAPIME <i>Enterobacteriaceae Paeruginosa, Acinetobacter Haemophilus infleznaze & haemophilus parainfluenzae Neisseria gonorrhoeae Streptococcus spp. Viridians group Streptococcus spp. beta haemolytic group</i>	CPM	30 mcg	25	19-24	18	31-37	-	25-31	-	-	-	-	-	-	-	-	26-31		
				18	15-17	14	-	23-29	-	-	-	-	-	-	-	-	-	-	-	
				26	-	-	-	-	-	-	-	25-31	-	-	-	-	-	-	-	
				31	-	-	-	-	-	-	-	-	-	-	-	-	37-46	-	-	
				24	22-23	21	-	-	-	-	-	-	-	-	-	-	-	28-35	-	-
TBD010	CEFAXIME <i>Enterobacteriaceae Haemophilus infleznaze & Haemophilus parainfluenzae Neisseria gonorrhoeae</i>	CFM	5 mcg	19	16-18	15	20-26	-	-	-	-	-	-	-	-	-	16-23	-		
				21	-	-	-	-	-	-	-	25-33	-	-	-	-	-	-	-	
TBD 011	CEFOPERAZONE <i>Enterobacteriaceae Staphylococcus spp.</i>	CS	75 mcg	21	16-20	15	28-34	24-33	23-29	-	-	-	-	-	-	-	-	-		
				21	16-20	15	28-34	24-33	23-29	-	-	-	-	-	-	-	-	-	-	
TBD 012	CEFOTAXIME (CEPHOTAXIME) <i>Enterobacteriaceae Acinetobacter & Staphylococcus spp. Haemophilus infleznaze & Haemophilus parainfluenzae Neisseria meningitidis Neisseria gonorrhoeae Streptococcus spp. Viridians group Streptococcus spp. Beta haemolytic group</i>	CTX	30 mcg	26	23-25	22	29-35	-	-	-	-	-	-	-	17-25	-	-	-		
				23	15-22	14	-	25-31	18-22	-	-	-	-	-	-	-	-	-	-	
				26	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	
				34	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	
				31	-	-	-	-	-	-	-	-	-	-	-	-	-	38-48	-	-
				28	26-27	25	-	-	-	-	-	-	-	-	-	-	-	-	31-39	-
TBD 050	CEFPDODXIME <i>Enterobacteriaceae Staphylococcus spp. Haemophilus infleznaze & Haemophilus parainfluenzae Neisseria gonohoeae</i>	CPD	10 mcg	21	18-20	17	23-28	19-25	-	-	-	-	-	-	9-16	-	28-34	-		
				21	-	-	-	-	-	-	-	25-31	-	-	-	-	-	-	-	
TBD 051	CEFROZIL <i>Enterobacteriaceae Staphylococcs Haemophilus infleznaze & Haemophilus parainfluenzae</i>	CPR	30 mcg	18	15-17	14	21-27	27-33	-	-	-	-	-	20-27	-	-	25-32	-		
				18	15-17	14	21-27	27-33	-	-	-	-	-	20-27	-	-	25-32	-	-	
TBD 014	CEFTAZIDIME <i>Enterobacteriaceae, B.cepacia Paeruginosa, Acinetobacter & Staphylococcus spp. Haemophilus infleznaze & Haemophilus parainfluenzae Neisseria gonorrhoeae</i>	CAZ	30 mcg	21	18-20	17	25-32	-	-	-	-	-	-	-	10-18	-	-	-		
				18	15-17	14	-	16-20	22-29	-	-	-	-	-	-	-	-	-	-	
				26	-	-	-	-	-	-	-	27-35	-	-	-	-	-	-	-	-
TBD 052	CEFTAZIDIME/CLAVULANIC ACID	CAC	30/10 mcg	22	19-21	19	23-29	-	-	-	-	-	-	-	-	-	-	-		

Antibiotic Sensitivity Discs*

Concentration of Antibiotics as per CLSI (Formerly known as NCCLS) Standards

PRODUCT CODE	ANTIMICROBIAL AGENT	SYMBOL	DISC CONTENT	INTERPRETATIVE CRITERIA			QUALITY CONTROL LIMITS (mm)											
				SENSITIVE mm or more	INTERMEDIATE mm	RESISTANT mm or less	E.coli ATCC 25922	S.aureus ATCC 25923	Paeruginosa ATCC 27853	E.coli ATCC 35218	S.aureus ATCC 25913	E.faecalis ATCC 29212	H.influenzae ATCC 49247	H.influenzae ATCC 49766	K.pneumoniae ATCC 700603	N.gonorrhoeae ATCC 49266	S.pneumoniae ATCC 49619	C.jejuni ATCC 33560
TBD 053	CEFTIZOXIME <i>Enterobacteriaceae</i> <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria gonarhoeae</i>	CZX	30 mcg	25 26 38	22-24 - -	21 - -	30-36 - -	27-35 - -	12-17 - -	- - -	- - -	- 29-39 -	- - -	- - -	- - 42-51	28-34 - -	- - -	
TBD 013	CEFTRIAXONE <i>Enterobacteriaceae</i> <i>Paeruginosa</i> , <i>Acientobacter</i> & <i>Staphylococcus</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria meningitidis</i> <i>Neisseria gonorrhoeae</i> <i>Streptococcus</i> spp. <i>Viridians</i> group <i>Streptococcus</i> spp. <i>Beta haemolytic</i> group	CTR	30 mcg	23 21 26 34 35 27 24	20-22 14-20 - - - 25-26 -	19 13 - - - 24 -	29-35 - - - - - -	- 22-28 - - - - -	- 17-23 - - - - -	- - - - - -	- - 31-39 - - -	- - - - - -	16-24 - - - - -	- - - - - -	- - - - - 39-51 -	- - - - - 30-35 -	- - - - - -	
TBD 015	CEFUROXIME <i>Enterobacteriaceae</i> (Parental) <i>Staphylococcus</i> spp. <i>Enterobacteriaceae</i> (Oral) <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria gonorrhoeae</i>	CXM	30 mcg	18 23 20 31	15-17 15-22 17-19 26-30	14 14 16 25	20-26 - - -	27-35 - - -	- - - -	- - - -	- - - -	- - - -	- - 28-36 -	- - - -	- - 33-41 -	- - - -	- - - -	
TBD 016	CEPHALOTHIN	CEP	30 mcg	-	-	-	15-21	29-37	-	-	-	-	-	-	-	26-32	-	
TBD 017	CHLORAMPHENICOL <i>Enterobacteriaceae</i> , <i>Staphylococcus</i> & <i>Enterococcus</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria meningitidis</i> <i>S. pneumoniae</i> <i>Streptococcus</i> spp. <i>Viridians</i> group & <i>Streptococcus</i> spp. <i>Beta haemolytic</i> group	C	30 mcg	18 29 26 21 21	13-17 26-28 20-25 - 18-20	12 25 19 20 17	21-27 - - - -	19-26 - - - -	- - - - -	- - - - -	- - 31-40 - -	- - - - -	- - - - -	- - - - -	- - - - -	- - - 23-27 -	- - - -	
TBD 018	CIPROFLOXACIN <i>Enterobacteriaceae</i> other than <i>S.typhi</i> & extraintestinal <i>Salmonella</i> spp., <i>Acientobacter</i> , <i>Staphylococcus</i> & <i>Enterococcus</i> spp. <i>Paeruginosa</i> For <i>S. typhi</i> and extraintestinal <i>Salmonella</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria meningitidis</i> <i>Neisseria gonorrhoeae</i>	CIP	5 mcg	26 21 25 31 21 35 41	22-25 16-20 19-24 21-30 - 33-34 28-40	21 15 18 20 - 32 27	29-37 - - - - - -	- 22-30 25-33 - - -	- - - - - -	- - - - - -	- - - - - -	- - - - 34-42 -	- - - - -	- - - - -	- - - - -	48-58	- - - -	
TBD 019	CLARITHROMYCIN <i>Staphylococcus</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>S. pneumoniae</i> , <i>Streptococcus</i> spp. <i>Viridians</i> group, <i>Streptococcus</i> spp. <i>Beta haemolytic</i> group	CLP	15 mcg	18 13 21	14-17 11-12 17-20	13 10 18	- - -	26-32 - -	- - -	- - -	- - -	- 11-17 -	- - -	- - -	- - -	- - 25-31	- -	
TBD 020	CLINDAMYCIN <i>Staphylococcus</i> spp. <i>S.pneumoniae</i> , <i>Streptococcus</i> spp. <i>Viridians</i> group, <i>Streptococcus</i> spp. <i>Beta haemolytic</i> group	CD	2 mcg	21 19	15-20 16-18	14 15	- -	24-30 -	- -	- -	- -	- -	- -	- -	- -	19-25	-	
TBD 054	COLISTIN (METHANE SULPHONSTE) <i>Paeruginosa</i>	CL	10 mcg	-	-	-	11-17	-	11-17	-	-	-	-	-	-	-	-	
TBD 022	Co-TRIMOXAZOLE (Sulphamethoxazole) <i>Enterobacteriaceae</i> , <i>Acientobacter</i> , <i>B. cepacia</i> , <i>S. maltophilia</i> , <i>staphylococcus</i> , <i>Haemophilus</i> <i>inflenzae</i> & <i>Haemophilus</i> <i>parainflenzae</i> <i>Neisseria meningitidis</i> <i>S. pneumoniae</i>	COT	25 mcg (23.75/ 1.25mcg)	16 30 19	11-15 26-29 16-18	10 25 15	23-29 - -	24-32 - -	- - -	- - -	- - -	>=20 24-32 -	- - -	- - -	- - -	- - 20-28	- -	
TBD 023	DOXYCYCLINE HCL <i>Enterobacteriaceae</i> <i>Acientobacter</i> <i>Staphylococcus</i> & <i>Enterococcus</i> spp. <i>S.pneumoniae</i>	DO	30 mcg	14 13 16 28	11-13 10-12 13-15 25-27	10 9 12 24	18-24 - - -	23-29 - - -	- - - -	- - - -	- - - -	- - - -	- - - -	- - - -	- - - -	- - - 25-34	- -	
TBD 055	ERTAPENEM <i>Enterobacteriaceae</i> <i>Staphylococcus</i> spp. <i>Neisseria gonorrhoeae</i>	ETP	10 mcg	22 19 19	19-21 16-18 -	18 15 -	29-36 - -	- 24-31 -	13-21 - -	- - -	- - -	- - 20-28	- - -	27-33 - -	- - -	28-35 - -	- -	
TBD 024	ERYTHROMYCIN <i>Staphylococcus</i> & <i>Enterococcus</i> spp. <i>S. pneumoniae</i> , <i>Streptococcus</i> spp. <i>Viridians</i> group, <i>Streptococcus</i> spp. <i>Beta haemolytic</i> group	E	15 mcg	23 21	14-22 16-20	13 15	- -	22-30 -	- -	- -	- -	- -	- -	- -	- -	- 25-30	- -	
TBD 056	FAROPENEM <i>Enterobacteriaceae</i>	FAR	5 mcg	-	-	-	20-26	27-34	-	-	-	-	15-22	-	-	27-35	-	
TBD 058	GATIFLOXACIN <i>Enterobacteriaceae</i> , <i>Paeruginosa</i> , <i>Acientobacter</i> & <i>Enterococcus</i> spp. <i>Staphylococcus</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria gonorrhoeae</i> <i>S. pneumoniae</i> , <i>Streptococcus</i> spp. <i>Viridians</i> group, <i>Streptococcus</i> spp. <i>Beta harmolytic</i> group	GAT	5 mcg	18 23 18 38 21	15-17 20-22 - 34-37 18-20	14 19 - 33 17	30-37 - - -	- 27-33 - -	20-28 - -	- - -	- - -	- - 33-41 -	- - -	- - -	- - 45-56 -	- - -	- 24-31 -	

Antibiotic Sensitivity Discs*

Concentration of Antibiotics as per CLSI (Formerly known as NCCLS) Standards

PRODUCT CODE	ANTIMICROBIAL AGENT	SYMBOL	DISC CONTENT	INTERPRETATIVE CRITERIA			QUALITY CONTROL LIMITS (mm)											
				SENSITIVE mm or more	INTERMEDIATE mm	RESISTANT mm or less	E.coli ATCC 25922	S.aureus ATCC 25923	P.aeruginosa ATCC 27853	E.coli ATCC 35218	S.aureus ATCC 259213	E.faecalis ATCC 29212	H.influenzae ATCC 49247	H.influenzae ATCC 49766	K.pneumoniae ATCC 700603	N.gonorrhoeae ATCC 49266	S.pneumoniae ATCC 49619	C.jejuni ATCC 33560
TBD 059	GEMIFLOXACIN <i>Enterobacteriaceae</i> <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>S.pneumoniae</i>	GEM	5 mcg	20	16-19	15	29-36	27-33	19-25	-	-	-	-	-	-	-	-	-
				18	-	-	-	-	-	-	30-37	-	-	-	-	-	-	
				23	20-22	19	-	-	-	-	-	-	-	-	-	-	28-34	
TBD 060	GENTAMICIN <i>Enterococcus</i> spp.	HLG	120 mcg	10	7-9	6	-	-	-	-	-	16-23	-	-	-	-	-	-
TBD 025	GENTAMICIN <i>Enterobacteriaceae</i> , <i>P.aeruginosa</i> , <i>Acinetobacter</i> & <i>Staphylococcus</i> spp.	GEN	10 mcg	15	13-14	12	19-26	19-27	17-23	-	-	-	-	-	-	-	-	-
TBD 061	IMIPENEM <i>Enterobacteriaceae</i> <i>P.aeruginosa</i> <i>Staphylococcus</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Acinetobacter</i> spp.	IPM	10 mcg	23	20-22	19	26-32	-	-	-	-	-	-	-	-	-	-	-
				19	16-18	15	-	-	20-28	-	-	-	-	-	-	-	-	
				16	-	-	-	-	-	-	-	21-29	-	-	-	-	-	
22	19-21	18	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	
TBD 026	KANAMYCIN <i>Enterobacteriaceae</i> , <i>Staphylococcus</i> spp.	K	30 mcg	18	14-17	13	17-25	19-26	-	-	-	-	-	-	-	-	-	-
TBD 062	LEVOFLOXACIN <i>Enterobacteriaceae</i> , <i>S.typhi</i> , <i>P.aeruginosa</i> , <i>Acinetobacter</i> spp <i>S.maltophilia</i> , <i>Enterococcus</i> spp., <i>S.pneumoniae</i> , <i>Streptococcus</i> spp. <i>Viridians</i> group, <i>Streptococcus</i> spp. Beta haemolytic group <i>Staphylococcus</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i>	LE	5 mcg	21	17-20	16	29-37	-	-	-	-	-	-	-	-	-	-	-
				22	15-21	14	-	-	19-26	-	-	-	-	-	-	-	-	-
				17	14-16	13	-	-	-	-	-	-	-	-	-	-	20-25	-
				19	16-18	15	25-30	25-30	-	-	-	-	-	-	-	-	-	-
				17	-	-	-	-	-	-	32-40	-	-	-	-	-	-	-
TBD 063	LINEZOLID <i>Staphylococcus</i> spp. <i>Enterococcus</i> spp. <i>Streptococcus</i> group A, B, C & G <i>Corynebacterium</i> spp. <i>S.pneumoniae</i>	LZ	30 mcg	21	-	20	-	25-32	-	-	-	-	-	-	-	-	-	-
				23	21-22	20	-	-	-	-	-	-	-	-	-	-	-	
				21	-	-	-	-	-	-	-	-	-	-	-	-	25-34	-
TBD 027	LOMEFLOXACIN <i>Enterobacteriaceae</i> , <i>P.aeruginosa</i> & <i>Staphylococcus</i> spp. <i>B.cepacia</i> <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria gonorrhoeae</i>	LOM	10 mcg	22	19-21	18	27-33	23-29	22-28	-	-	-	-	-	-	-	-	-
				22	-	-	-	-	-	-	-	33-41	-	-	-	-	-	-
				38	27-37	26	-	-	-	-	-	-	-	-	45-54	-	-	-
TBD 028	MEROPENEM <i>Enterobacteriaceae</i> <i>P.aeruginosa</i> <i>Staphylococcus</i> spp. <i>B.cepacia</i> <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria meningitidis</i> <i>Acinetobacter</i> spp.	MRP	10 mcg	23	20-22	19	28-35	-	-	-	-	-	-	-	-	-	-	-
				19	16-18	15	-	-	27-33	-	-	-	-	-	-	-	-	-
				16	14-15	13	-	29-37	-	-	-	-	-	-	-	-	-	-
				20	16-19	15	-	-	-	-	-	-	-	-	-	-	-	-
				20	-	-	-	-	-	-	-	20-28	-	-	-	-	-	-
30	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	28-35	-	
18	5-17	14	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	
TBD 064	MINOCYCLINE <i>Enterobacteriaceae</i> , <i>Acinetobacter</i> <i>B.cepacia</i> , <i>S.maltophilia</i> , <i>Staphylococcus</i> & <i>Enterococcus</i> <i>Neisseria meningitidis</i>	MI	30 mcg	16	13-15	12	19-25	25-30	-	-	-	-	-	-	-	-	-	-
				19	15-18	14	-	-	-	-	-	-	-	-	-	-	-	-
				26	-	-	-	-	-	-	-	-	-	-	-	-	-	
TBD 065	MOXIFLOXACIN <i>Staphylococcus</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>S.pneumoniae</i>	MO	5 mcg	24	21-23	20	28-35	28-35	17-25	-	-	-	-	-	-	-	-	-
				18	-	-	-	-	-	-	-	31-39	-	-	-	-	-	-
				18	15-17	14	-	-	-	-	-	-	-	-	-	-	25-31	-
TBD 029	NALIDIXIC ACID <i>Enterobacteriaceae</i>	NA	30 mcg	19	14-18	13	22-28	-	-	-	-	-	-	-	-	-	-	-
TBD 030	NETILLIN (Netimicin Sulphate) <i>Enterobacteriaceae</i> , <i>P.aeruginosa</i> <i>Staphylococcus</i> spp.	NET	30 mcg	15	13-14	12	22-30	22-31	17-23	-	-	-	-	-	-	-	-	-
TBD 031	NITROFURANTOIN <i>Enterobacteriaceae</i> , <i>Staphylococcus</i> & <i>Enterococcus</i> spp.	NIT	300 mcg	17	15-16	14	20-25	18-22	-	-	-	-	-	-	-	-	23-29	-
TBD 032	NORFLOXACIN <i>Enterobacteriaceae</i> , <i>P.aeruginosa</i> <i>Staphylococcus</i> & <i>Enterococcus</i>	NX	10 mcg	17	13-16	12	28-35	17-28	22-29	-	-	-	-	-	-	-	15-21	-
TBD 033	OFLOXACIN <i>Enterobacteriaceae</i> , <i>P.aeruginosa</i> , <i>S.pneumoniae</i> , <i>Streptococcus</i> spp. <i>Viridians</i> group, <i>Streptococcus</i> spp. Beta haemolytic group <i>Staphylococcus</i> spp. <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria gonorrhoeae</i>	OF	5 mcg	16	13-15	12	29-33	-	17-21	-	-	-	-	-	-	-	16-21	-
				18	15-17	14	-	24-48	-	-	-	-	-	-	-	-	-	
				16	-	-	-	-	-	-	-	-	-	-	-	-	-	
				31	25-30	24	-	-	-	-	-	31-40	-	-	-	-	43-51	-
TBD 066	OXLOXACIN <i>Staphylococcus (S.pseudintermedius)</i> <i>S.pneumoniae</i>	OX	1 mcg	18	-	17	-	18-24	-	-	-	-	-	-	-	-	-	-
20	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	<=12"	-	
TBD 067	PEFLOXACIN <i>Enterobacteriaceae (S.typhi)</i>	PF	5 mcg	24	-	23	25-33	-	-	-	-	-	-	-	-	-	-	-
TBD 034	PENICILLIN-G <i>Staphylococcus</i> spp. <i>Enterococcus</i> spp. <i>Neisseria gonorrhoeae</i> <i>Streptococcus</i> spp. Beta haemolytic group	P	10 units	29	-	28	-	26-37	-	-	-	-	-	-	-	-	-	-
				15	-	14	-	-	-	-	-	-	-	-	-	-	-	
				47	27-46	26	-	-	-	-	-	-	-	-	-	26-34	-	-
				24	-	-	-	-	-	-	-	-	-	-	-	-	24-30	-

ABBR.: Abbreviation | CON.: Concentration | Store between 2-8°C *
Pack Size: 5 ct = 5*50 Discs in Plastic Container | 1 Blister Set (B. set) = 5*50 Discs in Blister Pack
1 vl = 100 Discs in vial | 5 vl = 5*100 Discs in vial



Antibiotic Sensitivity Discs*

Concentration of Antibiotics as per CLSI (Formerly known as NCCLS) Standards

PRODUCT CODE	ANTIMICROBIAL AGENT	SYMBOL	DISC CONTENT	INTERPRETATIVE CRITERIA			QUALITY CONTROL LIMITS (mm)											
				SENSITIVE mm or more	INTERMEDIATE mm	RESISTANT mm or less	E.coli ATCC 25922	S.aureus ATCC 25923	Paeruginosa ATCC 27853	E.coli ATCC 35218	S.aureus ATCC 25913	E.faecalis ATCC 29212	H.influenzae ATCC 49247	H.influenzae ATCC 49766	K.pneumoniae ATCC 700603	N.gonorrhoeae ATCC 49266	S.pneumoniae ATCC 49619	C.jejuni ATCC 33560
TBD 035	PIPERACILLIN <i>Enterobacteriaceae</i> & <i>Acinetobacter</i> spp. <i>Paeruginosa</i>	PI	100 mcg	21 21	18-20 15-20	17 14	24-30 -	- -	25-33 -	12-18 -	- -	- -	- -	- -	- -	- -	- -	- -
TBD 068	PIPERACILLIN/TOZOBACTAM <i>Enterobacteriaceae</i> & <i>Acinetobacter</i> spp. <i>Paeruginosa</i> <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i>	PIT	100/10 mcg	21 21 21	18-20 15-20 -	17 14 -	24-30 - -	27-36 - -	25-33 - -	24-30 - -	- - -	- - 33-38	- - -	- - -	- - -	- - -	- - -	- - -
TBD 036	POLYMYXIN-B	PB	300 units	-	-	-	13-19	-	14-18	-	-	-	-	-	-	-	-	-
TBD 069	PRULIFLOXACIN (ULIFLOXACIN)	PRU	5 mcg	-	-	-	32-38	20-26	27-33	-	-	-	-	-	-	-	-	-
TBD 037	RIFAMPICIN <i>Staphylococcus</i> , <i>Enterococcus</i> spp., <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria meningitidis</i> <i>S.pneumoniae</i>	RIF	5 mcg	- 20 25 19	- 17-19 29-24 17-18	- 16 19 16	8-10 - -	- 26-34 -	- - -	- - -	- - -	- 22-30 -	- - -	- - -	- - -	- - -	- 25-30 -	- -
TBD 038	SPARFLOXACIN <i>Enterobacteriaceae</i> <i>Staphylococcus</i> , <i>S.pneumoniae</i>	SPX	5 mcg	- 19	- 16-18	- 15	30-38 -	- 27-33	21-29 -	- -	- -	- -	32-40 -	- -	- -	- -	43-51 -	- 21-27
TBD 070	SPECTINOMYCIN <i>Neisseria gonorrhoeae</i>	SPT	100 mcg	18	15-17	14	-	-	-	-	-	-	-	-	-	23-29	-	-
TBD 039	SREPTOMYCIN <i>Enterococcus</i> spp.	S	10 mcg	15	12-14	1	12-20	14-22	-	-	-	-	-	-	-	-	-	-
TBD 071	STREPTOMYCIN <i>Enterobacteriaceae</i>	HLS	300 mcg	15	12-14	11	-	14-22	-	-	-	-	-	-	-	-	-	-
TBD 072	TEICoplanin <i>Enterococcus</i> spp.	TEI	30 mcg	14	11-13	10	-	15-21	-	-	-	-	-	-	-	-	-	-
TBD 040	TETRACYCLINE <i>Enterobacteriaceae</i> , <i>Acinetobacter</i> <i>Staphylococcus</i> , <i>Enterococcus</i> & <i>Neisseria meningitidis</i> <i>Haemophilus influenzae</i> & <i>Haemophilus parainfluenzae</i> <i>Neisseria gonorrhoeae</i> <i>S.pneumoniae</i> , <i>Streptococcus</i> spp. Beta haemolytic group & <i>Viridians</i> group	TE	30 mcg	15 19 29 38 28 23	12-14 15-18 26-28 31-37 25-27 19-22	11 14 25 30 24 18	18-25 - -	- 24-30 -	- - -	- - -	- -	- 14-22 -	- -	- -	- -	- 30-42 -	- 27-31 -	- -
TBD 073	TICARICILLIN	TI	75 mcg	-	-	-	24-30	-	21-27	6	-	-	-	-	-	-	-	-
TBD 041	TICARICILLIN/CLAVULANIC ACID <i>Enterobacteriaceae</i> & <i>Acinetobacter</i> <i>Paeruginosa</i> <i>Staphylococcus</i> spp.	TCC	75/10 mcg	20 24 23	15-19 16-23 -	14 15 22	24-30 - -	- 20-28 29-37	21-25 -	- -	- -	- -	- -	- -	- -	- -	- -	- -
TBD 042	TOBRAMYCIN <i>Enterobacteriaceae</i> , <i>Paeruginosa</i> , <i>Acinetobacter</i> & <i>Staphylococcus</i> spp.	TOB	10 mcg	15	13-14	12	18-26	19-29	20-26	-	-	-	-	-	-	-	-	-
TBD 043	TRIMETHOPRIM <i>Enterobacteriaceae</i> , <i>Staphylococcus</i>	TR	5 mcg	16	11-15	10	19-26	-	-	-	-	-	-	-	-	-	-	-
TBD 044	VANCOMYCIN <i>Enterococcus</i> spp. <i>S.pneumoniae</i> , <i>Streptococcus</i> spp. Beta haemolytic group & <i>Streptococcus</i> spp. <i>Viridians</i>	VA	30 mcg	- 17 17	- 15-16 -	- 14 -	- -	17-21 -	- -	- -	- -	- -	- -	- -	- -	- -	- 20-27	- -

ZONE SIZE INTERPRETATIVE CHART FOR ANTIFUNGAL AGENT										
CODE	ANTIFUNGAL AGENT	SYMBOL	DISC CONTENT	ZONE DIAMETER, NEAREST WHOLE (mm)			QUALITY CONTROL LIMITS (mm)			
				RESISTANT mm or less	S-DD*	SUSCEPTIBLE MM OR MORE	<i>C.albicans</i> ATCC 90028	<i>C.parapsilosis</i> ATCC 22019	<i>C.tropicalis</i> ATCC 750	<i>C.krusei</i> ATCC 6258
TBD 057	FLUCONAZOLE	FLC	25 mcg	14	15-18	19	28-39	22-33	26-37	-
TBD 074	VORICONAZOLE	VRC	1 mcg	13	14-16	17	31-42	28-37	-	16-25

- Zone size interpretative criteria is per as CLSI standard.
- For *E. coli*, *S. aureus*, *P. aeruginosa*: Mueller Hinton Agar (MHA); For *Haemophilus* spp.: *Haemophilus* Test Medium.
- For *S. pneumoniae*: Mueller Hinton Agar with 5% Sheep Blood; For *N. gonorrhoeae*: GC Agar Base with 1% defined growth supplement.

Microbiological Testing & Analytical Solutions

Differentiation Discs

CODE	PRODUCT NAME	PACK SIZE	
TBL 046	BACITRACIN DISCS * for identification of <i>Streptococcus pyogenes</i>	50 Discs/vl	
TBL 021	BILE ESCULIN DISCS * for detection of esculin hydrolysis in the presence of bile	50 Discs/vl	
TBL 022	DMACA INDOLE DISCS * for indole testing	50 Discs/vl	
TBL 027	NITRATE DISCS * substrate for nitrate reduction	50 Discs/vl	

Differentiation Test Strips

CODE	PRODUCT NAME	PACK SIZE	
TBL 023	KOVAC'S REAGENT STRIPS* for indole testing	25 Strips /1 vl	

Staining and Indicator Solutions

CODE	PRODUCT NAME	PACK SIZE	
TBL 041	ACID FAST DECOLOURIZER	125 ml 500 ml	
TBL 090	BASIC FUCHSIN 0.1%	125 ml 500 ml	
TBL 058	EOSIN, 2% W/V	125 ml 500 ml	
TBL 025	GRAM'S CRYSTAL VIOLET Bacteriological Stain	125 ml 500 ml	
TBL 032	GRAM'S DECOLOURIZER	125 ml 500 ml	
TBL 034	GRAM'S IODINE Bacteriological Mordant	125 ml 500 ml	
TBL 026	GRAM'S SAFRANINE Counter Stain for Bacteria	125 ml 500 ml	
TBL 059	GENTIAN VIOLET	125 ml 500 ml	
TBL 016	GIEMSA'S STAIN SOLUTION	100 ml 500 ml	
TBL 008	HAEMATOXYLIN STAIN (Delafield) For nuclear staining in histology and cytology	125 ml 500 ml	
TBL 042	HARRIS HAEMATOXYLIN (Papanicolaou's Solution 1) Cytological Stain	125 ml 500 ml	
TBL 053	JSB STAIN SOLUTION NO. 1	125 ml	
TBL 054	JSB STAIN SOLUTION NO. 2	125 ml	
TBL 033	LACTOPHENOL COTTON BLUE	100 ml 500 ml	
TBL 092	LEISHMAN'S STAIN (Twin Pack)	250 ml 500 ml	
TBL 093	LEISHMAN STAIN SOLUTION	250 ml 500 ml	

CODE	PRODUCT NAME	PACK SIZE	
TBL 028	NITRATE REAGENT DISCS * (Double Pack) for detection of nitrate reduction Part I : Discs Part II (a) : Reagent A Part II (b) : Reagent B	50 Discs/vl 2.5 ml/vl 2.5 ml/vl	
TBL 029	ONPG DISCS * testing for ONPG	50 Discs/vl	
TBL 030	OXIDASE DISCS * for detection of oxidase production by microorganisms	50 Discs/vl	
TBL 056	STERILE DISCS *	50 Discs/vl	

CODE	PRODUCT NAME	PACK SIZE	
TBL 024	LEAD ACETATE PAPER STRIPS* for detection of H ₂ S production	25 Strips /1 vl	

CODE	PRODUCT NAME	PACK SIZE	
TBL 051	LITMUS BLUE SOLUTION	125 ml 500 ml	
TBL 052	LITMUS RED SOLUTION	125 ml 500 ml	
TBL 019	LOEFFLER'S METHYLENE BLUE	125 ml 500 ml	
TBL 035	LUGOL'S IODINE	125 ml 500 ml	
TBL 060	MALACHITE GREEN, 1% W/V	100 ml 500 ml	
TBL 063	METHYLENE BLUE (AQUEOUS)	125 ml 500 ml	
TBL 094	NEISSER'S METHYLENE BLUE	100 ml	
TBL 089	NEUTRAL RED	125 ml	
TBL 039	PAPANICOLAOU'S-EA-36	125 ml 500 ml	
TBL 055	PAPANICOLAOU'S-EA-50 (Papanicolaou's 3b)	125 ml 500 ml	
TBL 049	PAPANICOLAOU'S-EA-65	125 ml 500 ml	
TBL 038	PAPANICOLAOU'S-OG-6 (Papanicolaou's 2b Orange)	125 ml 500 ml	
TBL 070	PHENOLPHTHALEIN, 1% w/v	125 ml 500 ml	
TBL 062	POTASSIUM PERMANGANATE SOLUTION	125 ml 500 ml	
TBL 095	READYMADE LEISHMAN STAIN WITH BUFFER POUCH / TABLET	500 ml	
TBL 040	STRONG CARBOL FUCHSIN (ZN)	125 ml 500 ml	
TBL 013	UNIVERSAL INDICATOR SOLUTION PH 4-11	125 ml 500 ml	

Microbiological Testing & Analytical Solutions

Readymade Staining Kits

CODE	PRODUCT NAME	PACK SIZE	
TK 005	ALBERT'S METACHROMATIC STAIN-KIT i) ALBERT'S STAIN A-125 ml ii) ALBERT'S STAIN B-125 ml	1 Kit	
TK 003	CAPSULE STAIN-KIT for Bacterial Capsule Staining i) METHYLENE BLUE (AQUEOUS)-125 ml ii) NIGROSIN STAIN-100 ml	1 Kit	
TK 004	GRAM'S STAIN-KIT Bacteriological Stain i) GRAM'S CRYSTAL VIOLET Bacteriological Stain ii) GRAM'S DECOLOURIZER Bacteriological Stain iii) GRAM'S IODINE Bacteriological Stain iv) GRAM'S SAFRANINE 0.5% W/V Counter Stain for Bacteria	1 Kit (Each 125 ml) 1 Kit (Each 500 ml)	

CODE	PRODUCT NAME	PACK SIZE	
TK 002	FIELD'S STAIN-KIT for Malarial Parasite i) FIELD'S STAIN A-125 ml ii) FIELD'S STAIN B-125 ml	1 Kit	
TK 011	PAPANICOLAOU'S STAIN-KIT I) PAPANICOLAOU'S-EA-50-125 ml ii) PAPANICOLAOU'S-OG-6-125 ml iii) HARRIS HAEMATOXYLIN-125 ml (Papanicolaou's Solution 1) Cytological Stain	1 Kit	
TK 009	SCHAEFFER & FULTON'S SPORE STAIN-KIT i) SCHAEFFER & FULTON'S SPORE STAIN A-100 ml ii) SCHAEFFER & FULTON'S SPORE STAIN B-100 ml	1 Kit	
TK 001	ZN ACID FAST STAIN-KIT for Acid Fast Bacteria i) ACID FAST DECOLOURIZER-125 ml ii) STRONG CARBOL FUCHSIN (ZN)-125 ml iii) LOEFFLER'S METHYLENE BLUE-125 ml	1 Kit (Each 125 ml) 1 Kit (Each 500 ml)	

Analytical Reagent

CODE	PRODUCT NAME	PACK SIZE	
TR 002	BARIUM CHLORIDE SOLUTION, 10% W/V	500 ml	
TR 105	BARRITT REAGENT A (for VP test)	100 ml	
TR 106	BARRITT REAGENT B (for VP test)	100 ml	
TR 120	BUTTERFIELD'S PHOSPHATE BUFFERED DILUTION WATER	500 ml	
TR 044	C.S.F. DILUTING FLUID	125 ml	
TS 207	DMACA REAGENT (10 ml/vl)	1 vl	
TR 107	EDTA (di-Sodium) 5%	125 ml	
TR 108	EOSINE YELLOW STAIN SOLUTION 2% W/V	125 ml 500 ml	
TR 052	FEHLING SOLUTION A	500 ml	
TR 053	FEHLING SOLUTION B	500 ml	
TR 109	FIELD'S STAIN SOLUTION A	125 ml 500 ml	
TR 110	FIELD'S STAIN SOLUTION B	125 ml 500 ml	
TR 111	FOLIN & WU'S PHOSPHATE, MOLYBDATE SOLUTION	100 ml 500 ml	

CODE	PRODUCT NAME	PACK SIZE	
TR 008	KOVAC'S INDOLE REAGENT	100 ml	
TR 029	MILLON'S REAGENT (For Protein)	125 ml 500 ml	
TR 060	NESSLER'S REAGENT (For detection of Ammonia and Ammonium Salt)	100 ml	
TR 113	NESSLER'S REAGENT (For Serum Urea Nitrogen)	100 ml	
TR 020	OXIDASE REAGENT	100 ml	
TR 114	POTASSIUM CHROMATE 5% W/V	125 ml	
TR 115	POTASSIUM OXALATE 5% W/V	125 ml	
TS 170	PYR REAGENT (10 ml/vl)	1 vl	
TR 013	SODIUM CITRATE 3.8% W/V	500 ml	
TR 015	SULPHANILIC ACID 0.8%	100 ml	
TR 085	5-SULPHOSALICYLIC ACID 3%	125 ml	
TS 208	TDA REAGENT (10 ml/vl)	1 vl	
TR 009	O-TOLIDINE REAGENT (Free Chloride)	500 ml	

Microbiological Testing & Analytical Solutions

Standard Solutions

CODE	PRODUCT NAME	PACK SIZE	
TR 102	BUFFER Sol. pH 4.0	500 ml	
TR 103	BUFFER Sol. pH 7.0	500 ml	
TR 104	BUFFER Sol. pH 9.2	500 ml	
TR 048	EDTA 0.1M (0.2N) VOLUMETRIC SOLUTION	500 ml	
TR 049	EDTA 0.5M (1N) VOLUMETRIC SOLUTION	500 ml	
TR 118	HYDROCHLORIC ACID 0.1N	500 ml 2.5 Ltr	
TR 132	IODINE SOLUTION N/10 (0.1 N)	500 ml	

CODE	PRODUCT NAME	PACK SIZE	
TR 112	McFARLAND STANDARD KIT (Each set contains 1 tube of 0.5, 1, 2, 3, 4 Mcfarland Standard)	1 Nos	
TR 028	α-NAPHTHYLAMINE SOLN.	100 ml	
TR 070	SILVER NITRATE SOLUTION 0.1N	500 ml	
TR 071	SILVER NITRATE SOLUTION 0.05N	500 ml	
TR 117	SODIUM HYDROXIDE 1N	500 ml 2.5 Ltr	
TR 116	SODIUM HYDROXIDE SOLUTION (N/10)	500 ml 2.5 Ltr	
TR 074	SODIUM THIOSULPHATE 0.1N	500 ml 1 Ltr	

pH Indicators Paper Strips

CODE	PRODUCT NAME	PACK SIZE	
TBL 031	pH INDICATOR PAPER 1.0-14.0	200 lvs	
TBL 036	pH INDICATOR PAPER 2.0-10.5	200 lvs	
TBL 037	pH INDICATOR PAPER 2.0-4.5	200 lvs	
TBL 044	pH INDICATOR PAPER 3.8-5.3	200 lvs	
TBL 047	pH INDICATOR PAPER 3.5-6.0	200 lvs	
TBL 048	pH INDICATOR PAPER 5.0-7.5	200 lvs	

CODE	PRODUCT NAME	PACK SIZE	
TBL 057	pH INDICATOR PAPER 6.5-9.0	200 lvs	
TBL 061	pH INDICATOR PAPER 8.0-10.5	200 lvs	
1634	LITMUS PAPER-BLUE	200 lvs	
1636	LITMUS PAPER-RED	200 lvs	
1635	STARCH IODIDE PAPERS	200 lvs	

pH Buffer Capsules

CODE	PRODUCT NAME	PACK SIZE	
284	BUFFER CAPSULES pH 4.0	10 Cap. 10*10	
285	BUFFER CAPSULES pH 7.0	10 Cap. 10*10	

CODE	PRODUCT NAME	PACK SIZE	
286	BUFFER CAPSULES pH 9.2	10 Cap. 10*10	



Culture Strains

Quality, Convenience, and Reliability in Every Dish

Features

- **Accredited Strain:** ISO certified for unparalleled quality
- **Extensive Collection:** Over 90 strains to suit your specific needs
- **All-in-One Design:** Minimize contamination risks with our innovative packaging
- **Ready-to-Use Format:** Save time and money with pre-prepared discs
- **Online Available Certificate of Analysis:** Access detailed strain information instantly
- **Easy-to-Use:** Discs Liquify easily

Benefits

- **Guaranteed Purity:** First-generation derivatives from trusted sources
- **High Viability:** Discs contain at least 10^6 organisms at purchase
- **Versatility:** Suitable for both Solid and Liquid cultures
- **Long Shelf Life:** Store easily and maintain viability for extended periods
- **Convenient Packaging:** Available in packs of 10 & 25 discs per vial

For Agar Plates

Place disc on suitable growth medium

Leave disc for a few minutes to liquify, then spread plate and incubate to produce isolated colonies

For Broth Media

Place disc in a small volume of suitable broth medium such as Brain-heart Infusion

Allow disc a few minutes to dissolve, then spread aliquot onto a plate of suitable growth medium

Obtain a stock culture

Culture Strains

CODE	PRODUCT NAME	PACK SIZE	
TMCS 001	<i>Aspergillus brasiliensis</i> NCPF 2275	10	
TMCS 086	<i>Acinetobacter baumannii</i> NCTC 12156	10	
TMCS 002	<i>Bacillus cereus</i> (Recently renamed as <i>Bacillus toyonensis</i>) NCTC 10320	10	
TMCS 003	<i>Bacillus cereus</i> NCTC 7464	10 25	
TMCS 004	<i>Bacillus subtilis</i> NCTC 10400	10 25	
TMCS 005	<i>Bacteroides fragilis</i> NCTC 9343	10	
TMCS 006	<i>Burkholderia cepacia</i> NCTC 10661	10	
TMCS 007	<i>Camphylobacter jejuni</i> NCTC 11322	10	
TMCS 008	<i>Camphylobacter jejuni</i> NCTC 11351	10	
TMCS 009	<i>Candida albicans</i> NCPF 3255	10 25	
TMCS 010	<i>Candida albicans</i> NCPF 3179	10 25	
TMCS 091	<i>Candida guilliermondii</i> NCPF 3099	10	
TMCS 011	<i>Candida krusei</i> NCPF 3953	10	
TMCS 092	<i>Candida lusitanae</i> NCPF 3954	10	
TMCS 088	<i>Candida parapsilosis</i> NCPF 3104	10	
TMCS 093	<i>Candida tropicalis</i> NCPF 3111	10	
TMCS 013	<i>Citrobacter freundii</i> NCTC 9750	10	
TMCS 014	<i>Clostridium perfringens</i> NCTC 8237	10	
TMCS 015	<i>Clostridium sporogenes</i> NCTC 532	10	
TMCS 017	<i>Enterobacter cloacae</i> NCTC 13380	10	
TMCS 018	<i>Enterococcus faecalis</i> NCTC 775	10 25	
TMCS 019	<i>Enterococcus faecalis</i> NCTC 12697	10 25	
TMCS 020	<i>Enterococcus faecalis</i> NCTC 13379	10	
TMCS 085	<i>Enterococcus faecium</i> NCTC 7171	10	
TMCS 021	<i>Enterococcus hirae</i> NCTC 13383	10	
TMCS 022	<i>Enterobacter cloacae</i> NCTC 13406	10	
TMCS 023	<i>Escherichia coli</i> (CRE) NCTC 13476	10	
TMCS 024	<i>Escherichia coli</i> NCTC 12241	10 25	
TMCS 025	<i>Escherichia coli</i> NCTC 11954	10 25	
TMCS 026	<i>Escherichia coli</i> NCTC 10418	10 25	
TMCS 027	<i>Escherichia coli</i> NCTC 12923	10 25	

CODE	PRODUCT NAME	PACK SIZE	
TMCS 028	<i>Escherichia coli</i> NCTC 11560	10	
TMCS 029	<i>Escherichia coli</i> NCTC 9001	10	
TMCS 030	<i>Escherichia coli</i> (mcr-1) NCTC 13846	10	
TMCS 031	<i>Escherichia coli</i> 0157 non-toxigenic NCTC 12900	10	
TMCS 032	<i>Haemophilus influenzae</i> NCTC 12699	10	
TMCS 033	<i>Haemophilus influenzae</i> NCTC 11931	10	
TMCS 034	<i>Haemophilus influenzae</i> NCTC 8468	10	
TMCS 035	<i>Haemophilus influenzae</i> NCTC 12975	10	
TMCS 036	<i>Klebsiella aerogenes</i> / <i>Raoultella planticola</i> NCTC 9528	10	
TMCS 087	<i>Klebsiella aerogenes</i> NCTC 10006	10 25	
TMCS 037	<i>Klebsiella pneumoniae</i> (CRE) NCTC 13438	10	
TMCS 038	<i>Klebsiella pneumoniae</i> (CRE) NCTC 13440	10	
TMCS 039	<i>Klebsiella pneumoniae</i> (CRE) NCTC 13442	10	
TMCS 040	<i>Klebsiella pneumoniae</i> (CRE) NCTC 13443	10	
TMCS 041	<i>Klebsiella pneumoniae</i> NCTC 9633	10 25	
TMCS 042	<i>Klebsiella pneumoniae</i> NCTC 13368	10	
TMCS 043	<i>Lactobacillus brevis</i> NCTC 13386	10	
TMCS 044	<i>Legionella anisa</i> NCTC 11974	10	
TMCS 045	<i>Legionella pneumophila</i> serogroup 1 NCTC 11192	10	
TMCS 046	<i>Listeria innocua</i> NCTC 11288	10	
TMCS 089	<i>Listeria ivanovii</i> NCTC® 11007	10	
TMCS 047	<i>Listeria monocytogenes</i> NCTC 7973	10 25	
TMCS 048	<i>Listeria monocytogenes</i> NCTC 13372	10	
TMCS 049	<i>Listeria monocytogenes</i> NCTC 11994	10	
TMCS 050	<i>Listeria monocytogenes</i> NCTC 10527	10	
TMCS 090	<i>Mycobacterium smegmatis</i> NCTC® 8159	10	
TMCS 051	<i>Neisseria gonorrhoea</i> NCTC 8375	25	
TMCS 052	<i>Neisseria gonorrhoea</i> NCTC 12700	25	
TMCS 053	<i>Proteus mirabilis</i> NCTC 13376	10	
TMCS 054	<i>Proteus mirabilis</i> NCTC 10975	10	

CODE	PRODUCT NAME	PACK SIZE	
TMCS 055	<i>Proteus mirabilis</i> NCTC 11938	10	
TMCS 056	<i>Proteus vulgaris</i> NCTC 4175	10 25	
TMCS 057	<i>Pseudomonas aeruginosa</i> NCTC 12903	10 25	
TMCS 058	<i>Pseudomonas aeruginosa</i> NCTC 12924	10 25	
TMCS 059	<i>Pseudomonas aeruginosa</i> NCTC 13359	10	
TMCS 060	<i>Pseudomonas aeruginosa</i> NCTC 10662	10	
TMCS 061	<i>Rhodococcus equi</i> NCTC 1621	10	
TMCS 062	<i>Saccharomyces cerevisiae</i> NCPF 3178	10	
TMCS 063	<i>Saccharomyces cerevisiae</i> NCTC 10716	10	
TMCS 064	<i>Salmonella enteritidis</i> NCTC 12694	10	
TMCS 065	<i>Salmonella nottingham</i> NCTC 7832	10	
TMCS 066	<i>Salmonella poona</i> NCTC 4840	10 25	
TMCS 067	<i>Salmonella typhimurium</i> NCTC 12023	10 25	
TMCS 068	<i>Serratia marcescens</i> NCTC 13382	10	
TMCS 069	<i>Shigella sonnei</i> NCTC 12984	10	
TMCS 070	<i>Staphylococcus aureus</i> NCTC 12981	10 25	

CODE	PRODUCT NAME	PACK SIZE	
TMCS 071	<i>Staphylococcus aureus</i> NCTC 12973	10 25	
TMCS 072	<i>Staphylococcus aureus</i> NCTC 7447	10 25	
TMCS 073	<i>Staphylococcus aureus</i> NCTC 10788	10 25	
TMCS 074	<i>Staphylococcus aureus</i> (MRSA) NCTC 12493	10	
TMCS 075	<i>Staphylococcus aureus</i> NCTC 6571	10	
TMCS 076	<i>Staphylococcus aureus</i> (MRSA) NCTC 13373	10	
TMCS 077	<i>Staphylococcus epidermidis</i> NCTC 13360	10 25	
TMCS 078	<i>Streptococcus agalactiae</i> NCTC 8181	10	
TMCS 079	<i>Streptococcus dysgalactiae</i> subsp. <i>Equisimilis</i> NCTC 13762	10	
TMCS 080	<i>Streptococcus pneumoniae</i> NCTC 12695	10 25	
TMCS 081	<i>Streptococcus pneumoniae</i> NCTC 12977	10	
TMCS 082	<i>Streptococcus pyogenes</i> NCTC 12696	10 25	
TMCS 083	<i>Vibrio parahaemolyticus</i> NCTC 10885	10	
TMCS 084	<i>Yersinia enterocolitica</i> NCTC 12982	10 25	

Microbiological Lab Consumables

Anaerobic Box Systems



TMC 011

ANAEROBIC BOX-L

Transparent unbreakable poly box, Size: 28.5 cm x 22 cm x 11 cm | Capacity: 7 Ltr. Accessories required but not provided. Anaerobe Gas Pack 3.5 Ltr. (TMC 001D) 2 Packs and Anaerobe Indicator Tablet (TMC 001C) 1 Pack

Pack: Nos



TMC 012

ANAEROBIC BOX-S

Transparent unbreakable poly box, Size: 24 cm x 17.7 cm x 8 cm | Capacity: 3.5 Ltr. Accessories required but not provided. Anaerobe Gas Pack 3.5 Ltr. (TMC 001D) 1 Packs and Anaerobe Indicator Tablet (TMC 001C) 1 Pack

Pack: 1 Nos

Anaerobic Jar Systems



TMC 001

ANAEROBIC JAR SYSTEM

(3.5 Ltr) Poly-carbonate Jar with sturdy Aluminum Lid Clamp and Sealing Ring With built-in safety features, Pressure Valve With Safety Valve & Two Way Pressure Gauge, Petri Plate Carrier, Indicator Tablets & Gas Pack

Pack: 1 Nos



TMC 013

ANAEROBIC JAR SYSTEM

(3.5 Ltr) Poly-carbonate Jar with sturdy Aluminum Lid Clamp and Sealing Ring With built-in safety features, Petri Plate Carrier, Indicator Tablets & Gas Pack

Pack: 1 Nos

Anaerobic Jar System Accessories

CODE	PRODUCT NAME		PACK SIZE
TMC 001E	POLYCARBONATE JAR (WITHOUT PRESSURE GAUGE) Anaerobic system consists of plastic container	 	1 Nos
TMC 001A	POLYCARBONATE JAR (WITH PRESSURE GAUGE) Anaerobic system consists of plastic container		1 Nos
TMC 001B	PETRI PLATES CARRIER Stainless Steel Carrier for 10 petri plates		1 Nos
TMC 001C	ANAEROBE INDICATOR TABLET Indicator for detection of anaerobic condition in the jar		1*5 Nos 1*10 Nos
TMC 001D	ANAEROBE GAS PACK Carbon dioxide generating pack	 	1*5 Nos

Microbiological Lab Consumable

Spreader/Needle/Inoculating Loops



TMC 002

L - SHAPED SPREADERS (P.P.)

For spreading bacterial samples on agar plates

Pack: 5*10 Nos



TMC 003

STERILE NEEDLE (P. P.)

Straight loop, for stabbing into medium

Pack: 5*10 Nos



TMC 004

STERILE DISPOSABLE LOOP (P. P.)

Size: Inner Diameter 0.7 mm, Outer Diameter 2.1 mm, Loop calibrated to 1 µl

Pack: 5*10 Nos

Pack: 1*50 Nos



TMC 005

STERILE DISPOSABLE LOOP (P. P.)

Size: Inner Diameter 3.5 mm, Outer Diameter 6.0 mm, Loop calibrated to 10 µl

Pack: 5*10 Nos

Pack: 1*50 Nos



TMC 006

STERILE DISPOSABLE LOOP (ABS)

Size: Inner Diameter 0.7 mm, Outer Diameter 2.1 mm, Loop calibrated to 1 µl

Pack: 5*10 Nos



TMC 007

STERILE DISPOSABLE LOOP (ABS)

Size: Inner Diameter 3.5 mm, Outer Diameter 6.0 mm, Loop calibrated to 10 µl

Pack: 5*10 Nos



TMC 008

METAL LOOP HOLDER (WITHOUT LOOP)

Size: Length 200 mm

Pack: 10 Nos



TMC 009

NICHROME LOOP

Size: Replacement Dia 2.5 mm loop end, length 50 mm

Pack: 50 Nos




TMC 047

NICHROME LOOP

Size: Dia 1.3 double wound, calibrated to 1.0 µl

Pack: 50 Nos

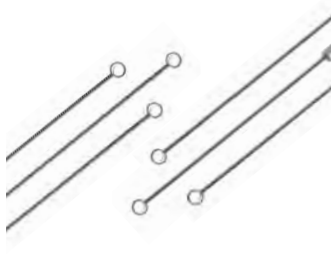

TMC 048

NICHROME LOOP

Size: Dia 3 mm, double wound calibrated to 0.06 ml

Pack: 50 Nos



TMC 049

NICHROME LOOP

Size: Dia 4 mm, double wound calibrated to 0.01 ml

Pack: 50 Nos



TMC 050

ASSORTED NICHROME LOOP

Each packet contains 3 pieces of each Nichrome loop having diameter 4 mm (0.01 ml), 2.5 mm (0.05 ml), 1.3 mm (1 µl) and straight wire, double wound


Pack: 50 Nos



Petri Plates

Reliable Surfaces for Reliable Results

- Crafted from Medical-grade optically clear polystyrene
- Available in standard (90 mm) and speciality (55 mm, 100 mm, 120 mm, and 150 mm) sizes
- Vented and non-vented lid options
- Pre-sterilized by gamma irradiation
- Triple-vented design for better air exchange
- Stackable design for optimized incubation space
- Available in single or compartment (Bi/Tri/Quad) formats

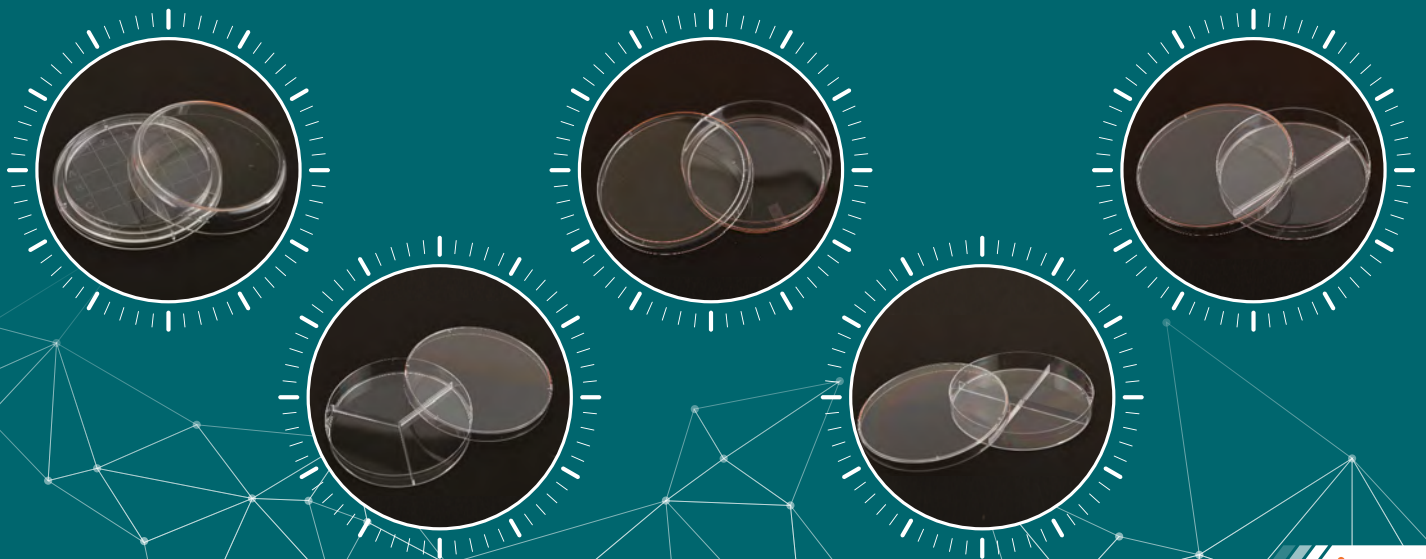

Microbiological Lab Consumables

Petri Plates (Culture Vessels & Accessories)

CODE	PRODUCT NAME	PACK SIZE
55 mm Petri Plates		
PP 008	STERILE DISPOSABLE SCORED TOP PETRI PLATES Polystyrene, Optically Clear, Flat Bottom, Clean Room Size : 65 mm x 12 mm Pack of 10 Plates	10*50 Nos
PP 007	STERILE DISPOSABLE SCORED TOP PETRI PLATES (γ-Irradiated) Polystyrene, Optically Clear, Flat Bottom (γ -Irradiated) Size : 65 mm x 12 mm Pack of 10 Plates	10*50 Nos
PP 010	STERILE DISPOSABLE SCORED TOP PETRI PLATES Polystyrene, Optically Clear, Domed Bottom, Clean Room Size : 65 mm x 12 mm Pack of 10 Plates	10*50 Nos
PP 009	STERILE DISPOSABLE SCORED TOP PETRI PLATES (γ-Irradiated) Polystyrene, Optically Clear, Domed Bottom (γ -Irradiated) Size : 65 mm x 12 mm Pack of 10 Plates	10*50 Nos
90 mm Petri Plates		
PP 004	STERILE DISPOSABLE PETRI PLATES Polystyrene, Optically Clear, Clean Room Size : 90 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 005	STERILE DISPOSABLE PETRI PLATES Polystyrene, Optically Clear, Clean Room Size : 90 mm Diameter x 15 mm Pack of 10 Plates	10*60 Nos
PP 001	STERILE DISPOSABLE PETRI PLATES (γ-Irradiated) Polystyrene, Optically Clear (γ -Irradiated) Size : 90 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 002	STERILE DISPOSABLE PETRI PLATES (γ-Irradiated) Polystyrene, Optically Clear (γ -Irradiated) Size : 90 mm Diameter x 15 mm Pack of 10 Plates	10*60 Nos
90 mm Petri Plates ECO		
PP 026	STERILE DISPOSABLE PETRI PLATES ECO Polystyrene, Optically Clear, Clean Room Size : 90 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 027	STERILE DISPOSABLE PETRI PLATES ECO Polystyrene, Optically Clear, Clean Room Size : 90 mm Diameter x 15 mm Pack of 10 Plates	10*50 Nos
PP 028	STERILE DISPOSABLE PETRI PLATES ECO (γ-Irradiated) Polystyrene, Optically Clear, Clean Room (γ -Irradiated) Size : 90 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 029	STERILE DISPOSABLE PETRI PLATES ECO (γ-Irradiated) Polystyrene, Optically Clear, Clean Room (γ -Irradiated) Size : 90 mm Diameter x 15 mm Pack of 10 Plates	10*50 Nos
90 mm Petri Plates Non-vented		
PP 032	STERILE DISPOSABLE PETRI PLATES NON-VENTED Polystyrene, Optically Clear, Clean Room Size : 90 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 033	STERILE DISPOSABLE PETRI PLATES NON-VENTED Polystyrene, Optically Clear, Clean Room Size : 90 mm Diameter x 15 mm Pack of 10 Plates	10*50 Nos
PP 030	STERILE DISPOSABLE PETRI PLATES NON-VENTED (γ-Irradiated) Polystyrene, Optically Clear, Clean Room (γ -Irradiated) Size : 90 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 031	STERILE DISPOSABLE PETRI PLATES NON-VENTED (γ-Irradiated) Polystyrene, Optically Clear, Clean Room (γ -Irradiated) Size : 90 mm Diameter x 15 mm Pack of 10 Plates	10*50 Nos
100 mm Petri Plates		
PP 011	STERILE DISPOSABLE PETRI PLATES Polystyrene, Optically Clear, Clean Room Size : 100 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 012	STERILE DISPOSABLE PETRI PLATES (γ-Irradiated) Polystyrene, Optically Clear, Clean Room (γ -Irradiated) Size : 100 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 013	STERILE DISPOSABLE PETRI PLATES Polystyrene, Optically Clear, Clean Room Size : 100 mm Diameter x 15 mm Pack of 10 Plates	10*50 Nos
PP 025	STERILE DISPOSABLE PETRI PLATES (γ-Irradiated) Polystyrene, Optically Clear, Clean Room (γ -Irradiated) Size : 100 mm Diameter x 15 mm Pack of 10 Plates	10*50 Nos

Microbiological Lab Consumables

Petri Plates (Culture Vessels & Accessories)

CODE	PRODUCT NAME	PACK SIZE
120 mm Petri Plates		
PP 014	STERILE DISPOSABLE SQUARE PETRI PLATES Polystyrene, Optically Clear, transparent useful for carrying out bioassay and antibacterial sensitivity test, Clean Room Size: 120 mm x 120 mm Individually Packed	1*100 Nos
150 mm Petri Plates		
PP 015	STERILE DISPOSABLE PETRI PLATES Polystyrene, Optically Clear, transparent useful for carrying out bioassay and antibacterial sensitivity test, Clean Room Size: 150 mm Diameter x 20 mm Individually Packed	1*100 Nos
PP 006	STERILE DISPOSABLE PETRI PLATES Polystyrene, Optically Clear, transparent useful for carrying out bioassay and antibacterial sensitivity test, Clean Room Size: 150 mm Diameter x 20 mm Set of 10 Plates	10*10 Nos
PP 018	STERILE DISPOSABLE PETRI PLATES (γ-Irradiated) Polystyrene, Optically Clear, transparent useful for carrying out bioassay and antibacterial sensitivity test, Clean Room (γ -Irradiated) Size: 150 mm Diameter x 20 mm Individually Packed	1*100 Nos
PP 003	STERILE DISPOSABLE PETRI PLATES (γ-Irradiated) Polystyrene, Optically Clear, transparent useful for carrying out bioassay and antibacterial sensitivity test, Clean Room (γ -Irradiated) Size: 150 mm Diameter x 20 mm Set of 10 Plates.	10*10 Nos
Bi-Petri Plates		
PP 016	STERILE DISPOSABLE PETRI PLATES (BI-PETRI PLATES) (Divided Petri Plates) Polystyrene, Optically clear, having partition, useful for culturing work Size: 90 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 017	STERILE DISPOSABLE PETRI PLATES (BI-PETRI PLATES) (Divided Petri Plates) Polystyrene, Optically clear, having partition, useful for culturing work Size: 90 mm Diameter x 15 mm Set of 10 Plates	10*50 Nos
PP 019	STERILE DISPOSABLE PETRI PLATES (BI-PETRI PLATES) (γ-Irradiated) (Divided Petri Plates) Polystyrene, Optically clear, having partition, useful for culturing work (γ -Irradiated) Size: 90 mm Diameter x 15 mm Individually Packed	1*500 Nos
PP 020	STERILE DISPOSABLE PETRI PLATES (BI-PETRI PLATES) (γ-Irradiated) (Divided Petri Plates) Polystyrene, Optically clear, having partition, useful for culturing work (γ -Irradiated) Size: 90 mm Diameter x 15 mm Set of 10 Plates	10*50 Nos
Tri-Petri Plates		
PP 021	STERILE DISPOSABLE PETRI PLATES (TRI-PETRI PLATES) (Divided Petri Plates) Polystyrene, Optically clear, having partition, useful for culturing work Size: 90 mm Diameter x 15 mm Set of 10 Plates	10*50 Nos
PP 022	STERILE DISPOSABLE PETRI PLATES (TRI-PETRI PLATES) (γ-Irradiated) (Divided Petri Plates) Polystyrene, Optically clear, having partition, useful for culturing work (γ -Irradiated) Size: 90 mm Diameter x 15 mm Set of 10 Plates	10*50 Nos
Quad-Petri Plates		
PP 023	STERILE DISPOSABLE PETRI PLATES (QUAD-PETRI PLATES) (Divided Petri Plates) Polystyrene, Optically clear, having partition, useful for culturing work Size: 90 mm Diameter x 15 mm Set of 10 Plates	10*50 Nos
PP 024	STERILE DISPOSABLE PETRI PLATES (QUAD-PETRI PLATES) (γ-Irradiated) (Divided Petri Plates) Polystyrene, Optically clear, having partition, useful for culturing work (γ -Irradiated) Size: 90 mm Diameter x 15 mm Set of 10 Plates	10*50 Nos
Mono-Petri Plates		
PP 034	STERILE DISPOSABLE PETRI PLATES (MONO PLATES) Polystyrene, Optically Clear, Clean Room Size : 60 mm x 15 mm Set of 10 Plates	10*50 Nos

Microbiological Lab Consumables

Stainless Steel Petri Plate Carriers



TMC 051
PETRI PLATE CARRIER W/ HANDLE
 Made out of S.S 304, **Size:** Diameter 90 mm petri plates for 12 plates
Pack: 1 Nos



TMC 052
PETRI PLATE CARRIER W/ HANDLE
 Made out of S.S 304, **Size:** Diameter 55 mm petri plates for 15 plates
Pack: 1 Nos

Durham Tubes



TMC 053
DURHAM TUBES
 Neutral glass, autoclavable, Length: 25 mm-27 mm, Diameter: 6-7 mm,
Pack: 1*100 Nos

Maccartney Bottle



TMC 032
MCCARTNEY BOTTLE W/ ALUMINIUM CAP
 Neutral glass, autoclavable
Pack: 1*100 Nos

Disposable Sterile Bottles



TMC 054
DISPOSABLE STERILE MEDIA BOTTLES 125 ML
 Alternative to Glass Media Bottles, Square (PETG), graduated bottles, sterilized by γ -Irradiation
Pack: 1*50 Nos
Pack: 1*100 Nos



TMC 055
DISPOSABLE STERILE MEDIA BOTTLES 250 ML
 Alternative to Glass Media Bottles, Square (PETG), graduated bottles, sterilized by γ -Irradiation
Pack: 1*50 Nos
Pack: *100 Nos



TMC 056
DISPOSABLE STERILE MEDIA BOTTLES 500 ML
 Alternative to Glass Media Bottles, Square (PETG), graduated bottles, sterilized by γ -Irradiation
Pack: 1*50 Nos
Pack: 1*100 Nos

Microbiological Lab Consumables


High Temperature Liquid Sample Bottles

CODE	PRODUCT NAME	PACK SIZE
TMC 077	STERILE HIGH TEMPERATURE LIQUID SAMPLE BOTTLES (SQUARE) 250 ml (γ -Irradiated)	1*50 Nos 1*100 Nos
TMC 078	STERILE HIGH TEMPERATURE LIQUID SAMPLE BOTTLES (SQUARE) 500 ml (γ -Irradiated)	1*50 Nos 1*100 Nos



Swabs

CODE	PRODUCT NAME	PACK SIZE
TMC 036	STERILE COTTON SWAB W/ TUBE in screw capped polypropylene tube, cotton bud w/ polypropylene stick, Size: 75 mm, Packed in 12 mm diameter tube, Individually packed	1*100 Nos
TMC 033	STERILE COTTON SWAB W/ TUBE in polypropylene tube, cotton bud w/ polypropylene stick, Size: 150 mm, Packed in 12 mm diameter tube, Individually packed	1*100 Nos
TMC 037	STERILE COTTON SWAB with polypropylene stick, Size: 150 mm, Individually packed	1*100 Nos
TMC 038	CLEAN SWAB with polypropylene stick, recommended for surface sampling	1*50 Nos 1*100 Nos
TMC 039	STERILE CLEAN SWAB W/ TUBE in polypropylene tube, polyester swab w/ polypropylene stick, recommended for surface sampling, Individually packed	1*100 Nos
TMC 040	CLEAN SWAB with polypropylene stick, sterile (γ -Irradiated), recommended for surface sampling, Individually packed	1*100 Nos
TMC 152	CLEAN SWAB W/ BREAKPOINT Polyester foam swab, Size: 125 mm, Non-sterile	1*100 Nos
TMC 153	CLEAN SWAB W/ BREAKPOINT Polyester foam swab, Size: 125 mm, Individually packed (γ -Irradiated)	1*100 Nos
TMC 154	STERILE DRY SPONGE SWAB STICK Size: 1.5" x 3", Sponge stick	10*5 Nos
TMC 041	PURE VISCOSE SWAB with polypropylene stick, Size: 150 x 2.5 mm	1*500 Nos
TMC 035	STERILE PURE VISCOSE SWAB with polypropylene stick (γ -Irradiated), Size: 150 x 2.5 mm, Individually packed	1*500 Nos
TMC 034	STERILE NASOPHARYNGEAL FLOCKED NYLON SWAB WITH BREAKPOINT (SCORED) for nasopharyngeal or urethral specimen collection and viral transport	1*50 Nos 1*100 Nos 1*500 Nos
TMC 042	STERILE NASOPHARYNGEAL FLOCKED NYLON SWAB WITH BREAKPOINT (SCORED) (RNase / DNase free) recommended for viral transport	1*50 Nos
TMC 043	STERILE OROPHARYNGEAL FLOCKED NYLON SWAB WITH BREAKPOINT recommended for clinical sampling	1*50 Nos 1*100 Nos 1*500 Nos
TMC 044	STERILE OROPHARYNGEAL FLOCKED NYLON SWAB WITH BREAKPOINT (RNase / DNase FREE) recommended for viral transport	1*50 Nos
TMC 045	STERILE POLYESTER TIPPED SWAB recommended for tracheal specimen collection	1*50 Nos 1*500 Nos
TMC 046	STERILE POLYESTER SWAB WITH BREAKPOINT (SCORED) recommended for tracheal specimen collection	1*50 Nos
TMC 071	STERILE NYLON SWAB W/ TUBE in polypropylene tube, nylon bud w/ polypropylene stick, Size: 150 mm, packed in 12 mm diameter tube, Individually packed	1*100 Nos
TMC 074	STERILE NYLON SWAB FOR HPV SPECIMEN COLLECTION recommended for HPV specimen collection	1*50 Nos
TMC 075	STERILE VISCOSE SWAB W/ TUBE in polypropylene tube, cotton bud w/ polypropylene stick, Size: 150 mm, packed in 12 mm diameter tube	1*100 Nos

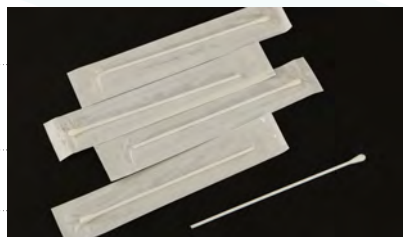


TMC 069 Filtration Funnel

Membrane Filtration Funnel is a ready-to-use, pre-sterilized, individually packed disposable filtration device to prevent cross-contamination of samples. These are easy to handle and speed up microbial analysis of up to 70% of water as well as products in the Pharmaceutical, Beverage, and Food processing industries.

Benefits

- Superior performance
- High total throughput
- Sterilized and ready-to-use
- Easy handling
- Time saving – up to 70%
- Converts into a Petri Dish
- Available with a Membrane Filter
- Gridded

How to use

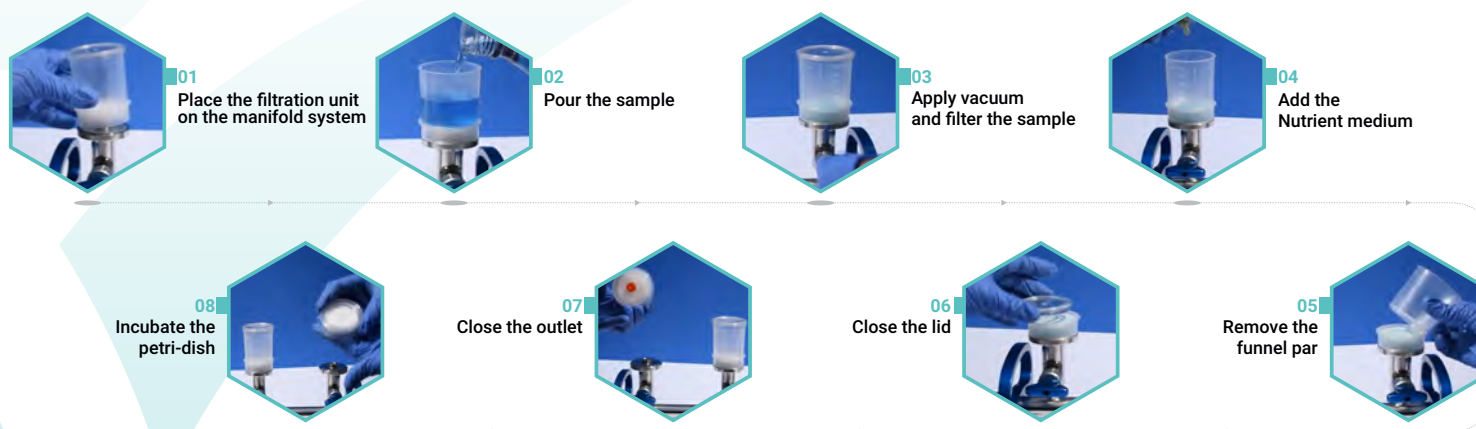

CODE	PRODUCT NAME	PACK SIZE
TMC 063	FILTRATION FUNNEL (PC) Dia: 47 mm, Funnel Capacity: 300 ml	1 Unit

CODE	PRODUCT NAME	PACK SIZE
TMC 064	FILTRATION FUNNEL (PP) Dia: 47 mm, Funnel Capacity: 300 ml	1 Unit

Membrane Filtration Funnels

CODE	PRODUCT NAME	PACK SIZE
TMC 065	MEMBRANE FILTRATION FUNNEL W/ CN GRIDDED MEMBRANE Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed, Funnel capacity: 100 ml, Pack size of 100 Nos	100 Nos
TMC 066	MEMBRANE FILTRATION FUNNEL W/ MCE GRIDDED MEMBRANE Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed, Funnel capacity: 100 ml, Pack size of 100 Nos	100 Nos
TMC 067	MEMBRANE FILTRATION FUNNEL W/ CN PLAIN MEMBRANE Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed, Funnel capacity: 100 ml, Pack size of 100 Nos	100 Nos

CODE	PRODUCT NAME	PACK SIZE
TMC 068	MEMBRANE FILTRATION FUNNEL W/ MCE PLAIN MEMBRANE Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed, Funnel capacity: 100 ml, Pack size of 100 Nos	100 Nos
TMC 069	MEMBRANE FILTRATION FUNNEL W/ CN GRIDDED MEMBRANE WITH ABSORBANT PAD Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed, Funnel capacity: 100 ml, Pack size of 100 Nos	100 Nos

Membrane Filters

CODE	PRODUCT NAME	PACK SIZE
TMC 125	STERILE MCE GRIDDED MEMBRANE Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed	100 Nos

CODE	PRODUCT NAME	PACK SIZE
TMC 127	STERILE MCE GRIDDED MEMBRANE Pore size: 0.22 µm, Dia of membrane: 47 mm, Individually packed	100 Nos

Microbiological Lab Consumables

Membrane Filters

CODE	PRODUCT NAME	PACK SIZE	
TMC 128	STERILE CN GRIDDED MEMBRANE Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed	100 Nos	
TMC 130	STERILE MCE PLAIN MEMBRANE Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed	100 Nos	
TMC 132	STERILE CN PLAIN MEMBRANE Pore size: 0.45 µm, Dia of membrane: 47 mm, Individually packed	100 Nos	

CODE	PRODUCT NAME	PACK SIZE	
TMC 129	STERILE CN GRIDDED MEMBRANE Pore size: 0.22 µm, Dia of membrane: 47 mm, Individually packed	100 Nos	
TMC 131	STERILE MCE PLAIN MEMBRANE Pore size: 0.22 µm, Dia of membrane: 47 mm, Individually packed	100 Nos	
TMC 133	STERILE CN PLAIN MEMBRANE Pore size: 0.22 µm, Dia of membrane: 47 mm, Individually packed	100 Nos	

Cold Pack



**TMC 060
ICE COLD PACK**
Keep colder for a long time
Content non-toxic and safe

Pack: 10*200 gm



**TMC 137
CLEANROOM PEN**
Black, γ-Irradiated

Pack: 10 Nos




**TMC 138
CLEANROOM MARKER**
Black, γ-Irradiated

Pack: 10 Nos

Marker Pens

Manifold Systems

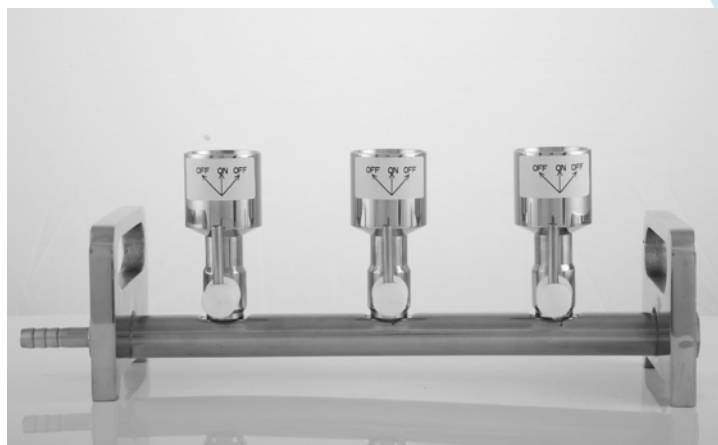

CODE	PRODUCT NAME	PACK SIZE	
TMC 070	SS. STERILITY TEST MANIFOLD SYSTEM Made of Stainless steel 316 Available with 3 stations with PTFE Valve	1 Unit	
TMC 126	SS. STERILITY TEST MANIFOLD SYSTEM Made of Stainless steel 316 Available with 6 stations with PTFE Valve	1 Unit	



Autoclave Bags

Designed for Biohazardous Waste Containment & Sterilization

Key Features

- High-grade polypropylene or polyethylene construction
- Withstands temperatures up to 135°C
- Available with biohazard symbols
- Leakproof and tear-resistant
- Self-standing bottom gusset design
- Multiple sizes for different waste volumes
- Color-coded for easy segregation
- Compatible with standard Autoclaves

How to Use

Step 01



Step 02



Step 03



Step 04



Microbiological Lab Consumables

Autoclave Accessories/Disposable Autoclave Bags

CODE	PRODUCT NAME	PACK SIZE	
	Clear Bags		
TMC 014	DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 11" X (B) 7", Clear, Transparent autoclavable disposable bag, Maximum weight holding capacity: 0.8 kg	250 Nos 500 Nos	
TMC 015	DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 10" X (B) 6", Clear, Transparent autoclavable disposable bag, Maximum weight holding capacity: 0.5 kg	250 Nos 500 Nos	
TMC 016	DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 12" X (B) 10", Clear, Transparent autoclavable disposable bag, Maximum weight holding capacity: 1 kg	250 Nos 500 Nos	
TMC 017	DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 20" X (B) 14", Clear, Transparent autoclavable disposable bag, Maximum weight holding capacity: 2.5 kg	250 Nos 500 Nos	
TMC 018	DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 32" X (B) 24", Clear, Transparent autoclavable disposable bag, Maximum weight holding capacity: 9.5 kg	250 Nos 500 Nos	
TMC 019	DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 40" X (B) 20", Clear, Transparent autoclavable disposable bag, Maximum weight holding capacity: 14 kg	250 Nos 500 Nos	
	Red Color Bags		
TMC 020	RED COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 11" X (B) 7", Red autoclavable disposable bag, Maximum weight holding capacity: 0.8 kg	250 Nos 500 Nos	
TMC 021	RED COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 10" X (B) 6", Red autoclavable disposable bag, Maximum weight holding capacity: 0.5 kg	250 Nos 500 Nos	
TMC 022	RED COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 12" X (B) 10", Red autoclavable disposable bag, Maximum weight holding capacity: 1 kg	250 Nos 500 Nos	
TMC 023	RED COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 20" X (B) 14", Red autoclavable disposable bag, Maximum weight holding capacity: 2.5 kg	250 Nos 500 Nos	
TMC 024	RED COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 32" X (B) 24", Red autoclavable disposable bag, Maximum weight holding capacity: 9.5 kg	250 Nos 500 Nos	
TMC 025	RED COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 40" X (B) 20", Red autoclavable disposable bag, Maximum weight holding capacity: 14 kg	250 Nos 500 Nos	
	Yellow Color Bags		
TMC 026	YELLOW COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 11" X (B) 7", Yellow autoclavable disposable bag, Maximum weight holding capacity: 0.8 kg	250 Nos 500 Nos	
TMC 027	YELLOW COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 10" X (B) 6", Yellow autoclavable disposable bag, Maximum weight holding capacity: 0.5 kg	250 Nos 500 Nos	
TMC 028	YELLOW COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 12" X (B) 10", Yellow autoclavable disposable bag, Maximum weight holding capacity: 1 kg	250 Nos 500 Nos	
TMC 029	YELLOW COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 20" X (B) 14", Yellow autoclavable disposable bag, Maximum weight holding capacity: 2.5 kg	250 Nos 500 Nos	
TMC 030	YELLOW COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 32" X (B) 24", Yellow autoclavable disposable bag, Maximum weight holding capacity: 9.5 kg	250 Nos 500 Nos	
TMC 031	YELLOW COLOR DISPOSABLE AUTOCLAVABLE BAGS Size: (H) 40" X (B) 20", Yellow autoclavable disposable bag, Maximum weight holding capacity: 14 kg	250 Nos 500 Nos	



Pack of 250 / 500 Nos in multiple packing of 50 bags



Pack of 250 / 500 Nos in multiple packing of 50 bags



Pack of 250 / 500 Nos in multiple packing of 50 bags



Precision at Every Drop

- Manufactured using high-quality, virgin polypropylene
- Universal compatibility (fits most Micropipettes)
- DNase, RNase, and pyrogen-free
- Crystal-clear for easy volume verification
- Available in filtered & non-filtered variants
- Sterile and non-sterile options
- Autoclavable and low-retention designs

CODE	PRODUCT NAME	PACK SIZE
Tips with Filter		
TMC 079	10 µl MICROPIPETTE TIPS, WITH FILTER Max capacity 10 µl	1*96 Nos 10*96 Nos
TMC 080	10 µl MICROPIPETTE TIPS, WITH FILTER Max capacity 20 µl	1*96 Nos 10*96 Nos
TMC 081	10 µl MICROPIPETTE LONG TIPS, WITH FILTER Max capacity 20 µl	1*96 Nos 10*96 Nos
TMC 082	20 µl MICROPIPETTE TIPS, WITH FILTER Max capacity 20 µl	1*96 Nos 10*96 Nos
TMC 083	100 µl MICROPIPETTE TIPS, WITH FILTER Max capacity 100 µl	1*96 Nos 10*96 Nos
TMC 084	200 µl MICROPIPETTE TIPS, WITH FILTER Max capacity 200 µl	1*96 Nos 10*96 Nos
TMC 085	300 µl MICROPIPETTE TIPS, WITH FILTER Max capacity 300 µl	1*96 Nos 10*96 Nos
TMC 086	1000 µl MICROPIPETTE TIPS, WITH FILTER Max capacity 1000 µl	1*96 Nos 10*96 Nos
TMC 087	1000 µl MICROPIPETTE TIPS, WITH FILTER Max capacity 1200 µl	1*96 Nos 10*96 Nos
Tips without Filter		
TMC 088	10 µl MICROPIPETTE TIPS Max capacity 10 µl	1*96 Nos 10*96 Nos
TMC 089	200 µl MICROPIPETTE TIPS Max capacity 200 µl	1*96 Nos 10*96 Nos
TMC 090	1000 µl MICROPIPETTE TIPS Max capacity 1000 µl	1*96 Nos 10*96 Nos
Tips in Bulk Pack		
TMC 091	200 µl MICROPIPETTE TIPS Max capacity 200 µl	1*1000 Nos
TMC 092	1000 µl MICROPIPETTE TIPS Max capacity 1000 µl	1*1000 Nos

Transfer Pipettes

CODE	PRODUCT NAME	PACK SIZE
TMC 093	TRANSFER PIPETTES 3 ML INDIVIDUALLY WRAP STERILE (Pasteur pipette / Dropper)	1*50 Nos 1*500 Nos



TMFlex™ Bags

Designed to Streamline Lab sampling

Available in various sizes and capacities to meet a wide range of laboratory sampling needs. From compact sizes for small samples to larger sizes for bulk collection, our range ensures strength and reliability.

TMFlex™ Bags are essential in laboratories for:

Sample Collection: Collect and store samples conveniently for further analysis.

Transportation: Safely move samples between lab areas or external facilities.

Storage: Securely store samples for extended periods, maintaining integrity and quality.

Why choose TMFlex™ Bags

Comprehensive Range: Various sizes and capacities accommodate your different sampling needs.

TMFlex™ Technology: Ensure durability, flexibility, and reliability in sample containment.

Specialized Variants: Write-On Bags, Stand-Up Bags, and Homogenizer Blender Filter Bags, providing tailored solutions for specific sampling requirements.

Sizes Available: Choose from sizes ranging from 60 ml to over 5000 ml, ensuring there's a bag for every sampling requirement.

How to Use

Step 01



Step 02



Step 03



Step 04



Microbiological Lab Consumables

TMFlex™ Sampling Bags/Whirl Bags

CODE	PRODUCT NAME	PACK SIZE	
	Homogenizer Blender Filter Sampling Bag		
TMC 115	TMFlex™ HOMOGENIZER BLENDER FILTER SAMPLING BAG Size: 15 x 23 cm, 710 ml capacity	1*100 Nos	
TMC 116	TMFlex™ HOMOGENIZER BLENDER FILTER SAMPLING BAG Size: 19 x 30 cm, 1627 ml capacity	1*100 Nos	
	Sampling Bag		
TMC 094	TMFlex™ SAMPLING BAG Size: 7.5 x 12.5 cm, 60 ml capacity	1*100 Nos	
TMC 095	TMFlex™ SAMPLING BAG Size: 7.5 x 18.5 cm, 120 ml capacity	1*100 Nos	
TMC 096	TMFlex™ SAMPLING BAG Size: 9.5 x 18 cm, 207 ml capacity	1*100 Nos	
TMC 097	TMFlex™ SAMPLING BAG Size: 13 x 19 cm, 384 ml capacity	1*100 Nos	
TMC 098	TMFlex™ SAMPLING BAG Size: 11.5 x 23 cm, 540 ml capacity	1*100 Nos	
TMC 099	TMFlex™ SAMPLING BAG Size: 15 x 23 cm, 720 ml capacity, Thickness: 0.07 mm	1*100 Nos	
TMC 100	TMFlex™ SAMPLING BAG Size: 15 x 23 cm, 720 ml capacity, Thickness: 0.1 mm	1*100 Nos	
TMC 101	TMFlex™ SAMPLING BAG Size: 15 x 38 cm, 1242 ml capacity	1*100 Nos	
TMC 102	TMFlex™ SAMPLING BAG Size: 19 x 30 cm, 1650 ml capacity	1*100 Nos	
TMC 103	TMFlex™ SAMPLING BAG Size: 19 x 38 cm, 2041 ml capacity	1*100 Nos	

CODE	PRODUCT NAME	PACK SIZE	
	Sampling Bag		
TMC 104	TMFlex™ SAMPLING BAG Size: 25.4 x 38 cm, 2721 ml capacity	1*100 Nos	
TMC 105	TMFlex™ SAMPLING BAG Size: 38 x 50.8 cm, 5441 ml capacity	1*100 Nos	
	Stand-up Sampling Bag		
TMC 112	TMFlex™ STAND-UP SAMPLING BAG Size: 11.5 x 23 cm, 532 ml capacity	1*100 Nos	
TMC 113	TMFlex™ STAND-UP SAMPLING BAG Size: 7.5 x 18.5 cm, 118 ml capacity	1*100 Nos	
TMC 114	TMFlex™ STAND-UP SAMPLING BAG Size: 15 x 23 cm, 710 ml capacity	1*100 Nos	
	Write-on Sampling Bag		
TMC 106	TMFlex™ WRITE-ON SAMPLING BAG Size: 11.5 x 23 cm, 532 ml capacity	1*100 Nos	
TMC 107	TMFlex™ WRITE-ON SAMPLING BAG, YELLOW TAPE Size: 15 x 23 cm, 710 ml capacity	1*100 Nos	
TMC 108	TMFlex™ WRITE-ON SAMPLING BAG Size: 15 x 23 cm, 710 ml capacity	1*100 Nos	
TMC 109	TMFlex™ WRITE-ON SAMPLING BAG Size: 7.5 x 12.5 cm, 58 ml capacity	1*100 Nos	
TMC 110	TMFlex™ WRITE-ON SAMPLING BAG Size: 7.5 x 18.5 cm, 118 ml capacity	1*100 Nos	
TMC 134	TMFlex™ WRITE-ON SAMPLING BAG Size: 9.5 x 18 cm, 207 ml capacity	1*100 Nos	
TMC 111	TMFlex™ WRITE-ON SAMPLING BAG Size: 19 x 30 cm, 1627 ml capacity	1*100 Nos	

Microbiological Lab Consumables

Centrifuge Tubes

CODE	PRODUCT NAME	PACK SIZE
15 ml Tubes		
TMC 120	CENTRIFUGE TUBE, 15 ML, NON STERILE, AUTOCLAVABLE, CONICAL BOTTOM, RCF VALUE: 17000G	4*25 Nos 20*25 Nos
TMC 119	CENTRIFUGE TUBE, 15 ML, STERILE, AUTOCLAVABLE, CONICAL BOTTOM, RCF VALUE: 17000G	4*25 Nos 1*100 Nos 20*25 Nos
50 ml Tubes		
TMC 122	CENTRIFUGE TUBE, 50 ML, NON STERILE, AUTOCLAVABLE, CONICAL BOTTOM, RCF VALUE: 20000G	4*25 Nos 20*25 Nos
TMC 124	CENTRIFUGE TUBE, 50 ML, SELF STANDING, NON STERILE, AUTOCLAVABLE, CONICAL BOTTOM, RCF VALUE: 5000G	4*25 Nos 20*25 Nos
TMC 123	CENTRIFUGE TUBE, 50 ML, SELF STANDING, STERILE, AUTOCLAVABLE, CONICAL BOTTOM, RCF VALUE: 5000G	4*25 Nos 1*100 Nos 20*25 Nos
TMC 121	CENTRIFUGE TUBE, 50 ML, STERILE, AUTOCLAVABLE, CONICAL BOTTOM, RCF VALUE: 20000G	4*25 Nos 1*100 Nos 20*25 Nos
Micro Centrifuge Tubes		
TMC 117	MICRO CENTRIFUGE TUBES, AUTOCLAVABLE, WITH GRADUATION, CONICAL BOTTOM, CAPACITY 1.5 ML	1*500 Nos
TMC 118	MICRO CENTRIFUGE TUBES, AUTOCLAVABLE, WITH GRADUATION, 'U' SHAPED BOTTOM, CAPACITY 2.0 ML	1*500 Nos

Cleanroom Tapes & Location Markers

CODE	PRODUCT NAME	PACK SIZE
TMC 139	CLEANROOM MASKING TAPE Size & Length: 18 mm x 50 Mtr, γ -Irradiated, Individually packed	10 Roll
TMC 140	CLEANROOM MASKING TAPE Size & Length: 18 mm x 50 Mtr	10 Roll
TMC 141	CLEANROOM FLOOR MARKING TAPES Size & Length: 24 mm x 20 Mtr-1 inch, γ -Irradiated, Individually packed, Yellow	10 Roll
TMC 142	CLEANROOM FLOOR MARKING TAPES Size & Length: 24 mm x 20 Mtr-1 inch, γ -Irradiated, Individually packed, Red	10 Roll
TMC 143	CLEANROOM FLOOR MARKING TAPES Size & Length: 24 mm x 20 Mtr-1 inch, γ -Irradiated, Individually packed, Green	10 Roll
TMC 144	CLEANROOM FLOOR MARKING TAPES Size & Length: 24 mm x 20 Mtr-1 inch, γ -Irradiated, Individually packed, Yellow & Black	10 Roll
TMC 145	CLEANROOM ENVIRONMENTAL MONITORING LOCATION MARKERS WITHOUT NUMBERS Size: 100 mm, γ -Irradiated, Individually packed, Yellow	50 Nos
TMC 146	CLEANROOM ENVIRONMENTAL MONITORING LOCATION MARKERS WITHOUT NUMBERS Size: 100 mm, γ -Irradiated, Individually packed, Red	50 Nos
TMC 147	CLEANROOM ENVIRONMENTAL MONITORING LOCATION MARKERS WITHOUT NUMBERS Size: 50 mm, γ -Irradiated, Individually packed, Yellow	50 Nos
TMC 148	CLEANROOM ENVIRONMENTAL MONITORING LOCATION MARKERS WITHOUT NUMBERS Size: 50 mm, γ -Irradiated, Individually packed, Red	50 Nos

Microbiological Lab Consumables

Receiver Flasks

CODE	PRODUCT NAME	PACK SIZE	
TMC 149	VACUUM FLASK Moc: Glass, Volume: 1000 ml	1 Nos	

CODE	PRODUCT NAME	PACK SIZE	
TMC 150	VACUUM FLASK Moc: Glass, Volume: 2000 ml	1 Nos	

Spatulas

CODE	PRODUCT NAME	PACK SIZE	
TMC 135	SPATULA PP, 200 mm, Non Sterile	50 Nos	

CODE	PRODUCT NAME	PACK SIZE	
TMC 136	SPATULA PP, 200 mm, Sterile	50 Nos	

Silicon Cork

CODE	PRODUCT NAME	PACK SIZE	
TMC 151	SILICON CORK Size: No #8	10 Nos	



Syringe Filters

CODE	PRODUCT NAME	PACK SIZE	
	Cellulose Acetate		
TMC 155	SYRINGE FILTER Cellulose Acetate, Pore size: 0.2 µm, 33 mm diameter	100 Nos	
TMC 156	SYRINGE FILTER Cellulose Acetate, Pore size: 0.45 µm, 33 mm diameter	100 Nos	
	Nylon-66		
TMC 157	SYRINGE FILTER Nylon-66, Hydrophilic membrane, Pore size: 0.2 µm, 13 mm diameter	100 Nos	
TMC 158	SYRINGE FILTER Nylon-66, Hydrophilic membrane, Pore size: 0.2 µm, 25 mm diameter	100 Nos	
TMC 159	SYRINGE FILTER Nylon-66, Hydrophilic membrane, Pore size: 0.2 µm, 33 mm diameter	100 Nos	
TMC 160	SYRINGE FILTER Nylon-66, Hydrophilic membrane, Pore size: 0.45 µm, 25 mm diameter	100 Nos	
TMC 161	SYRINGE FILTER Nylon-66, Hydrophilic membrane, Pore size: 0.45 µm, 33 mm diameter	100 Nos	

CODE	PRODUCT NAME	PACK SIZE	
	PVDF		
TMC 162	SYRINGE FILTER PVDF, Hydrophilic membrane, Pore size: 0.2 µm, 13 mm diameter	100 Nos	
TMC 163	SYRINGE FILTER PVDF, Hydrophilic membrane, Pore size: 0.2 µm, 33 mm diameter	100 Nos	
TMC 164	SYRINGE FILTER PVDF, Hydrophilic membrane, Pore size: 0.45 µm, 13 mm diameter	100 Nos	
TMC 165	SYRINGE FILTER PVDF, Hydrophilic membrane, Pore size: 0.45 µm, 33 mm diameter	100 Nos	
	PTFE		
TMC 166	SYRINGE FILTER PTFE, Hydrophilic membrane, Pore size: 0.2 µm, 33 mm diameter	100 Nos	
TMC 167	SYRINGE FILTER PTFE, Hydrophilic membrane, Pore size: 0.45 µm, 33 mm diameter	100 Nos	
	PES		
TMC 168	SYRINGE FILTER PES, Hydrophilic membrane, Pore size: 0.2 µm, 33 mm diameter	100 Nos	
TMC 169	SYRINGE FILTER PES, Hydrophilic membrane, Pore size: 0.45 µm, 33 mm diameter	100 Nos	

Sterile Zip Lock Bags

Secure Every Sample. Preserve Every Result

From collection to storage, **Sterile Zip Lock Bags** are designed to protect sample integrity with dependable contamination control. Manufactured using high-quality sterile LDPE, they combine durability, clarity, and secure sealing for reliable laboratory handling.

Key Advantages

- Preserves sterility throughout handling and transport
- Durable, transparent, and easy to handle
- Strong, tear-resistant material for safe handling
- Leak-resistant zip lock for secure containment

Microbiological Lab Consumables

Zip Lock Bags

CODE	PRODUCT NAME	PACK SIZE
TMC 170	STERILE ZIP LOCK BAGS Size: 6" X 8"	100 Nos
TMC 171	STERILE ZIP LOCK BAGS Size: 8" X 12"	100 Nos

Deodorising Pearls

NEW

CODE	PRODUCT NAME	PACK SIZE
TMC 172	DEODORISING PEARLS	5*50 Nos 10*50 Nos

TME 006-INFRA-RED LOOP STERILIZER

- High sterilization temperature (up to 825°C)
- Adjustable time range
- Fast heating time
- Intelligent temperature controller

- Stainless steel body
- Flameless & pollution-free
- Compact & space-saving
- Lightweight




CODE	PRODUCT NAME	PACK SIZE
TME 006	INFRA-RED LOOP STERILIZER	1 Nos

TME 007-INFRA-RED STERILIZER FOR INOCULATION LOOPS

- Reaches 900–1300°C in just a few seconds
- Instant sterilization of metal tools like inoculation loops, needles, scissors, tweezers
- Automatic start with sensor detection
- Stops heating automatically after 7 seconds
- Compact and portable weigh 2.6 Kg only

- Ideal for labs where flammable gases are restricted
- Low energy uses with tungsten halogen lamp heating
- Works well in anaerobic chambers
- Durable quartz glass heating zone

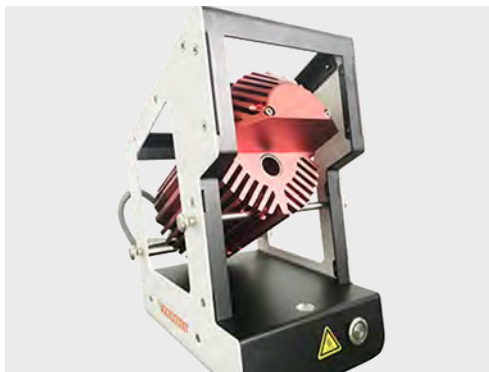

CODE	PRODUCT NAME	PACK SIZE
TME 007	INFRA-RED STERILIZER FOR INOCULATION LOOPS	1 Nos

Microbiological Lab Equipments

Microbial Air Monitoring Systems (Air Sampler) & Accessories

TME 003-MICROBIAL AIR MONITORING SYSTEM

- Validated as per BS EN ISO 14698-Part-I
- Special design prevents re-circulation of air for more accurate results
- Comes with built-in adaptors for 90 mm and 55 mm Petri plates
- Allows custom air volume settings as per your need
- Air Flow Rate: 100 litres per minute
- Easy to carry and handle
- Weigh only 2.68 Kg
- Strong and lightweight aviation-grade aluminium
- Battery runs up to 8 hours



CODE	PRODUCT NAME	PACK SIZE
TME 003	MICROBIAL AIR MONITORING SYSTEM (AIR SAMPLER)	1 Nos
	Add-on Accessories	
TME 003A	COMPRESSED GAS HEAD COMPONENT OF AN ENVIRONMENTAL MONITORING PROGRAM USED FOR EVALUATION OF COMPRESSED GAS	1 Nos
TME 003B	POROUS SAMPLING HEAD AND SAFETY LID	1 Nos
TME 003C	POWER CHARGER	1 Nos
TME 003D	BATTERY	1 Nos
TME 003E	COMPRESSED GAS HEAD COMPONENT WITH FLOW METER TO ACHIEVE THE REQUIRED AIR FLOW DURING SAMPLING	1 Nos

TME 004-MICROBIAL AIR MONITORING SYSTEM

- Validated as per BS EN ISO 14698-Part-I
- Connects to a PC to meet 21 CFR Part 11 guidelines
- Special design prevents recirculation of air for accurate results
- Comes with built-in adaptors for both 90 mm and 55 mm petri plates
- Nominal air flow rate: 100 liters per minute
- Easy to carry and handle weigh only 3 Kg
- Strong body made of anodized aluminium
- Battery runs up to 8 hours
- Comes in a high-strength carry case
- Infrared remote control for operating from a distance



CODE	PRODUCT NAME	PACK SIZE
TME 004	MICROBIAL AIR MONITORING SYSTEM (AIR SAMPLER)	1 Nos
	Add-on Accessories	
TME 004A	COMPRESSED GAS HEAD COMPONENT OF AN ENVIRONMENTAL MONITORING PROGRAM USED FOR EVALUATION OF COMPRESSED GAS	1 Nos
TME 004B	POROUS SAMPLING HEAD AND SAFETY LID	1 Nos
TME 004C	POWER CHARGER	1 Nos
TME 004D	BATTERY	1 Nos
TME 004E	REMOTE	1 Nos
TME 004F	CORD	1 Nos
TME 004G	COMPRESSED GAS HEAD COMPONENT WITH FLOW METER TO ACHIEVE THE REQUIRED AIR FLOW DURING SAMPLING	1 Nos

Microbiological Lab Equipments

Microbial Air Monitoring Systems (Air Sampler) & Accessories

TME 005-MICROBIAL AIR MONITORING SYSTEM

The equipment uses isokinetic sampling, which delivers uniform airflow across the surface of an agar plate at high velocity for quick analysis. A specified volume of air is passed through small holes or slits onto a Petri Plate containing Culture Medium for a specified period of time. After the cycle is complete, the plate is removed from the Air sampler and incubated, and visible colonies are counted to determine the level of contamination. This will result in the number of viable microorganisms present in each liter of air in the area being calculated.



Easy-to-Use



Lightweight
& Portable



Unidirectional
airflow



Suitable for
90 mm plates



Easily export
& print data



Individual user
credentials



Can store up to
10,000 entry logs



3-level authority
management



Validated as per
BS EN ISO:14698-Part-1



Easily customize sample
air volumes or choose
from pre-sets



PC connectivity to meet
21 CFR:Part-11 requirements



SPECIFICATIONS

Nominal Air Flow Rate	100 liters of air per minute
Air Sample Volume	1- 9999 liters can be selected
Recommended Volume of Air per Cycle	1000 liters
Recipe Storage	Upto 6 recipes can be created and stored
User Configuration	Multiple user support with individual username and password
Sample Data Records	10,000 data can be stored
Software Compliance	Comply with 21CFR Part 11
Software Communication	USB output, to wireless data export through Wi-Fi

Levels of Authority	3 Levels of authority management
Audit Trail	Audit trail function available
Alarm Logs	Completed cycle data can be reviewed
Plate Size Compatible	90±5 mm
Battery Pack	Rechargeable Li-Ion (without memory effect). After charging, the battery will work up to 8 hours
Weight of Instrument	2.5 Kg
Add-on Accessories	Bluetooth printer (Range within 10 m)
Accessories	Empty Petri-plate, Air sampler Case, Power Charger

Microbiological Lab Equipments

Petri Dish Revolving Table

TME 008-PETRI DISH REVOLVING TABLE-MANUAL



- Specially designed for even coating of Agar plates
- Ideal for use in Microbiology and Molecular biology labs
- Manual operation with a smooth and stable rotating base
- Comes with 2 surfaces to support Petri Dishes of different sizes (60-150 mm diameter)
- Made for both Glass and Plastic Petri Dishes
- Rubber non-slip mat ensures Petri Dishes stay secure during rotation

- Lightweight (2.4 Kg) and durable for regular lab use
- Allows easy and uniform spreading of samples using a spreader
- Compact and easy to use on any flat lab bench
- Supports quick and efficient plate preparation for Culture Media
- Double-sided design allows both sides of the surface to be used

CODE	PRODUCT NAME	PACK SIZE
TME 008	PETRI DISH REVOLVING TABLE-MANUAL	1 Nos

TME 009-PETRI DISH REVOLVING TABLE-ELECTRIC



- Uses infrared sensor technology for contact-free, hand-controlled operation
- Helps with uniform spreading of bacteria on agar media during inoculation
- Adjustable rotation timing from 30 to 180 seconds or continuous operation mode
- Rotates from 25 rpm to 150 rpm, ensuring even and consistent sample coating
- Compatible with Petri dishes sized 60–150 mm in diameter

- Can be autoclaved, making it safe for repeated use in sterile environments
- Ideal for Clinical, Research, and Industrial microbiology labs
- Comes with 2 work plates to support Petri dishes of different sizes
- Lightweight (1.2 Kg) and durable for regular lab use
- Equipped with a non-slip rubber lining for a secure grip of Petri dishes during spinning

CODE	PRODUCT NAME	PACK SIZE
TME 009	PETRI DISH REVOLVING TABLE-ELECTRIC	1 Nos

TME 010-COLONY COUNTER



- Digital Display Tool for accurate bacterial colony counting on Petri Dishes
- Clear LED Screen with 4-digit display and easy-to-read results
- Ring Fluorescent Light gives bright, even side lighting to highlight colonies clearly without glare
- Sensitive Probe Pen ensures precise colony marking and correction during counting
- Built-in Counting Pool with 118 mm diameter suitable for standard Petri Dishes
- Compact and light design, easy to place in labs with limited space
- Energy-saving lighting with brightness control for visual comfort and accuracy
- Magnifying Glass (90 mm) enhances view with 3x/6x magnification for better visibility
- Flexible counting options allows adding or removing colony numbers easily
- Sturdy base with anti-slip design keeps the counter stable during use
- Ideal for Clinical Labs, Pharmaceutical Companies, Food Testing, Water Testing, Biological Research, & Educational Institutes
- 100–240 V power supply

CODE	PRODUCT NAME	PACK SIZE
TME 010	COLONY COUNTER	1 Nos

Vacuum Pumps

CODE	PRODUCT NAME	PACK SIZE
TME 011	VACUUM PUMP 300 (DIAPHRAGM TYPE-OIL FREE)	1 Nos
TME 012	VACUUM PUMP 400 (DIAPHRAGM TYPE-OIL FREE)	1 Nos



Antibiotic Discs Dispenser

TME 014-TM OCTACCLICK DISPENSER (8 Position Antibiotic Disc Dispensing System)



Key Features

- 8-disc simultaneous dispensing for faster workflow
- Uniform and precise disc placement
- Single-hand ergonomic operation
- Self-tamping mechanism for proper agar contact
- Rotating turret for smooth alignment
- Clear cartridge labeling for easy identification
- Adjustable base ring for variable agar depth

Working Principle

- Load disc cartridges into slots until "click-fit"
- Place over inoculated Petri plate
- Press plunger in one smooth motion
- 8 discs are automatically embedded evenly
- Proceed for incubation as per AST protocol

Note:

Store in a protective container with desiccant at 2–8°C when not in use. Allow to reach room temperature before operation.

CODE	PRODUCT NAME	PACK SIZE
TME 014	TM OCTACCLICK DISPENSER	1 Nos



Sterilization Indicators

Evident Control of Sterilization Process

- Easy-to-Use & Interpret the result
- Reliable Sensitivity
- Accurate Performance
- Certified in accordance with ISO/IRAM/NS-EN/ANIMAT PM
- Trusted products

"Your Safety, Our Commitment"

(The most trusted partner for ensuring sterility)

Sterilization Indicators

Biological Indicators (Spore Strips)

CODE	PRODUCT NAME	PROCESS	PACK SIZE	
TBI 004	<i>Geobacillus stearothermophilus</i> 10 ⁶ Spore Strips 10 (ATCC 7953)	Steam Sterilization	1*100 Nos	

Biological Indicators (Spore Ampoules)

CODE	PRODUCT NAME	PROCESS	PACK SIZE	
TBI 006	<i>Geobacillus stearothermophilus</i> 10 ⁶ (Vol. 2 ml)	Steam Sterilization	1*50 Nos	

Biological Indicators (Self-Contained Biological Indicators-SCBI)

CODE	PRODUCT NAME	PROCESS	PACK SIZE	
TBI 010	<i>Geobacillus stearothermophilus</i> 10 ⁶	Steam Sterilization	1*100 Nos	

Sterilization Indicators

Chemical Indicators (Self Adhesive Sterilization Indicators for Radiation)

CODE	PRODUCT NAME	PROCESS	PACK SIZE
TCI 001	RADIATION CHEMICAL PROCESS INDICATOR (CLASS 1)	Plain Indicator Circular Dot (15 mm)	1*1000 Nos 5*1000 Nos

Chemical Indicators (Self Adhesive Sterilization Indicators for Ethylene Oxide)

CODE	PRODUCT NAME	PROCESS	PACK SIZE
TCI 002	ETHYLENE OXIDE (EO) CHEMICAL PROCESS INDICATOR (CLASS 1)	Plain Indicator Circular Dot (15 mm)	1*1000 Nos 5*1000 Nos

Chemical Indicators (Sterilization Indicators for Steam)

CODE	PRODUCT NAME	PROCESS	PACK SIZE
TCI 005	BOWIE & DICK TEST PACK (CLASS 2)	Steam Sterilization Process	1*1 Pack 1*20 Pack
TCI 006	STEAM CHEMICAL PROCESS INDICATOR (CLASS 1)	Plain Indicator Circular Dot (15 mm)	1*1000 Nos 5*1000 Nos
TCI 012	MOVING FRONT INTEGRATOR INDICATOR (CLASS 5)	Steam Sterilization Process	1*250 Nos
TCI 013	UNIQUE POINT CHEMICAL INDICATOR FOR STEAM STERILIZATION (TYPE 5, CLASS 1)	Steam Sterilization Process	200 Nos

Chemical Indicators Tapes

CODE	PRODUCT NAME	PROCESS	PACK SIZE
TCI 007	SELF ADHESIVE AUTOCLAVE TAPE (CLASS 1)	19 mm x 50 mm	1 Nos 1*5 Nos
TCI 008	SELF ADHESIVE ETO TAPE (CLASS 1)	19 mm x 50 mm	1 Nos

Chemical Indicators Test Strips

CODE	PRODUCT NAME	PROCESS	PACK SIZE
TCI 009	DRY HEAT INDICATOR STRIP (Class 4)	105 mm x 18 mm	1*250 Nos
TCI 014	STEAM INDICATOR STRIP (CLASS 4)	105 mm x 18 mm	2*250 Nos



Plant Tissue Culture Media & Ingredients

Fulfilling the Nutritional Needs of In-Vitro Plants

Plant Tissue Culture Media are specialized nutrient formulations that support the in vitro growth, development, and propagation of plant tissues. These formulations ensure optimal pH, nutrient balance, and sterility for successful Micro-propagation, Callus Formation, and Regeneration.

Features

- Reproducible results
- Batch-to-batch consistency
- In-house Manufacturing
- Fast solubility and high clarity
- Superior stability and shelf life
- Dust-free powder for easy handling
- Optimized powder flowability
- Tailored for different plant species and applications

Make Your Own PTCM

Macronutrients

Macronutrients are elements such as Nitrogen, Phosphorus, Potassium, Calcium, Magnesium, and Sulphur. These elements are required in large amounts (>0.5 mM) for robust plant structure and metabolism. They comprise at least 0.1% of the dry weight of plants.

Micronutrients

Micronutrients are elements such as Iron, Manganese, Zinc, Copper, Molybdenum, Boron, Iodine, and Cobalt. These trace elements (<0.5 mM) are essential for enzyme function and metabolic balance.

Vitamins

Vitamins are added to the plants in several forms and concentrations. Certainly, these compounds are essential for many biochemical reactions. In almost all media for plant cell and tissue culture, Thiamine (Vitamin B1) is included. Other frequently used Vitamins are Nicotinic Acid, Pyridoxine (B6) and Myo-inositol. Although Myo-inositol is a carbohydrate, not a vitamin, it has been shown to stimulate growth in certain cell cultures. It is hard to judge the virtual importance of Vitamins, as their effects on cells vary from species to species.

Gelling Agents

Gelling agents such as Agar and Gellan Gum provides the right consistency to the medium. Agar (0.5-1.0%) for firm gels Gellan Gum (1.5-2.5 g/L) for transparent, contamination-detectable media.


Plant Growth Regulators (PGRs)

Plant hormones are added to Plant Tissue Culture Media to regulate growth. In tissue culture, they are mainly used to stimulate adventitious regeneration of roots, shoots and embryos, outgrowth of axillary buds, and formation of callus. Moreover, cytokinin and auxin are often required to achieve quantitative growth (increase of cell number and volume). In tissue culture, usually only cytokinin and auxin are added. Plant hormones are typically added within the range 0.1-10 μ M.

Antibiotics

Antibiotics are used to prevent microbial contamination and eliminate Agrobacterium post-transformation and suppress bacterial infections and mold and yeast infections in cell cultures. They are suitable for sensitive research and genetic modification workflows.

Other Essential Ingredients

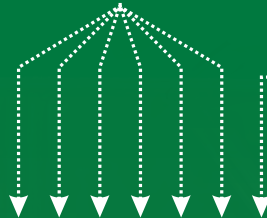

CODE	PRODUCT NAME	F. Wt. Gm./Ltr.	PACK SIZE	
TP 007	ANDERSON'S RHODODENDRON MEDIUM (W/O Vitamin, Sucrose, Agar) Micro & Macro Element	1.83	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 008	ANDERSON'S RHODODENDRON MEDIUM (W/ Vitamins, Sucrose, Agar & CaCl ₂ W/O IAA & 2iP)	42.17	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 082	ANDERSON'S RHODODENDRON MEDIUM (W/ Vitamins, Sucrose, Gelzan CM & CaCl ₂ W/O IAA & 2iP) for micro propagation of Musa species and Daucus species	35.17	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 095	AVOCADO-EMBRYO MEDIA	42.31	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 096	AVOCADO-PROLIFERATION MEDIUM	34.71	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 097	AVOCADO-PROPAGATION & MULTIPLICATION MEDIA	42.64	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 071	BANANA MICROPROPAGATION MEDIUM (W/ Vitamins, W/O NH ₄ NO ₃ , Sucrose, Agar)	2.58	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 072	BANANA MULTIPLICATION MEDIUM (W/ Vitamins, Glucose, Plant Growth Regulators & Agar) for in vitro multiplication of Musa species	41.5	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 073	BANANA MULTIPLICATION MEDIUM (W/ Vitamins, Tryptone, IAA W/O Sucrose, Agar)	4.77	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 075	BANANA MULTIPLICATION MEDIUM (W/ Vitamins, Sucrose, Tryptone, Agar & IAA)	42.77	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 076	BANANA MULTIPLICATION MEDIUM (W/ CaCl ₂ , Vitamins, Sucrose, Glucose, Ascorbic Acid, IAA, 6-BAP, Agar)	42.47	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 077	CARROT CALLUS INITIATION MEDIUM (W/ Vitamins, 2,4-D, W/O Sucrose, Agar)	3.21	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 027	CHEE AND POOL (C2D) VITIS MEDIUM (W/O Vitamins, Sucrose, Agar) Micro & Macro Element for regeneration and micro propagation of grapevine	4.44	5x1 Ltr 1x10 Ltr 1x50 Ltr	

CODE	PRODUCT NAME	F. Wt. Gm./Ltr.	PACK SIZE	
TP 028	CHEE AND POOL (C2D) VITIS MEDIUM (W/ Vitamins W/O Sucrose, Agar) Micro & Macro Element	4.49	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 029	CHU (N6) MEDIUM (W/O Vitamins, Sucrose, Agar) Micro & Macro Element	3.95	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 030	CHU (N) MEDIUM (W/ Vitamins, W/O Sucrose, Agar) Micro & Macro Element	4.03	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 031	CLC/ IPOMOEA CP MEDIUM (W/O Vitamins) Micro & Macro Element	6.61	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 032	CLC/ IPOMOEA EP MEDIUM (W/ Vitamins) Micro & Macro Element	3.55	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 083	CLC/ IPOMOEA BASAL MEDIUM (W/ Vitamins, W/O CaCl ₂ Sucrose & Agar)	6.40	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 078	DATE PALM CALLUS INITIATION MEDIUM (W/ Vitamins, Sucrose, 2,4-D, 2iP, Activated Charcoal, Agar) for in vitro micro propagation of date palm	45.76	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 033	DE GREEF & JACOBS MEDIUM (W/O Vitamins) Micro & Macro Element	3.77	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 034	DE GREEF & JACOBS MEDIUM (W/ Vitamins W/O Sucrose, Agar)	3.99	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 094	DKW BASAL MEDIUM (W/ Vitamins)	3.99	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 035	DKW MEDIUM (W/O Vitamins) Micro & Macro Element	5.48	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 036	DKW MEDIUM (W/ Vitamins)	5.36	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 009	ERIKSSON (ER) MEDIUM (W/O Vitamins, Sucrose & Agar)	4.00	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 010	ERIKSSON (ER) MEDIUM (W/ Vitamins, Sucrose W/O Agar)	34.14	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 011	GAMBORG MEDIUM (W/ Vitamins, Sucrose, CaCl ₂ , W/O IAA, Kinetin & Agar)	23.23	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 017	GAMBORG MEDIUM (W/O Vitamins, Sucrose & Agar)	3.05	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 012	GAMBORG MEDIUM (W/ Vitamins, Sugar & Agar)	31.23	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 074	GAMBORG'S (Vitamin 1000x Powder)	112.0	100 ml	

Plant Tissue Culture Media

CODE	PRODUCT NAME	F. Wt. Gm./Ltr.	PACK SIZE	
TP 085	GAMBORG B5 MEDIUM (W/ CaCl ₂ , Vitamins, Agar and Sucrose, W/O IAA & Kinetin)	31.23	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 079	GERBERA MULTIPLICATION MEDIUM (W/ CaCl ₂ Vitamins, Tyrosine, Sucrose, Adenine Sulphate, Agar, W/O IAA & Kinetin)	59.83	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 037	GRESSHOFF & DOY (DBM2) MEDIUM (W/O Vitamins) Micro & Macro Element	2.63	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 038	GRESSHOFF & DOY (DBM2) MEDIUM (W/ Vitamins) Micro & Macro Element	2.81	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 039	HELLER MEDIUM (W/O Vitamins) Micro & Macro Element	1.66	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 090	HOAGLAND NO.2 BASAL SALT MIXTURE	1.63	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 040	KAO & MICHAYLUK MEDIUM (W/O Vitamins) Micro & Macro Element	3.80	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 018	KNUDSON-C ORCHID MEDIUM MOREL MODIFICATION Micro & Macro Element	2.00	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 019	KNUDSON-C ORCHID MEDIUM (W/ Sucrose & Vitamins W/O Agar)	21.99	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 086	KNUDSON-C ORCHID MEDIUM, MOREL MODIFICATION (W/ Sucrose, Vitamins & Agar)	29.99	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 041	LINDEMANN ORCHID MEDIUM (W/O Vitamins) Micro & Macro Element	2.60	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 042	LINSMAIER & SKOOG MEDIUM (W/ Vitamins W/O Sucrose & Agar) Micro & Macro Element	4.54	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 043	LITVAY MEDIUM (W/O Vitamins) Micro & Macro with vitamins	4.95	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 044	LITVAY MEDIUM (W/ Vitamins) Micro & Macro Element	5.08	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 045	MCCOWN WOODY PLANT MEDIUM (W/O Vitamins) Micro & Macro Element	2.36	5x1 Ltr 1x10 Ltr 1x50 Ltr	

CODE	PRODUCT NAME	F. Wt. Gm./Ltr.	PACK SIZE	
TP 046	MCCOWN WOODY PLANT MEDIUM (W/ Vitamins) Micro & Macro Element	2.46	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 087	MURASHIGE & MILLER MEDIUM (W/ Na ₂ HPO ₄ and Vitamins W/O Sucrose & Agar)	4.69	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 093	MURASHIGE & SKOOG MICROELEMENT	9.79	1x10 Ltr 1x50 Ltr	
TP 047	MURASHIGE & SKOOG MEDIUM (W/ Vitamins & Folic Acid)	4.55	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 048	MURASHIGE & SKOOG MEDIUM (W/ Vitamins & MES Buffer)	5.04	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 001	MURASHIGE & SKOOG MEDIUM (W/ Vitamins, W/O Sucrose, Agar & Calcium Chloride)	4.10	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 001A	MURASHIGE & SKOOG MEDIUM (W/ Vitamins & Calcium Chloride W/O Sucrose, Agar)	4.54	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 002	MURASHIGE & SKOOG MEDIUM (W/ Sucrose, & Vitamins, W/O Agar & Calcium Chloride)	34.10	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 089	MURASHIGE & SKOOG MEDIUM W/ Calcium Chloride, Vitamins & Sucrose W/O Agar	34.54	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 060	MURASHIGE & SKOOG MEDIUM W/ Calcium Chloride, Vitamins, Sucrose & Agar	41.54	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 088	MURASHIGE & SKOOG MEDIUM W/ Calcium Chloride, Vitamins, Sucrose and Gellan Gum	37.54	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 003	MURASHIGE & SKOOG MEDIUM (W/ Sucrose, Agar, Vitamins, 6-BAP W/O Calcium Chloride, IAA & Kinetin)	42.10	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 004	MURASHIGE & SKOOG MEDIUM (S.B) (W/O Vitamins and with Calcium Chloride)	2.65	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 006	MURASHIGE & SKOOG MEDIUM (S.B) (W/O Calcium Chloride & Vitamins)	4.07	5x1 Ltr 1x10 Ltr 1x50 Ltr	

* Store between 2-8°C | Rare Chemicals



Plant Tissue Culture Media

CODE	PRODUCT NAME	F. Wt. Gm./Ltr.	PACK SIZE	
TP 020	MURASHIGE & SKOOG MEDIUM MODIFIED (SB) With 1/2 Macronutrient	2.18	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 021	MURASHIGE & SKOOG MEDIUM MODIFIED (SB) W/O NH_4NO_3	2.65	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 022	MURASHIGE & SKOOG MEDIUM MODIFIED (SB) $\frac{1}{2}$ x NH_4NO_3 and KNO_3	2.53	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 023	MURASHIGE & SKOOG MEDIUM MODIFIED (SB) W/ Micro and Macro nutrient W/O Ammonium Salt	4.40	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 024	MURASHIGE & SKOOG MEDIUM (SB) W/ minimal organics (MSMO)	4.39	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 025	MURASHIGE & SKOOG MODIFIED MEDIUM W/ Calcium Chloride & Gamborg's Vitamins	4.55	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 026	MURASHIGE & SKOOG MEDIUM Vitamin (1000x)	103.10	100 ml	
TP 013	NITSCH MEDIUM (W/ Vitamins, Sucrose & Agar) for anther culture	30.04	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 014	NITSCH MEDIUM (W/ Vitamins and Sucrose W/O Agar) Micro & Macro Element	22.04	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 049	NLN MEDIUM W/O $\text{CaNO}_3 \cdot 4\text{H}_2\text{O}$ Micro & Macro Element	0.39	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 050	NLN MEDIUM (With Vitamins) Micro & Macro Element	1.45	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 051	ORCHID MAINTENANCE MEDIUM (W/O Vitamins, MES, Sucrose, Tryptone, Charcoal & Agar)	2.17	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 052	ORCHID MAINTENANCE MEDIUM (W/ Vitamins, Buffer MES, Sucrose, Tryptone & Charcoal W/O Agar)	27.36	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 053	ORCHID MAINTENANCE MEDIUM (W/ Vitamins, Buffer MES, Sucrose and Tryptone, Charcoal & Agar)	35.36	5x1 Ltr 1x10 Ltr 1x50 Ltr	

CODE	PRODUCT NAME	F. Wt. Gm./Ltr.	PACK SIZE	
TP 080	ORCHID MAINTENANCE/ REPLATE MEDIUM (W/ CaCl_2 , Vitamins, Sucrose, Peptone, MES, Agar W/O Activated Charcoal)	33.36	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 081	POTATO MICROPROPAGATION MEDIUM (W/ CaCl_2 , Vitamins, IAA, Kinetin, W/O Sucrose & Agar)	4.54	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 054	QUOIRIN & LEPOIVRE MEDIUM (W/O Vitamins)	3.28	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 055	QUOIRIN & LEPOIVRE MEDIUM (W/ Vitamins W/O Sucrose & Agar)	3.41	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 057	RUGINI OLIVE MEDIUM (W/ Vitamins, W/O Sucrose & Agar) Micro & Macro Element	4.16	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 058	SCHENK & HILDEBRANDT MEDIUM (W/O Vitamins) Micro & Macro with Vitamins	3.18	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 059	SCHENK & HILDEBRANDT MEDIUM (W/ Vitamins W/O Sucrose & Agar) Micro & Macro with Vitamins	4.26	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 063	S-MEDIUM (Morel and Wetmore) (W/ Vitamins, MES Buffer & Sucrose)	13.03	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 070	VACIN & WENT MEDIUM (W/ Vitamins, Sucrose & Agar)	29.63	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 069	WEST VACO WU5 MEDIUM (for Embryogenic Culture)	5.22	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 015	WHITE MEDIUM Micro & Macro Element	0.94	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 016	WHITE MEDIUM (W/ Vitamins and Sucrose & W/O Agar)	21.06	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 084	WOODY PLANT MEDIUM (W/ CaCl_2 , Vitamins, W/O Sucrose, Agar)	2.43	5x1 Ltr 1x10 Ltr 1x50 Ltr	
TP 092	WOODY PLANT MEDIUM (W/ CaCl_2 , Vitamins & Sucrose W/O Agar)	22.40	5x1 Ltr 1x10 Ltr 1x50 Ltr	

Plant Tissue Culture Media

Ingredients

CODE	PRODUCT NAME	PACK SIZE	
	VITAMINS		
561	ADENINE 99% (6 Aminopurine, Vitamin B4) CAS 73-24-5 $C_5H_5N_5$ F. Wt. 135.13	5 gm 25 gm 1 Kg	
233	p-AMINOBENZOIC ACID (PABA) (Para Amino Benzoic Acid) CAS 150-13-0 $C_7H_7NO_2$ F. Wt. 137.14	100 gm 500 gm	
1706	L-ASCORBIC ACID (Vitamin C) CAS 50-81-7 $C_6H_8O_6$ F. Wt. 176.12	100 gm 500 gm	
1709	D(+)-BIOTIN (Vitamin H) * CAS 58-85-5 $C_{10}H_{16}N_2O_3$ S F. Wt. 244.31	1 gm 10 gm 25 gm	
1701	CALCIUM-D-PANTOTHENATE (Vitamin B5)* CAS 137-08-6 $C_{18}H_{32}CaN_4O_{10}$ F. Wt. 476.53	25 gm 100 gm 1 Kg	
335	FOLIC ACID (Vitamin B9)* (Pterogl-L-glutamic acid) CAS 59-30-3 $C_{19}H_{19}N_7O_6$ F. Wt. 441.40	5 gm 25 gm	
3151	HYPOXANTHINE CAS 68-94-0 $C_5H_4N_4O$ F. Wt. 136.112	5 gm 100 gm	
604	INOSITOL (Meso-Inositol) CAS 87-89-8 $C_6H_{12}O_6$ F. Wt. 180.2	25 gm 100 gm 1 Kg	
685	NICOTINAMIDE (Niacinamide) (Vitamin PP) CAS 98-92-0 $C_6H_6N_2O$ F. Wt. 122.1	25 gm 100 gm 1 Kg	
759	NICOTINIC ACID (Niacin) (Vitamin B) CAS 59-67-6 $C_6H_5NO_2$ F. Wt. 123.1	25 gm 100 gm 1 Kg	
4238	PYRROLOQUINOLINE QUINONE (PQQ) CAS 72909-34-3 $C_{11}H_6N_2O_3$ F. Wt. 330.21	100 mg 1 gm	
4239	PYRROLOQUINOLINE QUINONE DISODIUM CAS 122628-50-6 $C_{14}H_4Na_2O_8$ F. Wt. 347.17	100 mg 1 gm	
1705	PYRIDOXINE HCl (Vitamin B) CAS 58-56-0 $C_6H_{11}NO_3 \cdot HCl$ F. Wt. 205.64	10 gm 25 gm 100 gm	
1703	RIBOFLAVIN (Vitamin B2), Vitamin G2 CAS 83-88-5 $C_{17}H_{20}N_4O_6$ F. Wt. 376.36	10 gm 25 gm 100 gm	
1702	THIAMINE HCL (Vitamin B1) CAS 67-03-8 $C_{12}H_{17}ClN_4OS \cdot HC$ F. Wt. 337.27	25 gm 100 gm 500 gm	
1701	Vitamin B5 (Calcium-D-pantothenate) CAS 137-08-6 $C_{18}H_{32}CaN_4O_{10}$ F. Wt. 476.53	25 gm 100 gm 1 Kg	
1704	Vitamin B12 (Cyanocobalamin) CAS 68-19-9 $C_{63}H_{88}CoN_{14}P$ F. Wt. 1355.38	1 gm 10 gm	

CODE	PRODUCT NAME	PACK SIZE	
	AMINO ACID		
1103	L-ALANINE, 99% + PURITY CRYSTALLINE CAS 56-41-7 C_3H_7NO F. Wt. 89.10	25 gm 100 gm 500 gm	
1104	L- ARGININE 99% + PURITY (Free Base) CAS 74-79-3 $C_6H_{14}N_4O_2$ F. Wt. 174.20	25 gm 100 gm 1 Kg	
1106	L-ASPARAGINE, 99% + Purity (Mono) CAS 5794-13-8 $C_4H_8N_2O_3 \cdot H_2O_2$ F. Wt. 150.14	25 gm 100 gm 500 gm	
1108	L- ASPARTIC ACID 99% + PURITY, CRYSTALLINE, CAS 56-84-8 $C_4H_7NO_4$ F. Wt. 294.3	25 gm 100 gm 500 gm	
1138	L-CYSTEINE, (Free Base) 99% + CRYSTALLINE, CAS 52-90-4 $C_3H_7NO_2S$ F. Wt. 121.17 ₇₂	5 gm 25 gm 100 gm	
1139	L-CYSTEINE HCl (Mono) 99% + Purity CAS 7048-04-6 $C_3H_8ClN_2O_2 \cdot S \cdot HO$ F. Wt. 175.63	25 gm 100 gm 1 Kg 25 Kg	
1140	L-CYSTINE, 99% + CRYSTALLINE (Free Base) CAS 56-89-3 $C_6H_{12}N_2O_4S_2$ F. Wt. 240.3	25 gm 100 gm 1 Kg	
1111	L-GLUTAMIC ACID, 99% + CRYSTALLINE, AR GRADE, CAS 56-86-0 $C_5H_9NO_4$ F. Wt. 147.13	100 gm 500 gm 25 Kg	
1114	L-GLUTAMINE, 99% + Crystalline CAS 56-85-9 $C_5H_{10}N_2O_3$ F. Wt. 146.14	25 gm 100 gm 500 gm	
4340	GAMMA POLY GLUTAMIC ACID CAS 84960-48-5 $C_{15}H_{21}N_3O_7X_2$ F. Wt. 387.34	100 gm	
757	GLUTATHIONE (Reduced) CAS 70-18-8 $C_{10}H_{12}N_3O_6S$ F. Wt. 307.32	1 gm 5 gm 25 gm	
1109	GLYCINE 99% + CRYSTALLINE, AR GRADE CAS 56-40-6 $C_2H_5N_2O$ F. Wt. 75.07	100 gm 500 gm 5 Kg 25 Kg	
1115	L-HISTIDINE, 99% + Crystalline CAS 71-00-1 $C_6H_9N_3O_2$ F. Wt. 155.15	25 gm 500 gm	
1119	L-ISOLEUCINE, 99% + Crystalline CAS 73-32-5 $C_6H_{13}NO_2$ F. Wt. 131.17	5 gm 25 gm	
1121	L-LEUCINE, 99% + Crystalline CAS 61-90-5 $C_6H_{13}NO_2$ F. Wt. 131.17	25 gm 100 gm 500 gm	
1122	L-LYSINE MONO HCl, 99% + Crystalline CAS 657-27-2 $C_6H_{14}N_2O_2 \cdot HCl$ F. Wt. 182.65	100 gm 500 gm 25 Kg	
1123	L-METHIONINE, 99%+Crystalline CAS 63-68-3 $C_5H_{11}NO_2S$ F. Wt. 149.21	25 gm 100 gm 500 gm 25 Kg	
1125	L-ORNITHINE MONO HCl, 99% + Crystalline CAS 3184-13-2 $C_5H_{12}N_2O_2 \cdot HCl$ F. Wt. 168.62	25 gm 100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
1127	L-PHENYLALANINE, 99% + Crystalline CAS 63-91-2 $C_9H_9NO_2$ F. Wt. 165.2	25 gm 100 gm 500 gm	
4339	POLY GLUTAMIC ACID CAS 25513-46-6 $C_5H_9NO_4$ F. Wt. 147.13	100 gm	
1128	L-PROLINE, 99% + Crystalline CAS 147-85-3 $C_5H_9NO_2$ F. Wt. 115.13	5 gm 25 gm 500 gm	
1130	L-SERINE, 99% + Crystalline CAS 56-45-1 $C_3H_7NO_3$ F. Wt. 105.09	10 gm 25 gm 1 Kg	
956	TAURINE (2-Aminoethanesulphonic Acid) CAS 107-35-7 $C_2H_7NO_3S$ F. Wt. 125.15	25 gm 100 gm	
1131	L-THREONINE, 99% + Crystalline CAS 72-19-5 $C_4H_9NO_3$ F. Wt. 119.12	5 gm 25 gm 500 gm	
1134	L-TRYPTOPHAN, 99% + PURITY, CRYSTALLINE, CAS 73-22-3 $C_{11}H_{12}N_2O_2$ F. Wt. 204.23	5 gm 25 gm 500 gm	
1135	L-TYROSINE, 99 + Crystalline CAS 60-18-4 $C_9H_{11}NO_3$ F. Wt. 181.19	25 gm 100 gm 500 gm	
674	URACIL CAS 66-22-8 $C_4H_4N_2O_2$ F. Wt. 112.09	5 gm 25 gm	
1137	L-VALINE, 99% + PURITY, Crystalline, CAS 72-18-4 $C_6H_{11}NO_2$ F. Wt. 117.15	25 gm 100 gm	
	MICRONUTRIENTS		
245	AMMONIUM MOLYBDATE, (Tetra), EXTRA PURE CAS 12054-85-2 $(NH_4)_6Mo_7O_{24} \cdot 4H_2O$ F. Wt. 1235.86	100 gm 500 gm 25 Kg	
282	BORIC ACID, AR GRADE CAS 10043-35-3 H_3BO_3 F. Wt. 61.83	500 gm 5 Kg 50 Kg	
977	COBALT (II) NITRATE (Hexa), EXTRA PURE CAS 10026-22-9 $Co(NO_3)_2 \cdot 6H_2O$ F. Wt. 291.03	100 gm 500 gm	
978	COBALT (II) SULPHATE (Hepta), EXTRA PURE CAS 10026-24-1 $Co_2SO_3 \cdot 7H_2O$ F. Wt. 281.10	100 gm 500 gm	
320	CUPRIC SULPHATE (Penta), AR GRADE CAS 7758-99-8 $CuSO_4 \cdot 5H_2O$ F. Wt. 249.68	500 gm 25 Kg	
973	ETHYLENE DIAMINE TETRAACETIC ACID FERRIC-MONO SODIUM, EXTRA PURE CAS 15708-41-5 $C_{10}H_{12}N_2O_8FeNa$ F. Wt. 367	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
328	ETHYLENE DIAMINE TETRAACETIC ACID DISODIUM SALT, AR GRADE CAS 6381-92-6 $C_{10}H_{14}N_2O_8Fe \cdot Na_2 \cdot 2H_2O$ F. Wt. 372.2	100 gm 500 gm 5 Kg 25 Kg	
330	FERRIC AMMONIUM CITRATE (BROWN), EXTRA PURE CAS 1185-57-5 $C_{12}H_{18}FeN_2O_{14}$ F. Wt. 488.16	500 gm 5 Kg 25 Kg	
4223	FE-EDDHA (ETHYLENEDIAMINE-N,N'-BIS) (2-HYDROXYPHENYLACETATE ACID) FERRIC SODIUM COMPLEX CAS 16455-61-1 $C_{18}H_{18}FeN_2NaO_6$ F. Wt. 437.19	25 gm 100 gm	
756	FERRIC CITRATE, (Mono) for Bacteriology CAS 2338-05-8 $C_6H_5FeO_7$ F. Wt. 244.94	500 gm	
514	FERROUS SULPHATE (Hepta), AR GRADE CAS 7782-63-0 $FeSO_4 \cdot 7H_2O$ F. Wt. 278.01	500 gm 50 Kg	
704	MANGANESE (II) CHLORIDE, (Tetra), EXTRA PURE CAS 13446-34-9 $MnCl_2 \cdot 4H_2O$ F. Wt. 197.91	500 gm 50 Kg	
703	MANGANESE (II) SULPHATE (Mono), AR GRADE CAS 10034-96-5 $MnSO_4 \cdot H_2O$ F. Wt. 169.02	500 gm 50 Kg	
2131	NICKEL CHLORIDE (Hexa), EXTRA PURE CAS 7791-20-0 $NiCl_2 \cdot 6H_2O$ F. Wt. 237.69	500 gm	
394	POTASSIUM IODIDE (Confirming IP), EXTRA PURE CAS 7681-11-0 KI F. Wt. 166.00	100 gm 500 gm 25 Kg	
924	SODIUM CHLORIDE, MOLECULAR BIOLOGY GRADE, CAS 7647-14-5 NaCl F. Wt. 58.44	500 gm 5 Kg 50 Kg	
449	SODIUM MOLYBDATE, (Dihydrate), EXTRA PURE CAS 10102-40-6 $Na_2MoO_4 \cdot 2H_2O$ F. Wt. 241.95	100 gm 500 gm 25 Kg	
469	SODIUM THIOSULPHATE (Penta), EXTRA PURE CAS 10102-17-7 $Na_2S_2O_3 \cdot 5H_2O$ F. Wt. 248.18	500 gm 5 Kg 50 Kg	
2325	STRONTIUM NITRATE (ANH.), EXTRA PURE CAS 10042-76-9 $Sr(NO_3)_2$ F. Wt. 211.63	500 gm 50 Kg	
3508	ZINC SULPHATE (Hepta), AR GRADE CAS 7446-20-0 $ZnSO_4 \cdot 7H_2O$ F. Wt. 287.55	500 gm 5 Kg 50 Kg	

Plant Tissue Culture Media

Ingredients

CODE	PRODUCT NAME	PACK SIZE	
	MACRONUTRIENTS		
242	AMMONIUM CHLORIDE, AR GRADE CAS 12125-02-9 NH ₄ Cl F. Wt. 53.49	500 gm 5 Kg 50 Kg	
256	AMMONIUM DIHYDROGEN PHOSPHATE, AR GRADE CAS 7722-76-1 NH ₄ H ₂ PO ₄ F. Wt. 115.02	500 gm	
238	AMMONIUM HYDROGEN CARBONATE, AR GRADE CAS 1066-33-7 CH ₅ NO ₃ F. Wt. 79.056	500 gm	
246	AMMONIUM MOLYBDATE (Tetra), AR GRADE CAS 12054-85-2 (NH ₄) ₆ Mo ₇ O ₂₄ ·4H ₂ O F. Wt. 1235.86	100 gm 500 gm 25 Kg	
258	AMMONIUM SULPHATE, AR GRADE CAS 7783-20-2 (NH ₄) ₂ SO ₄ F. Wt. 132.16	500 gm 5 Kg 50 Kg	
289	CALCIUM CARBONATE (Precipitated), (ANH.), AR GRADE CAS 471-34-1 CaCO ₃ F. Wt. 100.09	500 gm 50 Kg	
292	CALCIUM CHLORIDE (Dihydrate), AR GRADE CAS 10035-04-8 CaCl ₂ ·2H ₂ O F. Wt. 147.02	500 gm 50 Kg	
828	CALCIUM CITRATE (Tetra), EXTRA PURE CAS 5785-44-4 C ₁₂ H ₁₀ Ca ₃ O ₁₄ ·4(H ₂ O) F. Wt. 570.50	500 gm 25 Kg	
4333	COPPER DISODIUM EDTA CAS No.14025-15-1 C ₁₀ H ₁₂ CuN ₂ NaO ₈ F. Wt. 374.75	100 gm	
160	CALCIUM NITRATE (Tetra), AR GRADE CAS 13477-34-4 CaH ₈ N ₂ O ₁₀ F. Wt. 236.15	500 gm 5 Kg	
163	tril-CALCIUM PHOSPHATE, EXTRA PURE CAS 7758-87-4 Ca ₃ (PO ₄) ₂ F. Wt. 310.2	500 gm 5 Kg 50 Kg	
614	MAGNESIUM CHLORIDE (Hexa), AR GRADE CAS 7791-18-6 MgCl ₂ ·6H ₂ O F. Wt. 203.3	500 gm 5 Kg 50 Kg	
4335	CALCIUM DISODIUM EDTA CAS No.62-33-9 C ₁₀ H ₁₂ CaN ₂ Na ₂ O ₈ ·xH ₂ O F. Wt. 374.27	500 gm	
4336	MAGNESIUM DISODIUM EDTA CAS No.14402-88-1 C ₁₀ H ₁₂ MgN ₂ NaO ₈ F. Wt. 335.51	500 gm	
4337	MANGANESE DISODIUM EDTA CAS No.15375-84-5 C ₁₀ H ₁₂ MnN ₂ NaO ₈ F. Wt. 366.14	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
357	MAGNESIUM NITRATE (Hexa), EXTRA PURE CAS 13446-18-9 Mg(NO ₃) ₂ ·6H ₂ O F. Wt. 256.41	500 gm	
359	MAGNESIUM SULPHATE (Hepta), AR GRADE CAS 10034-99-8 MgSO ₄ ·7H ₂ O F. Wt. 246.48	500 gm 5 Kg 50 Kg	
386	POTASSIUM CHLORIDE, AR GRADE CAS 7447-40-7 KCl F. Wt. 74.55	500 gm 5 Kg 50 Kg	
218	di-POTASSIUM HYDROGEN ORTHOPHOSPHATE (ANH.), AR GRADE CAS 7758-11-14 K ₂ HPO ₄ F. Wt. 174.17	500 gm 5 Kg 50 Kg	
220	POTASSIUM DIHYDROGEN ORTHOPHOSPHATE (ANH.), AR GRADE CAS 7778-77-0 KH ₂ PO ₄ F. Wt. 329.25	500 gm 5 Kg 50 Kg	
511	POTASSIUM HYDROXIDE PELLETS, AR GRADE CAS 1310-58-3 KOH F. Wt. 56.11	500 gm 5 Kg 50 Kg	
634	POTASSIUM IODIDE, AR GRADE CAS 7681-11-0 KI F. Wt. 166.0	100 gm 500 gm 25 Kg	
761	POTASSIUM NITRATE, EXTRA PURE CAS 7757-79-1 KNO ₃ F. Wt. 101.10	500 gm 5 Kg 50 Kg	
407	POTASSIUM SULPHATE, AR GRADE CAS 7778-80-5 K ₂ SO ₄ F. Wt. 176.26	500 gm 50 Kg	
763	SODIUM DIHYDROGEN ORTHOPHOSPHATE (Mono Sodium Phosphate Dihydrate), AR GRADE CAS 13472-35-0 NaH ₂ PO ₄ ·2H ₂ O F. Wt. 156.01	500 gm 5 Kg 50 Kg	
451	SODIUM NITRATE, AR GRADE CAS 7631-99-4 NaNO ₃ F. Wt. 85.0	500 gm 50 Kg	
4375	SODIUM NITROPHENOLATE (Sodium Salt) CAS 824-39-5 C ₆ H ₅ N NaO ₃ F. Wt. 161.09	25 gm 100 gm	
4334	ZINC DISODIUM EDTA CAS No.14025-21-9 C ₁₀ H ₁₂ N ₂ NaO ₈ Zn F. Wt. 376.58	500 gm	
	ANTIBIOTIC POWDER		
3300	ACETYLSALICYLIC ACID (Aspirin, 2-Acetoxybenzoic Acid) CAS 50-78-2 C ₉ H ₈ O ₄ F. Wt. 180.15	500 gm	
573	ACTIDIONE (Cycloheximide), AR GRADE CAS 66-81-9 C ₁₅ H ₂₃ NO ₄ F. Wt. 281.36	1 gm 5 gm	
2028	AMOXICILLIN (Trihydrate)* CAS 61336-70-7 C ₁₆ H ₂₅ N ₃ O ₅ S F. Wt. 419.5	1 gm	

* Store between 2-8°C | Rare Chemicals



CODE	PRODUCT NAME	PACK SIZE	
3078	AMPHOTERICIN B * CAS 1397-89-3 $C_{47}H_{79}NO_{17}$ F. Wt. 924.1	1 gm 5 gm	
2007	AMPICILLIN SODIUM * CAS 69-52-3 $C_{16}H_{18}N_3NaO_4S$ F. Wt. 371.4	1 gm 10 gm 25 gm 100 gm	
4002	AUGMENTIN (Amoxicillin : Clavulanic Acid Potassium 5:1)	1 gm	
2008	BACITRACIN (50,000 Units/vl) CAS 1405-87-4 $C_{66}H_{103}N_{17}O_{16}S$ F. Wt. 1422.69	1 vl	
2520	BACITRACIN ZINC CAS 1405-89-6 $C_{66}H_{101}N_{17}O_{16}SZn$ F. Wt. 1486.07	1 gm	
4228	CARBENDAZIM CAS 10605-21-7 $C_9H_8N_2O_2$ F. Wt. 191.19	1 gm 5 gm 25 gm	
3314	CARBENICILLIN DISODIUM * CAS 4800-94-6 $C_{17}H_{16}N_2Na_2O_6S$ F. Wt. 422.4	1 gm	
2016	CEFOTAXIME SODIUM * CAS 64485-93-4 $C_{16}H_{16}N_3NaO_7S_2$ F. Wt. 477.4	1 gm 5 gm	
2012	CEPHALEXIN (Mono) * CAS 23325-78-2 $C_{16}H_{17}N_3O_4S.H_2O$ F. Wt. 365.4	1 gm	
2017	CHLORAMPHENICOL * CAS 56-75-7 $C_{11}H_{12}Cl_2N_2O_5$ F. Wt. 323.1	5 gm 25 gm 100 gm	
4004	CLINDAMYCIN HCL CAS 18323-44-9 $C_{18}H_{33}ClN_2O_5S.HCl$ F. Wt. 424.98 g/mole	5 gm	
4207	D-CYCLOSERINE CAS 68-41-7 $C_3H_6N_2O_2$ F. Wt. 102.09	1 gm	
1805	COLCHICINE CAS 64-86-8 $C_{22}H_{25}NO_6$ F. Wt. 399.4 (Store at room temp.)	1 gm 10 gm	
2020	COLISTIN SULPHATE SODIUM (2,50,000/VL) * CAS 1264-72-8	1 gm	
4369	DA-6 (2-Diethylaminoethyl Hexanoate) CAS 10369-83-2 $C_{12}H_{25}NO_2$ F. Wt. 215.33	5 gm 25 gm	
4208	DICAMBA CAS 1918-00-9 $C_8H_6Cl_2O_3$ F. Wt. 221.04	100 mg 1 gm	
4370	2-(3,4-dichlorophenoxy) TRIETHYLAMINE (DCPTA) CAS 65202-07-05 $C_{12}H_{17}Cl_2NO$ F. Wt. 262.18	25 gm	
3539	DOXYCYCLINE HCl CAS 10592-13-9 $C_{22}H_{24}N_2O_8HCl$ F. Wt. 480.89	1 gm	

CODE	PRODUCT NAME	PACK SIZE	
2023	ERYTHROMYCIN CAS 114-07-8 $C_{37}H_{67}NO_{13}$ F. Wt. 733.9	1 gm	
3341	G-418-DISULPHATE (Geneticin Sulphate) CAS 108321-42-2 $C_{20}H_{40}N_4O_{10}$ F. Wt. 692.70	500 mg 1 gm	
2025	GENTAMYCIN SULPHATE * CAS 1405-41-0C $C_{21}H_{43}N_5O_7.H_2SO_4$ F. Wt. 575.67	1 gm 5 gm	
3357	HYGROMYCIN B CAS 31282-04-9 $C_{20}H_{37}N_3O_{13}$ F. Wt. 527.60	100 mg 1 gm	
2035	KANAMYCIN SULPHATE (Mono) * CAS 29701-07-3 $C_{18}H_{37}N_5O_{10}.H_2O_4S$ F. Wt. 581.68	1 gm 5 gm	
2026	LINCOMYCIN HCl (Mono) * (10,00000 Unit/vl) CAS 7179-49-9 $C_{18}H_{37}ClN_2O_7S$ F. Wt. 461.01	1 vl	
4209	METRONIDAZOLE CAS 443-48-1 $C_6H_9N_3O_3$ F. Wt. 171.2	1 gm	
4210	MINOCYCLINE HCL CAS 10118-90-8 $C_{23}H_{27}N_3O_7$ F. Wt. 457.47	100 mg 1 gm	
1713	NALIDIXIC ACID * CAS 389-08-2 $C_{12}H_{12}N_2O_3$ F. Wt. 232.2	1 gm 5 gm	
2003	NEOMYCIN SULPHATE * CAS 1405-10-3 $C_{23}H_{46}N_6O_{13}.3H_2SO_4$ F. Wt. 908.9	5 gm	
2463	NYSTATINE (10,00,000 u/vl) * CAS 1400-61-9 $C_{47}H_{75}NO_{17}$ F. Wt. 926.10	1 vl	
2004	OXYTETRACYCLINE HCl * CAS 2058-46-0 $C_{22}H_{24}N_2O_9.HCl$ F. Wt. 496.9	1 gm	
4229	PAROMOMYCIN SULPHATE CAS 1263-89-4 $C_{23}H_{47}N_5O_{18}S$ F. Wt. 713.71	1 gm 5 gm	
2005	POLYMYXIN B. SULPHATE (1 million units/vl) CAS 1405-20-5 $C_{55}H_{96}N_{16}O_{13}.2H_2SO_4$	1 vl 5 vl	
3393	PENICILLIN BENZYL SODIUM (PENCILLIN-G) CAS 69-57-8 $C_{16}H_{17}N_2NaO_4S(1MU/vl)$ F. Wt. 356.35	1 gm	
2200	RIFAMPICIN CAS 13292-46-1 $C_{43}H_{58}N_4O_{12}$	1 gm	

Plant Tissue Culture Media

Ingredients

CODE	PRODUCT NAME	PACK SIZE	
2030	STREPTOMYCIN SULPHATE CAS 3810-74-0 $C_{27}H_{39}O_{12}O_7 \cdot (H_2SO_4)_3$	5 gm 25 gm	
3439	SULPHAMETHOXAZOLE CAS 723-46-6 $C_{10}H_{11}N_3O_3S$	5 gm 25 gm	
070	TETRACYCLINE CAS 60-54-8 $C_{22}H_{24}N_2O_8$ F. Wt. 444.43	5 gm	
478	THIOMERSAL, EXTRA PURE CAS 54-64-8 $C_9H_9HgNaO_2S$ F. Wt. 404.81	25 gm 100 gm	
071	TICARCILLIN / CLAVULANIC ACID CAS 4697-14-7 $C_{15}H_{14}N_2Na_2O_6S_2$ F. Wt. 428.39	500 mg	
3441	TICARCILLIN SODIUM CAS 4697-14-7 $C_{15}H_{14}N_2Na_2O_6S_2$ F. Wt. 428.39	1gm 5 gm 25 gm	
2031	TRIMETHOPRIM CAS 738-70-5 $C_{11}H_{18}N_4O_3$ F. Wt. 290.3	5 gm	
072	TOBRAMYCIN SULPHATE CAS 49842-07-1 $C_{18}H_{37}N_5O_9$ F. Wt. 567.51	100 mg	
4211	VALIDAMYCINE CAS 37248-47-8 $C_{20}H_{35}NO_{13}$ F. Wt. 497.49	1 gm	
2032	VANCOMYCIN HCl CAS 1404-93-9 $C_{66}H_{75}Cl_2N_9O_{24} \cdot HCl$ F. Wt. 1485.7	500 mg	
4213	PLANT GROWTH REGULATORS 2, 8 HOMO BRASSINOLIDE CAS 80483-89-2 $C_{29}H_{50}O_6$ F. Wt. 494.70	10 mg 100 mg 1 gm	
4372	2, 8 HOMO BRASSINOLIDE 0.01% CAS 74174-44-0 $C_{29}H_{50}O_6$ F. Wt. 494.70	100 gm	
4371	2, 8 HOMO BRASSINOLIDE 90% CAS 74174-44-0 $C_{29}H_{50}O_6$ F. Wt. 494.70	1 gm	
4368	2-(3-Chlorophenoxy)-PROPIONIC ACID 99% CAS 101-10-0 $C_9H_9ClO_3$ F. Wt. 200.62	10 gm 25 gm	
3051	ABSCISIC ACID (ABA) CAS 21293-29-8 $C_{15}H_{20}O_4$ F. Wt. 264.3	100 mg 500 mg 1 gm	
3051A	ABSCISIC ACID (ABA) 10%* CAS 21293-29-8 $C_{15}H_{20}O_4$ F. Wt. 264.3	10 gm 25 gm 100 gm	
4240	ACETOSYRINGONE CAS 2478-38-8 $C_{10}H_{12}O_4$ F. Wt. 196.2	1 gm 5 gm	

CODE	PRODUCT NAME	PACK SIZE	
561	ADENINE 99% (6 Aminopurine, Vitamin B) CAS 73-24-5 $C_5H_5N_5$ F. Wt. 135.13	5 gm 25 gm 1 Kg	
562	ADENINE SULPHATE (Dihydrate) (6 Aminopurine Sulphate) CAS 106-60-5 $C_{10}H_{12}N_{10}O_4S$ F. Wt. 368.34	10 gm 100 gm	
4241	5-AMINOLEVULINIC ACID HYDROCHLORIDE * CAS 5451-09-2 $C_5H_7NO_3 \cdot HCl$ F. Wt. 167.59	1 gm 5 gm	
816	6-BENZYLADENINE (6BAP) (6-Benzylaminopurine) CAS 1214-39-7 $C_{12}H_{11}N_5$ F. Wt. 225.25	1 gm 5 gm 25 gm 1 Kg 50 Kg	
4322	CHLORMEQUAT CHLORIDE CAS 999-81-5 $C_5H_{13}Cl_2N$ F. Wt. 158.07	1 gm 5 gm 25 gm	
2565	p-CHLOROPHENOXY ACETIC ACID (4 CPA) CAS 122-88-3 $C_8H_7ClO_3$ F. Wt. 186.6	25 gm 100 gm 500gm	
4242	DAMINOZIDE (Alar) CAS 1596-84-5 $C_6H_{12}N_2O_3$ F. Wt. 160.17	5 gm	
2590	2,4, DICHLORO PHENOXY ACETIC ACID (2,4-D) CAS 94-75-7 $C_8H_6Cl_2O_3$ F. Wt. 221.04	500 gm	
4323	ETHEPHON CAS 16672-87-0 $ClCH_2CH_2P(O)(OH)_2$ F. Wt. 144.49	100 mg 500 mg 1 gm	
4212	FORCHLORFENURON (4-CPPU) N-(2-chloro-4-pyridinyl)-N'- phenyl CAS 68157-60-8 $C_{12}H_{10}ClN_3O$ F. Wt. 247.68	5 gm	
4243	GENISTEIN CAS 466-72-0 $C_{15}H_{10}O_5$ F. Wt. 270.24	100 mg	
3547	GIBBERELIC ACID (GA) ⁴⁺⁷ CAS 202467-69-4 $C_{20}H_{19}N_3O_5S$ F. Wt. 445.45	1 gm 5 gm 100 gm	
849	GIBBERELIC ACID A3 (GA) ³ CAS 77-06-5 $C_{19}H_{22}O_6$ F. Wt. 346.37	1 gm 10 gm 100 gm 1 Kg	
606	INDOLE-3-ACETIC ACID (3-Indole Acetic Acid, IAA) CAS 87-51-4 $C_{10}H_9NO_2$ F. Wt. 175.18 (Store at 2-15°C)	5 gm 25 gm 100 gm 1 Kg	
606A	INDOLE-3-ACETIC ACID POTASSIUM SALT (IAA-K) CAS 2338-19-4 $C_{10}H_8KNO_2$ F. Wt. 213.27	25 gm 100 gm 1 Kg	
607	INDOLE-3-BUTYRIC ACID * (3-Indole Butyric Acid, IBA) CAS 133-32-4 $C_{12}H_{13}NO_2$ F. Wt. 203.33	5 gm 25 gm 100 gm 1 Kg	

* Store between 2-8°C | Rare Chemicals



CODE	PRODUCT NAME	PACK SIZE	
607A	INDOLE-3-BUTYRIC ACID POTASSIUM SALT CAS 60096-23-3 $C_{12}H_{12}KNO_2$ F. Wt. 241.33	25 gm 100 gm 1 Kg	
4158	JASMONIC ACID CAS 77026-92-7 $C_{12}H_{18}O_3$ F. Wt. 210.27	100 mg 250 mg	
4215	N6-(2-ISOPENTENYL) ADENINE; (2IP) [6-(g-g-Dimethylallylamino) purine] CAS 2365-40-4 $C_{10}H_{13}N_5$ F. Wt. 203.24	1 gm 5 gm	
4324	N6-(2-ISOPENTENYL) ADENINE; (2IP) RIBOSIDE [6-(g-γγ-Dimethylallylamino)purine] CAS 7724-76-7 $C_{15}H_{21}N_5O_4$ F. Wt. 335.4	1 gm 5 gm	
863	KINETIN (6-Furfurylaminopurine), AR GRADE * CAS 525-79-1 $C_{10}H_9N_5O$ F. Wt. 215.21	1 gm 5 gm	
4214	META-TOPOLINE CAS 75737-38-1 $C_{12}H_{10}N_2OH$ F. Wt. 241.25	25 mg 100 mg	
2122	α-NAPHTHALENE ACETIC ACID (NAA) CAS 86-87-3 $C_{12}H_{10}O_2$ F. Wt. 186.21	25 gm 100 gm 25 Kg	
4374	1-NAPHTHYLACETAMIDE (98%) CAS 86-86-2 $C_{12}H_{11}NO$ F. Wt. 185.22	5 gm 25 gm	
2122A	α-NAPHTHALENE ACETIC ACID POTASSIUM SALT (NAA-K) CAS 15165-79-4 $C_{12}H_9O_2K$ F. Wt. 224.30	25 gm 100 gm 1 Kg	
3565	β-NAPHTHOXY ACETIC ACID (BONA) CAS 120-23-0 $C_{12}H_{10}O_3$ F. Wt. 202.21	25 gm 100 gm 500 gm	
3570	PACLOBUTRAZOL * CAS 76738-62-0 $C_{15}H_{20}ClN_3O$ F. Wt. 293.80	5 gm 25 gm	
906	PHENYLACETIC ACID (PAA; Benzeneacetic acid) CAS 103-82-2 $C_8H_8O_2$ F. Wt. 136.15	500 gm	
2457	PHLOROGLUCINOL, AR GRADE (1,3,5, Trihydroxybenzene) CAS 108-73-6 $C_6H_4O_3$ F. Wt. 126.11	25 gm 100 gm	
3576	PICLORAM CAS 1918-02-1 $C_6H_3Cl_3N_2O$ F. Wt. 241.46	1 gm 5 gm	
4245	PROHEXADIONE-CALCIUM CAS 127277-53-6 $C_{10}H_{10}CaO_5$ F. Wt. 250.263	100 mg 1 gm	
4325	PUTRESCINE DIHYDROCHLORIDE (1,4-Butanediamine dihydrochloride) CAS 333-39-7 (Store below 30° C) $C_4H_{12}H_2 \cdot 2HCl$ F. Wt. 161.07	1 gm 5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
680	SALICYCLIC ACID, AR GRADE (2-Hydroxybenzoic acid) CAS 69-72-7 $C_7H_6O_3$ F. Wt. 138.12	500 gm	
3436	SPERMIDINE [N-(3-Aminoopropyl)-1,4-diaminobutane] CAS 124-20-9 $C_7H_{19}N_3$ F. Wt. 145.2 *	1 gm 5 gm	
4270	THIABENDAZOLE CAS 148-79-8 $C_{10}H_7N_2S$ F. Wt : 201.25	50 gm 100 gm	
3584	THIDIAZURON CAS 51707-55-2 $C_9H_8N_4OS$ F. Wt. 220.25	250 mg 1 gm	
4294	2,3,5-TRIIODOBENZOIC ACID (TIBA)	5 gm	
3587	2,4,5, TRICHLOROPHENOXY ACETIC ACID CAS 93-76-5 $C_8H_5Cl_3O_3$ F. Wt. 255.49	25 gm 100 gm	
4364	TRIACONTANOL 98% CAS 593-50-0 $C_{30}H_{62}O$ F. Wt. 438.81 (Store at 2-4° C)	500 mg 1 gm	
516	ZEATIN CAS 1637-39-4 $C_{10}H_{13}N_5O$ F. Wt. 219.25 (Store at 2-4° C)	100 mg 500 mg 1 gm	
	CARBOHYDRATES		
4326	Trans-ZEATIN RIBOSIDE * CAS 6025-53-2 $C_{18}H_{21}N_5O_5$ F. Wt. 351.36	10 mg 25 mg	
3625	β-CYCLODEXTRIN CAS 7585-39-9 $C_{42}H_{70}O_{35}$ F. Wt. 1134.98	500 gm	
1315	D(-) FRUCTOSE (Bacteriological Grade) CAS 57-48-7 $C_6H_{12}O_6$ F. Wt. 180.16	100 gm 500 gm	
TC 010	D(+)-GALACTOSE CAS 59-23-4 $C_6H_{12}O_6$ F. Wt. 180.16	25 gm 100 gm 500 gm	
1303	DEXTROSE (Mono), EXTRA PURE CAS 5996-10-1 $C_6H_{12}O_6 \cdot H_2O$ F. Wt. 198.17	500 gm 5 Kg 25 Kg	
1304	DEXTROSE (ANH.), EXTRA PURE CAS 50-99-7 $C_6H_{12}O_6$ F. Wt. 180.16	500 gm 5 Kg 25 Kg	
1305	LACTOSE (Mono) CAS 10039-26-6 $C_{12}H_{22}O_{11} \cdot H_2O$ F. Wt. 360.31	500 gm 25 Kg	
1244	MALTOSE (Mono), EXTRA PURE CAS 6363-53-7 $C_{12}H_{24}O_{11} \cdot H_2O$ F. Wt. 360.31	100 gm 500 gm 25 Kg	

Plant Tissue Culture Media

Ingredients

CODE	PRODUCT NAME	PACK SIZE	
365	D+MANNOSE CAS 3458-28-4 $C_6H_{12}O_6$ F. Wt. 180.16	25 gm 100 gm	
1246	D-MANNITOL, AR GRADE CAS 69-65-8 $C_6H_{14}O_6$ F. Wt. 182.17	100 gm 500 gm 25 Kg	
TC 018	D-RAFFINOSE (Penta) CAS 17629-30-0 $C_{18}H_{32}O_{16} \cdot 5H_2O$ F. Wt. 594.52	10 gm 25 gm	
1306	D-RIBOSE CAS 50-69-1 $C_5H_{10}O_5$ F. Wt. 150.13	10 gm 25 gm 100 gm	
465	D-SORBITOL POWDER (D-Glucitol) CAS 50-70-4 $C_6H_{14}O_6$ F. Wt. 182.17	100 gm 500 gm 25 Kg	
1309	STARCH, POTATO (Insoluble) CAS 9005-25-8	500 gm	
TC 023	SUCROSE for Tissue Culture & Molecular Biology CAS 57-50-1 $C_{12}H_{22}O_{11}$ F. Wt. 342.3	500 gm 5 Kg	
TC 024	D-TREHALOSE (Dihydrate) CAS 6138-23-4 $C_{12}H_{22}O_{11} \cdot 2H_2O$ F. Wt. 378.33	5 gm 25 gm	
1312	D(+)-XYLOSE CAS 58-86-6 $C_5H_{10}O_5$ F. Wt. 150.13	25 gm 100 gm 500 gm	
	GELLING AGENT		
3500	AGAR AGAR Type I CAS 9002-18-0	100 gm 500 gm 5 Kg 25 Kg	
243 M	AGAR AGAR Type II CAS 9002-18-0	100 gm 500 gm 5 Kg 25 Kg	
1227	AGAROSE (High EEO) CAS 9012-36-6	5 gm 25 gm 100 gm 500 gm	
1270	AGAROSE SPL (Low EEO) (Counterelectrophoresis) CAS 9012-36-6	5 gm 25 gm 100 gm 500 gm	
4224	CHITOSAN CAS 9012-76-4 $(C_6H_{11}NO_2)_n$ F. Wt. 309.54	25 gm 100 gm 1 Kg	
4338	CHITOSAN OLIGOSACCHARIDE CAS No.148411-57-8 $C_{12}H_{24}N_2O_9$ F. Wt. 340.32	100 gm	

CODE	PRODUCT NAME	PACK SIZE	
3526	GELATIN CRYSTAL, Bloom Type B (Bacteriological Grade) CAS 9000-70-8	500 gm 5 Kg 50 Kg	
3318	GELLAN GUM Plant Tissue Culture Grade CAS 71010-52-1	100 gm 500 gm 1 Kg	
2115	POLYVINYL PYRROLIDINE K-30 (PVP) CAS 9003-39-8 $(C_6H_9O)_n$ F. Wt. 165.19	100 gm 500 gm 25 Kg	
	BUFFERING AGENT		
3086	BES, SODIUM, EXTRA PURE CAS 66992-27-6 $C_6H_4NNaO_3S$ F. Wt. 235.23	5 gm 25 gm	
3088	BIS-TRIS, EXTRA PURE (Bis (2-hydroxyethyl)amino-tris hydroxymethyl methane) * CAS 6976-37-0 $C_8H_{19}NO_5$ F. Wt. 209.24	25 gm 100 gm	
503	HEPES, FREE ACID CAS 7365-45-9 $C_8H_{19}N_2O_4S$ F. Wt. 238.30	25 gm 100 gm 500 gm 1 Kg	
2110	MES, (MONO)(2 (N-MORPHILINO) ETHANE SULFONIC ACID) CAS 4432-31-9 $C_6H_{13}NO_4SH_2O$ F. Wt. 213.2	25 gm 100 gm 1 Kg	
2219	PIPES (Piperazine-N,N-bis (2 Ethane Sulphonic Acid) CAS 5625-37-6 $C_8H_{18}N_2O_6S_2$ F. Wt. 302.4	5 gm 25 gm	
487	TRIS, 99%+Purity AR GRADE (Tris Hydroxymethyl Aminomethane) CAS 77-86-1 $C_4H_{11}NO_3$ F. Wt. 121.14	100 gm 500 gm 25 Kg	
488	TRIS HYDROCHLORIDE, AR GRADE CAS 1185-53-1 $C_4H_{11}NO_3 \cdot HCL$ F. Wt. 157.6	100 gm 500 gm 25 Kg	
	ADSORBING AGENT		
2903	CHAROCAL ACTIVATED, AR GRADE CAS 7740-44-0	500 gm	
	DISINFECTING AGENTS		
2632	8-HYDROXYQUINOLINE, AR GRADE, (Oxine) CAS 148-24-3 C_9H_7NO F. Wt. 145.2	100 gm 500 gm	
	OTHER ESSENTIAL INGREDIENTS		
3093	5-BROMO-4-CHLORO-3-INDOLYL-β-D-GALCTOPYRANOSIDE (X-Gal) CAS 7240-90-6 $C_{14}H_{14}BrClNO_6$ F. Wt. 408.63	100 mg 500 mg 1 gm	
3095	5-BROMO-4-CHLORO-3-INDOLYL-β-D-GLUCOPYRANOSIDE (X-Glucoside) CAS 15548-60-4 $C_{14}H_{14}BrClNO_6$ F. Wt. 408.63	100 mg	

* Store between 2-8°C | Rare Chemicals



CODE	PRODUCT NAME	PACK SIZE
3300	ACETYLSALICYLIC ACID (Aspirin)	500 gm
3509	ADENOSINE (9-β-Riboluranosylademine) * CAS 58-61-7 C ₁₀ H ₁₃ N ₅ O ₄ F. Wt. 267.25	5 gm 25 gm 100 gm
3063	ADENOSINE 5-DIPHOSPHATE DISODIUM (ADP-Na ₂) CAS 16178-48-6 C ₁₀ H ₁₃ N ₅ Na ₂ O ₁₀ P ₂ F. Wt. 391.19	1 gm 5 gm
3066	ALCIAN BLUE 8GX	5 gm 25 gm
1614	BROMO CRESOL PURPLE, AR GRADE (pH 5.2-6.8, Yellow to Purple)	5 gm 25 gm 100 gm
1657	BROMOTHYMOL BLUE (pH 5.8-7.6, Yellow to Blue)	5 gm 25 gm 100 gm
3523	CELLULASE * CAS 9012-54-8, Activity 10,000 U/gm	1 gm
827	CHITIN [Poly (N-acetyl-1,4-D-glucopyranosamine)]	100 gm 500 gm 25 Kg
4224	CHITOSAN CAS 9012-76-4 (C ₆ H ₁₁ NO ₄) _n F. Wt. 309.54	25 gm 100 gm 1 Kg
310	CITRIC ACID (Mono), AR GRADE CAS 5949-29-1 C ₆ H ₁₀ O ₈ F. Wt. 210.14	500 gm 5 Kg
1850	DITHIOERYTHREITOL (D.T.E.) * (1,4 Dithioerythreitol) CAS 6892-68-8 C ₄ H ₁₀ O ₂ S ₂ F. Wt. 154.25	1 gm 5 gm
1852	DITHIOTHREITOL (DTT)* (DL-DITHIOTHREITOL) CAS 3483-12-3 C ₄ H ₁₀ O ₂ S ₂ F. Wt. 154.25	1 gm 5 gm
304	ESCULIN (Aesculin) CAS 531-75-9 C ₁₅ H ₁₆ O ₉ F. Wt. 340.29	5 gm 25 gm
597	GLYCEROL (Glycerine), AR GRADE CAS 56-81-5 C ₃ H ₈ O ₃ F. Wt. 92.09	500 ml 2.5 Ltr 5 Ltr 25 Ltr
227	GUANIDINE HYDROCHLORIDE, AR GRADE CAS 50-01-1 CH ₅ N ₃ HCl F. Wt. 95.53	25 gm 100 gm
3156	ISOPROPYL-β-D-1-THIOGALACTOPYRANO SIDE (IPTG, DIOXAN FREE) * CAS 367-93-1 C ₉ H ₁₈ O ₅ S F. Wt. 238.3	1 gm 5 gm 25 gm
3556	MALEIC HYDRAZIDE CAS 123-33-1 C ₄ H ₂ N ₂ O ₂ F. Wt. 112.09	100 gm
363	(DL)-MALIC ACID CAS 6915-15-7 C ₄ H ₆ O ₅ F. Wt. 134.09	500 gm 5 Kg

CODE	PRODUCT NAME	PACK SIZE
4142	MALT EXTRACT POWDER Bacteriological Grade	500 gm 25 Kg
4373	MEPIQUAT CHLORIDE CAS 24307-26-4 C ₇ H ₁₆ N ₂ Cl F.Wt. 149.66	25 gm
369	MERCURIC CHLORIDE, EXTRA PURE (Mercury (II) Chloride) CAS 7487-94-7 HgCl ₂ F. Wt. 271.52	100 gm 250 gm 500 gm
2651	METHYLENE BLUE, AR GRADE	25 gm 100 gm 1 Kg
3564	MTT Methylthiazolyldiphenyl Tetrazolium Bromide * CAS 298-93-1 C ₁₈ H ₁₆ N ₆ SBBr F. Wt. 414.3	1 gm
2153	NITRO BLUE TETRAZOLIUM (Nitro Tetrazolium Blue, NBT), AR GRADE CAS 298-83-9 C ₄₀ H ₃₀ Cl ₂ N ₁₀ O ₆ F. Wt. 817.6	100 mg 250 mg 1 gm
3184	O-NITROPHENYL-b-D- GALACTOPYRANO SIDE * (ONPG) CAS 369-07-3 C ₁₂ H ₁₅ NO ₈ F. Wt. 301.25	500 mg 1 gm 5 gm
3518	PEPTONE (BACTO GRADE) CAS 73049-73-7	500 gm 5 Kg
908	POLYETHYLENE GLYCOL 4000 (P.E.G. 4000) CAS 25322-68-3	500 gm 25 Kg
1013	POLYETHYLENE GLYCOL 6000 (P.E.G. 6000) CAS 25322-68-3	500 gm 25 Kg
670	POLYOXYETHYLENE SORBITAN MONOLAUATE (Tween 20) CAS 9005-64-5 C ₅₈ H ₁₁₃ O ₂₆ F. Wt. 1226.5	500 ml 25 Ltr
511	POTASSIUM HYDROXIDE PELLETS, AR GRADE	500 gm 5 Kg 50 Kg
3258	POTASSIUM SILICATE POWDER	500 gm 50 Kg
673	POLYOXYETHYLENE SORBITAN MONOOLEATE (Tween 80) CAS 9005-65-6 C ₃₂ H ₆₀ O ₁₀ F. Wt. 604.8	100 ml 500 ml 25 Ltr
648	PROPYLENE GLYCOL CAS 57-55-6 C ₃ H ₈ O ₂ F. Wt. 76.09	500 ml 5 ltr 50 ltr
649	SILVER NITRATE, AR GRADE CAS 7761-88-8 AgNO ₃ F. Wt. 169.87	10 gm 25 gm
4247	SILICON DIOXIDE CAS-7631-86-9 SiO ₂ F. Wt. 60.08	100 gm
542	tri-SODIUM CITRATE (Dihydrate), AR GRADE CAS 6132-04-3 C ₆ H ₅ Na ₃ O ₇ ·2H ₂ O F. Wt. 294.1	500 gm 5 Kg 25 Kg

Plant Tissue Culture Media

Ingredients

CODE	PRODUCT NAME	PACK SIZE	
658	di-SODIUM HYDROGEN PHOSPHATE (Dihydrate), AR GRADE CAS 10028-24-7 $\text{Na}_2\text{HPO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 159.97	500 gm 5 Kg 50 Kg	
440	SODIUM HYDROXIDE PELLETS, AR GRADE CAS 1310-73-2 NaOH F. Wt. 39.99	500 gm 5 Kg 50 Kg	
659	SODIUM PYRUVATE, AR GRADE CAS 113-24-6 $\text{C}_3\text{H}_3\text{NaO}_3$ F. Wt. 110.04	100 gm 500 gm	
663	SODIUM THIOSULPHATE (ANH.), AR GRADE CAS 7772-98-2 $\text{Na}_2\text{S}_2\text{O}_3$ F. Wt. 158.1	500 gm 50 Kg	
3436	SPERMIDINE * [N-(3-Aminoopropyl)-1,4-diaminobutane] CAS 124-20-9 $\text{C}_7\text{H}_{19}\text{N}_3$ F. Wt. 145.2	1 gm 5 gm	

CODE	PRODUCT NAME	PACK SIZE	
1016	SUCCINIC ACID, EXTRA PURE CAS 110-15-6 $\text{C}_4\text{H}_6\text{O}_4$ F. Wt. 118.09	500 gm 25 Kg	
2742	TALCUM POWDER (400 mesh) CAS 14807-96-6	500 gm	
3226	TALCUM POWDER (700 mesh) CAS 14807-96-6	500 gm	
821	TETRAZOLIUM BLUE CHLORIDE, AR GRADE	1 gm 5 gm	
3514	TRYPTONE (REGULAR GRADE)	500 gm	
492	UREA, AR GRADE (Carbamide) CAS 57-13-6 $\text{CH}_4\text{N}_2\text{O}$ F. Wt. 60.06	500 gm 5 kg 50 kg	

Plant Tissue Culture Media

Pure Antibiotic Powder

CODE	PRODUCT NAME	PACK SIZE	
3525	ACTINOMYCIN D CAS 50-76-0 $C_{62}H_{86}N_2O_8S$ F. Wt. 1255.42	1 mg 5 mg	
2001	AMIKACIN SULPHATE CAS 39831-55-5 $C_{22}H_{43}N_5O_{13} \cdot 2H_2SO_4$ F. Wt. 781.76	1 gm	
4002	AMOXICLAV (AMOXICILLIN/CLAVULANIC ACID) CAS 61336-70-7 $C_{16}H_{25}N_3O_8S$ F. Wt. 419.45 CAS 58001-44-8 $C_8H_9NO_5$ F. Wt. 199.16	1 gm	
2028	AMOXICILLIN CAS 26787-78-0 $C_{16}H_{19}N_3O_5S$ F. Wt. 365.4	1 gm	
035	AMPICILLIN / SULBACTAM CAS 69388-84-7 CAS 69-52-3 $C_{16}H_{18}N_3NaO_4S$ F. Wt. 371.42 $C_8H_{10}NHaO_3S$ F. Wt. 255.22	1 gm	
2007	AMPICILLIN SODIUM CAS 69-52-3 $C_{16}H_{18}N_3O_4SNa$ F. Wt. 371.39	1 gm 10 gm 25 gm 100 gm	
3007	AMPICILLIN (Trihydrate) (Soluble in water) CAS 7177-48-2 $C_{16}H_{19}N_3O_5S \cdot 3H_2O$ F. Wt. 403.45	1 gm	
3078	AMPHOTERICIN B CAS 1397-89-3 $C_{27}H_{37}NO_{17}$ F. Wt. 924.10	1 gm 5 gm	
4002	AUGMENTIN (Amoxicillin : Clavulanic Acid Potassium 5:1)	1 gm	
4138	AZITHROMYCIN SULPHATE CAS 83905-01-5 $C_{38}H_{72}N_2O_{12}$ F. Wt. 748.99	25 mg	
2008	BACITRACIN (50,000 Units/vl) CAS 1405-87-4 $C_{66}H_{103}N_{17}O_{16}S$ F. Wt. 1422.69	1 vl	
2520	BACITRACIN ZINC CAS 1405-89-6 $C_{66}H_{101}N_{17}O_{16}SZn$ F. Wt. 1486.07	1 gm	
4228	CARBENDAZIM CAS 10605-21-7 $C_9H_9N_3O_2$ F. Wt. 191.19	1 gm 5 gm 25 gm	
3314	CARBENICILLIN DISODIUM CAS 4800-94-6 $C_{17}H_{16}N_2O_6SNa$ F. Wt. 422.36	1 gm	
2034	CEFACTOR CAS 53994-73-3 $C_{15}H_{14}ClN_3O_4S$ F. Wt. 367.81	1 gm	
2929	CEFADROXIL CAS 50370-12-2 $C_{16}H_{17}N_3O_5S$ F. Wt. 363.39	1 gm	

CODE	PRODUCT NAME	PACK SIZE	
2015	CEFAZOLIN SODIUM CAS 27164-46-1 $C_{14}H_{13}NaO_4S_2$ F. Wt. 477.49	1 gm	
4003	CEFIXIME CAS 79350-37-1 $C_{16}H_{15}N_5O_7S_2$ F. Wt. 453.452	1 gm	
2016	CEFOTAXIME SODIUM CAS 64485-93-4 $C_{16}H_{16}N_4NaO_7S_2$ F. Wt. 477.45	1 gm 5 gm	
3633	CEFOPERAZONE CAS 62893-19-0 $C_{22}H_{27}N_5O_7S_2$ F. Wt. 645.68	5 gm	
3636	CEFOXITIN CAS 35607-66-0 $C_{16}H_{17}N_5O_7S_2$ F. Wt. 427.45	100 mg	
3634	CEFSULODIN SODIUM SALT CAS 52152-93-9 $C_{22}H_{19}NaO_8S_2$ F. Wt. 554.52	10 mg 250 mg 1 gm	
2010	CEFTAZIDIME (Penta) CAS : 78439-06-2 $C_{22}H_{22}N_6O_7S_2 \cdot 5H_2O$ F.Wt : 636.65	1 gm	
2011	CEFTRIAZONE SODIUM CAS 104376-79-6 $C_{18}H_{16}N_8Na_2O_7S_3 \cdot \frac{1}{2}H_2O$ F. Wt. 598.54	1 gm	
3632	CEFUROXIME CAS 55268-75-2 $C_{16}H_{16}N_4O_8S$ F. Wt. 424.38	1 gm	
2012	CEPHALEXIN (MONO) CAS 23325-78-2 $C_{16}H_{19}N_3O_5S$ F. Wt. 365.4	1 gm	
3119	CEPHALOTHIN SODIUM CAS 58-71-9 $C_{16}H_{15}N_2O_6S_2Na$ F. Wt. 418.41	1 gm	
3120	CEPHALORIDINE CAS 50-59-9 $C_{19}H_{17}N_3O_4S_2$ F. Wt. 415.48	500 mg	
2017	CHLORAMPHENICOL CAS 56-75-7 $C_{11}H_{12}Cl_2N_2O_5$ F. Wt. 323.13	5 gm 25 gm 100 gm	
2019	CIPROFLOXACIN HCl CAS 86393-32-0 $C_{17}H_{18}FN_3O_3 \cdot HCl \cdot H_2O$ F. Wt. 385.82	1 gm	
4180	CLARITHROMYCIN CAS 81103-11-9 $C_{38}H_{69}NO_{13}$ F. Wt. 474.96	100 mg 500 mg	
4004	CLINDAMYCIN HCL CAS : 21462-39-5 $C_{18}H_{33}ClN_2O_5S \cdot HCl$ F.Wt : 461.44	5 gm	
4005	CLOXACILLIN CAS 61-72-3 $C_{15}H_{18}ClN_3O_5S$ F. Wt. 435.88	1 gm	

Plant Tissue Culture Media

Pure Antibiotic Powder

CODE	PRODUCT NAME	PACK SIZE	
2020	COLISTIN SULPHATE SODIUM (2,50,000 Units/vl) CAS 1264-72-8 $C_{58}H_{10}SN_{16}O_{28}S_5Na_5$ F. Wt. 1351.59	1 gm	
573	CYCLOHEXIMIDE (Actidione), AR GRADE CAS 66-81-9 $C_{15}H_{23}NO_4$ F. Wt. 281.36	1 gm 5 gm	
4207	D-CYCLOSERINE CAS 68-41-7 $C_3H_6N_2O_2$ F. Wt. 102.09	1 gm	
4208	DICAMBA CAS 1918-00-9 $C_6H_6Cl_2O_3$ F. Wt. 221.04	100 mg 1 gm	
3539	DOXYCYCLINE HCl CAS 10592-13-9 $C_{22}H_{24}N_2O_5HCl$ F. Wt. 480.89	1 gm	
2023	ERYTHROMYCIN CAS 114-07-8 $C_{37}H_{67}N_{13}$ F. Wt. 694.00	1 gm	
2024	ETHAMBUTOL Di-HCl CAS 1070-11-7 $C_{10}H_{24}N_2O_2 \cdot 2HCl$ F. Wt. 277.23	1 gm	
4006	FLUCONAZOLE CAS 86386-73-4 $C_{13}H_{12}F_2N_6O$ F. Wt. 306.2708	5 gm	
3341	G-418-DISULPHATE (Geneticin Sulphate) CAS 108321-42-2 $C_{20}H_{40}N_4O_{10}$ F. Wt. 692.70	500 mg 1 gm	
2025	GENTAMYCIN SULPHATE CAS 1405-41-0 $C_{21}H_{43}N_5O_7H_2SO_4$ F. Wt. 575.67	1 gm 5 gm	
3357	HYGROMYCIN B CAS 31282-04-9 $C_{20}H_{37}N_3O_{13}$ F. Wt. 527.60	100 mg 1 gm	
2948	ISONIAZIDE 99% CAS 54-85-3 $C_6H_7N_3O$ F. Wt. 137.14	5 gm	
2035	KANAMYCIN SULPHATE CAS 25389-94-0 $C_{18}H_{36}N_{16}N_4O_{11} \cdot H_2SO_4$ F. Wt. 582.58	1 gm 5 gm	
2026	LINCOMYCIN HCl (Mono) (10,00000 Unit/vl) * CAS 7179-49-9 $C_{18}H_{37}ClN_2O_5S$ F. Wt. 461.01	1 vl	
4168	LOMEFLOXACIN HCL CAS 98079-52-8 $C_{17}H_{20}ClF_2N_3O_3$ F. Wt. 387.81	1 gm	
4176	MEROPENEM CAS 96036-03-2 $C_{17}H_{25}N_5O_6S$ F. Wt. 387.81	10 mg	
4209	METRONIDAZOLE CAS 443-48-1 $C_6H_8N_4O_3$ F. Wt. 171.2	1 gm	

CODE	PRODUCT NAME	PACK SIZE	
4210	MINOCYCLINE HCL CAS 13614-98-7 $C_{23}H_{27}N_3O_7$ F.Wt : 493.94	100 mg 1 gm	
3637	MOXALACTAM CAS 64952-97-2 $C_{20}H_{20}N_6O_5S$ F. Wt. 520.47	500 mg	
2463	MYCOSTATINE (Nystatin) 10,00,000 unit/vl CAS 1400-61-9 $C_{47}H_{75}NO_{17}$ F. Wt. 926.09	1 vl	
1713	NALIDIXIC ACID CAS 389-08-2 $C_{12}H_{12}N_2O_3$ F. Wt. 232.24	1 gm 5 gm	
3528	NATAMYCIN (Pimaricin) CAS 7681-93-8 $C_{33}H_{47}NO_{13}$ F. Wt. 665.73	100 mg 1 gm	
2003	NEOMYCIN SULPHATE CAS 1405-10-3 $C_{23}H_{52}N_6O_{25}S_3$ F. Wt. 908.88	5 gm	
3384	NETILMICIN SULPHATE CAS 56391-57-2 $C_{42}H_{92}N_{10}O_{34}S_5$ F. Wt. 1441.55	100 mg	
3590	NISIN CAS 1414-45-5 $C_{143}H_{230}O_{37}N_{42}S_7$ F. Wt. 3354.07	1 gm 5 gm 25 gm	
2962	NORFLOXACIN CAS 70458-96-7 $C_{16}H_{18}FN_3O_3$ F. Wt. 319.33	1 gm	
3529	NOVOBIOCIN SODIUM CAS 1476-53-5 $C_{31}H_{35}N_2NaO_{11}$ F. Wt. 634.62	1 gm	
063	OFLOXACIN CAS 82419-36-1 $C_{18}H_{20}FN_3O_4$ F. Wt. 361.41	500 mg 5 gm	
3390	OXACILLIN SODIUM CAS 1173-88-2 $C_{19}H_{18}N_3NaO_5S$ F. Wt. 423.42	5 gm	
2004	OXYTETRACYCLINE HCL *	1 gm	
	CAS 2058-46-0 $C_{22}H_{24}N_2O_9 \cdot HCl$ F. Wt. 496.9		
4229	PAROMOMYCIN SULPHATE CAS 1263-89-4 $C_{23}H_{47}N_5O_{18}S$ F. Wt. 713.71	1 gm 5 gm	
3393	PENICILLIN BENZYL SODIUM (Pencillin-G) CAS 69-57-8 $C_{16}H_{17}N_2NaO_4S(1MU/vl)$ F. Wt. 356.35	1 gm	
3572	PENICILLIN V POTASSIUM CAS 132-98-9 $C_{16}H_{17}N_2O_5SK$ F. Wt. 388.48	500 mg	

* Store between 2-8°C | Rare Chemicals



CODE	PRODUCT NAME	PACK SIZE	
3397	PHOSPHOMYCIN DISODIUM CAS 26016-99-9 $C_3H_5O_4PNa_2$ F. Wt. 182.02	1 gm 5 gm	
3398	PIPERACILLIN SODIUM CAS 59703-84-3 $C_{23}H_{26}N_5NaO_7S$ F. Wt. 539.54	1 gm 5 gm	
2005	POLYMYXIN B. SULPHATE (1 million units/vl) CAS 1405-20-5 $C_{55}H_{96}N_{16}O_{13} \cdot 2H_2SO_4$	1 vl 5 vl	
3203	PYRAZINAMIDE CAS 98-96-4 $C_5H_5N_3O$	5 gm	
2200	RIFAMPICIN CAS 13292-46-1 $C_{43}H_{58}N_4O_{12}$	1 gm	
068	SPARFLOXACIN CAS 110871-86-8 $C_{19}H_{22}F_2N_4O_3$ F. Wt. 392.40	1 gm	
3435	SPECTINOMYCIN DIHYDROCHLORIDE, (Penta), CAS 22189-32-8 $C_{14}H_{24}N_2O_7 \cdot 2HCl \cdot 5H_2O$ F. Wt. 495.3	1 gm	
2030	STREPTOMYCIN SULPHATE CAS 3810-74-0 $C_{21}H_{39}O_{12}O_7 \cdot (H_2SO_4)_3$	5 gm 25 gm	
2033	SULFADIAZINE CAS 68-35-9 $C_{10}H_{10}N_4O_2S$	50 gm	
3244	SULPHACETAMIDE SODIUM CAS 127-56-0 $C_8H_9N_2NaO_3S \cdot aq$	5 gm	

CODE	PRODUCT NAME	PACK SIZE	
3439	SULPHAMETHOXAZOLE CAS 723-46-6 $C_{10}H_{11}N_3O_3S$	5 gm 25 gm	
3443	SULPHAMETHAZINE CAS 57-68-1 $C_{12}H_{14}N_4O_2S$ F. Wt. 278.33	25 gm	
070	TETRACYCLINE HYDROCHLORIDE CAS : 64-75-5 $C_{22}H_{24}N_2O_8$ F. Wt. 444.43	5 gm	
478	THIOMERSAL, EXTRA PURE CAS 54-64-8 $C_9H_9HgNaO_2S$ F. Wt. 404.81	25 gm 100 gm	
3441	TICARCILLIN DISODIUM CAS 4697-14-7 $C_{15}H_{14}N_2Na_2O_6S_2$ F. Wt. 428.39	1gm 5 gm 25 gm	
071	TICARCILLIN / CLAVULANIC ACID CAS 4697-14-7 $C_{15}H_{14}N_2Na_2O_6S_2$ F. Wt. 428.39	500 mg	
072	TOBRAMYCIN SULPHATE CAS 49842-07-1 $C_{18}H_{37}N_5O_9$ F. Wt. 567.51	100 mg	
2031	TRIMETHOPRIM CAS 738-70-5 $C_{14}H_{18}N_4O_3$ F. Wt. 290.32	5 gm	
4211	VALIDAMYCINE CAS 37248-47-8 $C_{20}H_{35}NO_{13}$ F. Wt. 497.49	1 gm	
2032	VANCOMYCIN HYDROCHLORIDE CAS 1404-93-9 $C_{66}H_{76}Cl_2N_9O_{24} \cdot HCl$ F. Wt. 1485.72	500 mg	



Chemicals

Pure & Precise

The lab deserves more than just **Chemicals**. It deserves the right grade for every breakthrough.

Available Grades

- Molecular Biology Grade
- Extra-Pure
- AR-Grade
- HPLC Grade
- ACS Grade
- Food Grade

High-purity Chemicals for



Quality Control



Diagnostics



Research



Drug Discovery



Cell Culture



Sensitive Analytical
Techniques



Reference
Standards



Nucleic Acid
Work



Enzyme Assays

1000 MT/Annum production capacity

CODE	PRODUCT NAME	PACK SIZE	
TMB 001	ACES (N-[2-Acetamido]-2-Aminoethane -Sulfonic Acid) CAS 7365-82-4 $C_4H_{10}N_2O_4S$ F. Wt. 182.20 Assay : $\geq 99\%$	5 gm 25 gm 100 gm	
TMB 002	ACETAMIDE (Amide C) CAS 60-35-5 C_5H_9NO F. Wt. 59.07 Assay : $\geq 99\%$	100 gm	
TMB 003	ACRYLAMIDE CAS 79-06-1 C_3H_5NO F. Wt. 71.08 Assay : $\geq 99\%$	25 gm 100 gm 500 gm 1 Kg	
TMB 004	ADENOSINE DIPHOSPHATE DISODIUM SALT (Dihydrate) * CAS 16178-48-6 $C_{10}H_{13}N_5Na_2O_{10}P_2 \cdot H_2O$ F. Wt. 507.2 Assay : $\geq 95\%$	1 gm 5 gm	
TMB 090	AGAR AGAR CAS 9002-18-0	100 gm 500 gm 1 kg	
TMB 005	AGAROSE DNA Grade, Low Melting Point CAS 9012-36-6	10 gm 100 gm	
TMB 006	AGAROSE, LOW EEO (Nuclease and Protease Free) CAS 9012-36-6	10 gm 100 gm 500 gm 1 Kg	
TMB 007	ALCIAN BLUE 8GX CAS 33864-99-2 $C_{56}H_{68}Cl_4CuN_{16}S_4$ F. Wt. 1298.86	25 gm 50 gm	
TMB 121	ALUMINIUM AMMONIUM SULPHATE CAS 7784-25-0 $AlNH_4(SO_4)_2 \cdot 12H_2O$ F. Wt. 453.32 Assay : $\geq 99.5\%$	500 gm	
TMB 008	AMMONIUM ACETATE CAS 631-61-8 $C_2H_7NO_2$ F. Wt. 77.08 Assay : $\geq 98\%$	100 gm 500 gm	
TMB 009	AMMONIUM CHLORIDE CAS 12125-02-9 NH_4Cl F. Wt. 53.49 Assay : $\geq 99.5\%$	500 gm	
TMB 010	DI-AMMONIUM HYDROGEN CITRATE CAS 3012-65-5 $C_6H_{14}N_2O_7$ F. Wt. 226.18 Assay : $\geq 99\%$	500 gm	
TMB 011	AMMONIUM MOLYBDATE (Tetra) CAS 12054-85-2 $(NH_4)_6Mo_7O_{24} \cdot 4H_2O$ F. Wt. 1235.86 Assay : $\geq 99\%$	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
TMB 012	AMMONIUM PERSULPHATE CAS 7727-54-0 $(NH_4)_2S_2O_8$ F. Wt. 228.2 Assay : $\geq 98\%$	25 gm 100 gm 500 gm	
TMB 013	AMMONIUM PHOSPHATE DIBASIC CAS 7783-28-0 $(NH_4)_2HPO_4$ F. Wt. 132.06 Assay : $\geq 98\%$	500 gm	
TMB 014	AMMONIUM PHOSPHATE MONOBASIC CAS 7722-76-1 $NH_4H_2PO_4$ F. Wt. 115.03 Assay : $\geq 99\%$	500 gm	
TMB 015	AMMONIUM SULPHATE CAS 7783-20-2 $(NH_4)_2SO_4$ F. Wt. 132.14 Assay : $\geq 99\%$	250 gm 500 gm	
TMB 016	BES-BUFFER CAS 10191-18-1 $C_6H_{15}NO_5S$ F. Wt. 213.25 Assay : $\geq 99.5\%$	25 gm 100 gm	
TMB 017	BICINE CAS 150-25-4 $C_6H_{13}NO_4$ F. Wt. 163.17 Assay : $\geq 99\%$	25 gm 100 gm 500 gm	
TMB 018	BIS-TRIS (Bis[2-hydroxyethyl Amino-Tris(hydroxymethyl)]Methane) CAS 6976-37-0 $C_8H_{19}NO_5$ F. Wt. 209.24 Assay : $\geq 98\%$	25 gm 100 gm	
TMB 019	BISACRYLAMIDE [N-N'-Methylene-Bis(Acrylamide)] CAS 110-26-9 $C_7H_{10}N_2O_2$ F. Wt. 154.17 Assay : $\geq 99\%$	25 gm 100 gm 250 gm 500 gm	
TMB 091	BORAX POWDER (Deca) (Sodium Tetraborate) CAS 1303-96-4 $Na_2B_4O_7 \cdot 10H_2O$ F. Wt. 381.37 Assay : $\geq 99.5\%$	500 gm	
TMB 092	BORIC ACID CAS 10043-35-3 H_3BO_3 F. Wt. 61.83 Assay : $\geq 99.5\%$	500 gm 1 Kg	
TMB 020	BOVINE SERUM ALBUMIN * CAS 9048-46-8 Assay : $\geq 98\%$ Protein	5 gm 25 gm 100 gm 500 gm	
TMB 021	BRILLIANT BLUE R-250 (Coomassie Brilliant Blue R-250) CAS 6104-59-2 C.I. No. 42660 $C_{45}H_{44}N_3NaO_7S_2$ F. Wt. 825.97	5 gm 25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
TMB 024	CALCIUM CHLORIDE (Dihydrate) CAS 10035-04-8 CaCl₂·2H₂O F. Wt. 147.01 Assay : ≥99%	100 gm 500 gm	
TMB 022	CAPS (3-[Cyclohexylamino]-1-Propanesulphonic Acid) CAS 1135-40-6 C₉H₁₉NO₃S F. Wt. 221.32 Assay : ≥99%	25 gm 100 gm 250 gm 1 Kg	
TMB 023	CHAPS CAS 75621-03-3 C₃₂H₅₈N₂O₇S F. Wt. 614.88 Assay : ≥98%	1 gm 5 gm 10 gm	
TMB 025	CITRIC ACID (ANH.) CAS 77-92-9 C₆H₈O₇ F. Wt. 192.12 Assay : ≥99.5%	500 gm	
TMB 093	CITRIC ACID (Mono) CAS 5949-29-1 C₆H₈O₇·H₂O F. Wt. 210.14 Assay : ≥99.5%	500 gm	
TMB 026	CITRIC ACID TRISODIUM SALT (Dihydrate), CAS 6132-04-3 C₆H₈O₇Na₃·2H₂O F. Wt. 294.1 Assay : ≥99%	500 gm	
TMB 027	COOMASSIE BRILLIANT BLUE G 250 CAS 6104-58-1 C.I. 42655 C₄₇H₄₈N₃NaO₇S₂ F. Wt. 854.02	5 gm 25 gm	
TMB 028	COPPER (II) SULPHATE (Penta) CAS 7758-99-8 CuSO₄·5H₂O F. Wt. 249.69 Assay : ≥99.5%	100 gm 500 gm	
TMB 029	CTAB (N-Cetyl-N,N,N-Trimethylammonium Bromide) CAS 57-09-0 CH₁₉BrN₄₂ F. Wt. 364.45 Assay : ≥99%	100 gm 500 gm	
TMB 030	DL-DTT (DL-Dithiothreitol) * CAS 3483-12-3 C₄H₁₀O₂S₂ F. Wt. 154.25 Assay : ≥99%	5 gm 25 gm	
TMB 031	1,4-Dithioerythritol (DTE) * CAS 6892-68-8 C₄H₁₀O₂S₂ F. Wt. 154.25 Assay : ≥99%	1 gm 5 gm	
TMB 032	EDTA SODIUM SALT (Dihydrate) (Ethylenediaminetetraacetic Acid Disodium Salt Dehydrate) CAS 6381-92-6 C₁₀H₁₄N₂Na₂O₈·2H₂O F. Wt. 372.24 Assay : ≥99%	500 gm 1 Kg	

CODE	PRODUCT NAME	PACK SIZE	
TMB 033	ETHIDIUM BROMIDE CAS 1239-45-8 C₂₂H₂₀BrN₃ F. Wt. 394.31 Assay : ≥95%	1 gm 5 gm	
TMB 094	D-FRUCTOSE (Levulose) CAS 57-48-7 C₆H₁₂O₆ F. Wt. 180.16 Assay : ≥99%	100 gm 500 gm 1 Kg	
TMB 095	D-GALACTOSE CAS 59-23-4 C₆H₁₂O₆ F. Wt. 180.16 Assay : ≥99%	100 gm 500 gm 1 Kg	
TMB 034	D(+)-GLUCOSE (ANH.) CAS 50-99-7 C₆H₁₂O₆ F. Wt. 180.16 Assay : ≥99.5%	100 gm 500 gm	
TMB 035	L-GLUTATHIONE REDUCED CAS 70-18-8 C₁₀H₁₇N₂O₆S F. Wt. 307.32 Assay : ≥98% Store Below 2-8° C	1 gm 5 gm 25 gm 100 gm	
TMB 036	GLYCINE CAS 56-40-6 C₂H₅NO₂ F. Wt. 75.07 Assay : ≥99%	100 gm 500 gm 1 Kg 5 Kg	
TMB 037	GUANIDINE HYDROCHLORIDE CAS 50-01-1 CH₅N₃·HCl F. Wt. 95.53 Assay : ≥99%	25 gm 100 gm 500 gm	
TMB 038	GUANINE CAS 73-40-5 C₅H₄N₆O F. Wt. 151.13 Assay : ≥98% Store Below 2-8° C	50 gm	
TMB 039	HEPES [N-(2-Hydroxyethyl) Piperazine-N'-(2-Ethanesulphonic Acid)] CAS 7365-45-9 C₈H₁₈N₂O₄S F. Wt. 238.3 Assay : ≥99%	25 gm 100 gm 500 gm	
TMB 040	HEPES SODIUM SALT CAS 75277-39-3 C₈H₁₇N₂NaO₄S F. Wt. 260.29 Assay : ≥99% Store Below 30° C	25 gm 100 gm 500 gm	
TMB 041	ISOPROPYL-β-D-THIOGALACTOPYRANOSIDE * CAS 367-93-1 C₉H₁₈O₅S F. Wt. 238.30 Assay : ≥99%	1 gm 5 gm 10 gm 25 gm	
TMB 042	LACTOSE (Monohydrate) CAS 64044-51-5 C₁₂H₂₂O₁₁·H₂O F. Wt. 360.31 Assay : ≥99.5%	500 gm 5 Kg	

CODE	PRODUCT NAME	PACK SIZE	
TMB 043	LAURYL SULPHATE SODIUM SALT (Dodecyl Sulphate Sodium Salt, SDS) CAS 151-21-3 $C_{12}H_{25}NaO_4S$ F. Wt. 288.38 Assay : $\geq 99\%$	25 gm 100 gm 500 gm 1 Kg	
TMB 096	LITHIUM CHLORIDE (ANH.) CAS 7447-41-8 LiCl F. Wt. 42.39 Assay : $\geq 99\%$	100 gm 500 gm	
TMB 044	MAGNESIUM ACETATE (Tetra) CAS 16674-78-5 $C_4H_6MgO_4 \cdot 4H_2O$ F. Wt. 214.45 Assay : $\geq 99\%$	100 gm	
TMB 045	MAGNESIUM CHLORIDE (Hexa) CAS 7791-18-6 $MgCl_2 \cdot 6H_2O$ F. Wt. 203.30 Assay : $\geq 99.5\%$	100 gm 500 gm	
TMB 046	MAGNESIUM SULPHATE (Hepta) CAS 10034-99-8 $MgSO_4 \cdot 7H_2O$ F. Wt. 246.47 Assay : $\geq 99.5\%$	500 gm	
TMB 047	DL-MALIC ACID CAS 6915-15-7 $C_4H_6O_5$ F. Wt. 134.09 Assay : $\geq 99\%$	500 gm	
TMB 048	D-(+)-MALTOSE (Mono) CAS 6363-53-7 $C_{12}H_{22}O_{11} \cdot H_2O$ F. Wt. 360.31 Assay : $\geq 95\%$	100 gm 500 gm	
TMB 049	MANGANESE (II) CHLORIDE (Tetra) CAS 13446-34-9 $MnCl_2 \cdot 4H_2O$ F. Wt. 197.91 Assay : $\geq 98\%$	100 gm 500 gm	
TMB 050	D-MANNITOL CAS 69-65-8 $C_6H_{14}O_6$ F. Wt. 182.17 Assay : $\geq 99\%$	500 gm 1 Kg	
TMB 051	2-MERCAPTOETHANOL CAS 60-24-2 C_2H_6OS F. Wt. 78.13 Assay : $\geq 98\%$	100 ml 500 ml	
TMB 052	MES (Monohydrate) [2-(N-Morpholino) Ethanesulphonic Acid Monohydrate] CAS 145224-94-8 $C_6H_{13}NO_4S \cdot H_2O$ F. Wt. 213.25 Assay : $\geq 90\%$	25 gm 100 gm	
TMB 053	MOPS (3-(N-Morpholino) Propanesulphonic Acid) CAS 1132-61-2 $C_7H_{15}NO_4S$ F. Wt. 209.26 Assay : $\geq 99.5\%$	25 gm 100 gm 500 gm	
TMB 054	MOPSO BUFFER CAS 68399-77-9 $C_7H_{15}NO_5S$ F. Wt. 225.26 Assay : $\geq 99\%$	25 gm	

CODE	PRODUCT NAME	PACK SIZE	
TMB 055	MTT * 3-(4,5-Dimethyl-2-Thiazolyl)-2,5 (Diphenyl -2H-Tetrazolium Bromide); Thiazolyl Blue CAS 298-93-1 $C_{18}H_{16}BrN_5S$ F. Wt. 414.32 Assay : $\geq 98\%$	100 mg 500 mg 1 gm	
TMB 056	NITROBLUE TETRAZOLIUM CHLORIDE (NBT) * CAS 298-83-9 $C_{40}H_{30}N_{10}O_6Cl_2$ F. Wt. 817.64 Assay : $\geq 98\%$	100 mg 250 mg 1 gm	
TMB 057	p-NITROPHENYL PHOSPHATE DISODIUM (Hexa) * CAS 333338-18-4 $C_6H_4NNa_2O P_6H_2O$ F. Wt. 371.14 Assay : $\geq 98\%$	5 gm 25 gm	
TMB 058	PHENOL, CRYSTALS (Hydroxybenzene) CAS 108-95-2 C_6H_6O F. Wt. 94.11 Assay : $\geq 99\%$ Store Below 2-8° C	100 gm 500 gm	
TMB 059	PIPES (Piperazine-N,N'-Bis) (2-Ethanesulphonic Acid) CAS 5625-37-6 $C_8H_{18}N_2O_6S_2$ F. Wt. 302.37 Assay : $\geq 99.5\%$	100 gm 500 gm	
TMB 060	POTASSIUM ACETATE CAS 127-08-2 $C_2H_3KO_2$ F. Wt. 98.14 Assay : $\geq 99\%$	100 gm 500 gm	
TMB 061	POTASSIUM CHLORIDE CAS 7447-40-7 KCl F. Wt. 74.55 Assay : $\geq 99\%$	500 gm	
TMB 097	POTASSIUM IODIDE CAS 7681-11-0 KI F. Wt. 166.00 Assay : $\geq 99.5\%$	100 gm 250 gm	
TMB 062	POTASSIUM NITRATE CAS 7757-79-1 KNO_3 F. Wt. 101.1 Assay : $\geq 99\%$	500 gm	
TMB 063	POTASSIUM PHOSPHATE DIBASIC (ANH.) CAS 7758-11-4 K_2HPO_4 F. Wt. 174.18 Assay : $\geq 99\%$	100 gm 500 gm	
TMB 064	POTASSIUM PHOSPHATE MONOBASIC (ANH.) CAS 7778-77-0 KH_2PO_4 F. Wt. 136.09 Assay : $\geq 99\%$	500 gm	
TMB 065	SILVER NITRATE CAS 7761-88-8 $AgNO_3$ F. Wt. 169.87 Assay : $\geq 99\%$	10 gm 25 gm 100 gm	
TMB 066	SODIUM ACETATE (Trihydrate) CAS 6131-90-4 $C_2H_3O_2Na \cdot 3H_2O$ F. Wt. 136.08 Assay : $\geq 99\%$	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
TMB 067	SODIUM BICARBONATE CAS 144-55-8 NaHCO₃ F. Wt. 84.01 Assay : ≥99.7%	100 gm 500 gm	
TMB 068	SODIUM CARBONATE (ANH.) (Soda Ash) CAS 497-19-8 Na₂CO₃ F. Wt. 105.99 Assay : ≥99.9%	100 gm 500 gm 1 Kg	
TMB 069	SODIUM CHLORIDE CAS 7647-14-5 NaCl F. Wt. 58.44 Assay : ≥99.5%	500 gm 1 Kg 5 Kg	
TMB 070	SODIUM HYDROXIDE PELLETS CAS 1310-73-2 NaOH F. Wt. 40.00 Assay : ≥98%	100 gm 500 gm 5 Kg	
TMB 071	SODIUM PHOSPHATE DIBASIC (Dihydrate) CAS 10028-24-7 Na₂HPO₄ .2H₂O F. Wt. 177.99 Assay : ≥99%	500 gm 5 Kg	
TMB 072	SODIUM PHOSPHATE MONOBASIC (ANH.) CAS 7558-80-7 NaH₂PO₄ F. Wt. 119.98 Assay : ≥98%	500 gm	
TMB 073	SODIUM SULPHATE (ANH.) CAS 7757-82-6 Na₂SO₄ F. Wt. 142.04 Assay : ≥99%	1 Kg	
TMB 074	SUCROSE CAS 57-50-1 C₁₂H₂₂O₁₁ F. Wt. 342.3 Assay : ≥99.5%	500 gm 5 Kg	
TMB 075	TEMED * (N,N,N',N'-Tetramethylethylenediamine) CAS 110-18-9 C₆H₁₆N₂ F. Wt. 116.20 Assay : ≥99%	100 ml	
TMB 076	TES (N-[Tris(Hydroxymethyl)Methyl]-2-Aminoethanesulphonic Acid) CAS 7365-44-8 C₆H₁₅NO₆S F. Wt. 229.25 Assay : ≥99%	100 gm	
TMB 077	TRICINE CAS 5704-04-1 C₆H₁₃NO₅ F. Wt. 179.17 Assay : ≥99%	100 gm	
TMB 078	2,3,5-TRIPHENYLTETRAZOLIUM CHLORIDE CAS 298-96-4 C₁₉H₁₅N₄Cl F. Wt. 334.8 Assay : ≥99%	10 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
TMB 079	TRIS, FREE BASE (Tris-[Hydroxymethyl] Aminomethane) CAS 77-86-1 C₄H₁₁NO₃ F. Wt. 121.14 Assay : ≥99%	100 gm 500 gm 1 Kg 5 Kg	
TMB 080	TRIS HYDROCHLORIDE CAS 1185-53-1 C₄H₁₁NO₃.HCl F. Wt. 157.6 Assay : ≥99%	100 gm 500 gm 1 Kg	
TMB 081	TRIS (HYDROXYMETHYL) AMINOMETHANE ACETATE CAS 6850-28-8 C₆H₁₅NO₅ F. Wt. 181.19 Assay : ≥99%	100 gm 500 gm	
TMB 082	TWEEN 20 CAS 9005-64-5 C₅₈H₁₁₂O₂₆ F. Wt. 1226.5	100 ml	
TMB 083	TWEEN 80 CAS 9005-65-6 C₃₂H₆₀O₁₀ F. Wt. 604.8	100 ml	
TMB 084	UREA CAS 57-13-6 CH₄N₂O F. Wt. 60.06 Assay : ≥99.5%	500 gm 1 Kg 5 Kg	
TMB 085	X-GAL * (5-Bromo-4-Chloro-3-Indolyl-β-D-Galactopyranoside) CAS 7240-90-6 C₁₄H₁₅BrClN₂O₆ F. Wt. 408.63 Assay : ≥98%	100 mg 500 mg 1 gm	
TMB 086	X-GLUCURONO CHA SALT * (5-Bromo-4-chloro-3-indolyl-β-D-glucuronide Cyclohexylammonium Salt) CAS 114162-64-0 C₂₀H₂₆BrClN₂O₇ F. Wt. 521.79 Assay : ≥98%	100 mg 500 mg 1 gm	
TMB 087	X-GLUCURONO SODIUM SALT * (5-Bromo-4-chloro-3-indolyl-β-D-glucuronide Sodium Salt) CAS 129541-41-9 C₁₄H₁₂BrClN₂NaO₇ F. Wt. 444.59 Assay : ≥98%	100 mg 500 mg 1 gm	
TMB 088	X-GLUCOSIDE * (5-Bromo-4-Chloro-3-Indolyl-β-D-Glucopyranoside) CAS 15548-60-4 C₁₄H₁₅BrClNO₆ F. Wt. 408.63 Assay : ≥97%	1 mg 5 mg 10 mg	
TMB 089	ZINC SULPHATE (Hepta) CAS 7446-20-0 ZnSO₄ .7H₂O F. Wt. 287.56 Assay : ≥99.5%	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
3051	ABSCISIC ACID (ABA, DORMIVE) * CAS 21293-29-8 $C_{15}H_{20}O_4$ F. Wt. 64.3	100 mg 500 mg 1 gm	
3051A	ABSCISIC ACID (ABA) 10%* CAS 21293-29-8 $C_{15}H_{20}O_4$ F. Wt. 264.3	10 gm 25 gm 100 gm	
210	ACACIA GUM POWDER (Spray Dried) (PHARMA GRADE) CAS 9000-01-5	500 gm 25 Kg	
3599	ACES BUFFER CAS 7365-82-4 $C_4H_{10}N_2O_4S$ F. WT. 182.20	5 gm 25 gm 100 gm	
681	ACETAMIDE (for Synthesis) CAS 60-35-5 C_2H_5NO F. Wt. 59.07	500 gm	
002	ACETANILIDE (N-Phenylacetamide), AR GRADE CAS 103-84-4 C_8H_9NO F. Wt. 135.17	500 gm	
982	ACETYL CHOLINE CHLORIDE, AR GRADE* CAS 60-31-1 $C_6H_{16}ClNO_2$ F. Wt. 181.68	5 gm 25 gm	
1143	N-ACETYL-L-CYSTEINE* CAS 616-91-1 $C_5H_9NO_3S$ F. Wt. 163.20	5 gm 25 gm 100 gm	
3058	N-ACETYL-L-TYROSINE* CAS 537-55-3 $C_{11}H_{13}NO_4$ F. Wt. 223.23	5 gm	
3300	ACETYLSALICYLIC ACID (Aspirin) CAS 50-78-2 $C_9H_8O_4$ F. Wt. 180.15	500 gm	
1681	ACID BLUE 9 (Erioglaucine) (Brilliant Blue FCF) CAS 3844-45-9 $C_{37}H_{34}Na_2N_2O_9S_3$ F. Wt. 792.86	10 gm	
1640	ACID FUCHSIN CAS 3244-88-0 C.I. 42685 $C_{20}H_{17}N_3Na_2O_5S_3$ F. Wt. 585.55	25 gm 100 gm	
1610	ACRIDINE ORANGE CAS 65-61-2 C.I. 46005 $C_{17}H_{20}ClN_3$ F. Wt. 301.81	10 gm 25 gm	
1695	ACRIFLAVINE (Euflavine) (For Biochemistry) CAS 8048-52-0 C.I. 46000 $C_{14}H_{14}ClN_3$ F. Wt. 259.73	5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
646	ACRIFLAVINE HYDROCHLORIDE CAS 8063-24-9 $C_{14}H_{13}N_3HCl$ F.Wt. 259.73	5 gm 25 gm	
201	ACRYLAMIDE, EXTRA PURE CAS 79-06-1 C_3H_5NO F.Wt. 71.80	500 gm	
573	ACTIDIONE (Cycloheximide), AR GRADE CAS 66-81-9 $C_{15}H_{23}NO_4$ F. Wt. 281.36	1 gm 5 gm	
561	ADENINE 99% (6 Aminopurine, Vitamin B) CAS 73-24-5 $C_5H_5N_5$ F. Wt. 135.13	5 gm 25 gm 1 Kg	
562	ADENINE SULPHATE (Adenine Hemi Sulphate) CAS 321-30-2 $C_{10}H_{12}N_6O_5S$ F. Wt. 368.34	10 gm 100 gm	
3509	ADENOSINE * CAS 58-61-7 $C_{10}H_{13}N_5O_4$ F. Wt. 267.25	5 gm 25 gm 100 gm	
3061	ADENOSINE 5-MONOPHOSPHATE* AMP Free Acid CAS 61-19-8 $C_{10}H_{14}N_5O_7P$ F. Wt. 347.22	1 gm 5 gm 25 gm	
3287	ADENOSINE-5-MONO PHOSPHATE DISODIUM (5'-ADP-Na) * CAS 4578-31-8 $C_{10}H_{13}N_5Na_2O_7P$ F. Wt. 391.18	1 gm 5 gm 25 gm	
3063	ADENOSINE 5-DIPHOSPHATE DISODIUM (5'-AMP-Na) * CAS 16178-48-6 $C_{10}H_{13}N_5Na_2O_{10}P_2$ F. Wt. 391.19	1 gm 5 gm	
3064	ADENOSINE-5-TRIPHOSPHATE DISODIUM (5'-ATP-Na) * CAS 987-65-5 $C_{10}H_{14}N_5Na_2O_{13}P_3$ F. Wt. 551.15	1 gm 5 gm 25 gm	
500	ADIPIC ACID, (Hexanedioic Acid) CAS 124-04-9 $C_6H_{10}O_4$ F. Wt. 146.14	500 gm	
2506	L-ADRENALINE, AR GRADE (Epinephrine) CAS 51-43-4 $C_9H_{13}NO_3$ F. Wt. 183.204	5 gm	
304	AESCULIN (Esculin) CAS 531-75-9 $C_{15}H_{16}O_9$ F. Wt. 340.29	5 gm 25 gm	
3500	AGAR AGAR Type I CAS 9002-18-0	100 gm 500 gm 5 Kg 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
3628	AGAR AGAR, FOOD GRADE CAS 9002-18-0	100 gm 500 gm	
1227	AGAROSE, HIGH EEO (0.23-0.26) CAS 9012-36-6 C ₂₄ H ₄₈ O ₁₉ F. Wt. 630.547	5 gm 25 gm 100 gm 500 gm	
1282	AGAROSE, MEDIUM EEO (0.16-0.19) CAS 9012-36-6 C ₂₄ H ₄₈ O ₁₉ F. Wt. 630.547	5 gm 25 gm 100 gm 500 gm	
1258	AGAROSE, LOW EEO (0.05-0.13) CAS 9012-36-6 C ₂₄ H ₄₈ O ₁₉ F. Wt. 630.547	5 gm 25 gm 100 gm 500 gm	
1270	AGAROSE, LOW EEO spl. for immunoelectrophoresis Counterelectrophoresis CAS 9012-36-6 C ₂₄ H ₄₈ O ₁₉ F. Wt. 630.547	5 gm 25 gm 100 gm 500 gm	
1283	AGAROSE, LOW EEO spl. for routine use CAS 9012-36-6 C ₂₄ H ₄₈ O ₁₉ F. Wt. 630.547	25 gm 100 gm 500 gm	
1101	β-ALANINE, 99% + PURITY Crystalline (L-Aminopropionic Acid) CAS 107-95-9 C ₃ H ₇ NO ₂ F. Wt. 89.10	25 gm 100 gm 500 gm	
1102	DL-ALANINE, 99% + PURITY Crystalline CAS 302-72-7 C ₃ H ₇ NO ₂ F. Wt. 89.10	25 gm 1 Kg	
1103	L-ALANINE, 99% + PURITY Crystalline CAS 56-41-7 C ₃ H ₇ NO ₂ F. Wt. 89.10	25 gm 100 gm 500 gm	
322	ALBUMIN EGG POWDER (Egg Albumin Powder) CAS 9006-59-1	500 gm	
3066	ALCIAN BLUE 8GX CAS 33864-99-2 C.I. 74240 C ₅₆ H ₆₈ Cl ₄ CuN ₁₆ S ₄ F. Wt. 1298.86	5 gm 25 gm	
684	ALGINIC ACID CAS 9005-32-7 (C ₆ H ₈ O ₆) _n F. Wt. 176.13	500 gm	
019	ALIZARIN CYANINE GREEN (Acid Green 25) CAS 4403-90-1 C.I. 61570 C ₂₈ H ₂₀ N ₆ O ₈ S ₂ Na ₂ F. Wt. 622.57	25 gm	
020	ALIZARIN RED S (Sodium Alizarin Sulphonate) CAS 130-22-3 C.I. 58005 C ₁₄ H ₇ NaO ₇ S F. Wt. 342.26	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
3069	ALLOXAN (Mono) CAS 2244-11-3 C ₄ H ₂ N ₂ O ₄ HO ₂ F. Wt. 160.10	25 gm 100 gm	
3288	ALLYLTHIOUREA (Thiosinamine) CAS 109-57-9 C ₄ H ₆ N ₂ S. F. Wt. 116.18	25 gm 100 gm	
4355	ALPHA OLEFIN SULFONATE (AOS), AR GRADE CAS 68439-57-6 C ₁₄ H ₂₇ NaO ₃ S F. Wt. 324.00	500 gm	
024	ALUMINIUM METAL POWDER, AR GRADE CAS 7429-90-5 Al F. Wt. 26.98	500 gm	
231	ALUMINIUM AMMONIUM SULPHATE (Ammonium Alum), EXTRA PURE CAS 7784-26-1 AlNH ₄ (SO ₄) ₂ .12H ₂ O F. Wt. 453.33	500 gm 50 Kg	
232	ALUMINIUM AMMONIUM SULPHATE (Ammonium Alum), AR GRADE CAS 7784-26-1 AlNH ₄ (SO ₄) ₂ .12H ₂ O F. Wt. 453.33	500 gm 50 Kg	
2404	ALUMINIUM HYDROXIDE GEL PASTE, AR GRADE CAS 21645-51-2 F. Wt. 78	500 gm 50 Kg	
2598	ALUMINIUM NITRATE (Nonahydrate), EXTRA PURE CAS 7784-27-2 Al(NO ₃) ₃ .9H ₂ O F. Wt. 375.134	500 gm	
753	ALUMINIUM POTASSIUM SULPHATE (Potash-Alum), EXTRA PURE CAS 7784-24-9 AlK ₂ O ₈ S ₂ .12H ₂ O F. Wt. 474.39	500 gm 5 Kg 50 Kg	
230	ALUMINIUM POTASSIUM SULPHATE (Potash-Alum), AR GRADE CAS 7784-24-9 AlK ₂ O ₈ S ₂ .12H ₂ O F. Wt. 474.39	500 gm 5 Kg 50 Kg	
791	ALUMINIUM SULPHATE, EXTRA PURE CAS 7784-31-8 Al ₂ (SO ₄) ₃ .18H ₂ O F. Wt. 666.42	500 gm 5 Kg 50 Kg	
1601	AMIDO BLACK 10 B (Naphthol Blue Black) CAS 1064-48-8 C.I. 20470 (Acid Black-1) C ₂₂ H ₁₄ N ₆ O ₉ S ₂ Na ₂ F. Wt. 616.50	25 gm 100 gm	
TCK 05	AMINO ACID KIT (Biochemistry) (20 Amino Acid 1 gm each & 4 Amino Acid 100 mg each) *	1 Kit	

CODE	PRODUCT NAME	PACK SIZE	
TCK 06	AMINO ACID KIT * (22 Essential Amino Acid 1 gm each)	1 Kit	
4020	1-AMINO-2-NAPHTHOL-4-SULPHONIC ACID, EXTRA PURE CAS 116-63-2 $C_{10}H_9NO_4S$ F. Wt. 239.25	25 gm	
036	4-AMINOANTIPYRINE (4-Amino Phenazone) CAS 83-07-8 $C_{11}H_{13}N_3O$ F. Wt. 203.25	25 gm 100 gm	
233	4-AMINOBENZOIC ACID (PABA) (Para Amino Benzoic Acid) CAS 150-13-0 $C_7H_7NO_2$ F. Wt. 137.14	100 gm 500 gm	
996	DL-2-AMINO BUTYRIC ACID CAS 2835-81-6 $C_4H_7NO_2$ F. Wt. 103.12	25 gm 100 gm	
796	p-AMINO-N, N-DIMETHYLANILINE OXALATE (N, N-Dimethyl-p-Phenylenediamine Oxalate) (for oxidase reagent) CAS 62778-12-5 $C_{10}H_{14}N_2O_2$ F. Wt. 226.23	1 gm 5 gm	
841	4-AMINO N,N-DIMETHYLANILINE SULPHATE (N, N-Dimethyl-p-Phenylenediamine Sulphate) CAS 536-47-0 $C_8H_{14}N_2SO_4$ F. Wt. 234.87	5 gm 25 gm	
042	4-AMINOPHENOL (p-Aminophenol) CAS 123-30-8 C_6H_7NO F. Wt. 109.13	100 gm 500 gm	
234	AMMONIUM ACETATE, EXTRA PURE CAS 631-61-8 $C_2H_3O_2NH_4$ F. Wt. 77.08	500 gm 50 Kg	
235	AMMONIUM ACETATE, AR GRADE CAS 631-61-8 $C_2H_3O_2NH_4$ F. Wt. 77.08	500 gm 50 Kg	
237	AMMONIUM BICARBONATE, EXTRA PURE (Ammonium Hydrogen Carbonate) CAS 1066-33-7 CH_5NO_3 F. Wt. 79.056	500 gm 5 Kg 50 Kg	
238	AMMONIUM BICARBONATE, AR GRADE (Ammonium Hydrogen Carbonate) CAS 1066-33-7 CH_5NO_3 F. Wt. 79.056	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
4023	AMMONIUM BIFLUORIDE, EXTRA PURE CAS 1341-49-7 F_2H_5N F. Wt. 57.04	500 gm	
240	AMMONIUM CARBONATE, EXTRA PURE CAS 506-87-6 $(NH_4)_2CO_3$ F. Wt. 96.09	500 gm	
053	AMMONIUM CARBONATE, AR GRADE CAS 506-87-6 $(NH_4)_2CO_3$ F. Wt. 96.09	500 gm	
054	AMMONIUM CERIC NITRATE, EXTRA PURE CAS 16774-21-3 $H_8N_6CeO_{18}$ F. Wt. 548.26	100 gm	
4025	AMMONIUM CERIC NITRATE, AR GRADE CAS 16774-21-3 $CeH_8N_6O_{18}$ F. Wt. 548.22	100 gm	
800	AMMONIUM CERIC SULPHATE (Dihydrate), EXTRA PURE CAS 10378-47-9 $CeH_{16}N_4O_{16}S_4 \cdot 2H_2O$ F. Wt. 632.55	100 gm 500 gm	
241	AMMONIUM CHLORIDE, EXTRA PURE CAS 12125-02-9 NH_4Cl F. Wt. 53.49	500 gm 5 Kg 50 Kg	
242	AMMONIUM CHLORIDE, AR GRADE CAS 12125-02-9 NH_4Cl F. Wt. 53.49	500 gm 5 Kg 50 Kg	
2514	AMMONIUM CITRATE DIBASIC (di-Ammonium Hydrogen Citrate) CAS 3012-65-5 $C_6H_{14}N_2O_7$ F. Wt. 226.19	500 gm	
802	tri-AMMONIUM CITRATE, EXTRA PURE CAS 3458-72-8 $C_6H_{17}N_3O_7$ F. Wt. 243.22	500 gm 50 Kg	
686	tri-AMMONIUM CITRATE, AR GRADE CAS 3458-72-8 $C_6H_{17}N_3O_7$ F. Wt. 243.22	500 gm	
551	AMMONIUM DICHROMATE, EXTRA PURE CAS 7789-09-5 $(NH_4)_2Cr_2O_7$ F. Wt. 252.07	500 gm	
255	AMMONIUM DIHYDROGEN ORTHOPHOSPHATE, EXTRA PURE (Ammonium Phosphate, Monobasic) CAS 7722-76-1 $NH_4H_2PO_4$ F. Wt. 115.03	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
256	AMMONIUM DIHYDROGEN ORTHOPHOSPHATE, AR GRADE CAS 7722-76-1 $\text{NH}_4\text{H}_2\text{PO}_4$ F. Wt. 115.03	500 gm	
330	AMMONIUM FERRIC CITRATE, (Brown) CAS 1185-57-5 $\text{C}_{12}\text{H}_{22}\text{FeN}_3\text{O}_{14}$ F. Wt. 488.16	500 gm 5 Kg 25 Kg	
243	AMMONIUM FERROUS SULPHATE, HEXAHYDRATE (Mohr's Salt), EXTRA PURE CAS 7783-85-9 $(\text{NH}_4)_2\text{Fe}(\text{SO}_4)_2 \cdot 6\text{H}_2\text{O}$ F. Wt. 392.16	500 gm 50 Kg	
244	AMMONIUM FERROUS SULPHATE, HEXAHYDRATE (Mohr's Salt), AR GRADE CAS 7783-85-9 $(\text{NH}_4)_2\text{Fe}(\text{SO}_4)_2 \cdot 6\text{H}_2\text{O}$ F. Wt. 392.16	500 gm 50 Kg	
804	AMMONIUM FLUORIDE, EXTRA PURE CAS 12125-1-8 FH_4N F. Wt. 37.04	500 gm	
055	AMMONIUM FORMATE, EXTRA PURE CAS 540-69-2 CH_3NO_2 F. Wt. 63.06	500 gm	
056	AMMONIUM FORMATE, AR GRADE CAS 540-69-2 CH_3NO_2 F. Wt. 63.06	500 gm	
2514	di-AMMONIUM HYDROGEN CITRATE, EXTRA PURE CAS 3012-65-5 $\text{C}_6\text{H}_{14}\text{N}_2\text{O}_7$ F. Wt. 226.19	500 gm	
253	di-AMMONIUM HYDROGEN ORTHOPHOSPHATE, EXTRA PURE (DiAmmonium phosphate dibasic) CAS 7783-28-0 $(\text{NH}_4)_2\text{HPO}_4$ F. Wt. 132.07	500 gm	
254	di-AMMONIUM HYDROGEN ORTHOPHOSPHATE, AR GRADE CAS 7783-28-0 $(\text{NH}_4)_2\text{HPO}_4$ F. Wt. 132.07	500 gm	
806	AMMONIUM META VANADATE, (Ammonium Vanadate), EXTRA PURE CAS 7803-55-6 NH_4VO_3 F. Wt. 116.98	100 gm 500 gm	
807	AMMONIUM META VANADATE, AR GRADE CAS 7803-55-6 NH_4VO_3 F. Wt. 116.98	100 gm 500 gm	
245	AMMONIUM MOLYBDATE (Tetra), EXTRA PURE CAS 12054-85-2 $(\text{NH}_4)_6\text{MO}_7\text{O}_{24} \cdot 4\text{H}_2\text{O}$ F. Wt. 1235.86	100 gm 500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
246	AMMONIUM MOLYBDATE (Tetra), AR GRADE CAS 12054-85-2 $(\text{NH}_4)_6\text{MO}_7\text{O}_{24} \cdot 4\text{H}_2\text{O}$ F. Wt. 1235.86	100 gm 500 gm 25 Kg	
057	AMMONIUM NICKEL SULPHATE (Hexa), (Nickel Ammonium Sulphate), EXTRA PURE CAS 7785-20-8 $\text{H}_8\text{N}_2\text{NiO}_8\text{S}_2 \cdot 6\text{H}_2\text{O}$ F. Wt. 395.00	500 gm	
249	AMMONIUM OXALATE (Mono), EXTRA PURE CAS 6009-70-7 $\text{C}_2\text{H}_8\text{N}_2\text{O}_4\text{HO}_2$ F. Wt. 142.12	500 gm 5 Kg	
250	AMMONIUM OXALATE (Mono), AR GRADE CAS 6009-70-7 $\text{C}_2\text{H}_8\text{N}_2\text{O}_4\text{HO}_2$ F. Wt. 142.12	500 gm	
251	AMMONIUM PERSULPHATE, EXTRA PURE CAS 7727-54-0 $\text{H}_8\text{N}_2\text{O}_8\text{S}_2$ F. Wt. 228.20	500 gm 25 Kg	
252	AMMONIUM PERSULPHATE, AR GRADE CAS 7727-54-0 $\text{H}_8\text{N}_2\text{O}_8\text{S}_2$ F. Wt. 228.20	500 gm 25 Kg	
060	AMMONIUM PURPURATE (Murexide) CAS 3051-09-0 C.I. 56085 $\text{C}_8\text{H}_8\text{N}_6\text{O}_6$ F. Wt. 284.19	5 gm 25 gm	
236	AMMONIUM SULPHAMATE, EXTRA PURE (Ammonium Amido Sulphonate) CAS 7773-06-0 $\text{H}_6\text{N}_2\text{O}_3\text{S}$ F. Wt. 114.125	100 gm 500 gm	
257	AMMONIUM SULPHATE, EXTRA PURE CAS 7783-20-2 $(\text{NH}_4)_2\text{SO}_4$ F. Wt. 132.14	500 gm 5 Kg 50 Kg	
258	AMMONIUM SULPHATE, AR GRADE CAS 7783-20-2 $(\text{NH}_4)_2\text{SO}_4$ F. Wt. 132.16	500 gm 5 Kg 50 Kg	
576	AMMONIUM(+)TARTRATE, AR GRADE CAS 3164-29-2 $\text{C}_4\text{H}_{12}\text{N}_2\text{O}_6$ F. Wt. 184.15	500 gm	
921	AMMONIUM THIOCYANATE, EXTRA PURE (Ammonium Sulphocyanate) (Amonium Rhodanide) CAS 1762-95-4 $\text{CH}_4\text{N}_2\text{S}$ F. Wt. 76.12	500 gm 50 Kg	
1417	α -AMYLASE (Diastase) 2000 Anson Unit/gm CAS 9000-90-2	100 gm 500 gm 5 Kg	

CODE	PRODUCT NAME	PACK SIZE	
1687	ANILINE BLUE (China Blue) (Water Soluble) CAS 28983-56-4 C.I. 42780 $C_{37}H_{27}N_3Na_2O_9S_3$ F. Wt. 799.80	25 gm 100 gm	
3295	ANILINE HYDROCHLORIDE CAS 142-04-1 $C_6H_5NH_2 \cdot HCl$ F. Wt. 129.59	500 gm	
3296	ANILINE SULPHATE CAS 542-16-5 $C_{12}H_{14}N_2 \cdot H_2SO_4$ F. Wt. 284.33	500 gm	
076	ANTHRACENE 99%, EXTRA PURE CAS 120-12-7 C.I. 10790 $C_{14}H_{10}$ F. Wt. 178.23	100 gm	
079	ANTHRAQUINONE CAS 84-65-1 $C_{14}H_8O_2$ F. Wt. 208.22	500 gm	
080	ANTHRONE (9,10, Dihydro-9-oxoanthracene) CAS 90-44-8 $C_{14}H_{10}O$ F. Wt. 194.24	25 gm	
2516	ANTIMONY METAL POWDER 99% CAS 7440-36-0 Sb F. Wt. 121.76	500 gm	
2517	ANTIMONY POTASSIUM TARTRATE, EXTRA PURE CAS 28300-74-5 $C_8H_4K_2O_{12}Sb_2 \cdot 3H_2O$ F. Wt. 667.87	100 gm 500 gm 25 Kg	
811	ANTIMONY POTASSIUM TARTRATE, (Potassium Antimony (III) Oxide Tartrate), AR GRADE CAS 28300-74-5 $C_8H_4K_2O_{12}Sb_2 \cdot 3H_2O$ F. Wt. 667.87	100 gm 500 gm	
809	ANTIMONY TRICHLORIDE, EXTRA PURE CAS 10025-91-9 $SbCl_3$ F. Wt. 228.13	100 gm 500 gm	
810	ANTIMONY TRICHLORIDE (Antimony (III) Chloride), AR GRADE CAS 10025-91-9 $SbCl_3$ F. Wt. 228.13	100 gm 500 gm	
085	ANTIMONY TRIOXIDE, EXTRA PURE CAS 1309-64-4 Sb_2O_3 F. Wt. 291.52	100 gm 500 gm	
086	ANTIMONY TRIOXIDE, AR GRADE CAS 1309-64-4 Sb_2O_3 F. Wt. 291.52	100 gm 500 gm	
2519	ANTIPYRINE, (Phenazone) CAS 60-80-0 $C_{11}H_{12}N_2O$ F. Wt. 188.23	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
087	D (-) ARABINOSE, AR GRADE CAS 10323-20-3 $C_5H_{10}O_5$ F. Wt. 105.13	5 gm 25 gm	
259	L (+) ARABINOSE CAS 87-72-9 $C_5H_{10}O_5$ F. Wt. 105.13	25 gm 100 gm	
1104	L- ARGININE 99% + PURITY (Free Base) CAS 74-79-3 $C_6H_{14}N_4O_2$ F. Wt. 174.20	25 gm 100 gm 1 Kg	
1105	L-ARGININE HYDROCHLORIDE (Mono) 99% + PURITY, CRYSTALLINE CAS 1119-34-2 $C_6H_{14}N_4O_2 \cdot HCl$ F. Wt. 210.67	25 gm 100 gm 1 Kg	
088	ARSENAZO I (Neothorin), AR GRADE CAS 520-10-5 $C_{16}H_{10}AsN_2Na_3O_{11}S_2$ F. Wt. 614.27	1 gm 5 gm	
090	ARSENIC TRIOXIDE, EXTRA PURE CAS 1327-53-3 As_2O_3 F. Wt. 197.84	500 gm	
1706	L-ASCORBIC ACID (Vitamin C) CAS 50-81-7 $C_6H_8O_6$ F. Wt. 176.12	100 gm 500 gm	
1715	L-ASCORBIC ACID (Vitamin C), AR GRADE CAS 50-81-7 $C_6H_8O_6$ F. Wt. 176.13	25 gm 100 gm 500 gm	
1714	L-ASCORBIC ACID SODIUM (Sodium Ascorbate) CAS 134-03-2 $C_6H_7NaO_6$ F. Wt. 198.11	100 gm 500 gm	
1148	DL-ASPARAGINE (Mono) 98.5% + Purity, Crystalline CAS 3130-87-8 $C_4H_8N_2O_3 \cdot H_2O$ F. Wt. 150.14	25 gm 100 gm	
1106	L- ASPARAGINE (Mono) 99% + Purity Crystalline CAS 5794-13-8 $C_4H_8N_2O_3 \cdot H_2O$ F. Wt. 150.14	25 gm 100 gm 500 gm	
1107	DL-ASPARTIC ACID 99% + Purity, Crystalline CAS 617-45-8 $C_4H_7NO_4$ F. Wt. 133.11	25 gm 100 gm 500 gm	
1108	L- ASPARTIC ACID 99% + Purity, Crystalline CAS 56-84-8 $C_4H_7NO_4$ F. Wt. 133.10	25 gm 100 gm 500 gm	
1638	AURAMINE CAS 2465-27-2 C.I. 41000 $C_{17}H_{21}N_3ClH$ F. Wt. 303.84	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
1022	AZOMETHIN-H-MONOSODIUM SALT CAS 206752-32-1 $C_{17}H_{12}NH_4O_8S_2 \cdot xH_2O$ F. Wt. 445.40	1 gm 5 gm	
1699	AZUR A CAS 531-53-3 C.I. 52005 $C_{14}H_{14}ClN_3S$ F. WT. 291.80	5 gm 25 gm 100 gm	
1689	AZUR I (AZUR B) CAS 531-55-5 $C_{15}H_{16}ClNS_3$ F. Wt. 305.83	5 gm 25 gm	
1697	AZUR II (Azur I and Methylene Blue) CAS 37247-10-2 C.I. 52010/52015 $C_{40}H_{52}Cl_2N_6S_2$ F. Wt. 751.94	25 gm	
1694	AZUR II EOSIN CAS 53092-85-6	25 gm	
094	BALSAM CANADA (Synthetic) CAS 8007-47-4	500 ml	
3301	BARBITURIC ACID CAS 67-52-7 $C_4H_4N_2O_3$ F. Wt. 128.09	100 gm 500 gm	
099	BARIUM ACETATE, EXTRA PURE CAS 543-80-6 $C_4H_6BaO_4$ F. Wt. 255.43	500 gm	
262	BARIUM ACETATE, AR GRADE CAS 543-80-6 $C_4H_6BaO_4$ F. Wt. 255.43	500 gm	
264	BARIUM CARBONATE, EXTRA PURE CAS 513-77-9 $BaCO_3$ F. Wt. 197.34	500 gm	
265	BARIUM CARBONATE, AR GRADE CAS 513-77-9 $BaCO_3$ F. Wt. 197.34	500 gm	
266	BARIUM CHLORIDE (Dihydrate), EXTRA PURE CAS 10326-27-9 $BaCl_2 \cdot 2H_2O$ F. Wt. 244.27	500 gm 5 Kg	
267	BARIUM CHLORIDE (Dihydrate), AR GRADE CAS 10326-27-9 $BaCl_2 \cdot 2H_2O$ F. Wt. 244.27	500 gm	
4026	BARIUM DIPHENYLAMINE SULPHONATE, AR GRADE CAS 6211-24-1 $C_{24}H_{20}BaN_2O_6S_2$ F. Wt. 633.88	5 gm 25 gm	
268	BARIUM HYDROXIDE (Octa), EXTRA PURE CAS 12230-71-6 $Ba(OH)_2 \cdot 8H_2O$ F. Wt. 315.48	500 gm 5 Kg	
270	BARIUM NITRATE, EXTRA PURE CAS 10022-31-8 $Ba(NO_3)_2$ F. Wt. 261.37	500 gm 5 Kg	

CODE	PRODUCT NAME	PACK SIZE	
4027	BARIUM NITRATE, AR GRADE CAS 10022-31-8 $Ba(NO_3)_2$ F. Wt. 261.37	500 gm	
273	BARIUM SULPHATE, EXTRA PURE CAS 7727-43-7 $BaSO_4$ F. Wt. 233.40	500 gm 5 Kg	
274	BARIUM SULPHATE, AR GRADE CAS 7727-43-7 $BaSO_4$ F. Wt. 233.40	500 gm	
1603	BASIC FUCHSIN (Fuchsin Basic, Margneta Basic) CAS 632-99-5 C.I. 42510 $C_{20}H_{20}ClN_3$ F. Wt. 337.86	25 gm 100 gm 1 Kg	
4028	BATHOPHENANTHROLINE, AR GRADE CAS 1662-01-7 $C_{24}H_{16}N_2$ F. Wt. 332.4	250 mg	
3501	B. MEAT EXTRACT PASTE (for general purpose, Bacteriological Grade)	500 gm 50 Kg	
3502	B. MEAT EXTRACT POWDER (for general purpose, Bacteriological Grade)	500 gm 25 Kg	
4029	BENTONITE POWDER CAS 1302-78-9 $Al_2O_3 \cdot 4(SiO_2) \cdot H_2O$ F. Wt. 360.31	500 gm	
4356	BENZALKONIUMCHLORIDE 50% CAS 68424-85-1 $C_{17}H_{30}ClN$ F. Wt. 283.88	500 gm	
4357	BENZALKONIUMCHLORIDE 80% CAS 68424-85-1 $C_{17}H_{30}ClN$ F. Wt. 283.88	500 gm	
108	BENZAMIDE CAS 55-21-0 C_7H_7NO F. Wt. 121.14	500 gm	
109	BENZANILIDE CAS 93-98-1 $C_{13}H_{11}NO$ F. Wt. 197.24	100 gm 500 gm	
4030	BENZENE SULPHONIC ACID, EXTRA PURE CAS 98-11-3 $C_6H_5O_3S$ F. Wt. 158.18	500 gm	
4031	BENZETHONIUM CHLORIDE (Hyamine 1622) CAS 121-54-0 $C_{27}H_{42}ClNO_2$ F. Wt. 448.08	25 gm	
3085	BENZIL CAS 134-81-6 $C_{14}H_{10}O_2$ F. Wt. 210.24	500 gm	
4352	1,2-BENZISOTHIAZOL- (2H)-ONE CAS 2634-33-5 C_7H_5NOS F. Wt. 151.19	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
4353	1,2-BENZISOTHAZOLE CAS 271-61-4 C ₇ H ₅ NS F. Wt. 135.19	100 gm 500 gm	
276	BENZOIC ACID, EXTRA PURE CAS 65-85-0 C ₆ H ₅ COOH F. Wt. 122.12	500 gm 50 Kg	
4032	BENZOIC ACID, AR GRADE (Meets Analytical Standards of IP, BP, USP, Ph.Eur) CAS 65-85-0 C ₇ H ₆ O ₂ F. Wt. 122.12	500 gm 25 Kg	
2526	BENZOIN CAS 119-53-9 C ₁₄ H ₁₂ O ₂ F. Wt. 212.25	100 gm 500 gm	
817	BENZOIN-α-OXIME (Cupron), AR GRADE CAS 441-38-3 C ₁₄ H ₁₃ NO ₂ F. Wt. 227.27	25 gm 100 gm	
3469	BENZOPHENONE (Diphenylketene) CAS 119-61-9 C ₁₃ H ₁₀ O F. Wt. 182.22	500 gm	
117	1,2,3, BENZOTRIAZOLE (Benzotriazole) CAS 95-14-7 C ₆ H ₅ N ₃ F. Wt. 119.12	100 gm	
816	6-BENZYLADENINE (6BAP) (6-Benzylaminopurine), Plant Growth Regulator CAS 1214-39-7 C ₁₂ H ₁₁ N ₅ F. Wt. 225.25	1 gm 5 gm 25 gm 1 Kg 25 Kg	
3086	BES, SODIUM, EXTRA PURE (N,N-Bis(2-hydroxyethyl)-2-amino ethane sulphonic acid sodium CAS 66992-27-6 C ₆ H ₁₄ NNaO ₅ S F. Wt. 235.23	5 gm 25 gm	
2388	BETAINE HYDROCHLORIDE CAS 590-46-5 C ₅ H ₁₂ ClNO ₂ F. Wt. 153.61	100 gm	
121	BILIRUBIN, AR GRADE CAS 635-65-4 C ₃₃ H ₃₆ N ₄ O ₆ F. Wt. 584.68	500 mg 1 gm	
1709	D-BIOTIN (Vitamin H) * CAS 58-85-5 C ₁₀ H ₁₆ N ₂ O ₃ S F. Wt. 244.31	1 gm 10 gm 25 gm	
2533	BIPHENYL (Diphenyl) CAS 92-52-4 C ₆ H ₅ C ₆ H ₅ F. Wt. 154.21	100 gm 500 gm	
2534	2-2 BIPYRIDYL (2,2-Dipyridyl), AR GRADE CAS 366-18-7 C ₁₀ H ₈ N ₂ F. Wt. 156.18	5 gm 25 gm	
3088	BIS-TRIS, EXTRA PURE (Bis(2-hydroxyethyl)amino-trishydroxy methy methane) CAS 6976-37-0 C ₈ H ₁₉ NO ₅ F. Wt. 209.24	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
371	BISACRYLAMIDE, AR GRADE (N-N-Meethylene-bisacrylamide) for Electrophoresis CAS 110-26-9 C ₇ H ₁₀ N ₂ O ₂ F. Wt. 154.17	25 gm 100 gm 500 gm	
4033	BISMUTH CARBONATE, EXTRA PURE CAS 5892-10-4 (BiO) ₂ CO ₃ F. Wt. 509.97	100 gm 500 gm	
4034	BISMUTH CHLORIDE, EXTRA PURE (Bismuth Trichloride) CAS 7787-60-2 BiCl ₃ F. Wt. 315.34	100 gm 500 gm	
4035	BISMUTH NITRATE, EXTRA PURE CAS 10035-06-0 BiN ₃ O ₉ F. Wt. 394.99	100 gm 500 gm	
4036	BISMUTH SUBNITRATE, EXTR APURE (Bismuth Oxynitrate) CAS 10361-46-3 BiNO ₄ F. Wt. 286.99	100 gm 500 gm	
4037	BISMUTH OXIDE, EXTRA PURE CAS 1304-76-3 Bi ₂ O ₃ F. Wt. 465.65	100 gm 500 gm	
682	BISMUTH SULPHATE CAS 7787-68-0 Bi ₂ O ₁₂ S ₃ F. Wt. 706.15	100 gm 500 gm	
2537	BISPHENOL A CAS 80-05-7 C ₁₅ H ₁₆ O ₂ F. Wt. 228.29	500 gm	
4358	BIT-20 (1,2-BENZISOTHAZOLI -3-ONE 20% SOLUTION) (BONBON 20) CAS 2634-33-5 C ₇ H ₅ NOS F. Wt. 151.19	500 ml 1 Ltr	
130	BIURET, AR GRADE (Allophanic Acid Amide) CAS 108-19-0 C ₂ H ₂ N ₃ O ₂ F. Wt. 103.08	25 gm	
821	BLUE TETRAZOLIUM CHLORIDE, AR GRADE CAS 1871-22-3 C ₄₀ H ₃₆ Cl ₂ N ₈ O ₂ F. Wt. 731.67	1 gm 5 gm	
279	BORAX POWDER (Deca) (Sodium Tetraborate), EXTRA PURE CAS 1303-96-4 Na ₂ B ₄ O ₇ ·10H ₂ O F. Wt. 381.37	500 gm 5 Kg 50 Kg	
280	BORAX POWDER (Deca) (Sodium Tetraborate), AR GRADE CAS 1303-96-4 Na ₂ B ₄ O ₇ ·10H ₂ O F. Wt. 381.37	500 gm 5 Kg 50 Kg	
281	BORIC ACID, EXTRA PURE CAS 10043-35-3 H ₃ BO ₃ F. Wt. 61.83	500 gm 5 Kg 50 Kg	
282	BORIC ACID, AR GRADE CAS 10043-35-3 H ₃ BO ₃ F. Wt. 61.83	500 gm 5 Kg 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
4170	di-BORON TRIOXIDE, EXTRA PURE CAS 1303-86-2 B ₂ O ₃ F. Wt. 69.62	250 gm 500 gm	
1226	BOVINE ALBUMIN FRACTION-V * (Albumin Bovine Fraction-V) CAS 9048-46-8	5 gm 25 gm 100 gm	
133	BRIJ 35 (Polyoxyethylene Lauryl ether) CAS 9002-92-0 (C ₁₂ H ₂₅ O. C ₂ H ₄ O) _n	500 ml	
1681	BRILLIANT BLUE FCF (Erioglaucine A) CAS 3844-45-9 C ₃₇ H ₃₄ N ₂ Na ₂ O ₅ S ₃ F. Wt. 792.86	10 gm	
1696	BRILLIANT BLUE G-250 (Coomassie Brilliant Blue G) for Electrophoresis CAS 6104-58-1 C.I. 42655 C ₄₇ H ₄₈ N ₃ NaO ₇ S ₂ F. Wt. 854.02	5 gm 25 gm	
1653	BRILLIANT CRESYL BLUE, AR GRADE CAS 81029-05-2 (C ₁₇ H ₂₀ CIN ₄) ₂ ZnCl ₂ F. Wt. 771.92	25 gm 100 gm	
1611	BRILLIANT GREEN CAS 633-03-4 C.I. 42040 C ₂₇ H ₃₃ N ₂ .HO ₄ S F. Wt. 482.65	25 gm 100 gm 1 Kg	
1630	BRILLIANT GREEN FCF CAS 2353-45-9 C.I. 42053 C ₃₇ H ₃₄ N ₂ O ₁₀ S ₃ Na ₂ F. Wt. 808.85	5 gm 25 gm	
4256	5-BROMO-4-CHLORO-3-INDOLYL- α-D-GALACTOPYRANOSIDE * CAS 107021-38-5 C ₁₄ H ₁₅ BrClNO ₆ F. Wt. 408.63	100 mg	
3093	5-BROMO-4-CHLORO-3-INDOLYL- β-D-GALACTOPYRANOSIDE (X-Gal) * CAS 7240-90-6 C ₁₄ H ₁₄ BrClNO ₆ F. Wt. 408.63	100 mg 500 mg 1 gm	
4161	5-BROMO-4-CHLORO-3-INDOLYL-β-D- GLUCURONIDE CYCLOHEXYL AMMONIUM SALT * CAS 114162-64-0 C ₂₀ H ₂₆ BrCIN ₂ O ₇ F. Wt. 521.79	1 gm	
3096	5-BROMO-4-CHLORO-3-INDOLYL- PHOSPHATE di-SODIUM (BCIP)* CAS 102185-33-1 C ₈ H ₇ BrClNO ₄ P.Na ₂ F. Wt. 370.43	500 mg 1 gm	
275	BROMO CRESOL GREEN SODIUM, AR GRADE CAS 62625-32-5 C ₂₁ H ₁₄ Br ₄ NaO ₅ S F. Wt. 720.0	5 gm 25 gm	
1612	BROMOCRESOL GREEN, AR GRADE (Bromo Cresol Blue) CAS 76-60-8 C ₂₁ H ₁₄ Br ₄ O ₅ S F. Wt. 698.01	5 gm 25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
1614	BROMOCRESOL PURPLE, AR GRADE (pH 5.2-6.8, Yellow to Purple) CAS 115-40-2 C ₂₁ H ₁₆ Br ₂ O ₅ S F. Wt. 540.24	5 gm 25 gm 100 gm	
1656	BROMOPHENOL BLUE INDICATOR, AR GRADE CAS 115-39-9 C ₁₉ H ₁₀ Br ₄ O ₅ S F. Wt. 669.96	5 gm 25 gm	
4039	BROMOPYROGALLOL RED, AR GRADE CAS 16574-43-9 C ₁₉ H ₁₀ Br ₂ O ₈ S F. Wt. 558.15	1 gm 5 gm	
1657	BROMOTHYMOL BLUE (pH 5.8-7.6, Yellow to Blue) CAS 76-59-5 C ₂₇ H ₂₈ Br ₂ O ₅ S F. Wt. 624.38	5 gm 25 gm 100 gm	
3327	BROMOTHYMOL BLUE SODIUM CAS 34722-90-2 C ₂₇ H ₂₇ Br ₂ O ₅ Na F. Wt. 646.35	5 gm 25 gm	
3308	BRUCINE CAS 357-57-3 C ₂₃ H ₂₆ N ₂ O ₄ F. Wt. 394.46	25 gm	
3309	BRUCINE SULPHATE (Hepta) CAS 5787-00-8 (C ₂₃ H ₂₆ N ₂ O ₄) ₂ .H ₂ SO ₄ .7H ₂ O F. Wt. 1013.13	25 gm	
284	BUFFER CAPSULES pH 4.0	10 Cap. 10x10	
285	BUFFER CAPSULES pH 7.0	10 Cap. 10x10	
286	BUFFER CAPSULES pH 9.2	10 Cap. 10x10	
4041	1-BUTANE SULPHONIC ACID SODIUM SALT (Anhydrous), (For HPLC) CAS 2386-54-1 C ₄ H ₉ NaO ₃ S F. Wt. 160.16	25 gm	
4042	1-BUTANE SULPHONIC ACID SODIUM SALT (Monohydrate), (For HPLC) CAS 2386-54-1 C ₄ H ₉ NaO ₃ .S.H ₂ O F. Wt. 178.16	25 gm	
581	BUTYLATED HYDROXY ANISOLE (B.H.A.)* CAS 25013-16-5 C ₁₁ H ₁₆ O ₂ F. Wt. 180.24	100 gm 500 gm 25 Kg	
582	BUTYLATED HYDROXY TOLUENE (B.H.T.), (2,6-di-tert-butyl-p-cresol) CAS 128-37-0 C ₁₅ H ₂₄ O F. Wt. 220.35	500 gm 5 Kg 25 Kg	
144	Tert-BUTYL HYDROQUINONE (TBHQ) CAS 1948-33-0 C ₁₀ H ₁₄ O ₂ F. Wt. 166.22	100 gm 500 gm 25 gm	
3109	CADMIUM SULPHATE (Octa), EXTRA PURE CAS 7790-84-3 3CdSO ₄ .8H ₂ O F. Wt. 769.50	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
3110	CADMIUM SULPHATE (Octa), AR GRADE CAS 7790-84-3 $3\text{CdSO}_4 \cdot 8\text{H}_2\text{O}$ F. Wt. 769.50	100 gm 500 gm	
3111	CAFFEINE (ANH) (as per IP) CAS 58-08-2 $\text{C}_8\text{H}_{10}\text{N}_4\text{O}_2$ F. Wt. 194.19	100 gm 500 gm 25 Kg	
148	CALAMINE (as per IP) CAS 8011-96-9	500 gm 25 Kg	
149	CALCEIN (Fluorescein Complexone), AR GRADE CAS 1461-15-0 $\text{C}_{30}\text{H}_{26}\text{N}_2\text{O}_{13}\text{Na}_2$ F. Wt. 622.55	1 gm 5 gm	
287	CALCIUM ACETATE CAS 62-54-4 $\text{C}_4\text{H}_6\text{CaO}_4$ F. Wt. 158.17	500 gm 25 Kg	
4043	CALCIUM ACETATE, AR GRADE CAS 62-54-4 $\text{C}_4\text{H}_6\text{CaO}_4$ F. Wt. 158.17	500 gm 25 Kg	
4007	CALCIUM BUTYRATE CAS 5743-36-2 $\text{C}_8\text{H}_{14}\text{CaO}_4$ F. Wt. 214.27	500 gm 25 Kg	
288	CALCIUM CARBONATE (Precipitated), (ANH.), EXTRA PURE CAS 471-34-1 CaCO_3 F. Wt. 100.09	500 gm 5 Kg 50 Kg	
289	CALCIUM CARBONATE (Precipitated), (ANH.), AR GRADE CAS 471-34-1 CaCO_3 F. Wt. 100.09	500 gm 50 Kg	
290	CALCIUM CHLORIDE (Fused), EXTRA PURE CAS 10043-52-4 CaCl_2 F. Wt. 110.98	500 gm 5 Kg 25 Kg	
291	CALCIUM CHLORIDE (Dihydrate), EXTRA PURE CAS 10035-04-8 $\text{CaCl}_2 \cdot 2\text{H}_2\text{O}$ F. Wt. 147.02	500 gm 5 Kg 50 Kg	
292	CALCIUM CHLORIDE (Dihydrate), AR GRADE CAS 10035-04-8 $\text{CaCl}_2 \cdot 2\text{H}_2\text{O}$ F. Wt. 147.02	500 gm 50 Kg	
587	CALCIUM CHLORIDE (ANH.), EXTRA PURE CAS 10043-52-4 CaCl_2 F. Wt. 110.98	500 gm	
828	CALCIUM CITRATE (Tetra), EXTRA PURE CAS 5785-44-4 $\text{Ca}_3(\text{C}_6\text{H}_5\text{O}_7)_2 \cdot 4(\text{H}_2\text{O})$ F. Wt. 570.50	500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
1701	CALCIUM-D-PANTOTHENATE (VITAMIN B) * CAS 137-08-6 $\text{C}_{18}\text{H}_{32}\text{CaN}_2\text{O}_{10}$ F. Wt. 476.53	25 gm 100 gm 1 Kg	
4044	CALCIUM FLUORIDE CAS 7789-75-5 CaF_2 F. Wt. 78.07	500 gm	
823	CALCIUM GLUCONATE (Mono) CAS 18016-24-5 $\text{C}_{12}\text{H}_{22}\text{CaO}_{14} \cdot \text{H}_2\text{O}$ F. Wt. 448.39	500 gm	
3113	CALCIUM GLYCEROPHOSPHATE CAS 17603-42-8 $\text{C}_3\text{H}_8\text{PO}_6\text{Na}$ F. Wt. 194.06	500 gm	
298	CALCIUM HYDROXIDE (Purified) CAS 1305-62-0 $\text{Ca}(\text{OH})_2$ F. Wt. 74.09	500 gm 5 Kg 50 Kg	
299	CALCIUM HYDROXIDE, AR GRADE CAS 1305-62-0 $\text{Ca}(\text{OH})_2$ F. Wt. 74.09	500 gm 50 Kg	
2915	CALCIUM HYPOPHOSPHITE	500 gm	
157	CALCIUM LACTATE, EXTRA PURE CAS 814-80-2 $\text{C}_6\text{H}_{10}\text{CaO}_6$ F. Wt. 218.22	500 gm	
159	CALCIUM NITRATE (Tetra), EXTRA PURE CAS 13477-34-4 $\text{Ca}(\text{NO}_3)_2 \cdot \text{H}_2\text{O}$ F. Wt. 236.15	500 gm 5 Kg 50 Kg	
160	CALCIUM NITRATE (Tetra), AR GRADE CAS 13477-34-4 $\text{Ca}(\text{NO}_3)_2 \cdot \text{H}_2\text{O}$ F. Wt. 236.15	500 gm 5 Kg	
695	CALCIUM OXIDE POWDER, EXTRA PURE CAS 1305-78-8 CaO F. Wt. 56.077	500 gm 25 Kg	
4046	di-CALCIUM PHOSPHATE (Dihydrate), EXTRA PURE CAS 7789-77-7 $\text{CaHPO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 172.1	500 gm 25 Kg	
296	CALCIUM PHOSPHATE DIBASIC (ANH.), EXTRA PURE CAS 7757-93-9 CaHPO_4 F. Wt. 136.06	500 gm 5 Kg	
163	CALCIUM PHOSPHATE TRIBASIC (tri-Calcium Phosphate), EXTRA PURE CAS 7758-87-4 $\text{Ca}_3(\text{PO}_4)_2$ F. Wt. 310.20	500 gm 5 Kg 50 Kg	
933	CALCIUM PROPIONATE, EXTRA PURE CAS 4075-81-4 $\text{C}_5\text{H}_{10}\text{CaO}_4$ F. Wt. 186.22	500 gm 5 Kg 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
164	CALCIUM SULPHATE (Dihydrate), EXTRA PURE CAS 10101-41-4 $\text{CaSO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 172.17	500 gm	
165	CALCIUM SULPHATE (Dihydrate), AR GRADE CAS 10101-41-4 $\text{CaSO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 172.17	500 gm	
5086	CALCIUM SULPHATE (Dihydrate), (FOOD GRADE) CAS 7778-18-9 $\text{CaSO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 172.17	500 gm 50 Kg	
2995	CALCON (Solo Chrome Dark Blue/ Eriochrome blue dark) CAS 2538-85-4 C.I. 15705 $\text{C}_{20}\text{H}_{13}\text{N}_2\text{NaO}_5\text{S}$ F. Wt. 416.38	5 gm 25 gm 100 gm	
1025	CALCONCARBOXYLIC ACID, AR GRADE (Patton & Reeder's Reagent) CAS 3737-95-9 $\text{C}_{21}\text{H}_{14}\text{N}_2\text{O}_7\text{S}$ F. Wt. 438.42	5 gm 25 gm	
4048	CALMAGITE, AR GRADE CAS 3147-14-6 $\text{C}_{17}\text{H}_{14}\text{N}_2\text{O}_5\text{S}$ F. Wt. 358.37	1 gm 5 gm	
094	CANADA BALSAM (Synthetic) CAS 8007-47-4	500 ml	
TCK 002	CARBOHYDRATE KIT (21 items of 1 gm each) Adonitol, L-Arabinose, Cellobiose, Dextrose, Dulcitol, Galactose, Fructose, Inositol, Inulin, Lactose, Maltose, Mannitol, Mannose, Melibiose, Raffinose, Rhamnose, Salicin, Sorbitol, Sucrose Trehalose, Xylose	1 Kit 5 Kit	
TCK 003	CARBOHYDRATE KIT, TYPE A (4 Amino Sugars of 1 gm each) N-Acetylgalactosamine, N-Acetylglucosamine D-Galactosamine HCl, D-Glucosamine Hcl	1 Kit	
1658	CARBOL FUCHSIN (Para Fuchsin & Phenol) CAS 4197-24-4 $\text{C}_{26}\text{H}_{26}\text{ClN}_3\text{O}$ F. Wt. 431.96	25 gm 100 gm	
293	CARBOLIC ACID, CRYSTAL (Phenol Crystal), EXTRA PURE CAS 108-95-2 $\text{C}_6\text{H}_5\text{OH}$ F. Wt. 94.11	500 gm	
170	CARBOPOL 934® CAS 9003-01-4	500 gm	
2546	CARBOPOL 940® CAS 9003-01-4	500 gm	
826	CARBOXYMETHYL CELLULOSE Sodium High Viscosity (CMC) (15 mpas) CAS 9000-11-7	500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
1659	CARMINE (Carminic Acid) CAS 1390-65-4 C.I. 75470 $\text{C}_{22}\text{H}_{20}\text{O}_{13}$ F. Wt. 492.38	5 gm 25 gm	
1721	L-CARNITINE (VITAMIN-B)* CAS 541-15-1 $\text{C}_7\text{H}_{15}\text{NO}_3$ F. Wt. 161.2	5 gm 25 gm	
209	CASEIN (Soluble in Alkali) (Fat Free, Vitamin Free) CAS 9000-71-9	500 gm	
510	CATALASE (2000-5000 Unit/Mg) (from Bovine Liver) CAS 9001-05-2	1 gm	
4231	CEDARWOOD OIL Optically clear, Nonfluorescent for microscopy CAS 8000-27-9	25 ml 30 ml 100 ml 500 ml	
4232	CEDARWOOD OIL, AR GRADE Optically clear, Nonfluorescent for microscopy CAS 8000-27-9	25 ml 30 ml 100 ml 500 ml	
939	D-CELLOBIOSE CAS 528-50-7 $\text{C}_{12}\text{H}_{22}\text{O}_{11}$ F. Wt. 342.29	5 gm 25 gm	
3523	CELLULASE * Activity 10,000 U/gm CAS 9012-54-8	1 gm	
TC 005	CELLULOSE POWDER CAS 9004-34-6 $(\text{C}_6\text{H}_{10}\text{O}_5)_n$	500 gm 25 Kg	
2550	CETOSTEARYL ALCOHOL, EXTRA PURE CAS 8005-44-5	500 gm	
584	CETRIMIDE (N-Cetyl N,N,N- Trimethylammonium Bromide), EXTRA PURE CAS 57-09-0 $\text{C}_{19}\text{H}_{42}\text{NBr}$ F. Wt. 364.62	100 gm 500 gm 25 Kg	
2556	CETRIMIDE (N-Cetyl N,N,N- Trimethylammonium Bromide), AR GRADE CAS 57-09-0 $\text{C}_{19}\text{H}_{42}\text{NBr}$ F. Wt. 364.45	100 gm 500 gm	
2553	CETYL ALCOHOL (1-Hexadecanol), EXTRA PURE CAS 36653-82-4 $\text{C}_{16}\text{H}_{34}\text{O}$ F. Wt. 242.44	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
3316	CHAPS BUFFER CAS 75621-03-3 $C_{35}H_{58}N_2O_7S$ F. Wt. 614.9	1 gm 5 gm	
203	CHARCOAL ACTIVATED POWDER, EXTRA PURE CAS 7440-44-0 F. Wt. 12.01	500 gm 25 Kg	
2903	CHARCOAL ACTIVATED, AR GRADE CAS 7440-44-0 F. Wt. 12.01	500 gm	
2931	CHARCOAL ACTIVATED, GRANULAR SPECIAL GRADE 2.0-5.0 mm CAS 7440-44-0 F. Wt. 12.01	500 gm	
4255	CHES BUFFER CAS 103-47-9	25 gm 100 gm	
827	CHITIN [Poly (N-acetyl-1,4-D-glucopyranosamine)] CAS 1398-61-4 ($C_8H_{13}NO_5$) _n	100 gm 500 gm 25 Kg	
4224	CHITOSAN CAS 9012-76-4 ($C_6H_{11}NO_4$) _n F. Wt. 309.54	25 gm 100 gm 1 Kg	
183	CHLORAMINE-T (Trihydrate), AR GRADE CAS 7080-50-4 $C_7H_{13}N_3ClNNaO_5S$ F. Wt. 281.69	500 gm	
2559	CHLOROACETIC ACID, EXTRA PURE (Monochloroacetic Acid) CAS 79-11-8 $C_2H_3ClO_2$ F. Wt. 94.50	500 gm	
2560	CHLOROAUIC ACID, Min. 49% Au (Gold Chloride) CAS 16961-25-4 $HAuCl_4 \cdot 3H_2O$ F. Wt. 339.79	1 gm	
4354	2-CHLOROETHYL PHOSPHONIC ACID CAS 16672-87-0 $C_2H_5ClO_3P$ F. Wt. 144.49	100 ml 500 ml	
182	CHLOROPHENOL RED, AR GRADE CAS 4430-20-0 $C_{19}H_{12}Cl_2O_5S$ F. Wt. 423.27	5 gm 25 gm	
2565	p-CHLOROPHENOXY ACETIC ACID (4 CPA) CAS 122-88-3 $C_8H_7ClO_3$ F. Wt. 186.59	25 gm 100gm 500gm	
1707	CHOLECALCIFEROL (Vitamin D) * CAS 67-97-0 $C_{27}H_{44}O$ F. Wt. 384.64	1 gm	

CODE	PRODUCT NAME	PACK SIZE	
305	CHOLESTEROL (as per USP), EXTRA PURE * CAS 57-88-5 $C_{27}H_{46}O$ F. Wt. 386.65	100 gm 500 gm	
4053	CHOLESTEROL, AR GRADE CAS 57-88-5 $C_{27}H_{46}O$ F. Wt. 386.65	5 gm 25 gm 100 gm	
307	CHOLIC ACID Purity 99%+ CAS 81-25-4 $C_{26}H_{40}O_5$ F. Wt. 408.57	25 gm 100 gm	
4178	CHOLIC ACID SODIUM CAS 361-09-1 $C_{26}H_{39}NaO_5$ F. Wt. 430.56	25 gm 100 gm 5 Kg	
3125	CHROMAZUROL-S CAS 1667-99-8 $C_{23}H_{13}Cl_2Na_3O_9S$ F. Wt. 605.28	10 gm	
3126	CHROMIUM ACETATE, EXTRA PURE CAS 1066-30-4 $Cr_6H_{16}Cr_2O_{10}$ F. Wt. 376.2	500 gm	
185	CHROMIUM (III) CHLORIDE (Hexa), EXTRA PURE CAS 10060-12-5 $Cl_3Cr \cdot 6H_2O$ F. Wt. 266.45	500 gm 25 Kg	
187	CHROMIUM (III) NITRATE (Nona), EXTRA PURE CAS 7789-02-8 $CrN_3O_9 \cdot 9H_2O$ F. Wt. 400.15	500 gm	
189	CHROMIUM PICOLINATE, EXTRA PURE CAS 14639-25-9 $C_{18}H_{12}N_3O_6Cr$ F. Wt. 418.33	25 gm 100 gm	
2570	CHROMIUM (III) SULPHATE BASIC, EXTRA PURE CAS 39380-78-4 $Cr_4(SO_4)_5(OH)_2$ F. Wt. 722.32	500 gm	
4054	CHROMIUM OXIDE GREEN, EXTRA PURE CAS 1308-38-9 Cr_2O_3 F. Wt. 151.99	500 gm	
192	CHROMIUM TRIOXIDE, EXTRA PURE (Chromium (VI) Oxide/Chromic acid) CAS 1333-82-0 CrO_3 F. Wt. 99.99	500 gm	
196	CHROMOTROPIC ACID DISODIUM (Dihydrate), AR GRADE CAS 5808-22-0 $C_{10}H_6Na_2O_8S_2 \cdot 2H_2O$ F. Wt. 400.29	10 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
199	CINNAMIC ACID, (β -phenyl acrylic acid) EXTRA PURE CAS 140-10-3 $C_9H_8O_2$ F. Wt. 148.05	100 gm 500 gm	
309	CITRIC ACID (Mono), EXTRA PURE CAS 5949-29-1 $C_6H_8O_7 \cdot H_2O$ F. Wt. 210.14	500 gm 5 Kg 25 Kg	
310	CITRIC ACID (Mono), AR GRADE CAS 5949-29-1 $C_6H_{10}O_8$ F. Wt. 210.14	500 gm 5 Kg 25 Kg	
311	CITRIC ACID (ANH.), EXTRA PURE CAS 77-92-9 $C_6H_8O_7$ F. Wt. 192.12	500 gm 5 Kg 25 Kg	
312	CITRIC ACID (ANH.), AR GRADE CAS 77-92-9 $C_6H_8O_7$ F. Wt. 192.12	500 gm 5 Kg 25 Kg	
975	COBALT (II) ACETATE (Tetra), EXTRA PURE CAS 6147-53-1 $Co(CH_3COO)_4 \cdot 4H_2O$ F. Wt. 249.08	100 gm 500 gm	
1801	COBALT (II) CARBONATE BASIC, EXTRA PURE CAS 12602-23-2 $CoCO_3 \cdot x H_2O$ F. Wt. 516.73	100 gm 500 gm	
976	COBALT (II) CHLORIDE (Hexa), EXTRA PURE CAS 7791-13-1 $CoCl_2 \cdot 6(H_2O)$ F. Wt. 237.93	100 gm 500 gm	
977	COBALT (II) NITRATE (Hexa), EXTRA PURE CAS 10026-22-9 $Co(NO_3)_2 \cdot 6H_2O$ F. Wt. 291.03	100 gm 500 gm	
1803	COBALT (II) NITRATE (Hexa), AR GRADE CAS 10026-22-9 $Co(NO_3)_2 \cdot 6H_2O$ F. Wt. 291.03	100 gm 500 gm	
1804	COBALT (II) OXIDE, EXTRA PURE CAS 1308-06-1 Co_3O_4 F. Wt. 240.80	100 gm 500 gm	
978	COBALT (II) SULPHATE (Hepta), EXTRA PURE CAS 10026-24-1 $Co_3SO_4 \cdot 7H_2O$ F. Wt. 281.10	100 gm 500 gm	
1805	COLCHICINE CAS 64-86-8 $C_{22}H_{25}NO_6$ F. Wt. 399.44	1 gm 10 gm	

CODE	PRODUCT NAME	PACK SIZE	
1615	CONGO RED CAS 573-58-0 C.I. 22120 $C_{32}H_{22}N_6Na_2O_6S_2$ F. Wt. 696.67	25 gm 100 gm	
1696	COOMASSIE BRILLIANT BLUE G 250 (Brilliant Blue G-250) For electrophoresis CAS 6104-58-1 C.I. 42655 $C_{47}H_{48}NaO_7S_2$ F. Wt. 854.02	5 gm 25 gm	
1616	COOMASSIE BRILLIANT BLUE R-250 (Brilliant Blue R-250) CAS 6104-59-2 C.I. 42660 $C_{45}H_{44}N_3NaO_7S_2$ F. Wt. 825.97	5 gm 25 gm	
2409	COPPER METAL POWDER 99.5% (325 mesh) Electrolytic Grade CAS 7440-50-8 Cu F. Wt. 63.55	500 gm	
1687	COTTON BLUE, (Water Blue) CAS 28983-56-4 C.I. 42780 $C_{37}H_{27}N_3Na_2O_9S_3$ F. Wt. 799.82	25 gm 100 gm	
2574	COUMARIN [1-Benzopyran-2-One] CAS 91-64-5 $C_9H_6O_2$ F. Wt. 146.143	100 gm 500 gm	
1807	CREATINE (Mono) CAS 6020-87-7 $C_4H_9N_3O_2 \cdot H_2O$ F. Wt. 149.16	25 gm 100 gm	
1808	CREATININE, AR GRADE CAS 60-27-5 $C_4H_7N_3O$ F. Wt. 113.12	25 gm 100 gm 1 Kg	
1811	o-CRESOLPHTHALEIN (pH Indicator) CAS 596-27-0 $C_{22}H_{18}O_4$ F. Wt. 346.38	5 gm 25 gm	
4137	o-CRESOLPHTHALEIN COMPLEXONE CAS 2411-89-4 $C_{32}H_{32}N_2O_{12}$ F. Wt. 636.6	1 gm 5 gm	
1618	m-CRESOL PURPLE CAS 2303-01-7 $C_{21}H_{18}O_5S$ F. Wt. 382.43	1 gm 5 gm	
1619	CRESOL RED, AR GRADE CAS 1733-12-6 $C_{21}H_{18}O_5S$ F. Wt. 382.44	5 gm 25 gm	
1606	CRYSTAL VIOLET POWDER, EXTRA PURE (Gentian Violet or Methyl Violet) CAS 548-62-9 C.I. 42535 $C_{25}H_{30}ClN_3$ F. Wt. 407.99	25 gm 100 gm 1 Kg	

CODE	PRODUCT NAME	PACK SIZE	
1620	CRYSTAL VIOLET, AR GRADE (Gentian Violet or Methyl Violet) CAS 548-62-9 C.I. 42555 $C_{25}H_{30}ClN_3$ F. Wt. 407.98	25 gm 100 gm 1 Kg	
314	CUPRIC ACETATE (Mono), EXTRA PURE CAS 6046-93-1 $C_4H_6CuO_4 \cdot H_2O$ F. Wt. 199.64	100 gm 500 gm	
1866	CUPRIC CARBONATE BASIC, EXTRA PURE CAS 12069-69-1 $CH_2Cu_2O_5$ F. Wt. 221.12	500 gm	
315	CUPRIC CHLORIDE (Dihydrate), EXTRA PURE CAS 10125-13-0 $CuCl_2 \cdot 2H_2O$ F. Wt. 170.48	500 gm	
316	CUPRIC CHLORIDE (Dihydrate), AR GRADE CAS 10125-13-0 $CuCl_2 \cdot 2H_2O$ F. Wt. 170.48	500 gm	
317	CUPRIC NITRATE (Trihydrate), EXTRA PURE CAS 10031-43-3 $Cu(NO_3)_2 \cdot 3H_2O$ F. Wt. 241.60	500 gm	
1818	CUPRIC OXIDE, EXTRA PURE CAS 1317-38-0 CuO F. Wt. 79.55	100 gm 500 gm	
319	CUPRIC SULPHATE (Penta), EXTRA PURE CAS 7758-99-8 $CuSO_4 \cdot 5H_2O$ F. Wt. 249.68	500 gm 25 Kg	
320	CUPRIC SULPHATE (Penta), AR GARDE CAS 7758-99-8 $CuSO_4 \cdot 5H_2O$ F. Wt. 249.68	500 gm 25 Kg	
1819	CUPRIC SULPHATE (ANH.), EXTRA PURE CAS 7758-98-7 $CuSO_4$ F. Wt. 159.61	500 gm 25 Kg	
837	CUPRIC SULPHATE (ANH.), AR GRADE CAS 7758-98-7 $CuSO_4$ F. Wt. 159.61	500 gm 25 Kg	
2893	CUPROUS CHLORIDE, EXTRA PURE CAS 7758-89-6 $CuCl$ F. Wt. 98.999	500 gm	
817	CUPRON (Benzoin α -Oxime), AR GRADE CAS 441-38-3 $C_{14}H_{13}NO_2$ F. Wt. 227.27	25 gm 100 gm	
2936	CURCUMINE (Turmeric) CAS 458-37-7 $C_{21}H_{20}O_6$ F. Wt. 368.39	5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
1704	CYANOCOBALAMIN (Vitamin B) CAS 68-19-9 $C_{63}H_{88}CON_{14}O_{14}P$ F. Wt. 1355.38	1 gm 10 gm	
3625	β -CYCLODEXTRIN CAS 7585-39-9 $C_{42}H_{70}O_{35}$ F. Wt. 1134.98	500 gm	
573	CYCLOHEXIMIDE (Actidione), AR GRADE CAS 66-81-9 $C_{15}H_{23}NO_4$ F. Wt. 281.36	1 gm 5 gm	
1138	L-CYSTEINE (Free Base) 99% + Crystalline CAS 52-90-4 $C_3H_7NO_2S$ F. Wt. 121.17	5 gm 25 gm 100 gm	
1139	L-CYSTEINE HYDROCHLORIDE (Mono) 99% + Crystalline CAS 7048-04-6 $C_3H_8NO_2SCI \cdot H_2O$ F. Wt. 175.63	25 gm 100 gm 1 Kg 25 Kg	
1140	L-CYSTINE, 99% + CRYSTALLINE (Free Base) CAS 56-89-3 $C_6H_{12}N_2O_4S_2$ F. Wt. 240.3	25 gm 100 gm 1 Kg	
4056	CYTIDINE (for Biochemistry) CAS 65-46-3 $C_9H_{13}N_3O_5$ F. Wt. 243.22	1 gm 5 gm 25 gm	
4057	CYTOSINE (For Biochemistry) CAS 71-30-7 $C_4H_5N_3O$ F. Wt. 111.1	25 gm 100 gm 500 gm	
4058	1-DECANE SULPHONIC ACID SODIUM SALT (ANH), (For HPLC) CAS 13419-61-9 $CH_3(CH_2)_9SO_3Na$ F. Wt. 244.33	25 gm	
4059	1-DECANE SULPHONIC ACID SODIUM SALT (Mono), (For HPLC) CAS 13419-61-9 $C_{10}H_{21}NaO_3S \cdot H_2O$ F. Wt. 178.16	25 gm	
666	DETERGENT TABOLENE (Neutral pH, W/O Phosphate) (For glassware cleaning)	500 ml 5 Ltr	
1823	DEVARDA'S ALLOY POWDER, EXTRA PURE CAS 8049-11-4	100 gm 500 gm	
1313	DEXTRIN WHITE CAS 9004-53-9 $(C_6H_{10}O_5)_n \cdot xH_2O$ F. Wt. 414.55	500 gm	
1303	DEXTROSE (Mono), EXTRA PURE CAS 5996-10-1 $C_6H_{12}O_6 \cdot H_2O$ F. Wt. 198.17	500 gm 5 Kg 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
1317	DEXTROSE (Mono), AR GRADE CAS 5996-10-1 $C_6H_{12}O_6 \cdot H_2O$ F. Wt. 198.17	500 gm 5 Kg 25 Kg	
1304	DEXTROSE (ANH.), EXTRA PURE CAS 50-99-7 $C_6H_{12}O_6$ F. Wt. 180.16	500 gm 5 Kg 25 Kg	
1314	DEXTROSE (ANH.), AR GRADE CAS 50-99-7 $C_6H_{12}O_6$ F. Wt. 180.16	500 gm 5 Kg 25 Kg	
1825	DIACETYL MONOXIME, AR GRADE (2,3, Butanedione Monoxime) CAS 57-71-6 $C_4H_7NO_2$ F. Wt. 101.10	25 gm 100 gm	
1417	DIASTASE * (α -Amylase, Activity 2000 A U/g) CAS 9000-90-2	100 gm 500 gm 5 Kg	
1827	trans-1,2-DIAMINO CYCLOHEXANE N,N,N',N'-TETRA ACETIC ACID (CDTA), EXTRA PURE CAS 13291-61-7 $C_{14}H_{22}N_2O_8 \cdot H_2O$ F. Wt. 364.36	25 gm 100 gm	
3136	2,7 DICHLORO FLUORESCIN, AR GRADE CAS 76-54-0 $C_{20}H_{10}Cl_2O_5$ F. Wt. 401.2	10 gm	
2589	2,6-DICHLOROPHENOL INDOPHENOL SODIUM, AR GRADE CAS 620-45-1 $C_{12}H_6Cl_2NNaO_2$ F. Wt. 290.08	1 gm 5 gm 25 gm	
2590	2,4-DICHLOROPHENOXY ACETIC ACID (2,4-D) CAS 94-75-7 $C_8H_6Cl_2O_3$ F. Wt. 221.04	500 gm	
2591	2,6-DICHLOROQUINONE-4-CHLORIMIDE * (Gibbs Reagent) CAS 101-38-2 $C_6H_2Cl_3NO$ F. Wt. 210.45	5 gm 25 gm	
2942	DIETHYLENE GLYCOL, EXTRA PURE CAS 111-46-6 $C_4H_{10}O_3$ F. Wt. 106.12	500 ml 25 Kg	
2944	DIETHYLENE GLYCOL, AR GRADE CAS 111-46-6 $C_4H_{10}O_3$ F. Wt. 106.12	500 ml	

CODE	PRODUCT NAME	PACK SIZE	
754	DIGITONIN, AR GRADE CAS 11024-24-1 $C_{56}H_{92}O_{29}$ F. Wt. 1229.31	1 gm	
839	L-3, 4-DIHYDROXYPHENYL L-ALANINE (L-DOPA) CAS 59-92-7 $C_9H_{11}NO_4$ F. Wt. 197.19	5 gm 25 gm	
596	p -DIMETHYLAMINO BENZALDEHYDE (Ehrlich Reagent), EXTRA PURE CAS 100-10-7 $C_9H_{11}NO$ F. Wt. 149.19	25 gm 100 gm 500 gm 25 Kg	
795	p-DIMETHYL AMINO BENZALDEHYDE (Ehrlich Reagent), AR GRADE CAS 100-10-7 $C_9H_{11}NO$ F. Wt. 149.19	25 gm 100 gm 500 gm	
1836	DIMETHYL GLYOXIME, EXTRA PURE (2,3, butanedionedioxime) CAS 95-45-4 $C_4H_8N_2O_2$ F. Wt. 116.12	100 gm 500 gm	
1840	N,N-DIMETHYL p-PHENYLENEDIAMINE di-HCL, AR GRADE CAS 536-46-9 $C_8H_{14}Cl_2N_2$ F. Wt. 209.12	25 gm	
796	N,N-DIMETHYL-p-PHENYLENE DIAMINE OXALATE * (p-Amino-N,N dimethylaniline oxalate) CAS 62778-12-5 $C_{10}H_{14}N_2O_4$ F. Wt. 226.23	1 gm 5 gm	
841	N,N-DIMETHYL-p-PHENYLENE DIAMINE SULPHATE (4-Amino-N,N dimethylaniline sulphate) CAS 536-47-0 $C_8H_{14}N_2SO_4$ F. Wt. 234.87	5 gm 25 gm	
1842	DIMETHYL YELLOW (pH Indicator) CAS 60-11-7 C.I. 11020 $C_{14}H_{15}N_3$ F. Wt. 225.29	25 gm 100 gm	
4060	DIMIDIUM BROMIDE 98% (For Tenside Testing) CAS 518-67-2 $C_{20}H_{18}BrN_3$ F. Wt. 380.28	100 mg 1 gm	
2941	3,5 DINITROBENZOIC ACID, EXTRA PURE CAS 99-34-3 $C_7H_4N_2O_6$ F. Wt. 212.12	100 gm 500 gm	
1844	3,5- DINITROBENZOIC ACID, AR GRADE CAS 99-34-3 $C_7H_4N_2O_6$ F. Wt. 212.12	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
1845	2,4-DINITROPHENYL HYDRAZINE, EXTRA PURE CAS 119-26-6 $C_6H_6N_4O_4$ F. Wt. 198.14	25 gm 100 gm	
2599	3,5-DINITROSALICYLIC ACID, AR GRADE CAS 609-99-4 $C_7H_4N_2O_7$ F. Wt. 228.12	25 gm 100 gm	
4061	DIOCTYL SODIUM SULPHOSUCCINATE 98% (DOSS) (AEROSOL OT) CAS 577-11-7 $C_{20}H_{37}NaO_7S$ F. Wt. 444.56	500 gm	
3138	DIPHENYLAMINE, EXTRA PURE CAS 122-39-4 $C_{12}H_{11}N$ F. Wt. 169.22	100 gm 500 gm	
1849	DIPHENYLAMINE, AR GARDE CAS 122-39-4 $C_{12}H_{11}N$ F. Wt. 169.22	100 gm 500 gm	
2601	N,N'-DIPHENYL BENZIDINE, AR GRADE CAS 531-91-9 $C_{24}H_{20}N_2$ F. Wt. 336.43	5 gm	
2602	1,5-DIPHENYL CARBAZIDE, AR GRADE * CAS 140-22-7 $C_{13}H_{14}N_2O$ F. Wt. 242.28	25 gm	
4062	1,5-DIPHENYL CARBAZONE, AR GRADE (Reagent For Mercury) CAS 538-62-5 $C_{13}H_{14}N_2O$ F. Wt. 242.26	5 gm 25 gm	
4359	DISODIUM OCTABORATE (BORAX 20) CAS 12008-41-2 $B_8Na_2O_{24}$ F. Wt. 340.50	500 gm	
1850	1,4-DITHIOERYTHRITOL (D.T.E.) * (Cleland Reagent) CAS 6892-68-8 $C_4H_{10}O_2S_2$ F. Wt. 154.25	1 gm 5 gm	
1852	DL-DITHIOTHREITOL (D.T.T.) * CAS 3483-12-3 $C_4H_{10}O_2S_2$ F. Wt. 154.25	1 gm 5 gm	
2998	DITHIZONE (Diphenylthiocarbazone), EXTRA PURE CAS 60-10-6 $C_{13}H_{12}N_4S$ F. Wt. 256.33	5 gm 25 gm	
1853	DITHIZONE, AR GRADE CAS 60-10-6 $C_{13}H_{12}N_4S$ F. Wt. 256.33	5 gm 25 gm	
755	DNA * (Deoxyribonucleic Acid) CAS 9007-49-2	5 gm 25 gm	
1023	DODECYL BENZENE SULPHONIC ACID SODIUM SALT CAS 25155-30-0 $C_{18}H_{29}NaO_3S$ F. Wt. 326.49	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
3522	DPX MOUNTANT for Histology	250 ml	
TC 008	DULCITOL (Galactitol, Dulcite) CAS 608-66-2 $C_6H_{14}O_6$ F. Wt. 182.17	25 gm 100 gm	
349	EDTA DIPOTASSIUM, EXTRA PURE CAS 25102-12-9 $C_{10}H_{14}K_2N_2O_8 \cdot 2H_2O$ F. Wt. 404.46	100 gm 500 gm	
326	EDTA DIPOTASSIUM, AR GRADE CAS 25102-12-9 $C_{10}H_{14}K_2N_2O_8 \cdot 2H_2O$ F. Wt. 404.46	100 gm 500 gm 25 Kg	
327	EDTA DISODIUM, EXTRA PURE CAS 6381-92-6 $C_{10}H_{14}K_2Na_2O_8 \cdot 2H_2O$ F. Wt. 372.24	100 gm 500 gm 5 Kg 25 Kg	
328	EDTA DISODIUM, AR GRADE CAS 6381-92-6 $C_{10}H_{14}N_2Na_2O_8 \cdot 2H_2O$ F. Wt. 372.24	100 gm 500 gm 5 Kg 25 Kg	
973	EDTA FERRIC MONOSODIUM, EXTRA PURE CAS 15708-41-5 $C_{10}H_{12}FeN_2NaO_8$ F. Wt. 367.05	100 gm 500 gm	
2608	EDTA MAGNESIUM DISODIUM (Tetra), EXTRA PURE CAS 14402-88-1 $C_{10}H_{12}N_2O_8MgNa_2$ F. Wt. 358.50	500 gm	
596	EHRlich's REAGENT CAS 100-10-7 $C_9H_{11}NO$ F. Wt. 149.19	25 gm 100 gm 500 gm 25 Kg	
1661	EOSINE (Spirit Soluble) CAS 6359-05-3 C.I. 45386 $C_{22}H_{11}Br_4KO_5$ F. Wt. 714.03	25 gm 100 gm	
1669	EOSINE YELLOW (Water Soluble) CAS 548-26-5 C.I. 45380 $C_{20}H_6Br_4Na_2O_5$ F. Wt. 691.85	25 gm 100 gm 1 Kg	
4065	EOSIN YELLOW (Water Soluble) CAS 17372-87-1 C.I. 45380 $C_{20}H_6Br_4Na_2O_5$ F. Wt. 691.88	25 gm 100 gm	
1623	ERIOCHROME BLACK T, EXTRA PURE CAS 1787-61-7 C.I. 14645 $C_{20}H_{12}N_3NaO_7S$ F. Wt. 461.38	25 gm 100 gm	
1855	ERIOCHROME BLACK-T, AR GRADE CAS 1787-61-7 C.I. 14645 $C_{20}H_{12}N_3NaO_7S$ F. Wt. 461.38	25 gm 100 gm	
1856	ERIOCHROME CYANINE-R, AR GRADE CAS 3564-18-9 C.I. 43820 $C_{23}H_{15}Na_3O_9S$ F. Wt. 536.39	5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
1681	ERIOGLAUCINE (Acid Blue 9) (Brilliant Blue FCF) CAS 3844-45-9 C.I. 42090 $C_{37}H_{34}N_2Na_2O_9S_3$ F. Wt. 792.86	10 gm	
TC 009	ERYTHRITOL (Meso Erythritol) * CAS 149-32-6 $C_4H_{10}O_4$ F. Wt. 122.12	5 gm 25 gm 100 gm	
1859	ERYTHROSIN B SODIUM CAS 16423-68-0 C.I. 45430 $C_{20}H_{16}I_4Na_2O_5$ F. Wt. 879.86	25 gm	
304	ESCULIN (Aesculin) CAS 531-75-9 $C_{15}H_{16}O_9$ F. Wt. 367.30	5 gm 25 gm	
3141	ETHIDIUM BROMIDE, EXTRA PURE CAS 1239-45-8 $C_{21}H_{20}BrN_3$ F. Wt. 394.31	1 gm 5 gm	
4067	ETHIDIUM BROMIDE, AR GRADE (For Electrophoresis)	1 gm 5 gm	
324	ETHYLENEDIAMINE TETRAACETIC ACID, (EDTA), EXTRA PURE CAS 60-00-4 $C_{10}H_{16}N_2O_8$ F. Wt. 292.24	100 gm 500 gm	
325	ETHYLENE DIAMINE TETRA ACETIC ACID, (EDTA), AR GRADE CAS 60-00-4 $C_{10}H_{16}N_2O_8$ F. Wt. 292.24	100 gm 500 gm	
2610	ETHYL-P-HYDROXY BENZOATE (Ethyl Paraben) CAS 120-47-8 $C_9H_{10}O_3$ F. Wt. 166.17	500 gm 5 Kg	
2940	ETHYLENE GLYCOL, EXTRA PURE CAS 107-21-1 $C_2H_6O_2$ F. Wt. 62.07	500 ml	
2943	ETHYLENE GLYCOL, AR GRADE CAS 107-21-1 $C_2H_6O_2$ F. Wt. 62.07	500 ml	
1625	EVANS BLUE CAS 314-13-6 C.I. 23860 $C_{34}H_{24}N_6Na_4O_{14}S_4$ F. Wt. 960.80	5 gm 25 gm	
2613	FAST BLUE B SALT (for Microscopy) CAS 14263-94-6 C.I. 37235 $C_{14}H_{12}Cl_2N_4O_2 \cdot ZnCl_2$ F. Wt. 475.49	25 gm 100 gm	
1627	FAST BLUE BB SALT, AR GRADE CAS 5486-84-0 C.I. 37175 $C_{34}H_{36}Cl_4N_6O_6Zn$ F. Wt. 831.89	25 gm	

CODE	PRODUCT NAME	PACK SIZE	
1630	FAST GREEN FCF CAS 2353-45-9 C.I. 42053 $C_{37}H_{34}N_2Na_2O_{10}S_3$ F. Wt. 808.85	5 gm 25 gm	
4223	FE-EDDHA (ETHYLENEDIAMINE-N,N'-BIS (2-HYDROXYPHENYLACETATE ACID) FERRIC-SODIUM COMPLEX CAS 16455-61-1 $C_{18}H_{18}FeN_2NaO_6$ F. Wt. 437.19	25 gm 100 gm	
TR 052	FEHLING SOLUTION A	500 ml	
TR 053	FEHLING SOLUTION B	500 ml	
330	FERRIC AMMONIUM CITRATE (BROWN), EXTRA PURE CAS 1185-57-5 $C_{12}H_{22}FeN_3O_{14}$ F. Wt. 488.16	500 gm 5 Kg 25 Kg	
331	FERRIC CHLORIDE (ANH.), EXTRA PURE CAS 7705-08-0 $FeCl_3$ F. Wt. 162.20	500 gm 5 Kg 50 Kg	
332	FERRIC CHLORIDE (ANH.), AR GRADE CAS 7705-08-0 $FeCl_3$ F. Wt. 162.20	500 gm 50 Kg	
756	FERRIC CITRATE (Mono) (for Bacteriology) CAS 2338-05-8 $C_6H_5FeO_7$ F. Wt. 244.94	500 gm	
4068	FERRIC NITRATE, EXTRA PURE CAS 7782-61-8 FeN_3O_9 F. Wt. 241.86	500 gm	
228	FERRIC ORTHOPHOSPHATE (Tetra), EXTRA PURE CAS 10045-86-0 $FePO_4$ F. Wt. 150.82	500 gm	
990	FERRIC OXIDE RED (Iron Oxide Red), EXTRA PURE CAS 1309-37-1 Fe_2O_3 F. Wt. 159.69	500 gm	
845	FERRIC SULPHATE (Hydrate), EXTRA PURE CAS 10028-22-5 $Fe_2(SO_4)_3 \cdot xH_2O$ F. Wt. 399.87	500 gm	
243	FERROUS AMMONIUM SULPHATE, EXTRA PURE CAS 7783-85-9 $(NH_4)_2Fe(SO_4)_2 \cdot 6H_2O$ F. Wt. 392.16	500 gm 50 Kg	
244	FERROUS AMMONIUM SULPHATE, AR GRADE CAS 7783-85-9 $FeH_{20}N_2O_{14}S_2$ F. Wt. 392.14	500 gm 50 Kg	
513	FERROUS SULPHATE (Hepta), EXTRA PURE CAS 7782-63-0 $FeSO_4 \cdot 7H_2O$ F. Wt. 278.01	500 gm 5 Kg 50 Kg	
514	FERROUS SULPHATE (Hepta), AR GRADE CAS 7782-63-0 $FeSO_4 \cdot 7H_2O$ F. Wt. 278.01	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
1879	FERROUS SULPHIDE STICKS (Iron sulphide sticks) CAS 1317-37-9 FeS F. Wt. 87.91	1 Kg	
1682	FIELD'S STAIN-A, POWDER	25 gm 1 Kg	
1683	FIELD'S STAIN-B, POWDER	25 gm 1 Kg	
1880	FLUORESCIN CAS 2321-07-5 C.I. 45350 C ₂₀ H ₁₂ O ₅ F. Wt. 332.31	25 gm 100 gm	
149	FLUORESCIN COMPLEXONE (Calcein), AR GRADE CAS 1461-15-0 C ₃₀ H ₂₆ N ₂ O ₁₃ F. Wt. 622.55	1 gm 5 gm	
2616	FLUORESCIN SODIUM CAS 518-47-8 C ₂₀ H ₁₀ Na ₂ O ₅ F. Wt. 376.27	25 gm 100 gm	
3147	5-FLUOROURACIL CAS 51-21-8 C ₄ H ₅ FN ₂ O ₂ F. Wt. 130.08	1 gm 5 gm	
335	FOLIC ACID (Vitamin B) * CAS 59-30-3 C ₁₉ H ₁₉ N ₇ O ₆ F. Wt. 441.40	5 gm 25 gm	
4212	FORCHLORFENURON (4-CPPU) N-(2-chloro-4-pyridinyl)-N'- phenyl CAS 68157-60-8 C ₁₂ H ₁₀ ClN ₃ O F. Wt. 247.68	5 gm	
1301	D-FRUCTOSE (Levulose), EXTRA PURE CAS 57-48-7 C ₆ H ₁₂ O ₆ F. Wt. 180.16	100 gm 500 gm 25 Kg	
1315	D-FRUCTOSE (Levulose) (Bacteriological Grade) CAS 57-48-7 C ₆ H ₁₂ O ₆ F. Wt. 180.16	100 gm 500 gm	
1640	FUCHSIN ACID (Acid Fuchsin) CAS 3244-88-0 C.I. 42685 C ₂₀ H ₁₇ N ₃ Na ₂ O ₉ S ₃ F. Wt. 585.55	25 gm 100 gm	
1603	FUCHSIN BASIC (Rosaniline HCl) (Basic Fuchsin) CAS 632-99-5 C.I. 42510 C ₂₀ H ₂₀ N ₃ .HCl F. Wt. 337.86	25 gm 100 gm 1 Kg	
771	FUMARIC ACID, EXTRA PURE CAS 110-17-8 C ₄ H ₄ O ₄ F. Wt. 116.07	500 gm	
4071	FUSION MIXTURE	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
TC 010	D-GALACTOSE CAS 59-23-4 C ₆ H ₁₂ O ₄ F. Wt. 180.16	25 gm 100 gm 500 gm	
337	GALLIC ACID (Mono) (3,4,5, Trihydroxybenzoic Acid) CAS 5995-86-8 C ₇ H ₆ O ₅ .H ₂ O F. Wt. 188.14	100 gm 500 gm	
3318	GELLAN GUM (Plant Tissue Culture Grade) CAS 71010-52-1	100 gm 500 gm 1 Kg	
4072	GELATIN POWDER (for Bacteriology) CAS 9000-70-8	500 gm	
1606	GENTIAN VIOLET, EXTRA PURE (Crystal violet or Methyl Violet) CAS 548-62-9 C.I. 42555 C ₂₅ H ₃₀ ClN ₃ F. Wt. 407.99	25 gm 100 gm 1 Kg	
1620	GENTIAN VIOLET, AR GRADE (Crystal Violet or Methyl Violet) CAS 548-62-9 C.I. 42555 C ₂₅ H ₃₀ ClN ₃ F. Wt. 407.98	25 gm 100 gm 1 Kg	
849	GIBBERELIC ACID (GA) (Plant Growth Regulator) CAS 77-06-5 C ₁₉ H ₂₂ O ₆ F. Wt. 346.37	1 gm 10 gm 100 gm 1 Kg	
3547	GIBBERELIC ACID (GA) CAS 202467-69-4 C ₂₀ H ₁₉ N ₃ O ₇ S F. Wt. 445.45	1 gm 5 gm 100 gm	
2591	GIBBS REAGENT * (2,6, Dichloroquinone-4-chlorimide) CAS 101-38-2 C ₆ H ₂ Cl ₃ NO F. Wt. 210.45	5 gm 25 gm	
1672	GIEMSA'S STAIN (Azur Eosin Methylene Blue) CAS 51811-82-6	10 gm 25 gm 100 gm 1 Kg	
1303	D-GLUCOSE (Mono) (Dextrose) CAS 5996-10-1 C ₆ H ₁₂ O ₆ .H ₂ O F. Wt. 198.17	500 gm 5 Kg 25 Kg	
1304	D-GLUCOSE (ANH.) (Dextrose) CAS 50-99-7 C ₆ H ₁₂ O ₆ F. Wt. 180.15	500 gm 5 Kg 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
1145	L-GLUTAMIC ACID 99%, EXTRA PURE (2- Amino Glutaric Acid) CAS 56-86-0 $C_5H_9NO_4$ F. Wt. 147.13	100 gm 500 gm	
1111	L-GLUTAMIC ACID, 99% + Crystalline, AR GRADE CAS 56-86-0 $C_5H_9NO_4$ F. Wt. 147.13	100 gm 500 gm 25 Kg	
1112	L-GLUTAMIC MONO SODIUM, AR GRADE (Sodium Glutamate) 99% + Crystalline CAS 142-47-2 $C_5H_8NNaO_4$ F. Wt. 169.11	100 gm 500 gm 25 Kg	
1114	L-GLUTAMINE, 99% + Crystalline CAS 56-85-9 $C_5H_{10}N_2O_3$ F. Wt. 146.14	25 gm 100 gm 500 gm	
793	GLUTARIC ACID * CAS 110-94-1 $C_5H_8O_4$ F. Wt. 132.12	100 gm 500 gm	
850	GLUTATHIONE (Oxidized) * CAS 27025-41-8 $C_{20}H_{32}N_6O_{12}S_2$ F. Wt. 612.60	500 mg	
757	GLUTATHIONE (Reduced) * CAS 70-18-8 $C_{10}H_{17}N_5O_6S$ F. Wt. 307.32	1 gm 5 gm 25 gm	
598	GLYCEROL (Glycerine), EXTRA PURE CAS 56-81-5 $C_3H_8O_3$ F. Wt. 92.09	500 ml 2.5 Ltr 5 Ltr 25 Ltr	
597	GLYCEROL (Glycerine), AR GRADE CAS 56-81-5 $C_3H_8O_3$ F. Wt. 92.09	500 ml 2.5 Ltr 5 Ltr 25 Ltr	
851	GLYCEROL TRIBUTYRATE (Tributyrin) CAS 60-01-5 $C_{15}H_{26}O_3$ F. Wt. 302.36	100 ml 500 ml	
1110	GLYCINE, EXTRA PURE (Amino Acetic Acid) CAS 56-40-6 $C_2H_5NO_2$ F. Wt. 75.07	100 gm 500 gm 5 Kg 25 Kg	
1109	GLYCINE 99% + CRYSTALLINE, AR GRADE CAS 56-40-6 $C_2H_5NO_2$ F. Wt. 75.07	100 gm 500 gm 5 Kg 25 Kg	
2560	CHLOROauric ACID, Min. 49% Au (Gold Chloride)	1 gm	
4701	GOLD CHLORIDE, ACS GRADE CAS 16903-35-8 $HAuCl_4 \cdot 3H_2O$ F. Wt. 339.79	1 gm	
1672 A	GRAM'S IODINE CAS 12298-68-9 I_2 F. Wt. 419.81	25 gm	
227	GUANIDINE HYDROCHLORIDE, AR GRADE (2-Amino-6-hydroxypurine HCl) CAS 50-01-1 $CH_5N_3 \cdot HCl$ F. Wt. 95.33	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
TK 004	GRAM'S STAIN KIT I Crystal Violet Gram 15-125 ml II Decolourizer Gram 15-125 ml III Iodine Gram 15-125 ml IV Safranin 0.5% w/v OR Basic Fuchsin 0.1%-125 ml	1 Kit (125 ml) 1 Kit (500 ml)	
697	GUANINE (2-Amino 6-Hydroxypurine) CAS 73-40-5 $C_5H_5N_5O$ F. Wt. 151.13	5 gm 25 gm	
600	GUANINE HCL * (2-Amino 6-Hydroxy Purine HCl) CAS 635-39-2 $C_5H_5N_5O \cdot HCl$ F. Wt. 187.59	5 gm 25 gm	
3346	GUANOSINE CAS 118-00-3 $C_{10}H_{13}N_5O_5$ F. Wt. 283.241	10 gm 25 gm	
210	GUM ACACIA, PHARMA GRADE Spray Dried Powder CAS 9000-01-5	500 gm 25 Kg	
1890	GUM XANTHAN (Xanthan Gum), AR GRADE CAS 11138-66-2	100 gm 500 gm 25 Kg	
1673	HAEMATOTOXYLIN CAS 517-28-2 C.I. 75290 $C_{16}H_{14}O_6$ F. Wt. 302.28	5 gm 25 gm 1 Kg	
603	HEMIN (From Bovine) CAS 16009-13-5 $C_{34}H_{32}ClFeN_4O_4$ F. Wt. 651.95	1 gm	
341	HEPARIN SODIUM SALT 20,000 I.U./Vial CAS 9041-08-1	1 vl	
1891	HEPARIN SODIUM SALT * 1,00,000 I.U./Vial CAS 9041-08-1	1 vl	
503	HEPES, FREE ACID CAS 7365-45-9 $C_8H_{18}N_2O_4S$ F. Wt. 238.30	25 gm 100 gm 500 gm 1 Kg	
599	HEPES SODIUM (for Tissue Culture) CAS 75277-39-3 $C_8H_{17}O_4N_2Na$ F. Wt. 260.29	25 gm 100 gm 1 Kg	
1894	1-HEPTANE SULPHONIC ACID SODIUM (ANH.), (For HPLC), AR GRADE CAS 22767-50-6 $C_7H_{13}NaO_3S$ F. Wt. 202.24	25 gm 100 gm	
4073	1-HEPTANE SULPHONIC ACID SODIUM SALT (MONO), (for HPLC) CAS 207300-90-1 $C_7H_{13}NaO_3S \cdot H_2O$ F. Wt. 220.24	25 gm 100 gm	
3351	HEXAMINE CAS 100-97-0 $C_6H_{12}N_4$ F. Wt. 140.19	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
4074	1-HEXANE SULPHONIC ACID SODIUM SALT (ANH), (for HPLC) CAS 2832-45-3 $C_6H_{13}NaO_3S$ F. Wt. 188.22	25 gm 100 gm	
4075	1-HEXANE SULPHONIC ACID SODIUM SALT (Monohydrate), (for HPLC) CAS 207300-91-2 $C_6H_{13}NaO_3S$ F. Wt. 188.22	25 gm 100 gm	
792	HIPPURIC ACID, 99% + Crystalline CAS 495-69-2 $C_9H_9NO_3$ F. Wt. 179.17	100 gm 500 gm	
1115	L-HISTIDINE, 99% + Crystalline CAS 71-00-1 $C_6H_9N_3O_2$ F. Wt. 155.15	25 gm 500 gm	
1116	L-HISTIDINE HCl (Mono), 99% + Crystalline CAS 5934-29-2 $C_6H_9N_3O_2 \cdot HCl \cdot H_2O$ F. Wt. 209.63	25 gm 100 gm 1 Kg 25 Kg	
4213	2, 8 HOMOBRASSINOLIDE CAS 80483-89-2 $C_{29}H_{50}O_6$ F. Wt. 494.70	10 mg 100 mg	
1899	HYDRAZINE SULPHATE, EXTRA PURE CAS 10034-93-2 $H_4N_2 \cdot H_2SO_4$ F. Wt. 130.12	100 gm 500 gm 25 Kg	
2631	HYDRAZINE SULPHATE, AR GRADE CAS 10034-93-2 $H_4N_2 \cdot H_2SO_4$ F. Wt. 130.12	100 gm 500 gm 25 Kg	
1907	HYDROXYLAMINE HYDROCHLORIDE, EXTRA PURE CAS 5470-11-1 $NH_2OH \cdot HCl$ F. Wt. 69.49	100 gm 500 gm 25 Kg	
1908	HYDROXYLAMINE HYDROCHLORIDE, AR GRADE CAS 5470-11-1 ClH_4NO F.Wt. 69.49	100 gm 500 gm	
4077	HYDROXYLAMINE SULPHATE, EXTRA PURE (Hydroxylammonium Sulphate) CAS 10039-54-0 $H_8N_2O_6S$ F. Wt. 164.14	500 gm	
4078	HYDROXYLAMINE SULPHATE, AR GRADE (Hydroxylammonium Sulphate) CAS 10039-54-0 $H_8N_2O_6S$ F. Wt. 164.14	500 gm	
1911	HYDROXY NAPHTHOL BLUE (SODIUM), AR GRADE CAS 63451-35-4 $C_{20}H_{12}N_2Na_2O_{11}S_3$ F. Wt. 598.50	5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
1117	L-HYDROXY-L-PROLINE, 99% + Crystalline CAS 51-35-4 $C_5H_9NO_3$ F. Wt. 131.13	5 gm 25 gm	
1912	HYDROXY PROPYL METHY CELLULOSE (HPMC) E 15 LV CAS 9004-65-3 $C_{32}H_{60}O_{19}$ F. Wt. 748.81	500 gm	
2632	8-HYDROXYQUINOLINE (Oxine), AR GRADE CAS 148-24-3 C_9H_7NO F. Wt. 145.16	100 gm 500 gm	
1015	HYDROQUINONE, EXTRA PURE (Quinol) CAS 123-31-9 $C_6H_6O_2$ F. Wt. 110.11	500 gm	
1904	HYDROQUINONE, AR GRADE (Quinol) CAS 123-31-9 $C_6H_6O_2$ F. Wt. 110.11	500 gm	
3151	HYPOXANTHINE CAS 68-94-0 $C_5H_4N_4O$ F. Wt. 136.112	5 gm 100 gm	
342	IMIDAZOLE, AR GRADE (Buffer Substance for pH 6.0-7.8) CAS 288-32-4 $C_3H_4N_2$ F. Wt. 68.08	100 gm 500 gm	
345	IMMERSION OIL Optically clear non fluorescent for microscopy CAS 8000-27-9	25 ml 100 ml 500 ml	
346	IMMERSION OIL, AR GRADE Optically clear non fluorescent for microscopy CAS 8000-27-9	25 ml 100 ml 500 ml	
1915	INDIGO CARMINE, AR GRADE CAS 860-22-0 C.I. 73015 $C_{16}H_8N_2Na_2O_8S_2$ F. Wt. 466.36	25 gm 100 gm	
1916	INDOLE, AR GRADE (1-Benzo-Pyrrole) CAS 120-72-9 C_8H_7N F. Wt. 117.15	10 gm 100 gm	
606	INDOLE-3-ACETIC ACID (IAA) * (3-Indole acetic acid) CAS 87-51-4 $C_{10}H_9NO_2$ F. Wt. 175.18	5 gm 25 gm 100 gm 1 Kg	
606A	INDOLE-3-ACETIC ACID POTASSIUM SALT CAS 2338-19-4 $C_{10}H_8KNO_2$ F. Wt. 213.27	25 gm 100 gm 1 Kg	
607	INDOLE-3-BUTYRIC ACID (IBA) * (3-Indole butyric acid) CAS 133-32-4 $C_{12}H_{13}NO_2$ F. Wt. 203.33	5 gm 25 gm 100 gm 1 Kg	

CODE	PRODUCT NAME	PACK SIZE	
607A	INDOLE-3-BUTYRIC ACID POTASSIUM SALT CAS 60096-23-3 C ₁₂ H ₁₂ KNO ₂ F. Wt. 241.33	25 gm 100 gm 1 Kg	
604	INOSITOL (Meso Inositol) CAS 87-89-8 C ₆ H ₁₂ O ₆ F. Wt. 180.16	25 gm 100 gm 1 Kg	
TC 014	INULIN, AR GRADE CAS 9005-80-5 (C ₆ H ₁₀ O ₅) _n F. Wt. 5,000 (Approx)	25 gm 100 gm	
347	IODINE (Resublimed), EXTRA PURE CAS 7553-56-2 I ₂ F. Wt. 253.81	100 gm 500 gm 25 Kg	
348	IODINE (Resublimed), AR GRADE CAS 7553-56-2 I ₂ F. Wt. 253.81	100 gm 500 gm	
4079	2-IODOBENZOIC ACID, Min. 98% CAS 88-67-5 C ₇ H ₅ IO ₂ F. Wt. 248.02	100 gm	
1927	IRON METAL POWDER (300 Mesh) (ELECTROLYTIC GRADE) CAS 7439-89-6 Fe F. Wt. 55.85	1 Kg	
3154	ISATIN (Indole-2,3-dione), AR GRADE CAS 91-56-5 C ₈ H ₅ NO ₂ F. Wt. 147.13	25 gm 100 gm	
1119	L-ISOLEUCINE, 99% + Crystalline CAS 73-32-5 C ₆ H ₁₃ NO ₂ F. Wt. 131.17	5 gm 25 gm	
4215	N ⁶ -(2-ISOPENTENYL)ADENINE; (2IP) [6-(g-g-Dimethylallylamino)purine] CAS 2365-40-4 C ₁₀ H ₁₃ N ₅ F. Wt. 203.24	1 gm 5 gm	
3156	ISOPROPYL β-D-THIOGALACTOPYRANOSIDE * (IPTG Dioxan Free) CAS 367-93-1 C ₉ H ₁₈ O ₅ S F. Wt. 238.30	1 gm 5 gm 25 gm	
1933	JACK BEAN MEAL (Uric Acid activator) CAS 9002-13-5	100 gm 500 gm	
1934	JANUS GREEN-B CAS 2869-83-2 C ₃₀ H ₃₁ CIN ₆ F. Wt. 511.06	5 gm 25 gm	
1935	JENNER'S STAIN (Eosin Methylene Blue) CAS 62851-42-7	25 gm 100 gm	
4158	JASMONIC ACID CAS 77026-92-7 C ₁₂ H ₁₈ O ₃ F. Wt. 210.27	100 mg 250 mg	
1937	KAOLIN LIGHT, EXTRA PURE (as per IP) CAS 1332-58-7 H ₂ Al ₂ Si ₂ O ₈ ·H ₂ O F. Wt. 258.16	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
1120	α-KETOGLUTARIC ACID * (2-Oxo-Glutaric Acid) 99% + Purity, Crystalline CAS 328-50-7 C ₅ H ₆ O ₅ F. Wt. 146.10	25 gm 100 gm	
863	KINETIN (6-Furfurylaminopurine) AR GRADE * CAS 525-79-1 C ₁₀ H ₉ N ₅ O F. Wt. 215.21	1 gm 5 gm	
1938	KOJIC ACID CAS 501-30-4 C ₆ H ₆ O ₄ F. Wt. 142.11	5 gm 25 gm	
TR 008	KOVAC'S INDOLE REAGENT Indole Test for Microbiology	100 ml	
1940	LACTIC ACID (Hydroxy Propionic Acid), EXTRA PURE CAS 50-21-5 C ₃ H ₆ O ₃ F. Wt. 90.08	500 ml 2.5 Ltr 25 Ltr	
1305	LACTOSE (Mono) CAS 10039-26-6 C ₁₂ H ₂₂ O ₁₁ ·H ₂ O F. Wt. 360.31	500 gm 25 Kg	
1316	LACTOSE (Mono) 99% + Crystalline Bacteriological Grade CAS 10039-26-6 C ₁₂ H ₂₂ O ₁₁ ·H ₂ O F. Wt. 360.31	500 gm 25 Kg	
3158	LAURIC ACID CAS 143-07-7 C ₁₂ H ₂₄ O ₂ F. Wt. 200.32	500 gm	
444	LAURYL SULPHATE SODIUM, EXTRA PURE (Sodium Dodecyl Sulphate) CAS 151-21-3 C ₁₂ H ₂₅ NaO ₄ S F. Wt. 288.38	500 gm 25 Kg	
445	LAURYL SULPHATE SODIUM, AR GRADE CAS 151-21-3 C ₁₂ H ₂₅ NaO ₄ S F. Wt. 288.38	100 gm 500 gm	
1943	LEAD METAL POWDER CAS 7439-92-1 Pb F. Wt. 207.19	500 gm	
351	LEAD ACETATE (Trihydrate), EXTRA PURE CAS 6080-56-4 C ₆ H ₆ O ₄ Pb ₃ (H ₂ O) F. Wt. 379.33	500 gm	
352	LEAD ACETATE BASIC (Lead Subacetate Basic) for Sugar Analysis by Horne CAS 51404-69-4 C ₂ H ₃ O ₃ Pb F. Wt. 283.25	500 gm 2.5 Kg	
1944	LEAD CARBONATE BASIC, Min 98% CAS 1344-36-1 C ₂ H ₃ O ₃ Pb F. Wt. 329.23	500 gm	
2639	LEAD (II) CHLORIDE (ANH.) (for Synthesis) CAS 7758-95-4 PbCl ₂ F. Wt. 278.11	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
1946	LEAD DIOXIDE, EXTRA PURE (Lead Peroxide) CAS 1309-60-0 O ₂ Pb F. Wt. 239.19	500 gm	
1947	LEAD MONOXIDE (Litharge, Lead Oxide Yellow) CAS 1317-36-8 PbO F. Wt. 223.20	500 gm	
864	LEAD NITRATE, EXTRA PURE CAS 10099-74-8 N ₂ O ₆ Pb F. Wt. 331.2	500 gm	
1674	LEISHMAN'S STAIN C.I. 881038 (Eosin Methylene Blue) CAS 12627-53-1	25 gm 100 gm 1 Kg	
1121	L-LEUCINE, 99% + Crystalline CAS 61-90-5 C ₆ H ₁₃ NO ₂ F. Wt. 131.17	25 gm 100 gm 500 gm	
1950	LIGHT GREEN CAS 7114-03-06 C.I. 42040 C ₂₇ H ₃₄ N ₂ O ₄ S F. Wt. 482.63	25 gm 100 gm	
1686	LIGHT GREEN-SF, 80% CAS 5141-20-8 C.I. 42095 C ₃₇ H ₃₄ N ₂ Na ₂ O ₅ S ₃ F. Wt. 792.85	25 gm 100 gm	
1686A	LIGHT GREEN-SF, 65% CAS 5141-20-8 C.I. 42095 C ₃₇ H ₃₄ N ₂ Na ₂ O ₅ S ₃ F. Wt. 792.85	25 gm 100 gm 500 gm	
2640	LIPASE * CAS 9001-62-1	25 gm 100 gm	
4085	LITHIUM BROMIDE (ANH.), EXTRA PURE CAS 7550-35-8 LiBr F. Wt. 86.85	500 gm	
668	LITHIUM CARBONATE, EXTRA PURE CAS 554-13-2 Li ₂ CO ₃ F. Wt. 73.89	100 gm 500 gm 25 Kg	
4140	LITHIUM CARBONATE, AR GRADE CAS 554-13-2 Li ₂ CO ₃ F. Wt. 73.89	100 gm 500 gm	
3322	LITHIUM CHLORIDE (ANH.), EXTRA PURE CAS 7447-41-8 LiCl F. Wt. 42.39	500 gm	
1955	LITHIUM HYDROXIDE (Mono), EXTRA PURE CAS 1310-66-3 LiOH.H ₂ O F. Wt. 41.96	500 gm 25 Kg	
4086	LITHIUM NITRATE (ANH.), EXTRA PURE CAS 7790-69-4 LiNO ₃ F. Wt. 68.95	500 gm	
613	LITHIUM SULPHATE (Mono), EXTRA PURE CAS 10102-25-7 LiSO ₄ .H ₂ O F. Wt. 121.01	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
4087	di-LITHIUM TETRABORATE, AR GRADE CAS 12007-60-2 B ₄ Li ₂ O ₇ F. Wt. 169.12	100 gm 500 gm	
1634	LITMUS BLUE, Indicator paper (10 books of 20 leaves in Box)	200 lvs	
1636	LITMUS RED, Indicator paper (10 books of 20 leaves in box)	200 lvs	
1688	LITMUS INDICATOR, EXTRA PURE CAS 1393-92-6	10 gm 25 gm 100 gm	
4377	LIVER EXTRACT POWDER, BACTO GRADE  CAS 1393-92-6	500 gm	
1122	L-LYSINE, MONO Hcl, 99% + Crystalline CAS 657-27-2 C ₆ H ₁₄ N ₂ O ₂ .HCl F. Wt. 182.65	100 gm 500 gm 25 Kg	
4233	LYSOL (Cresol and Soap Solution) CAS 12772-68-8	500 ml 5 Ltr	
2951	LYSOZYME * CAS 12650-88-3	1 gm 5 gm	
2950	LYSOZYME CHLORIDE * CAS 52219-07-5	1 gm 5 gm	
4088	MAGNESIUM (Metal Turning) CAS 7439-95-4 Mg F. Wt. 24.31	500 gm	
1964	MAGNESIUM ACETATE (Tetra), EXTRA PURE CAS 16674-78-5 C ₄ H ₁₄ MgO ₈ F. Wt. 214.45	500 gm	
355	MAGNESIUM CARBONATE BASIC (Light), EXTRA PURE CAS 546-93-0 MgCO ₃ F. Wt. 84.31	500 gm	
2393	MAGNESIUM CARBONATE (Heavy), EXTRA PURE CAS 546-93-0 MgCO ₃ F. Wt. 84.313	500 gm	
2394	MAGNESIUM CHLORIDE (ANH.), EXTRA PURE CAS 7786-30-3 MgCl ₂ F. Wt. 95.21	500 gm	
356	MAGNESIUM CHLORIDE (Hexa), EXTRA PURE CAS 7791-18-6 MgCl ₂ .6H ₂ O F. Wt. 203.30	500 gm 5 Kg 50 Kg	
614	MAGNESIUM CHLORIDE (Hexa), AR GRADE CAS 7791-18-6 MgCl ₂ .6H ₂ O F. Wt. 203.30	500 gm 5 Kg 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
758	MAGNESIUM CITRATE (Nonahydrate) CAS 153531-96-5 $Mg_3(C_6H_5O_7)_2 \cdot 9H_2O$ F. Wt. 613.25	500 gm	
357	MAGNESIUM NITRATE (Hexa), EXTRA PURE CAS 13446-18-9 $Mg(NO_3)_2 \cdot 6H_2O$ F. Wt. 256.41	500 gm	
1968	MAGNESIUM OXIDE (Light), EXTRA PURE CAS 1309-48-4 MgO F. Wt. 40.30	500 gm	
2953	MAGNESIUM PHOSPHATE DIBASIC, EXTRA PURE CAS 7782-75-4 $MgHPO_4 \cdot 3H_2O$ F. Wt. 174.33	500 gm	
358	MAGNESIUM SULPHATE (Hepta), EXTRA PURE CAS 10034-99-8 $MgSO_4 \cdot 7H_2O$ F. Wt. 246.48	500 gm 5 Kg 50 Kg	
359	MAGNESIUM SULPHATE (Hepta), AR GRADE CAS 10034-99-8 $MgSO_4 \cdot 7H_2O$ F. Wt. 246.48	500 gm 5 Kg 50 Kg	
871	MAGNESIUM SULPHATE (ANH.), EXTRA PURE CAS 7487-88-9 $MgSO_4$ F. Wt. 120.37	500 gm 50 Kg	
1973	MAGNESIUM TRISILICATE (Hydrate) CAS 39365-87-2 $2MgO \cdot 3SiO_2 \cdot xH_2O$ F. Wt. 260.86+ xH_2O	500 gm	
1974	MAGNESON-I, AR GRADE (Reagent for Mg) CAS 74-39-5 $C_{12}H_{19}N_3O_3$ F. Wt. 259.22	25 gm 100 gm	
1975	MAGNESON-II, AR GRADE (Reagent for Mg) CAS 5290-62-0 $C_{16}H_{11}N_3O_3$ F. Wt. 293.28	25 gm 100 gm	
4360	METHYLSULFONYLMETHANE CAS 67-71-0 $(CH_3)_2SO_2$ F. Wt. 94.13	25 gm	
1641	MALACHITE GREEN (Malachite Green Oxalate) CAS 2437-29-8 C.I. 42000 $C_{25}H_{26}N_2O_4$ F. Wt. 927.03	25 gm 100 gm 500 gm	
360	MALEIC ACID, EXTRA PURE CAS 110-16-7 $C_4H_2O_4$ F. Wt. 116.07	500 gm 25 Kg	
362	MALEIC ANHYDRIDE [2,5-Furandione] CAS 108-31-6 $C_4H_2O_4$ F. Wt. 98.06	500 gm 5 Kg 25 Kg	
3556	MALEIC HYDRAZIDE CAS 123-33-1 $C_4H_4N_2O_2$ F. Wt. 112.09	100 gm	

CODE	PRODUCT NAME	PACK SIZE	
363	DL-MALIC ACID (DL-Hydroxysuccinic Acid) CAS 6915-15-7 $C_4H_6O_5$ F. Wt. 134.09	500 gm 5 Kg	
2955	MALONIC ACID, EXTRA PURE CAS 141-82-2 $C_3H_4O_4$ F. Wt. 104.06	100 gm 500 gm	
446	MALONIC ACID SODIUM CAS 141-95-7 $C_3H_2Na_2O_4$ F. Wt. 148.03	25 gm 100 gm	
3527	MALT EXTRACT PASTE (Food Grade) CAS 8002-48-0	500 gm 50 Kg	
3506	MALT EXTRACT POWDER (Food Grade) CAS 8002-48-0	500 gm 5 Kg 25 Kg	
4142	MALT EXTRACT POWDER Bacteriological Grade CAS 8002-48-0	500 gm 25 Kg	
1321	MALTODEXTRINE POWDER, (D/E 20) (Maltodextrin) CAS 9050-36-6 $C_4H_{12}O_6$ F. Wt. 180.16	500 gm 5 Kg 25 Kg	
1244	MALTOSE (Mono), EXTRA PURE (for Bacteriology) CAS 6363-53-7 $C_{12}H_{22}O_{11} \cdot H_2O$ F. Wt. 360.31	100 gm 500 gm 25 Kg	
766	MALTOSE (Mono), AR GRADE (for Bacteriology) CAS 6363-53-7 $C_{12}H_{22}O_{11} \cdot H_2O$ F. Wt. 360.31	100 gm 500 gm	
1979	DL-MANDELIC ACID CAS 90-64-2 $C_8H_8O_3$ F. Wt. 152.15	100 gm 500 gm	
2396	MANGANESE CARBONATE, EXTRA PURE CAS 598-62-9 $MnCO_3$ F. Wt. 114.95	500 gm 25 Kg	
704	MANGANESE (II) CHLORIDE (Tetra), EXTRA PURE CAS 13446-34-9 $MnCl_2 \cdot 4H_2O$ F. Wt. 197.91	500 gm 50 Kg	
705	MANGANESE (II) CHLORIDE (Tetra), AR GRADE CAS 13446-34-9 $MnCl_2 \cdot 4H_2O$ F. Wt. 197.91	500 gm 50 Kg	
1980	MANGANESE DIOXIDE (Black), TECH. CAS 1313-13-9 MnO_2 F. Wt. 86.94	500 gm	
702	MANGANESE (II) SULPHATE (Mono), EXTRA PURE CAS 10034-96-5 $MnSO_4 \cdot H_2O$ F. Wt. 169.02	500 gm 5 Kg 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
703	MANGANESE (II) SULPHATE (Mono), AR GRADE CAS 10034-96-5 MnSO ₄ ·H ₂ O F. Wt. 169.02	500 gm 50 Kg	
1245	D-MANNITOL, EXTRA PURE CAS 69-65-8 C ₆ H ₁₄ O ₆ F. Wt. 182.17	100 gm 500 gm 5 Kg 25 Kg	
1246	D-MANNITOL, AR GRADE CAS 69-65-8 C ₆ H ₁₄ O ₆ F. Wt. 182.17	100 gm 500 gm 25 Kg	
365	D + MANNOSE CAS 3458-28-4 C ₆ H ₁₂ O ₆ F. Wt. 180.16	25 gm 100 gm	
1982	D-MELIBIOSE 99% (Mono) CAS 585-99-9 C ₁₂ H ₂₂ O ₁₁ F. Wt. 342.30	5 gm	
1716	MENADIONE SODIUM BISULPHITE* (Vitamin K3 Sodium Sulphate) Water Soluble CAS 130-37-0 C ₁₁ H ₉ NaO ₅ S F. Wt. 276.24	25 gm 100 gm	
3513	MERCAPTOACETIC ACID SODIUM SALT (Sodium Thioglycollate) CAS 367-51-1 C ₂ H ₃ NaO ₂ S F. Wt. 114.10	100 gm 500 gm	
1988	2-MERCAPTOETHANOL, AR GRADE CAS 60-24-2 C ₂ H ₆ OS F. Wt. 78.13	100 ml 500 ml	
1989	MERCURIC ACETATE, EXTRA PURE CAS 1600-27-7 C ₆ H ₆ HgO ₄ F. Wt. 318.70	100 gm 500 gm	
1990	MERCURIC ACETATE, AR GRADE CAS 1600-27-7 C ₆ H ₆ HgO ₄ F. Wt. 318.70	100 gm	
369	MERCURIC CHLORIDE, EXTRA PURE (Mercury (II) Chloride) CAS 7487-94-7 HgCl ₂ F. Wt. 271.52	100 gm 250 gm 500 gm	
370	MERCURIC CHLORIDE, AR GRADE (Mercury (II) Chloride) CAS 7487-94-7 HgCl ₂ F. Wt. 271.52	100 gm 250 gm	
789	MERCURIC IODIDE RED, EXTRA PURE CAS 7774-29-0 HgI ₂ F. Wt. 454.40	100 gm	

CODE	PRODUCT NAME	PACK SIZE	
880	MERCURIC NITRATE (Mono), EXTRA PURE CAS 7783-34-8 Hg(NO ₃) ₂ ·H ₂ O F. Wt. 342.62	100 gm 250 gm	
616	MERCURIC OXIDE RED, EXTRA PURE CAS 21908-53-2 HgO F. Wt. 216.59	100 gm 500 gm	
618	MERCURIC OXIDE YELLOW, EXTRA PURE CAS 21908-53-2 HgO F. Wt. 216.59	100 gm	
882	MERCURIC SULPHATE, EXTRA PURE (Mercury (II) Sulphate) CAS 7783-35-9 HgSO ₄ F. Wt. 296.65	100 gm 250 gm 500 gm	
883	MERCURIC SULPHATE, AR GRADE CAS 7783-35-9 HgO ₄ S F. Wt. 296.65	250 gm	
4090	MERCURIC THIOCYANATE, EXTRA PURE CAS 592-85-8 C ₂ HgN ₂ S ₂ F. Wt. 316.75	100 gm	
4091	MERCUROUS CHLORIDE, AR GRADE CAS 10112-91-1 Cl ₂ Hg ₂ F. Wt. 472.08	100 gm	
1996	MERCUROUS NITRATE (Dihydrate), EXTRA PURE CAS 14836-60-3 Hg ₂ (NO ₃) ₂ ·2H ₂ O F. Wt. 561.22	100 gm	
2414	MERCURY (Metal), EXTRA PURE CAS 7439-97-6 Hg F. Wt. 200.59	100 gm	
4214	META-TOPOLINE CAS 75737-38-1 C ₁₂ H ₁₀ N ₅ OH F. Wt. 241.25	25 mg 100 mg	
1124	DL-METHIONINE 99% + Purity, Crystalline CAS 59-51-8 C ₆ H ₁₁ NO ₂ S F. Wt. 149.21	100 gm 500 gm 1 Kg 25 Kg	
1123	L-METHIONINE 99% + Purity, Crystalline CAS 63-68-3 C ₆ H ₁₁ NO ₂ S F. Wt. 149.21	25 gm 100 gm 500 gm 25 Kg	
4260	4-METHYLUMBELLIFERYL-β-D-GALACTOPYRANOSIDE (MUGA:MU-GAL) CAS 6160-78-7 C ₁₆ H ₁₈ O ₈ F. Wt. 338.31	250 mg	
1687	METHYL BLUE (Water Blue) CAS 28983-56-4 C.I. 42780 C ₃₇ H ₂₇ N ₃ Na ₂ O ₉ S ₃ F. Wt. 799.82	25 gm 100 gm	
2433	METHYL CELLULOSE (High Viscosity) 4000 CPS CAS 9004-67-5	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
4257	METHYL- α -D-GALACTOPYRANOSIDE (METHYL- α -D-GALACTOSIDE) * CAS 3396-99-4 $C_7H_{14}O_6$ F. Wt. 194.18	1 gm 5 gm	
4258	METHYL- β -D-GALACTOPYRANOSIDE (METHYL- β -D-GALACTOSIDE) * CAS 1824-94-8 $C_7H_{14}O_6$ F. Wt. 194.18	100 gm	
4259	METHYL- α -D-GLUCOPYRANOSIDE (METHYL- α -D-GLUCOSIDE) * CAS 1824-94-8 $C_7H_{14}O_6$ F. Wt. 194.18	100 gm	
4225	METHYL EUGENOL CAS 93-15-2 $C_{11}H_{14}O_2$ F. Wt. 178.23	500 ml	
2103	METHYL GREEN, AR GRADE CAS 7114-03-6 C.I. 42590 $C_{27}H_{35}Cl_4N_3Zn$ F. Wt. 608.79	10 gm 25 gm	
1604	METHYL ORANGE (ORANGE-III) CAS 547-58-0 C.I. 13025 $C_{14}H_{14}N_3NaO_3S$ F. Wt. 327.33	25 gm 100 gm 1 Kg	
884	METHYL PARABEN (Methyl 4-Hydroxy Benzoate) CAS 99-76-3 $C_8H_8O_3$ F. Wt. 152.15	500 gm	
619	METHYL PARABEN SODIUM (Methyl-P-Hydroxy Benzoate Sodium) CAS 5026-62-0 $C_8H_7NaO_3$ F. Wt. 174.13	500 gm 25 Kg	
1605	METHYL RED CAS 493-52-7 C.I. 13020 $C_{15}H_{15}N_2O_2$ F. Wt. 269.30	25 gm 100 gm	
1633	METHYL THYMOL BLUE SODIUM (Methyl Red Thymol Blue) CAS 1945-77-3 $C_{37}H_{43}N_2O_{13}NaS$ F. Wt. 778.82	5 gm 25 gm	
1606	METHYL VIOLET, EXTRA PURE (Crystal Violet or Gentian Violet) CAS 548-62-9 C.I. 42555 $C_{25}H_{30}ClN_3$ F. Wt. 407.99	25 gm 100 gm 1 Kg	
1620	METHYL VIOLET (Crystal Violet or Gentian Violet 10B) CAS 548-62-9 C.I. 42555 $C_{25}H_{30}ClN_3$ F. Wt. 407.98	25 gm 100 gm 1 Kg	
371	N,N-METHYLENE BISACRYLAMIDE (for Electrophoresis), AR GRADE CAS 110-26-9 $C_7H_{10}N_2O_2$ F. Wt. 154.17	25 gm 100 gm 500 gm	
1607	METHYLENE BLUE, EXTRA PURE CAS 61-73-4 $C_{16}H_{18}N_3SCl$ F. Wt. 319.85	25 gm 100 gm 1 Kg	
2651	METHYLENE BLUE, AR GRADE CAS 61-73-4 $C_{16}H_{18}N_3SCl$ F. Wt. 319.85	25 gm 100 gm 1 Kg	

CODE	PRODUCT NAME	PACK SIZE	
4092	METHYL GREEN (for Microscopical Staining) (C.I.42590) CAS 22383-16-0 $C_{26}H_{33}Cl_2N_3$ F. Wt. 458.47	10 gm 25 gm	
4093	MOLYBDENUM DISULPHIDE CAS 1317-33-5 MoS_2 F. Wt. 160.09	100 gm	
4094	MOLYBDENUM TRIOXIDE CAS 1313-27-5 MoO_3 F. Wt. 143.95	100 gm 500 gm	
948	MOLYBDIC ACID (Mono), EXTRA PURE CAS 7782-91-4 $MoO_3 \cdot H_2O$ F. Wt. 1161.95	100 gm 25 Kg	
4095	MOLYBDIC ACID (Hydrate), AR GRADE CAS 7782-91-4 H_2MoO_4 F. Wt. 161.95	100 gm 500 gm	
1112	MONO SODIUM GLUTAMATE CAS 142-47-2 $C_5H_9NNaO_4$ F. Wt. 169.11	100 gm 500 gm 25 Kg	
1542	MORIN DIHYDRATE, AR GRADE CAS 6472-38-4 C.I. 75660 $C_{15}H_{10}O_7 \cdot 2H_2O$ F. Wt. 338.27	1 gm 5 gm	
2110	2-MORPHOLINO ETHANE SULPHONIC ACID (MES) BUFFER CAS 4432-31-9 $C_6H_{13}NO_3S$ F. Wt. 195.2	25 gm 100 gm 1 Kg	
3596	3-MORPHOLINOPROPANE SULPHONIC ACID (MOPS-BUFFER) CAS 1132-61-2 $C_7H_{13}NO_3S$ F. Wt. 209.26	25 gm 1 Kg	
3564	MTT (Methyl Thiazolyl Diphenyl Tetrazolium Bromide) * CAS 298-93-1 $C_{18}H_{16}N_5SBr$ F. Wt. 414.32	1 gm	
060	MUREXIDE (Ammonium Purpurate) CAS 3051-09-0 C.I. 56085 $C_8H_8N_6O_6$ F. Wt. 284.19	5 gm 25 gm	
1602	MUREXIDE (Ammonium Purpurate), AR GRADE CAS 3051-09-0 C.I. 56085 $C_8H_8N_6O_6$ F. Wt. 284.19	5 gm 25 gm	
2379	NAPHTHALENE CAS 91-20-3 $C_{10}H_8$ F. Wt. 128.17	500 gm	
2122	α -NAPHTHALENE ACETIC ACID CAS 86-87-3 $C_{12}H_{10}O_2$ F. Wt. 128.17	25 gm 100 gm 500 gm 25 Kg	
2122A	α -NAPHTHALENE ACETIC ACID POTASSIUM SALT CAS 15165-79-4 $C_{12}H_9O_2K$ F. Wt. 224.30	25 gm 100 gm 1 Kg	

CODE	PRODUCT NAME	PACK SIZE	
893	1-NAPHTHOL (1-Hydroxynapthalene) CAS 90-15-3 $C_{10}H_8O$ F. Wt. 144.17	100 gm 500 gm	
2381	2-NAPHTHOL (β -Naphthol) CAS 135-19-3 $C_{10}H_8O$ F. Wt. 144.17	500 gm	
2465	1-NAPHTHOLBENZEIN (pH Indicator) CAS 145-50-6 $C_{27}H_{18}O_2$ F. Wt. 374.43	5 gm 25 gm	
1632	α -NAPHTHOLPHTHALEIN CAS 596-01-0 $C_{27}H_{18}O_2$ F. Wt. 418.44	1 gm 5 gm	
2656	1,2-NAPHTHOLQUINONE-4 SULPHONIC ACID SODIUM SALT CAS 521-24-4 $C_{10}H_5NaO_5S$ F. Wt. 260.19	5 gm 25 gm	
3565	β -NAPHTHOXY ACETIC ACID (BONA) CAS 120-23-0 $C_{12}H_{10}O_3$ F. Wt. 202.21	25 gm 100 gm 500 gm	
2121	1-NAPHTHYLACETATE, AR GRADE (Naphthol free) CAS 830-81-9 $C_{12}H_{10}O_2$ F. Wt. 186.21	5 gm 25 gm	
2122	1-NAPHTHYLACETIC ACID (α -Naphthalene Acetic Acid) CAS 86-87-3 $C_{12}H_{10}O_2$ F. Wt. 186.21	25 gm 100 gm 500 gm 25 Kg	
2124	1-NAPHTHYLAMINE, AR GRADE CAS 134-32-7 $C_{10}H_9N$ F. Wt. 143.19	100 gm 500 gm	
TR 028	α -NAPHTHYLAMINE SOLUTION CAS 134-32-7 $C_{10}H_9N$ F. Wt. 143.19	100 ml	
1913	N-(1-NAPHTHYL) ETHYLENE DIAMINE DIHYDROCHLORIDE (2-NEDA) CAS 1465-25-4 $C_{12}H_{14}N_2$ F. Wt. 259	5 gm 25 gm	
2126	NEOCUPROINE, AR GRADE (2,9, Dimethyl-1-10 Phenanthroline) CAS 484-11-7 $C_{14}H_{12}N_2$ F. Wt. 208.26	1 gm 5 gm	
2127	NEOCUPROINE HYDROCHLORIDE (Mono) (2,9, Dimethyl-1-10 Phenanthroline Hcl), AR GRADE CAS 7296-20-0 $C_{14}H_{12}N_2 \cdot HCl \cdot H_2O$ F. Wt. 262.72	1 gm 5 gm	
1643	NEUTRAL RED (pH Indicator) CAS 553-24-2 C.I.50040 $C_{15}H_{17}ClN_4$ F. Wt. 288.78	5 gm 25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
1643 A	NEUTRAL RED, AR GRADE CAS 553-24-2 $C_{15}H_{17}ClN_4$ F. Wt. 288.78	5 gm 25 gm 100 gm	
3178	NEW METHYLENE BLUE (Basic Blue 24) CAS 6586-05-0 C.I.52030 $C_{18}H_{22}ClN_3S \cdot ZnCl_2$ F. Wt. 416.05	5 gm	
759	NIACIN (Nicotinic Acid) CAS 59-67-6 $C_6H_5NO_2$ F. Wt. 123.11	25 gm 100 gm 1 Kg	
685	NIACINAMIDE (Nicotinamide) CAS 98-92-0 $C_6H_6N_2O$ F. Wt. 122.12	25 gm 100 gm 1 Kg	
2130	NICKEL CARBONATE BASIC PURIFIED CAS 12607-70-4 $CH_4Ni_3O_7$ F. Wt. 304.18	500 gm	
2131	NICKEL CHLORIDE (Hexa), EXTRA PURE CAS 7791-20-0 $NiCl_2 \cdot 6H_2O$ F. Wt. 237.69	500 gm	
3323	NICKEL (II) NITRATE (Hexa), EXTRA PURE CAS 13478-00-7 $N_2NiO_6 \cdot 6H_2O$ F. Wt. 290.79	500 gm	
2135	NICKEL OXIDE BLACK, EXTRA PURE CAS 1314-06-3 NiO F. Wt. 74.69	100 gm 500 gm	
3324	NICKEL SULPHATE (Hexa) CAS 10101-97-0 $NiSO_4 \cdot 6H_2O$ F. Wt. 262.85	500 gm	
3179	NICOTINAMIDE ADENINE DINUCLEOTIDE (DPN;NAD) * CAS 53-84-9 $C_{21}H_{27}N_7O_{14}P_2$ F. Wt. 663.43	1 gm 5 gm	
3180	NICOTINAMIDE ADENINE DINUCLEOTIDE REDUCED SODIUM (NADH Disodium) * CAS 606-68-8 $C_{21}H_{27}N_7Na_2O_{14}P_2$ F. Wt. 709.41	1 gm 5 gm	
3181	NICOTINAMIDE ADENINE DINUCLEOTIDE PHOSPHATE DISODIUM (NADP Disodium) * CAS 24292-60-2 $C_{21}H_{26}N_7Na_2O_{17}P_3$ F. Wt. 787.37	100 mg 1 gm	
1678	NIGROSINE B (A/S) CAS 8005-02-5 C.I. 50415	25 gm 100 gm 1 Kg	
1645	NIGROSINE (Water Soluble) CAS 8005-02-5 C.I. 50420	25 gm 100 gm	
2143	NILE BLUE CHLORIDE CAS 2381-85-3 C.I. 51180 $C_{20}H_{20}ClN_3O$ F. Wt. 353.85	25 gm	

CODE	PRODUCT NAME	PACK SIZE	
2144	NILE BLUE SULPHATE (Nile Blue A), AR GRADE CAS 3625-57-8 C.I. 51180 $C_{40}H_{40}N_6O_6S$ F. Wt. 732.84	25 gm	
2145	NINHYDRIN 99%, EXTRA PURE CAS 485-47-2 $C_9H_4O_3 \cdot H_2O$ F. Wt. 178.14	10 gm 25 gm 100 gm	
1608	NINHYDRIN, AR GRADE for detection and assay of Amino Acid CAS 485-47-2 $C_9H_4O_3 \cdot H_2O$ F. Wt. 178.14	10 gm 25 gm 100 gm	
2146	NIOBIUM PENTOXIDE 99.9%, AR GRADE CAS 1313-96-8 Nb_2O_5 F. Wt. 265.81	25 gm	
2892	2-NITROBENZALDEHYDE, AR GRADE CAS 552-89-6 $C_7H_5NO_3$ F. Wt. 151.11	10 gm 25 gm	
2702	4-(4-NITROBENZYL) PYRIDINE, AR GRADE CAS 1083-48-3 $C_{12}H_{10}N_2O_2$ F. Wt. 214.22	5 gm 25 gm	
3385	4-NITROBENZOIC ACID CAS 62-23-7 $C_7H_5NO_4$ F. Wt. 167.12	500 gm	
2153	NITRO B.T., AR GRADE * (Nitro Blue Tetrazolium Chloride) CAS 298-83-9 $C_{40}H_{30}Cl_2N_{10}O_6$ F. Wt. 817.64	100 mg 250 mg 1 gm	
2155	P-NITROPHENOL, EXTRA PURE CAS 100-02-7 $C_6H_5NO_3$ F. Wt. 139.12	500 gm 5 Kg 25 Kg	
3184	O-NITROPHENYL-β-D GALACOPYRANOSIDE (ONPG) * CAS 369-07-3 $C_{12}H_{15}NO_8$ F. Wt. 301.25	500 mg 1 gm 5 gm	
1974	4-(4-NITRO PHENYL AZO)-1-NAPHTHOL (Magneson-I), AR GRADE CAS 74-39-5 $C_{12}H_9N_3O_4$ F. Wt. 259.22	25 gm 100 gm	
1975	4-(4-NITROPHENYL AZO) RESORCINOL (Magneson-II), AR GRADE CAS 5290-62-0 $C_{16}H_{11}N_3O_3$ F. Wt. 293.28	25 gm 100 gm	
2158	4-NITROPHENYL PHOSPHATE DISODIUM (Hexa), AR GRADE CAS 333338-18-4 $C_6H_4NO_6PNa_2 \cdot 6H_2O$ F. Wt. 371.14	5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
2159	1-NITROSO-2-NAPHTHOL, AR GRADE CAS 131-91-9 $C_{10}H_7NO_2$ F. Wt. 173.17	25 gm 100 gm	
2161	NITROSO-R-SALT, AR GRADE CAS 525-05-3 $C_{10}H_5NNaO_5S_2$ F. Wt. 377.26	25 gm 100 gm	
2463	NYSTATIN (10,00,000 Units/vl) * CAS 1400-61-9 $C_{47}H_{75}NO_{17}$ F. Wt. 926.10	1 vl	
3386	1-OCTANE SULPHONIC ACID SODIUM SALT (MONO), (For HPLC) CAS 207596-29-0 $C_8H_{17}NaO_3S \cdot H_2O$ F. Wt. 234.29	25 gm	
3387	1-OCTANE SULPHONIC ACID SODIUM SALT (ANH), (For HPLC) CAS 5324-84-5 $C_8H_{17}NaO_3S$ F. Wt. 216.28	25 gm	
4098	OCTYL SULPHATE SODIUM SALT, AR GRADE, (for HPLC) CAS 142-31-4 $C_8H_{17}NaO_4S$ F. Wt. 232.3	5 gm	
1693	OIL RED-IV (Sudan-IV) (Seavlet-R) CAS 85-83-6 C.I. 26105 $C_{24}H_{20}N_4O$ F. Wt. 380.44	25 gm 100 gm	
2163	OLEIC ACID (Elainic Acid), EXTRA PURE * (for Biochemistry) CAS 112-80-1 $C_{18}H_{34}O_2$ F. Wt. 282.46	500 ml	
3520	OLIVE OIL CAS 8001-25-0	100 ml 500 ml	
1639	ORANGE G (Acid Orange 10) CAS 1936-15-8 C.I. 16230 $C_{16}H_{10}N_2O_7S_2Na_2$ F. Wt. 452.38	25 gm 100 gm 1 Kg	
2165	ORCEIN (for Microscopy) CAS 1400-62-0 $C_{28}H_{24}N_2O_7$ F. Wt. 500.49	5 gm 10 gm	
1125	L-ORNITHINE MONOHYDROCHLORIDE 99% + PURITY, CRYSTALLINE CAS 3184-13-2 $C_5H_{12}N_2O_2 \cdot HCl$ F. Wt. 168.62	25 gm 100 gm 500 gm	
2166	ORCINOL (Mono) (3,5-Dihydroxytoluence Mono) CAS 6153-39-5 $C_7H_8O_2 \cdot H_2O$ F. Wt. 142.15	10 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
372	OXALIC ACID (Dihydrate), EXTRA PURE CAS 6153-56-6 C ₂ H ₂ O ₄ ·2H ₂ O F. Wt. 126.07	500 gm 5 Kg 50 Kg	
TR 020	OXIDASE REAGENT (Gorden-Mcleod Reagent)	100 ml	
3570	PACLOBUTRAZOL CAS 76738-62-0 C ₁₅ H ₂₀ ClN ₃ O F. Wt. 293.80	5 gm 25 gm	
3186	PAN, INDICATOR 1-(2-Pyridylazo)-2-Naphthol CAS 85-85-8 C ₁₅ H ₁₁ N ₃ O F. Wt. 249	1 gm 5 gm	
4100	PAR INDICATOR, AR GRADE CAS 16593-81-0 C ₁₁ H ₈ N ₃ NaO ₂ ·H ₂ O F. Wt. 255.21	1 gm 5 gm	
1401	PANCREATIN 1NF CAS 8049-47-6	100 gm 500 gm	
1403	PANCREATIN 4 NF * CAS 8049-47-6	100 gm	
1404	PAPAIN (SPRAY DRIED) (Activity 100 T.U./gm or 6,000 USP/gm) CAS 9001-73-4	100 gm 500 gm	
609	PARAFFIN LIQUID (HEAVY) CAS 8012-95-1	500 ml 2.5 Ltr 25 Ltr	
622	PARAFFIN LIQUID (LIGHT) CAS 8012-95-1	500 ml 2.5 Ltr 25 Ltr	
375	PARAFFIN WAX (Block Form) Congealing point about 58-60°C (Packed in SS Box) CAS 8002-74-2	1 Kg	
375A	PARAFFIN WAX (Pellet Form) Congealing point about 58-60°C CAS 8002-74-2	1 Kg	
4173	PARAFFIN WAX (Solid) Congealing point about 58-60°C (Packed in SS Box) CAS 8002-74-2	1 Kg	
4174	PARAFFIN WAX (Block Form) Congealing point about 60-62°C (Packed in SS Box) CAS 8002-74-2	1 Kg	
4174A	PARAFFIN WAX (Pellet Form) Congealing point about 60-62°C CAS 8002-74-2	1 Kg	
2175	PARAFFIN WAX (Solid) Congealing point about 60-62°C (Packed in SS Box) CAS 8002-74-2	1 Kg 2 Kg	

CODE	PRODUCT NAME	PACK SIZE	
4175	PARAFFIN WAX (With Ceresin) (Block Form) Congealing point about 60°-62°C (Packed in SS Box) CAS 8002-74-2	1 Kg	
4350	PARAFFIN WAX (With Ceresin) (Block Form) Congealing point about 58-60°C CAS 8002-74-2	1 Kg	
4350A	PARAFFIN WAX (With Ceresin)  (Pallet Form) Congealing point about 58-60°C CAS 8002-74-2	1 Kg	
376	PARAFFIN WAX (With Ceresin) (Solid) Congealing point about 60°-62°C (Packed in SS Box) CAS 8002-74-2	500 gm 1 Kg 2 Kg	
623	PARAFORMALDEHYDE 96% (For Synthesis) CAS 30525-89-4 (CH ₂ O) ₄ F. Wt. 90.08	500 gm	
2178	PARAROSANILINE (BASE) (Para Fuchsin) CAS 25620-78-4 C.I. 42500 C ₁₉ H ₁₉ N ₃ O F. Wt. 305.37	25 gm	
2179	PARAROSANILINE HYDROCHLORIDE (Para Fuchsin HCl) CAS 569-61-9 C.I. 42500 C ₁₉ H ₁₇ N ₃ ·HCl F. Wt. 323.82	25 gm 100 gm	
1025	PATTON & REEDER'S REAGENT, AR GRADE (Calcon Carboxylic Acid) CAS 3737-95-9 C ₂₁ H ₁₄ N ₂ O ₇ S F. Wt. 438.42	5 gm 25 gm	
1405	PECTIN, PURE GRADE (Poly-D-galaeturonic acid methyl ester) CAS 9000-69-5	100 gm 500 gm	
3392	PECTINASE * CAS 9032-75-1	25 gm	
3188	1-PENTANE SULPHONIC ACID SODIUM SALT (ANH.), (For HPLC) CAS 22767-49-3 C ₅ H ₁₁ NaO ₃ S F. Wt. 192.21	25 gm	
3395	1-PENTANE SULPHONIC ACID SODIUM SALT (MONO), (For HPLC) CAS 207605-40-1 C ₅ H ₁₁ NaO ₃ S·H ₂ O F. Wt. 192.2	25 gm	
1406	PEPSIN 1:3,000 * CAS 9001-75-6	100 gm 500 gm 5 Kg	
2185	PEPSIN 1:10,000 * CAS 9001-75-6	25 gm 100 gm	
1266	PEPTONE, PHARMA GRADE (See more Peptones in Biological Media Base)	500 gm 25 Kg	
3518	PEPTONE, BACTO GRADE (See more Peptones in Biological Media Base)	500 gm 5 Kg 25 Kg	
3521	PEPTONE PASTE Bacteriological Grade	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
3519	PEPTONE GR (Granular), BACTO GRADE	500 gm 5 Kg 25 Kg	
2197	1,10, PHENANTHROLINE (Mono), AR GRADE CAS 5144-89-8 C ₁₂ H ₈ N ₂ ·H ₂ O F. Wt. 198.22	5 gm 25 gm	
2198	1,10, PHENANTHROLINE HYDROCHLORIDE, AR GRADE CAS 3829-86-5 C ₁₂ H ₈ N ₂ ·H ₂ O F. Wt. 234.69	5 gm 25 gm	
2519	PHENAZONE (Antipyrine) CAS 60-80-0 C ₁₁ H ₁₂ N ₂ O F. Wt. 188.23	100 gm 500 gm	
293	PHENOL CRYSTAL, EXTRA PURE CAS 108-95-2 C ₆ H ₆ O F. Wt. 94.11	500 gm	
1670	PHENOLPHTHALEIN (White) (pH 8.2-10.0) CAS 77-09-8 C ₂₀ H ₁₄ O ₄ F. Wt. 318.32	25 gm 100 gm 500 gm 25 Kg	
TBL 070	PHENOLPHTHALEIN, 1% W/V CAS 77-09-8 C ₂₀ H ₁₄ O ₄ F. Wt. 318.32	125 ml 500 ml	
1671	PHENOL RED (pH Indicator) CAS 143-74-8 C ₁₉ H ₁₄ O ₅ S F. Wt. 354.38	5 gm 25 gm 100 gm	
1671A	PHENOL RED SODIUM (Water Soluble), AR GRADE CAS 34487-61-1 C ₁₉ H ₁₃ NaO ₅ S F. Wt. 376.36	25 gm	
906	PHENYLACETIC ACID (α-Tolylic Acid) CAS 103-82-2 C ₈ H ₈ O ₂ F. Wt. 136.15	500 gm	
1126	DL-PHENYLALANINE 99% + Crystalline CAS 150-30-1 C ₉ H ₁₁ NO ₂ F. Wt. 165.19	25 gm 100 gm	
1127	L-PHENYLALANINE 99%+Crystalline CAS 63-91-2 C ₉ H ₁₁ NO ₂ F. Wt. 165.19	25 gm 100 gm 500 gm	
2669	N-PHENYLANTHRANILIC ACID, AR GRADE (Redox Indicator) CAS 91-40-7 C ₁₃ H ₁₁ NO ₂ F. Wt. 213.23	25 gm 100 gm	
2208	PHENYLHYDRAZINE HYDROCHLORIDE, EXTRA PURE CAS 59-88-1 C ₆ H ₈ N ₂ ·HCl F. Wt. 144.60	100 gm 250 gm	
903	PHENYLMERCURIC ACETATE CAS 62-38-4 C ₈ H ₈ HgO ₂ F. Wt. 336.74	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
904	PHENYLMERCURIC NITRATE, AR GRADE CAS 8003-05-02 C ₆ H ₅ HgNO ₃ F. Wt. 634.40	25 gm 100 gm	
2457	PHLOROGLUCINOL, AR GRADE (1,3,5, Trihydroxybenzene) CAS 108-73-6 C ₆ H ₆ O ₃ F. Wt. 126.11	25 gm 100 gm	
2211	PHLOXIN B (Cyanosine, Eosin 10 B) CAS 18472-87-2 C ₂₀ H ₂ Br ₄ Cl ₄ Na ₂ O ₅ F. Wt. 829.64	25 gm 100 gm	
2972	PHOSPHOMOLYBDIC ACID, EXTRA PURE CAS 51429-74-4 H ₃ PO ₄ ·12MoO ₃ ·xOH ₂ O F. Wt. 1825.25	25 gm	
773	PHOSPHOMOLYBDIC ACID, AR GRADE CAS 51429-74-4 H ₃ PO ₄ ·12MoO ₃ ·xOH ₂ O F. Wt. 1825.25	25 gm	
2000	PHOSPHOROUS PENTOXIDE, AR GRADE CAS 1314-56-3 P ₂ O ₅ F. Wt. 283.89	500 gm	
981	PHOSPHOTUNGSTIC ACID, EXTRA PURE CAS 12501-23-4 H ₃ PO ₄ ·12WO ₃ ·XH ₂ O F. Wt. 2880.05	100 gm 500 gm	
2673	PTHALIC ACID, EXTRA PURE CAS 88-99-3 C ₈ H ₆ O ₄ F. Wt. 166.13	500 gm	
907	PTHALIC ANHYDRIDE, EXTRA PURE CAS 85-44-9 C ₈ H ₄ O ₃ F. Wt. 148.1	500 gm 25 Kg	
2450	PTHALIMIDE CAS 85-41-6 C ₈ H ₅ NO ₂ F. Wt. 147.13	500 gm	
3576	PICLORAM CAS 1918-02-1 C ₆ H ₃ Cl ₃ N ₂ O ₂ F. Wt. 241.46	1 gm 5 gm	
4101	PICRIC ACID, EXTRA PURE (Moistened With Water) CAS 88-89-1 C ₆ H ₃ N ₃ O ₇ F. Wt. 229.10	500 gm	
770	PIPERAZINE (ANH.) CAS 110-85-0 C ₄ H ₁₀ N ₂ F. Wt. 86.14	100 gm 500 gm 25 Kg	
2219	PIPES (Buffer) (Piperazine-N, N-bis (2-ethane Sulphonic Acid) CAS 5625-37-6 C ₆ H ₁₆ N ₂ O ₆ S ₂ F. Wt. 302.4	5 gm 25 gm	
4367	γ-PGA (poly-γ-glutamic acid) CAS 49717-32-0 (C ₅ H ₉ NO ₄) _x F. Wt. 147.13	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
908	POLYETHYLENE GLYCOL-4000 (PEG-4000) CAS 25322-68-3	500 gm 25 Kg	
1013	POLYETHYLENE GLYCOL-6000 (PEG-6000) CAS 25322-68-3	500 gm 25 Kg	
4361	POLYSORBATE 20 (Commercial Grade) CAS 9005-64-5	500 ml 25 Ltr	
4362	POLYSORBATE 80 (Commercial Grade) CAS 9005-65-6	500 ml 25 Ltr	
4363	POLYSORBATE 85 (Commercial Grade) CAS 9005-70-3	500 ml 25 Ltr	
2224	POLYVINYL ALCOHOL M.W. 125,000 CAS 9002-89-5 (C ₂ H ₄ O) _n F. Wt. 165.19	500 gm	
2115	POLYVINYL PYRROLIDINE K-30 (PVP) CAS 9003-39-8 (C ₆ H ₉ O) _n F. Wt. 165.19	100 gm 500 gm 25 Kg	
1675	PONCEAU-S, Sodium Salt CAS 6226-79-5 C.I 27195 C ₂₂ H ₁₂ N ₄ Na ₄ O ₁₃ S ₄ F. Wt. 760.57	25 gm	
380	POTASSIUM ACETATE, EXTRA PURE CAS 127-08-2 C ₂ H ₃ KO ₂ F. Wt. 98.15	500 gm 50 Kg	
381	POTASSIUM ACETATE, AR GRADE CAS 127-08-2 C ₂ H ₃ KO ₂ F. Wt. 98.15	500 gm 50 Kg	
753	POTASSIUM ALUM, EXTRA PURE CAS 7784-24-9 ALK ₂ O ₃ ·12H ₂ O F. Wt. 474.39	500 gm 5 Kg 50 Kg	
230	POTASSIUM ALUM, AR GARDE CAS 7784-24-9 ALK ₂ O ₃ ·12H ₂ O F. Wt. 474.39	500 gm 5 Kg 50 Kg	
2517	POTASSIUM ANTIMONY TARTRATE (Trihydrate), EXTRA PURE CAS 28300-74-5 C ₈ H ₄ K ₂ O ₁₂ Sb ₂ ·3H ₂ O F. Wt. 667.87	100 gm 500 gm 50 Kg	
811	POTASSIUM ANTIMONY TARTRATE (Trihydrate), AR GRADE CAS 28300-74-5 C ₈ H ₄ K ₂ O ₁₂ Sb ₂ ·3H ₂ O F. Wt. 667.87	100 gm 500 gm	
382	POTASSIUM BICARBONATE, EXTRA PURE (Potassium Hydrogen Carbonate) CAS 298-14-6 KHCO ₃ F. Wt. 100.115	500 gm 50 Kg	
4226	POTASSIUM BICARBONATE, AR GRADE CAS 298-14-6 KHCO ₃ F. Wt. 100.12	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
3511	POTASSIUM BICARBONATE (R.O. SPL), For R.O. Water CAS 298-14-6 KHCO ₃ F. Wt. 100.12	500 gm 5 Kg 50 Kg	
2686	POTASSIUM BISULPHATE, EXTRA PURE (Potassium Hydrogen Sulphate) CAS 7646-93-7 KHSO ₄ F. Wt. 136.17	500 gm	
2405	POTASSIUM BISULPHITE, EXTRA PURE (Potassium Hydrogen Sulphite) CAS 7773-03-7 KHSO ₃ F. Wt. 120.16	500 gm 5 Kg	
919	POTASSIUM BITARTRATE (Potassium Hydrogen Tartrate) CAS 868-14-4 Kc ₄ H ₆ O ₆ F. Wt. 188.18	500 gm	
711	POTASSIUM BROMATE, EXTRA PURE CAS 7758-01-2 KBrO ₃ F. Wt. 167.00	500 gm 50 Kg	
628	POTASSIUM BROMIDE, EXTRA PURE CAS 7758-02-3 Kbr F. Wt. 119.002	500 gm 50 Kg	
2680	POTASSIUM BROMIDE, AR GRADE CAS 7758-02-3 Kbr F. Wt. 119.002	500 gm	
4102	POTASSIUM BROMIDE IR (for Spectroscopy) CAS 7758-02-3 Kbr F. Wt. 119.00	100 gm	
383	POTASSIUM CARBONATE (ANH.), EXTRA PURE CAS 584-08-7 K ₂ CO ₃ F. Wt. 138.21	500 gm 50 Kg	
384	POTASSIUM CARBONATE (ANH.), AR GRADE CAS 584-08-7 K ₂ CO ₃ F. Wt. 138.21	500 gm 50 Kg	
385	POTASSIUM CHLORIDE, EXTRA PURE CAS 7447-40-7 KCl F. Wt. 74.55	500 gm 5 Kg 50 Kg	
386	POTASSIUM CHLORIDE, AR GRADE CAS 7447-40-7 KCl F. Wt. 74.55	500 gm 5 Kg 50 Kg	
2226	POTASSIUM CHLOROPLATINATE (Hexa) (40% Platinum) CAS 16921-30-5 K ₂ PtCl ₆ F. Wt. 485.98	1 gm	
911	POTASSIUM CHROMATE, EXTRA PURE CAS 7789-00-6 K ₂ CrO ₄ F. Wt. 194.19	500 gm 50 Kg	
912	POTASSIUM CHROMATE, AR GRADE CAS 7789-00-6 K ₂ CrO ₄ F. Wt. 194.19	500 gm	
913	tri-POTASSIUM CITRATE (Mono) CAS 6100-05-6 C ₆ H ₈ K ₃ O ₇ ·H ₂ O F. Wt. 324.41	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
387	POTASSIUM DICHROMATE, EXTRA PURE CAS 7778-50-9 $K_2Cr_2O_7$ F. Wt. 294.18	500 gm 5 Kg 50 Kg	
388	POTASSIUM DICHROMATE, AR GRADE CAS 7778-50-9 $K_2Cr_2O_7$ F. Wt. 294.18	500 gm 50 Kg	
389	POTASSIUM DIHYDROGENORTHOPH PHOSPHATE (Mono Potassium Phosphate), (ANH.), EXTRA PURE CAS 7778-77-0 KH_2PO_4 F. Wt. 136.09	500 gm 5 Kg 50 Kg	
220	POTASSIUM DIHYDROGEN ORTHOPHOSPHATE (Mono Potassium Phosphate), (ANH.), AR GRADE CAS 7778-77-0 KH_2PO_4 F. Wt. 136.09	500 gm 5 Kg 50 Kg	
737	POTASSIUM FERRICYANIDE, EXTRA PURE CAS 13746-66-2 $K_3Fe(CN)_6$ F. Wt. 329.25	100 gm 500 gm 50 Kg	
631	POTASSIUM FERROCYANIDE, EXTRA PURE CAS 13943-58-3 $K_4Fe(CN)_6$ F. Wt. 368.35	500 gm 50 Kg	
4327	POTASSIUM FLUORIDE (ANH.) 97.0%, AR GRADE CAS 7789-23-3 F.Wt. 58.10	500 gm 25 Kg	
4328	POTASSIUM FLUORIDE (ANH) 97.0%, EXTRA PURE CAS 7789-23-3 F.Wt. 58.10	500 gm 25 Kg	
916	POTASSIUM FLUORIDE (ANH.), EXTRA PURE CAS 7789-23-3 KF F. Wt. 58.10	500 gm	
2229	POTASSIUM FORMATE, EXTRA PURE CAS 590-29-4 $CHKO_2$ F. Wt. 84.12	100 gm 500 gm	
636	di-POTASSIUM HYDROGEN ORTHOPHOSPHATE (ANH), EXTRA PURE CAS 7758-11-4 K_2HPO_4 F. Wt. 174.17	500 gm 5 Kg 50 Kg	
218	di-POTASSIUM HYDROGEN ORTHOPHOSPHATE (ANH.), AR GRADE CAS 7758-11-14 K_2HPO_4 F. Wt. 174.17	500 gm 5 Kg 50 Kg	
509	POTASSIUM HYDROGEN PHTHALATE (Potassium Biphthalate), EXTRA PURE CAS 877-24-7 $C_8H_4KO_4$ F. Wt. 204.22	500 gm 50 Kg	
877	POTASSIUM HYDROGEN PHTHALATE, AR GRADE CAS 877-24-7 $C_8H_4KO_4$ F. Wt. 204.22	500 gm 50 Kg	
2686	POTASSIUM HYDROGEN SULPHATE, EXTRA PURE CAS 7646-93-7 $KHSO_4$ F. Wt. 136.17	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
392	POTASSIUM HYDROGEN SULPHATE, AR GRADE CAS 7646-93-7 $KHSO_4$ F. Wt. 136.17	500 gm 50 Kg	
919	POTASSIUM HYDROGEN TARTRATE, (Potassium Bitartrate) CAS 868-14-4 $C_4H_5KO_6$ F. Wt. 188.18	500 gm 50 Kg	
2230	POTASSIUM HYDROXIDE FLAKES, EXTRA PURE CAS 1310-58-3 KOH F. Wt. 56.11	500 gm 5 Kg 50 Kg	
504	POTASSIUM HYDROXIDE PELLETS, EXTRA PURE CAS 1310-58-3 KOH F. Wt. 56.11	500 gm 5 Kg 50 Kg	
511	POTASSIUM HYDROXIDE PELLETS, AR GRADE CAS 1310-58-3 KOH F. Wt. 56.11	500 gm 5 Kg 50 Kg	
393	POTASSIUM IODATE, EXTRA PURE CAS 7758-05-6 KIO_3 F. Wt. 214.001	100 gm 500 gm 25 Kg	
2231	POTASSIUM IODATE, AR GRADE CAS 7758-05-6 KIO_3 F. Wt. 214.001	100 gm 500 gm 25 Kg	
394	POTASSIUM IODIDE (Confirming IP), EXTRA PURE CAS 7681-11-0 KI F. Wt. 166.00	100 gm 500 gm 25 Kg	
634	POTASSIUM IODIDE, AR GRADE CAS 7681-11-0 KI F. Wt. 166.00	100 gm 500 gm 25 Kg	
3200	POTASSIUM LACTATE 50% CAS 996-31-6 $C_3H_5KO_3$ F. Wt. 128.17	500 ml	
3128	POTASSIUM METABISULPHITE, (Potassium Pyrosulphite), EXTRA PURE CAS 16731-55-8 $K_2S_2O_5$ F. Wt. 222.32	500 gm 5 Kg 50 Kg	
397	POTASSIUM METAPERIODATE, AR GRADE CAS 7790-21-8 KIO_4 F. Wt. 230	100 gm 500 gm	
761	POTASSIUM NITRATE, EXTRA PURE CAS 7757-79-1 KNO_3 F. Wt. 101.10	500 gm 5 Kg 50 Kg	
399	POTASSIUM NITRATE, AR GRADE CAS 7757-79-1 KNO_3 F. Wt. 101.1	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
505	POTASSIUM OXALATE (Mono), (Dipotassium Oxalate), EXTRA PURE CAS 6487-48-5 $C_2K_2O_4 \cdot H_2O$ F. Wt. 85.10	500 gm 50 Kg	
920	POTASSIUM OXALATE, AR GRADE CAS 6487-48-5 $C_2K_2O_4 \cdot H_2O$ F. Wt. 85.10	500 gm 50 Kg	
400	POTASSIUM PERMANGANATE, EXTRA PURE CAS 7722-64-7 KMnO ₄ F. Wt. 158.03	500 gm 5 Kg 50 Kg	
2236	POTASSIUM PERMANGANATE, AR GRADE CAS 7722-64-7 KMnO ₄ F. Wt. 158.03	500 gm	
2237	POTASSIUM PERSULPHATE, EXTRA PURE (Potassium Peroxodisulphate) CAS 7727-21-1 $K_2S_2O_8$ F. Wt. 270.32	500 gm 50 Kg	
2238	POTASSIUM PERSULPHATE, AR GRADE CAS 7727-21-1 $K_2S_2O_8$ F. Wt. 270.32	500 gm	
636	POTASSIUM PHOSPHATE DIBASIC (ANH.) (Dipotassium Hydrogen Orthophosphate), EXTRA PURE CAS 7758-11-4 K_2HPO_4 F. Wt. 174.18	500 gm 5 Kg 50 Kg	
218	POTASSIUM PHOSPHATE DIBASIC (ANH.) (Dipotassium Hydrogen Orthophosphate), AR GRADE CAS 7758-11-4 K_2HPO_4 F. Wt. 174.18	500 gm 5 Kg 50 Kg	
389	POTASSIUM PHOSPHATE MONOBASIC, EXTRA PURE CAS 7778-77-0 KH ₂ PO ₄ F. Wt. 174.18	500 gm 5 Kg 50 Kg	
3259	tri-POTASSIUM PHOSPHATE CAS 7778-53-2 K ₃ PO ₄ F. Wt. 212.27	500 gm	
3258	POTASSIUM SILICATE POWDER CAS 1312-76-1 K_2O_3Si F. Wt. 154.28	500 gm 50 Kg	
403	POTASSIUM SODIUM TARTRATE (Tetra), EXTRA PURE CAS 6381-59-5 $C_4H_4KNaO_6 \cdot 4H_2O$ F. Wt. 282.22	100 gm 500 gm 25 Kg	
404	POTASSIUM SODIUM TARTRATE (Tetra), AR GRADE CAS 6381-59-5 $C_4H_4KNaO_6 \cdot 4H_2O$ F. Wt. 282.22	100 gm 500 gm 25 Kg	
405	POTASSIUM SORBATE (Granular), AR GRADE CAS 24634-61-5 $C_6H_7KO_2$ F. Wt. 150.22	500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
2687	POTASSIUM NITRITE, EXTRA PURE CAS 7758-09-0 KNO ₂ F. Wt. 85.10	500 gm 50 Kg	
4018	POTASSIUM STEARATE CAS 593-29-3 $C_{18}H_{35}KO_2$ F. Wt. 322.57	500 gm	
406	POTASSIUM SULPHATE, EXTRA PURE CAS 7778-80-5 K ₂ SO ₄ F. Wt. 174.26	500 gm 50 Kg	
407	POTASSIUM SULPHATE, AR GRADE CAS 7778-80-5 K ₂ SO ₄ F. Wt. 176.26	500 gm 50 Kg	
408	POTASSIUM TELLURITE, EXTRA PURE (for Bacteriology) CAS 7790-58-1 K_2TeO_3 F. Wt. 253.81	25 gm 100 gm	
2225	POTASSIUM TERT-BUTOXIDE (for Synthesis) CAS 865-47-4 C_4H_9KO F. Wt. 112.21	100 gm 500 gm 25 Kg	
699	POTASSIUM TETRAOXALATE (Dihydrate) CAS 6100-20-5 $C_4H_3KO_8 \cdot 2H_2O$ F. Wt. 254.19	500 gm	
409	POTASSIUM THIOCYANATE, EXTRA PURE (Potassium Sulphocyanide) CAS 333-20-0 KSCN F. Wt. 97.18	500 gm 50 Kg	
410	POTASSIUM THIOCYANATE, AR GRADE CAS 333-20-0 KSCN F. Wt. 97.18	500 gm	
2694	POTASSIUM TITANIUM OXALATE, AR GRADE CAS 14402-67-6 $C_4K_2O_5Ti \cdot 2H_2O$ F. Wt. 354.13	100 gm 500 gm	
1309	POTATO STARCH (Insoluble) CAS 9005-25-8	500 gm	
1128	L-PROLINE, 99% + Purity, Crystalline CAS 147-85-3 $C_5H_9NO_2$ F. Wt. 115.13	5 gm 25 gm 500 gm	
211	n-PROPYL GALLATE (Antioxidant & Preservative) CAS 121-79-9 $C_{10}H_{12}O_5$ F. Wt. 212.20	100 gm 500 gm	
922	PROPYL PARABEN (Propyl-P-Hydroxybenzoate) CAS 94-13-3 $C_{10}H_{12}O_3$ F. Wt. 180.2	500 gm	
637	PROPYL PARABEN SODIUM (Propyl-P-Hydroxy Benzoate Sodium) CAS 35285-69-9 $C_{10}H_{11}NaO_3$ F. Wt. 202.00	500 gm	
648	PROPYLENE GLYCOL (1,2-Propanediol) CAS 57-55-6 $C_3H_8O_2$ F. Wt. 76.09	500 ml 5 Ltr 50 Ltr	

CODE	PRODUCT NAME	PACK SIZE	
639	PYRIDOXAL HYDROCHLORIDE CAS 65-22-5 $C_8H_9NO_3 \cdot HCl$ F. Wt. 203.62	1 gm 5 gm	
1705	PYRIDOXINE HYDROCHLORIDE (VITAMIN B) CAS 58-56-0 $C_8H_{11}NO_3 \cdot HCl$ F. Wt. 205.64	10 gm 25 gm 100 gm	
171	PYROCATECHOL (1,2 Dihydroxybenzene) CAS 120-80-9 $C_6H_6O_3$ F. Wt. 110.10	100 gm 500 gm	
2148	PYROGALLOL, EXTRA PURE CAS 87-66-1 $C_6H_6O_3$ F. Wt. 126.11	100 gm 500 gm	
925	PYROGALLOL, AR GRADE (Pyrogalllic Acid, 1,2,3, Trihydroxybenzene) CAS 87-66-1 $C_6H_6O_3$ F. Wt. 126.11	100 gm 500 gm	
3204	PYRONINE-G (PYRONINE) CAS 92-32-0 C.I. 45005 $C_{17}H_{19}ClN_2O$ F. Wt. 302.8	5 gm	
659	PYRUVIC ACID SODIUM, AR GRADE (Sodium Pyruvate) CAS 113-24-6 $C_3H_3NaO_3$ F. Wt. 110.0	100 gm 500 gm	
3207	QUININE HYDROCHLORIDE (Dihydrate) CAS 6119-47-7 $C_{20}H_{24}N_2O_2 \cdot HCl \cdot 2H_2O$ F. Wt. 396.92	25 gm	
2632	8-QUINOLINOL (8-Hydroxyquinoline) (Oxine) CAS 148-24-3 C_8H_7NO F. Wt. 145.16	100 gm 500 gm	
2252	QUINOLINE YELLOW (Water Soluble) CAS 8004-92-0 C.I. 47005 $C_{18}H_{15}NNa_2O_8S_8$ F. Wt. 477.38	5 gm 25 gm	
698	RNA (From Torula Yeast) CAS 63231-63-0	25 gm 100 gm 500 gm	
TC 018	D-RAFFINOSE (Penta) CAS 17629-30-0 $C_{18}H_{32}O_{16} \cdot 5H_2O$ F. Wt. 594.52	10 gm 25 gm	
507	RESAZURIN SODIUM CAS 62758-13-8 $C_{12}H_6NO_4Na$ F. Wt. 251.18	1 gm 5 gm	
2981	RESORCINOL (1,3 Dihydroxybenzene), EXTRA PURE CAS 108-46-3 $C_6H_6O_2$ F. Wt. 110.1	250 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
411	L-RHAMNOSE 99% PURITY CAS 10030-85-0 $C_6H_{12}O_5 \cdot H_2O$ F. Wt. 182.17	5 gm 25 gm 1 Kg	
1609	RHODAMINE-B CAS 81-88-9 C.I. 45170 $C_{28}H_{31}ClN_2O_3$ F. Wt. 479.02	25 gm 100 gm	
1703	RIBOFLAVIN (Vitamin B), Vitamin G CAS 83-88-5 $C_{17}H_{20}N_4O_6$ F. Wt. 376.36	10 gm 25 gm 100 gm	
929	RIBOFLAVIN-5-PHOSPHATE, EXTRA PURE * CAS 130-40-5 $C_{17}H_{20}N_4NaO_9P_2H_2O$ F. Wt. 514.36	5 gm 25 gm	
1306	D-RIBOSE CAS 50-69-1 $C_5H_{10}O_5$ F. Wt. 150.13	10 gm 25 gm 100 gm	
1603	ROSANILINE HCl CAS 632-99-5 $C_{20}H_{20}N_3 \cdot HCl$ F. Wt. 337.86	25 gm 100 gm 1 Kg	
1677	ROSE BENGAL (Acid red 94) CAS 632-69-9 C.I. 45440 $C_{20}H_2Cl_2I_4Na_2O_5$ F. Wt. 1017.64	25 gm 100 gm	
1664	p-ROSOLIC ACID (Aurin) CAS 603-45-2 C.I. 43800 $C_{19}H_{14}O_3$ F. Wt. 290.32	25 gm	
2435	RUTIN (Trihydrate) (Vitamin P) (for Biochemistry) CAS 153-18-4 $C_{27}H_{30}O_{16} \cdot 3H_2O$ F. Wt. 628.50	25 gm 100 gm 1 Kg	
1646	SAFRANINE-T (Safranin-O) CAS 477-73-6 C.I. 50240 $C_{20}H_{19}ClN_4$ F. Wt. 350.85	25 gm 100 gm 500 gm	
TC 019	D-SALICIN CAS 138-52-3 $C_{13}H_{18}O_7$ F. Wt. 286.28	5 gm 10 gm 25 gm	
691	SALICYLIC ACID, EXTRA PURE (2-Hydroxy Benzoic Acid) CAS 69-72-7 $C_7H_6O_3$ F. Wt. 138.12	500 gm 25 Kg	
680	SALICYCLIC ACID, AR GRADE (2-Hydroxy Benzoic Acid) CAS 69-72-7 $C_7H_6O_3$ F. Wt. 138.12	500 gm	
4227	SAPONIN CAS 8047-15-2	100 gm 500 gm	
553	SELENIUM DIOXIDE CAS 7446-08-4 SeO_2 F. Wt. 110.96	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
3620	L-SELENOMETHIONINE CAS 3211-76-5 C ₆ H ₁₁ NO ₂ Se F. Wt. 196.11	1 gm	
1130	L-SERINE 99% + Purity, Crystalline CAS 56-45-1 C ₃ H ₇ NO ₃ F. Wt. 105.09	10 gm 25 gm 1 Kg	
412	SILICA GEL (Self Indicating Blue) 5-8 mm CAS 63231-67-4	500 gm 5 Kg 50 Kg	
2257	SILICA GEL (Self Indicating White) 5-8 mm CAS 63231-67-4	500 gm 5 Kg	
414	SILICA GEL, 60-120 Mesh (for Column Chromatography) for CC CAS 63231-67-4	500 gm	
2258	SILICA GEL, 100-200 Mesh CAS 63231-67-4	500 gm	
2259	SILICA GEL, 60-200 Mesh CAS 63231-67-4	500 gm	
2260	SILICA GEL, 200-400 Mesh For TLC CAS 63231-67-4	500 gm	
415	SILICA GEL-G For TLC CAS 112926-00-8	500 gm	
416	SILICA GEL HF 254 For TLC CAS 112926-00-8	500 gm	
417	SILICA GEL-H for TLC W/O Binder CAS 63231-67-4	500 gm	
418	SILICA GEL GF 254 for TLC W/O Binder CAS 63231-67-4	500 gm	
642	SILICONE, High Vacuum Grease (Vacuum Grease)	500 gm	
932	SILICONE ANTIFOAMING AGENT (NON IONIC)	100 ml 500 ml 35 Ltr	
644	SILICONE OIL (for Oil Bath, upto 250°C) CAS 63148-62-9 C ₆ H ₁₈ OSi ₂ F. Wt. 162.38	100 ml 500 ml	
645	SILICOTUNGSTIC ACID, EXTRA PURE CAS 12027-38-2 H ₄ O ₄₀ SiW ₁₂ F. Wt. 2878.17	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
2708	SILICOTUNGSTIC ACID, 99.5%, AR GRADE CAS 12027-38-2 H ₄ O ₄₀ SiW ₁₂ F. Wt. 2878.17	25 gm 100 gm	
2709	SILVER ACETATE (Acetic Acid Silver Salt), EXTRA PURE CAS 563-63-3 C ₂ H ₃ AgO ₂ F. Wt. 166.91	25 gm	
954	SILVER CARBONATE, EXTRA PURE CAS 534-16-7 CAg ₂ O ₃ F. Wt. 275.75	25 gm	
512	SILVER CHLORIDE, EXTRA PURE CAS 7783-90-6 AgCl F. Wt. 143.32	25 gm	
2982	SILVER IODIDE CAS 7783-96-2 AgI F. Wt. 234.77	25 gm	
419	SILVER NITRATE, EXTRA PURE CAS 7761-88-8 AgNO ₃ F. Wt. 169.87	10 gm 25 gm	
649	SILVER NITRATE, AR GRADE CAS 7761-88-8 AgNO ₃ F. Wt. 169.87	10 gm 25 gm	
2266	SILVER OXIDE, EXTRA PURE CAS 20667-12-3 Ag ₂ O F. Wt. 123.87	25 gm	
2267	SILVER SULPHATE, EXTRA PURE CAS 10294-26-5 Ag ₂ SO ₄ F. Wt. 311.80	25 gm	
650	SILVER SULPHATE, AR GRADE CAS 10294-26-5 Ag ₂ SO ₄ F. Wt. 311.80	25 gm	
651	SKIM MILK POWDER (Fat Free)	500 gm	
422	SODIUM ACETATE (ANH.), EXTRA PURE CAS 127-09-3 C ₂ H ₃ NaO ₂ F. Wt. 82.03	500 gm 5 Kg 50 Kg	
221	SODIUM ACETATE (ANH.), AR GRADE CAS 127-09-3 C ₂ H ₃ NaO ₂ F. Wt. 82.03	500 gm 5 Kg 50 Kg	
420	SODIUM ACETATE (Trihydrate), EXTRA PURE CAS 6131-90-4 C ₂ H ₃ NaO ₂ .3H ₂ O F. Wt. 100.04	500 gm 5 Kg 50 Kg	
2270	SODIUM ACETATE (Trihydrate), AR GRADE CAS 6131-90-4 C ₂ H ₃ NaO ₂ .3H ₂ O F. Wt. 136.08	500 gm 5 Kg 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
435	SODIUM ACID PHOSPHATE (Sodium dihydrogen orthophosphate) CAS 13472-35-0 $\text{NaH}_2\text{PO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 156.01	500 gm 5 Kg 50 Kg	
3623	SODIUM ACID PYROPHOSPHATE (SAPP), (Food Grade) CAS 7758-16-9 $\text{Na}_2\text{H}_2\text{P}_2\text{O}_7$ F. Wt. 221.94	500 gm 5 Kg	
423	SODIUM ALGINATE, EXTRA PURE (Sodium Polymannuronate) CAS 9005-38-3	500 gm 5 Kg 25 Kg	
2983	SODIUM ALGINATE (Alginic Acid Sodium Salt), AR GRADE CAS 9005-38-3 $\text{C}_{14}\text{H}_{22}\text{O}_{13}$ F. Wt. 101.1	500 gm	
3424	SODIUM ARSENATE (Hepta), EXTRA PURE CAS 10048-95-0 $\text{Na}_2\text{HAsO}_4 \cdot 7\text{H}_2\text{O}$ F. Wt. 312.01	25 gm 100 gm	
4104	SODIUM ARSENATE (Hepta), AR GRADE CAS 10048-95-0 $\text{Na}_2\text{HAsO}_4 \cdot 7\text{H}_2\text{O}$ F. Wt. 312.01	25 gm 100 gm	
2273	SODIUM ARSENITE, AR GRADE CAS 7784-46-5 NaAsO_2 F. Wt. 129.91	100 gm 500 gm	
1714	SODIUM ASCORBATE (Ascorbic Acid Sodium) CAS 134-03-2 $\text{C}_2\text{H}_3\text{NaO}_6$ F. Wt. 198.11	100 gm 500 gm	
424	SODIUM AZIDE, EXTRA PURE CAS 26628-22-8 NaN_3 F. Wt. 65.01	100 gm 500 gm	
425	SODIUM AZIDE, AR GRADE CAS 26628-22-8 NaN_3 F. Wt. 65.01	100 gm 500 gm	
426	SODIUM BENZOATE, EXTRA PURE CAS 532-32-1 $\text{C}_7\text{H}_5\text{NaO}_2$ F. Wt. 144.11	500 gm 25 Kg	
427	SODIUM BENZOATE, AR GRADE CAS 532-32-1 $\text{C}_7\text{H}_5\text{NaO}_2$ F. Wt. 144.11	500 gm 50 Kg	
429	SODIUM BICARBONATE, EXTRA PURE (Sodium Hydrogen Carbonate) CAS 144-55-8 NaHCO_3 F. Wt. 84.01	500 gm 5 Kg 50 Kg	
430	SODIUM BICARBONATE, AR GRADE CAS 144-55-8 NaHCO_3 F. Wt. 84.01	500 gm 5 Kg 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
2984	SODIUM BICARBONATE 'SPL' (for R.O. Water) CAS 144-55-8 NaHCO_3 F. Wt. 84.01	500 gm 5 Kg 50 Kg	
428	SODIUM BISELENITE (Sodium Hydrogen Selenite) CAS 7782-82-3 NaHSeO_3 F. Wt. 150.95	100 gm 1 Kg	
780	SODIUM BISULPHATE, EXTRA PURE (Sodium Hydrogen Sulphate) CAS 10034-88-5 $\text{NaHSO}_4 \cdot \text{H}_2\text{O}$ F. Wt. 138.07	500 gm	
2292	SODIUM BISULPHATE, AR GRADE (Sodium Hydrogen Sulphate) CAS 10034-88-5 $\text{NaHSO}_4 \cdot \text{H}_2\text{O}$ F. Wt. 138.07	500 gm	
946	SODIUM BISULPHITE, EXTRA PURE (Sodium Hydrogen Sulphite) CAS 7631-90-5 NaHSO_3 F. Wt. 104.06	500 gm	
655	SODIUM BISULPHITE, AR GRADE (Sodium Hydrogen Sulphite) CAS 7631-90-5 NaHSO_3 F. Wt. 104.06	500 gm	
2275	SODIUM BROMATE, EXTRA PURE CAS 7789-38-0 NaBrO_3 F. Wt. 150.89	500 gm 50 Kg	
936	SODIUM BROMIDE, EXTRA PURE CAS 7647-15-6 NaBr F. Wt. 102.89	500 gm 50 Kg	
2276	SODIUM BROMIDE, AR GRADE CAS 7647-15-6 NaBr F. Wt. 102.89	500 gm 50 Kg	
4008	SODIUM BUTYRATE CAS 156-54-7 $\text{C}_4\text{H}_7\text{NaO}_2$ F. Wt. 110.09	100 gm 500 gm 25 Kg	
2278	SODIUM CACODYLATE (Trihydrate) CAS 6131-99-3 $\text{C}_2\text{H}_6\text{AsNaO}_2 \cdot 3\text{H}_2\text{O}$ F. Wt. 214.03	25 gm 100 gm	
431	SODIUM CARBONATE (ANH.), EXTRA PURE CAS 497-19-8 Na_2CO_3 F. Wt. 105.99	500 gm 5 Kg 50 Kg	
432	SODIUM CARBONATE (ANH.), AR GRADE CAS 497-19-8 Na_2CO_3 F. Wt. 105.99	500 gm 5 Kg 50 Kg	
938	SODIUM CARBONATE (Mono), EXTRA PURE CAS 5968-11-6 $\text{Na}_2\text{CO}_3 \cdot \text{H}_2\text{O}$ F. Wt. 124.00	500 gm 50 Kg	
2280	SODIUM CARBONATE (Mono), AR GRADE CAS 5968-11-6 $\text{Na}_2\text{CO}_3 \cdot \text{H}_2\text{O}$ F. Wt. 124.00	500 gm 50 Kg	
4252	SODIUM CASEINATE, AR GRADE CAS 9005-46-3	500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
215	SODIUM CHLORIDE, EXTRA PURE CAS 7647-14-5 NaCl F. Wt. 58.44	500 gm 5 Kg 50 Kg	
216	SODIUM CHLORIDE, AR GRADE CAS 7647-14-5 NaCl F. Wt. 58.44	500 gm 5 Kg 50 Kg	
924	SODIUM CHLORIDE (for Tissue Culture), MOLECULAR BIOLOGY GRADE CAS 7647-14-5 NaCl F. Wt. 58.44	500 gm 5 Kg 50 Kg	
4178	SODIUM CHOLATE (Cholic Acid Sodium) CAS 361-09-1 C ₂₄ H ₃₉ NaO ₅ F. Wt. 430.56	25 gm 100 gm 25 Kg	
2281	SODIUM CHROMATE (Tetra), EXTRA PURE CAS 10034-82-9 Na ₂ CrO ₄ ·4H ₂ O F. Wt. 234.03	500 gm 25 Kg	
2882	SODIUM CHROMATE (Tetra), AR GRADE CAS 10034-82-9 Na ₂ CrO ₄ ·4H ₂ O F. Wt. 234.03	500 gm	
541	tri-SODIUM CITRATE (Dihydrate), EXTRA PURE CAS 6132-04-3 C ₆ H ₅ Na ₃ O ₇ ·2H ₂ O F. Wt. 294.10	500 gm 5 Kg 25 Kg	
542	tri-SODIUM CITRATE (Dihydrate), AR GRADE CAS 6132-04-3 C ₆ H ₅ Na ₃ O ₇ ·2H ₂ O F. Wt. 294.10	500 gm 5 Kg 25 Kg	
653	tri-SODIUM CITRATE (ANH.), EXTRA PURE CAS 68-04-2 C ₆ H ₅ Na ₃ O ₇ F. Wt. 258.07	500 gm 5 Kg 25 Kg	
2291	SODIUM CITRATE DIBASIC, EXTRA PURE (di-Sodium Hydrogen Citrate) CAS 6132-05-4 C ₆ H ₅ Na ₂ O ₇ F. Wt. 236.09	500 gm 25 Kg	
2283	SODIUM COBALTINITRIDE, EXTRA PURE CAS 13600-98-1 Na ₃ Co(NO ₂) ₆ F. Wt. 403.94	100 gm	
434	SODIUM DICHROMATE (Dihydrate), EXTRA PURE CAS 7789-12-0 Na ₂ Cr ₂ O ₇ ·2H ₂ O F. Wt. 298.00	500 gm 50 Kg	
2286	SODIUM DIETHYLDITHIOCARBAMATE, AR GRADE CAS 20624-25-3 (C ₂ H ₅) ₂ NCS ₂ Na·3H ₂ O F. Wt. 225.31	100 gm 500 gm	
435	SODIUM DIHYDROGEN ORTHOPHOSPHATE (Mono Sodium Phosphate Dihydrate), EXTRA PURE CAS 13472-35-0 NaH ₂ PO ₄ ·2H ₂ O F. Wt. 156.01	500 gm 5 Kg 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
763	SODIUM DIHYDROGEN ORTHOPHOSPHATE (Mono Sodium Phosphate Dihydrate), AR GRADE CAS 13472-35-0 NaH ₂ PO ₄ ·2H ₂ O F. Wt. 156.01	500 gm 5 Kg 50 Kg	
688	SODIUM DIPHENYLAMINE SULPHONATE, AR GRADE CAS 6152-67-6 C ₁₂ H ₁₀ NNaO ₃ S F. Wt. 271.27	5 gm 25 gm	
436	SODIUM DITHIONITE, EXTRA PURE CAS 7775-14-6 Na ₂ S ₂ O ₄ F. Wt. 174.11	500 gm	
444	SODIUM DODECYL SULPHATE (Sodium Lauryl Sulphate), EXTRA PURE CAS 151-21-3 C ₁₂ H ₂₅ NaO ₄ S F. Wt. 288.38	500 gm 25 Kg	
445	SODIUM DODECYL SULPHATE, AR GRADE CAS 151-21-3 C ₁₂ H ₂₅ NaO ₄ S F. Wt. 298.38	100 gm 500 gm	
777	SODIUM FLUORIDE, EXTRA PURE CAS 7681-49-4 NaF F. Wt. 41.99	500 gm 50 Kg	
706	SODIUM FLUORIDE, AR GRADE CAS 7681-49-4 NaF F. Wt. 41.99	500 gm	
2289	SODIUM FORMATE, EXTRA PURE CAS 141-53-7 CHNaO ₂ F. Wt. 68.01	500 gm	
4105	SODIUM FORMATE, AR GRADE CAS 141-53-7 CHNaO ₂ F. Wt. 68.01	500 gm	
941	SODIUM GLUCONATE (D-Gluconic Acid Sodium Salt) CAS 527-07-1 C ₆ H ₁₁ NaO ₇ F. Wt. 218.14	500 gm 50 Kg	
943	SODIUM β-GLYCEROPHOSPHATE, AR GRADE CAS 819-83-0 C ₃ H ₅ O ₆ P ₂ Na F. Wt. 306.11	25 gm 100 gm 500 gm	
944	SODIUM HEXAMETAPHOSPHATE, 62% (Sodium Polyphosphate, Graham's Salts), EXTRA PURE CAS 68915-31-1 (NaPO ₃) ₁₂₋₁₃ NaO ₂ F. Wt. 611.78	500 gm 5 Kg 50 Kg	
429	SODIUM HYDROGEN CARBONATE, EXTRA PURE CAS 144-55-8 NaHCO ₃ F. Wt. 84.01	500 gm 5 Kg 50 Kg	
430	SODIUM HYDROGEN CARBONATE, AR GRADE CAS 144-55-8 NaHCO ₃ F. Wt. 84.01	500 gm 5 Kg 50 Kg	
2984	SODIUM HYDROGEN CARBONATE 'SPL' for R.O Water CAS 144-55-8 NaHCO ₃ F. Wt. 84.01	500 gm 5 Kg 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
2291	di-SODIUM HYDROGEN CITRATE, EXTRA PURE CAS 6132-05-4 $C_6H_5Na_2O_7$ F. Wt. 236.09	500 gm 25 Kg	
460	di-SODIUM HYDROGEN ORTHOPHOSPHATE (Dihydrate), EXTRA PURE CAS 10028-24-7 $Na_2HPO_4 \cdot 2H_2O$ F. Wt. 177.99	500 gm 5 Kg 50 Kg	
658	di-SODIUM HYDROGEN ORTHOPHOSPHATE (Dihydrate), AR GRADE CAS 10028-24-7 $Na_2HPO_4 \cdot 2H_2O$ F. Wt. 177.99	500 gm 5 Kg 50 Kg	
459	di-SODIUM HYDROGEN ORTHOPHOSPHATE (ANH.), EXTRA PURE CAS 7558-79-4 Na_2HPO_4 F. Wt. 141.96	500 gm 5 Kg 50 Kg	
217	di-SODIUM HYDROGEN ORTHOPHOSPHATE (ANH.), AR GRADE CAS 7558-79-4 Na_2HPO_4 F. Wt. 141.96	500 gm 5 Kg 50 Kg	
446	SODIUM MALONATE CAS 141-95-7 $C_3H_2Na_2O_4$ F. Wt. 148.03	25 gm 100 gm	
4221	di-SODIUM OCTABORATE TETRAHYDRATE CAS 12280-03-4 $Na_2B_8O_{13} \cdot 4H_2O$ F. Wt. 412.53	500 gm	
2309	di-SODIUM TARTRATE (Dihydrate) CAS 6106-24-7 $C_4H_4Na_2O_6 \cdot 2H_2O$ F. Wt. 230.08	500 gm 50 Kg	
428	SODIUM HYDROGEN SELENITE for Bacteriology CAS 7782-82-3 $NaHSeO_3$ F. Wt. 150.95	100 gm 1 Kg	
780	SODIUM HYDROGEN SULPHATE (Mono), EXTRA PURE CAS 10034-88-5 $NaHSO_4 \cdot H_2O$ F. Wt. 138.07	500 gm	
2292	SODIUM HYDROGEN SULPHATE (Mono), AR GRADE CAS 10034-88-5 $NaHSO_4 \cdot H_2O$ F. Wt. 138.07	500 gm	
946	SODIUM HYDROGEN SULPHITE, EXTRA PURE (SODIUM BISULPHITE) CAS 7631-90-5 $NaHSO_3$ F. Wt. 104.061	500 gm	
655	SODIUM HYDROGEN SULPHITE, AR GRADE (SODIUM BISULPHITE) CAS 7631-90-5 $NaHSO_3$ F. Wt. 104.061	500 gm	
779	SODIUM HYDROXIDE FLAKES CAS 1310-73-2 NaOH F. Wt. 40.00	500 gm 5 Kg 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
439	SODIUM HYDROXIDE PELLETS, EXTRA PURE CAS 1310-73-2 NaOH F. Wt. 40.00	500 gm 5 Kg 50 Kg	
440	SODIUM HYDROXIDE PELLETS, AR GRADE CAS 1310-73-2 NaOH F. Wt. 40.00	500 gm 5 Kg 50 Kg	
2718	SODIUM HYPOPHOSPHITE (Mono), CAS 10039-56-2 $NaH_2PO_2 \cdot H_2O$ F. Wt. 105.99	500 gm	
544	SODIUM IODIDE, EXTRA PURE CAS 7681-82-5 NaI F. Wt. 149.89	100 gm 250 gm	
2294	SODIUM LACTATE 60% CAS 72-17-3 $C_3H_5NaO_3$ F. Wt. 112.06	500 ml	
444	SODIUM LAURYL SULPHATE, EXTRA PURE (Sodium Dodecyl Sulphate) CAS 151-21-3 $C_{12}H_{25}NaO_4S$ F. Wt. 288.38	500 gm 25 Kg	
445	SODIUM LAURYL SULPHATE 99%, AR GRADE CAS 151-21-3 $NaC_{12}H_{25}SO_4$ F. Wt. 288.38	100 gm 500 gm	
447	SODIUM METABISULPHITE, EXTRA PURE (Sodium Pyrosulfite) CAS 7681-57-4 $Na_2S_2O_5$ F. Wt. 190.11	500 gm 5 Kg 50 Kg	
448	SODIUM METABISULPHITE, AR GRADE CAS 7681-57-4 $Na_2S_2O_5$ F. Wt. 190.11	500 gm 5 Kg 50 Kg	
2296	SODIUM METABORATE (Tetra) CAS 10555-76-7 $BNaO_2 \cdot 4H_2O$ F. Wt. 137.86	500 gm	
2722	SODIUM METAPERIODATE, EXTRA PURE CAS 7790-28-5 $NaIO_4$ F. Wt. 213.89	100 gm 500 gm 25 Kg	
4106	SODIUM METAPERIODATE, AR GRADE (Sodium Periodate) CAS 7790-28-5 $INaO_4$ F. Wt. 213.89	100 gm 500 gm	
2298	SODIUM METAVANADATE (ANH.), AR GRADE CAS 13718-26-8 $NaVO_3$ F. Wt. 121.93	100 gm 500 gm	
2299	SODIUM METHOXIDE, EXTRA PURE CAS 124-41-4 $NaOCH_3$ F. Wt. 54.02	500 gm 5 Kg 25 Kg	
449	SODIUM MOLYBDATE (Dihydrate), EXTRA PURE CAS 10102-40-6 $Na_2MoO_4 \cdot 2H_2O$ F. Wt. 241.95	100 gm 500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
656	SODIUM MOLYBDATE (Dihydrate), AR GRADE CAS 10102-40-6 $\text{Na}_2\text{MoO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 241.9	100 gm 500 gm 25 Kg	
450	SODIUM NITRATE, EXTRA PURE CAS 7631-99-4 NaNO_3 F. Wt. 84.99	500 gm 5 Kg 50 Kg	
3430	SODIUM NITRATE, FOOD GRADE CAS 7631-99-4 NaNO_3 F. Wt. 84.99	500 gm 5 Kg 50 Kg	
451	SODIUM NITRATE, AR GRADE CAS 7631-99-4 NaNO_3 F. Wt. 85.0	500 gm 50 Kg	
452	SODIUM NITRITE, EXTRA PURE CAS 7632-00-0 NaNO_2 F. Wt. 68.99	500 gm 5 Kg	
772	SODIUM NITROPRUSSIDE, EXTRA PURE (Sodium Nitroferricyanide dihydrate) CAS 13755-38-9 $\text{Na}_2[\text{Fe}(\text{CN})_5\text{NO}] \cdot 2\text{H}_2\text{O}$ F. Wt. 297.95	100 gm 500 gm 10 Kg	
2720	SODIUM OLEATE (Pure), EXTRA PURE CAS 143-19-1 $\text{C}_{18}\text{H}_{33}\text{NaO}_2$ F. Wt. 304.44	500 gm 50 Kg	
2300	tri-SODIUM ORTHOPHOSPHATE (ANH.), EXTRA PURE CAS 7601-54-9 Na_3PO_4 F. Wt. 163.94	500 gm 50 Kg	
4107	tri-SODIUM ORTHO-PHOSPHATE (Dodecahydrate), EXTRA PURE CAS 10101-89-0 $\text{H}_{24}\text{Na}_3\text{O}_{16}\text{P}$ F. Wt. 380.12	500 gm 5 Kg 50 Kg	
4108	tri-SODIUM ORTHO-PHOSPHATE (Dodecahydrate), AR GRADE CAS 10101-89-0 $\text{H}_{24}\text{Na}_3\text{O}_{16}\text{P}$ F. Wt. 380.12	500 gm 5 Kg 50 Kg	
455	SODIUM OXALATE (di-Sodium Oxalate), EXTRA PURE CAS 62-76-0 $\text{Na}_2\text{C}_2\text{O}_4$ F. Wt. 133.999	500 gm	
456	SODIUM OXALATE, AR GRADE CAS 62-76-0 $\text{C}_2\text{Na}_2\text{O}_4$ F. Wt. 133.999	500 gm	
2721	SODIUM PENTACHLOROPHENATE CAS 131-52-2 $\text{C}_6\text{Cl}_5\text{NaO}$ F. Wt. 288.32	500 gm 50 Kg	
2301	SODIUM PERBORATE (Tetra), EXTRA PURE CAS 10486-00-7 $\text{NaBO}_3 \cdot 4\text{H}_2\text{O}$ F. Wt. 153.86	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
4154	SODIUM PERCARBONATE (Granules) CAS 15630-89-4 $\text{Na}_2\text{CO}_3 \cdot 1.5\text{H}_2\text{O}$ F. Wt. 156.98	1 Kg	
4097	SODIUM PERCARBONATE (Tablet) CAS 15630-89-4 $\text{Na}_2\text{CO}_3 \cdot 1.5\text{H}_2\text{O}$ F. Wt. 157.01	1 Kg	
2304	SODIUM PERSULPHATE, EXTRA PURE CAS 7775-27-1 $\text{Na}_2\text{S}_2\text{O}_8$ F. Wt. 238.03	500 gm 50 Kg	
459	SODIUM PHOSPHATE DIBASIC (ANH.) (Disodium Hydrogen Orthophosphate), EXTRA PURE CAS 7558-79-4 Na_2HPO_4 F. Wt. 141.96	500 gm 5 Kg 50 Kg	
217	SODIUM PHOSPHATE DIBASIC (ANH.) (Disodium Hydrogen Orthophosphate), AR GRADE CAS 7558-79-4 Na_2HPO_4 F. Wt. 141.96	500 gm 5 Kg 50 Kg	
460	SODIUM PHOSPHATE DIBASIC (Dihydrate) (Disodium Hydrogen Orthophosphate), EXTRA PURE CAS 10028-24-7 $\text{Na}_2\text{HPO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 177.99	500 gm 5 Kg 50 Kg	
658	SODIUM PHOSPHATE DIBASIC (Dihydrate) (Disodium Hydrogen Orthophosphate), AR GRADE CAS 10028-24-7 $\text{Na}_2\text{HPO}_4 \cdot 2\text{H}_2\text{O}$ F. Wt. 177.99	500 gm 5 Kg 50 Kg	
692	SODIUM POLYANETHOL SULPHONATE (SPS), AR GRADE CAS 55963-78-5	1 gm 5 gm	
403	SODIUM POTASSIUM TARTRATE (Tetra), EXTRA PURE CAS 6381-59-5 $\text{C}_4\text{H}_4\text{KNaO}_6 \cdot 4\text{H}_2\text{O}$ F. Wt. 282.22	100 gm 500 gm 25 Kg	
404	SODIUM POTASSIUM TARTRATE (Tetra), AR GRADE CAS 6381-59-5 $\text{C}_4\text{H}_4\text{KNaO}_6 \cdot 4\text{H}_2\text{O}$ F. Wt. 282.22	100 gm 500 gm 25 Kg	
955	SODIUM PROPIONATE, EXTRA PURE CAS 137-40-6 $\text{C}_3\text{H}_5\text{NaO}_2$ F. Wt. 96.07	500 gm 50 Kg	
2724	tetra-SODIUM PYROPHOSPHATE (ANH.), (TSPP), EXTRA PURE CAS 7722-88-5 $\text{Na}_4\text{P}_2\text{O}_7$ F. Wt. 265.90	500 gm 50 Kg	
659	SODIUM PYRUVATE, AR GRADE CAS 113-24-6 $\text{C}_3\text{H}_3\text{NaO}_3$ F. Wt. 110.0	100 gm 500 gm	
4110	SODIUM RHODIZONATE, AR GRADE (Sensitivity To BA 1:200000)	1gm 5gm	

CODE	PRODUCT NAME	PACK SIZE	
2727	SODIUM SALICYLATE, EXTRA PURE CAS 54-21-7 $C_7H_5NaO_3$ F. Wt. 160.11	500 gm 25 Kg	
2728	SODIUM SELENATE (ANH.), AR GRADE CAS 13410-01-0 Na_2SeO_4 F. Wt. 188.94	100 gm	
683	SODIUM SELENITE (ANH.) (Culture Media Ingredient) CAS 10102-18-8 Na_2SeO_3 F. Wt. 172.96	25 gm 100 gm 500 gm	
2240	SODIUM SUCCINATE (Hexa), EXTRA PURE CAS 6106-21-4 $C_4H_4Na_2O_4 \cdot 6H_2O$ F. Wt. 270.15	100 gm 500 gm	
701	SODIUM SULPHATE (ANH.), EXTRA PURE CAS 7757-82-6 Na_2SO_4 F. Wt. 142.04	500 gm 5 Kg 50 Kg	
543	SODIUM SULPHATE (ANH.), AR GRADE CAS 7757-82-6 Na_2SO_4 F. Wt. 142.04	500 gm 5 Kg 50 Kg	
661	SODIUM SULPHIDE FLAKES (Iron Free) CAS 1313-84-4 Na_2S F. Wt. 78.05	500 gm	
952	SODIUM SULPHITE (ANH.), EXTRA PURE CAS 7757-83-7 Na_2SO_3 F. Wt. 126.04	500 gm 5 Kg 50 Kg	
662	SODIUM SULPHITE (ANH.), AR GRADE CAS 7757-83-7 Na_2SO_3 F. Wt. 126.04	500 gm 5 Kg 50 Kg	
3260	SODIUM STEARATE (Stearic Acid Sodium Salt) CAS 822-16-2 $C_{18}H_{35}NaO_2$ F. Wt. 306.46	500 gm	
279	SODIUM TETRABORATE (Decahydrate) (BORAX), EXTRA PURE CAS 1303-96-4 $Na_2B_4O_7 \cdot 10H_2O$ F. Wt. 381.37	500 gm 5 Kg 50 Kg	
280	SODIUM TETRABORATE (Decahydrate) (BORAX), AR GRADE CAS 1303-96-4 $Na_2B_4O_7 \cdot 10H_2O$ F. Wt. 381.37	500 gm 5 Kg 50 Kg	
4083	SODIUM TETRAPHENYL BORON, AR GRADE (Kalignost) CAS 143-66-8 $C_{24}H_{20}BNa$ F. Wt. 342.22	10 gm 25 gm	
464	SODIUM THIOCYANATE, EXTRA PURE CAS 540-72-7 $NaSCN$ F. Wt. 81.07	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
2730	SODIUM THIOCYANATE, AR GRADE CAS 540-72-7 $NaSCN$ F. Wt. 81.07	500 gm 50 Kg	
3513	SODIUM THIOLYCOLLATE (for Bacteriology) * CAS 367-51-1 $C_2H_3NaO_2S$ F. Wt. 114.10	100 gm 500 gm 10 Kg	
469	SODIUM THIOSULPHATE (Penta), EXTRA PURE CAS 10102-17-7 $Na_2S_2O_3 \cdot 5H_2O$ F. Wt. 248.18	500 gm 5 Kg 50 Kg	
4111	SODIUM THIOSULPHATE (Penta), AR GRADE CAS 10102-17-7 $Na_2S_2O_3 \cdot 5H_2O$ F. Wt. 248.18	500 gm	
2235	SODIUM THIOSULPHATE (ANH.), EXTRA PURE CAS 7772-98-2 $Na_2S_2O_3$ F. Wt. 158.11	500 gm 50 Kg	
663	SODIUM THIOSULPHATE (ANH.), AR GRADE CAS 7772-98-2 $Na_2S_2O_3$ F. Wt. 158.11	500 gm 50 Kg	
2895	SODIUM TRIPOLYPHOSPHATE (STPP) (ANH.) CAS 7758-29-4 $Na_5S_3O_{10}$ F. Wt. 367.86	500 gm 50 Kg	
953	SODIUM TUNGSTATE (Dihydrate), EXTRA PURE CAS 10213-10-2 $Na_2O_4W \cdot 2H_2O$ F. Wt. 329.86	100 gm 500 gm 10 Kg	
1623	SOLOCHROME BLACK-T, EXTRA PURE (ERICHROME BLACK-T) CAS 1787-61-7 $C_{20}H_{12}N_3NaO_7S$ F. Wt. 461.38	25 gm 100 gm	
1855	SOLOCHROME BLACK-T, AR GRADE CAS 1787-61-7 C.I. 14645 $C_{20}H_{12}N_3NaO_7S$ F. Wt. 461.38	25 gm 100 gm	
1856	SOLOCHROME CYANINE-R, AR GRADE CAS 3564-18-9 C.I. 43820 $C_{23}H_{15}Na_3O_9S$ F. Wt. 536.40	5 gm 25 gm	
2995	SOLOCHROME DARK BLUE, (Calcon) CAS 2538-85-4 C.I. 15705 $C_{20}H_{13}N_2NaO_5S$ F. Wt. 416.38	5 gm 25 gm 100 gm	
2732	SORBITIC ACID, EXTRA PURE * CAS 110-44-1 $C_6H_{14}O_6$ F. Wt. 112.12	500 gm	
465	D-SORBITOL POWDER (D-Glucitol) CAS 50-70-4 $C_6H_{14}O_6$ F. Wt. 182.17	100 gm 500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
3065	L-SORBOSE CAS 87-79-6 $C_6H_{12}O_6$ F. Wt. 180.16	25 gm 100 gm	
466	SOYABEAN MEAL (Defatted) (Non GMO Powder)	500 gm 25 Kg	
467	SOYA LECITHIN POWDER * (L- α -phosphatidyl Choline 22%) CAS 8002-43-5 $C_{36}H_{72}NO_8P$ F. Wt. 677.93	100 gm 500 gm	
4351	SOYA LECITHIN LIQUID CAS 8002-43-5 $C_{36}H_{72}NO_8P$ F. Wt. 677.93	500 ml	
4205	SOYA PEPTONE (REGULAR GRADE) (for general purpose for bactriological work) CAS 91079-38-8	500 gm 50 Kg	
2733	SPADNS, AR GRADE CAS 23647-14-5 $C_{16}H_{19}N_2O_{11}S_3Na_3$ F. Wt. 570.41	1 gm 5 gm	
4329	SPAN 20 CAS 1338-39-2 $C_{18}H_{34}O_6$ F. Wt. 346.46	500 ml	
4330	SPAN 80 CAS 1338-43-8 $C_{24}H_{44}O_6$ F. Wt. 428.6	500 ml	
3436	SPERMIDINE * CAS 124-20-9 $C_7H_{19}N_3$ F. Wt. 145.25	1 gm 5 gm	
2317	STANNOUS CHLORIDE (Dihydrate), EXTRA PURE CAS 10025-69-1 $SnCl_2 \cdot 2H_2O$ F. Wt. 225.63	100 gm 250 gm	
2318	STANNOUS CHLORIDE (Dihydrate), AR GRADE CAS 10025-69-1 $SnCl_2 \cdot 2H_2O$ F. Wt. 225.63	100 gm 250 gm	
3006	STANNOUS FLOURIDE (ANH.) CAS 7783-47-3 SnF_2 F. Wt. 156.69	500 gm	
2734	STANNOUS OXIDE, EXTRA PURE CAS 21651-19-4 SnO F. Wt. 134.71	500 gm	
1310	STARCH CORN (Maize) (Insoluble) CAS 9005-84-9 $(C_6H_{10}O_5)_n$ F. Wt. 162.14	500 gm 25 Kg	
1309	STARCH POTATO (Insoluble) CAS 9005-25-8 $(C_6H_{10}O_5)_n$ F. Wt. 162.14	500 gm	
1311	STARCH SOLUBLE, EXTRA PURE CAS 9000-84-9 $(C_6H_{10}O_5)_n$	500 gm 25 Kg	
1318	STARCH SOLUBLE, AR GRADE CAS 9004-84-9 $(C_6H_{10}O_5)_n$	500 gm	
2321	STEARIC ACID, FLAKES (Carboxylic Acid) CAS 57-11-4 $C_{18}H_{36}O_2$ F. Wt. 284.48	500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
2322	STRONTIUM CARBONATE CAS 1633-05-2 $SrCO_3$ F. Wt. 147.63	500 gm 50 Kg	
2323	STRONTIUM CHLORIDE (Hexa), EXTRA PURE CAS 10025-70-4	500 gm 50 Kg	
2325	STRONTIUM NITRATE (ANH.), EXTRA PURE CAS 10042-76-9 $Sr(NO_3)_2$ F. Wt. 211.63	500 gm 50 Kg	
2735	STRONTIUM SULPHATE (ANH.) CAS 7759-02-6 $SrSO_4$ F. Wt. 183.68	500 gm	
1016	SUCCINIC ACID, EXTRA PURE CAS 110-15-6 $C_4H_6O_4$ F. Wt. 118.09	500 gm 25 Kg	
2328	SUCCINIC ANHYDRIDE CAS 108-30-5 $C_4H_4O_3$ F. Wt. 100.07	100 gm 500 gm 25 Kg	
4702	SUCROSE, ACS GRADE CAS 57-50-1 $C_{12}H_{22}O_{11}$ F. Wt. 342.30	500 gm	
1307	SUCROSE, EXTRA PURE CAS 57-50-1 $C_{12}H_{22}O_{11}$ F. Wt. 342.30	500 gm 5 Kg 50 Kg	
1308	SUCROSE, AR GRADE (Bacto Grade, Reducing Sugar Absent) CAS 57-50-1 $C_{12}H_{22}O_{11}$ F. Wt. 342.30	500 gm 5 Kg 50 Kg	
TC 023	SUCROSE for Tissue Culture & Molecular Biology CAS 57-50-1 $C_{12}H_{22}O_{11}$ F. Wt. 342.30	500 gm 5 Kg	
3503	SUCRALOSE (600 Times Sweeter than Sugar) CAS 56038-13-2 $C_{12}H_{22}O_{11}$ F. Wt. 397.64	25 gm 100 gm	
2737	SUDAN (I) CAS 842-07-9 $C_{16}H_{12}N_2O$ F. Wt. 248.28	25 gm 100 gm	
1692	SUDAN (III) CAS 85-86-9 C.I. 26100 $C_{22}H_{16}N_4O$ F. Wt. 352.39	25 gm 100 gm	
1693	SUDAN (IV) (Oil Red IV) (Seavlet-R) CAS 85-83-6 C.I. 26105 $C_{22}H_{20}N_4O$ F. Wt. 380.44	25 gm 100 gm	
1647	SUDAN BLACK-B (for Microscopy) CAS 4197-25-5 C.I. 26150 $C_{29}H_{24}N_6$ F. Wt. 456.54	10 gm 25 gm	
2333	SULPHAMIC ACID (Amido Sulfonic Acid), EXTRA PURE CAS 5329-14-6 NH_2SO_3H F. Wt. 97.10	500 gm 50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
2407	SULPHANILAMIDE, EXTRA PURE CAS 63-74-1 $C_6H_5N_2O_2S$ F. Wt. 172.20	500 gm	
2741	SULPHANILAMIDE, AR GRADE CAS 63-74-1 $C_6H_5N_2O_2S$ F. Wt. 172.20	100 gm 500 gm	
468	SULPHANILIC ACID, EXTRA PURE CAS 121-57-3 $C_6H_5NO_3S$ F. Wt. 173.19	100 gm 500 gm 50 Kg	
470	SULPHANILIC ACID, AR GRADE (4-Aminobenzene Sulphonic Acid) CAS 121-57-3 $C_6H_5NO_3S$ F. Wt. 173.19	100 gm 500 gm 50 Kg	
776	5-SULPHOSALICYLIC ACID DIHYDRATE, EXTRA PURE CAS 5965-83-3 $C_7H_6O_6S \cdot 2H_2O$ F. Wt. 254.22	100 gm 500 gm 50 Kg	
2412	SULPHUR POWDER, EXTRA PURE CAS 7704-34-9 S F. Wt. 32.07	500 gm	
665	TABOLENE GL (With Dispenser) for Glass & Computer Equipment Cleaning	500 ml 5 Ltr 50 Ltr	
666	TABOLENE GLW (DETERGENT TABOLENE) (With Neutral pH Phosphate Free Detergent) for Glassware Cleaning	500 ml 5 Ltr	
667	TABOLENE GP General Purpose Detergent Lab Floor Cleaning, Like Teepol	500 ml 5 Ltr 50 Ltr	
2742	TALCUM POWDER (400 mesh) (Magnesium Silicate) CAS 14807-96-6 $3MgO \cdot 4SiO_2 \cdot H_2O$ F. Wt. 379.29	500 gm	
3226	TALCUM POWDER (700 mesh) CAS 14807-96-6 $3MgO \cdot 4SiO_2 \cdot H_2O$ F. Wt. 379.29	500 gm	
2782	TANNIC ACID, (Gallotannin), EXTRA PURE CAS 1401-55-4 $C_{76}H_{52}O_{46}$ F. Wt. 1701.20	100 gm 500 gm 25 Kg	
4112	TANNIC ACID, AR GRADE CAS 1401-55-4 $C_{76}H_{52}O_{46}$ F. Wt. 1701.2	100 gm 250 gm	
4261	TAPS BUFFER (N-[TRIS(HYDROXYMETHYL) METHYL] -3-AMINOPROPANESULPHONIC ACID) CAS 29915-38-6 $C_7H_{17}NO_6S$ F. Wt. 243.28	25 gm 100 gm	
472	L (+) TARTARIC ACID, EXTRA PURE CAS 87-69-4 $C_4H_6O_6$ F. Wt. 150.09	100 gm 500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
4262	L-(+)-TARTARIC ACID, AR GRADE CAS 87-69-4 $C_4H_6O_6$ F. Wt. 150.09	500 gm	
474	DL-TARTARIC ACID (SYNTHETIC) CAS 133-37-9 $C_4H_6O_6$ F. Wt. 150.09	500 gm 50 Kg	
1648	TARTRAZINE, (for Microscopy) CAS 1934-21-0 C.I. 19140 $C_{16}H_9N_4Na_3O_9S_2$ F. Wt. 534.3	25 gm 100 gm	
956	TAURINE (2-Aminoethanesulphonic Acid) CAS 107-35-7 $C_2H_7NO_3S$ F. Wt. 125.15	25 gm 100 gm	
4113	TBTU (2-(1H-Benzotriazole-1-yl)-1,1,3,3 Tetramethyluronium Tetrafluoroborate) CAS 125700-67-6 $C_{11}H_{16}N_5O \cdot BF_4$ F. Wt. 321.08	5 gm 25 gm	
4263	TEICOPLANIN * CAS 61036-62-2	1 gm	
4264	TELLURIUM DIOXIDE CAS 7446-07-3 TeO_2 F. Wt. 159.60	50 gm	
4265	TELLURIUM (IV) CHLORIDE CAS 10026-07-0 $TeCl_4$ F. Wt. 269.41	25 gm	
775	TEMED, AR GRADE (N,N,N,N-Tetra-Methylethylene Diamine) CAS 110-18-9 $C_6H_{16}N_2$ F. Wt. 116.21	100 ml 500 ml	
4266	TERGITOL-7 CAS 68131-40-8	100gm 500gm	
2743	TERPINEOL (Mixture of Isomers), PHARMA GRADE CAS 8006-39-1 $C_{10}H_{18}O$ F. Wt. 154.25	500 ml 2.5 Ltr	
4114	TETRABUTYL AMMONIUM BROMIDE, AR GRADE CAS 1643-19-2	500 gm	
4117	TETRABUTYL AMMONIUM FLUORIDE CAS 87749-50-6 $C_{16}H_{36}FN$ F. Wt. 261.46	25 gm 100 gm	
4118	TETRABUTYL AMMONIUM IODIDE, AR GRADE CAS 311-28-4 $C_{16}H_{36}IN$ F. Wt. 369.37	25 gm 100 gm	
4119	TETRABUTYL-R (for HPLC) CAS 32503-27-8 $(CH_3CH_2CH_2CH_2)_4N(HSO_4)$ F. Wt. 339.53	100 gm	
144	Tert-BUTYL HYDROQUINONE (TBHQ) CAS 1948-33-0 $C_{10}H_{14}O_2$ F. Wt. 166.22	100 gm 500 gm 25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
775	N,N,N,N-TETRA METHYLETHYLENEDIAMINE (TEMED), AR GRADE CAS 110-18-9 $C_6H_{16}N_4$ F. Wt. 116.21	100 ml 500 ml	
4267	N,N,N,N'-TETRA ACETYLETHYLENEDIAMINE (TAED) CAS 10543-57-4 $C_{10}H_{16}N_2O_4$ F. Wt. 228.25	50 gm 100 gm	
654	N,N,N,N-TETRAMETHYL-p-PHENYLENEDIAMINE DIHYDROCHLORIDE, AR GRADE CAS 637-01-4 $C_{10}H_{16}N_2 \cdot 2HCl$ F. Wt. 237.17	5 gm 25 gm	
4268	TETRABUTYL AMMONIUM CHLORIDE CAS 1112-67-0 $C_{16}H_{36}ClN$ F. Wt. 261.46	25 gm	
821	TETRAZOLIUM BLUE CHLORIDE, AR GRADE CAS 1871-22-3 $C_{40}H_{32}Cl_2N_8O_8$ F. Wt. 731.67	1 gm 5 gm	
475	TETRAZOLIUM RED SALT (TTC), AR GRADE (2,3,5-Triphenyl Tetrazolium Chloride) CAS 298-96-4 $C_{19}H_{15}ClN_4$ F. Wt. 334.81	10 gm 25 gm	
4269	TETRAZOLIUM VIOLET CAS 1719-71-7 $C_{23}H_{17}ClN_4$ F. Wt. 384.86	1 gm	
4270	THIABENDAZOLE CAS 148-79-8 $C_{10}H_7N_3S$ F. Wt. 201.25	50 gm 100 gm	
1702	THIAMINE HYDROCHLORIDE * (Vitamin B1 Hydrochloride) CAS 67-03-8 $C_{12}H_{17}ClN_4OS \cdot HCl$ F. Wt. 337.27	25 gm 100 gm 500 gm	
3584	THIDIAZURON * CAS 51707-55-2 $C_9H_8N_4OS$ F. Wt. 220.25	250 mg 1 gm	
4272	THIOACETAMIDE (Ethanethioamide), AR GRADE CAS 62-55-5 C_2H_5NS F. Wt. 75.13	100gm	
959	THIOACETAMIDE (Ethanethioamide), EXTRA PURE * CAS 62-55-5 C_2H_5NS F. Wt. 75.13	100 gm	
4273	THIOACETIC ACID CAS 507-09-5 C_2H_3OS F. Wt. 76.12	100gm	
2750	2-THIOBARBITURIC ACID, AR GRADE CAS 504-17-6 $C_4H_4N_2O_2S$ F. Wt. 144.15	25 gm 100 gm	
4274	THIOGLYCOLIC ACID, 80% SOLUTION IN WATER CAS 68-11-1 $C_2H_2O_2S$ F. Wt. 92.12	100 ml 500 ml	

CODE	PRODUCT NAME	PACK SIZE	
2341	1-THIOGLYCEROL 98% (Mono Thioglycerol) CAS 96-27-5 $C_3H_8O_2S$ F. Wt. 108.15	100 ml 500 ml	
4275	2-THIOURACIL (4-HYDROXY-2-MERCAPTOPYRIMIDINE) CAS 141-90-2 $C_4H_4N_2OS$ F. Wt. 128.15	25 gm 100 gm	
3513	THIOGLYCOLLIC ACID SODIUM SALT (for Bacteriology) CAS 367-51-1 $C_2H_3NaO_2S$ F. Wt. 114.1	100 gm 500 gm 10 Kg	
478	THIOMERSAL, EXTRA PURE CAS 54-64-8 $C_9H_5HgNaO_2S$ F. Wt. 404.81	25 gm 100 gm	
2751	THIOSEMICARBAZIDE, AR GRADE CAS 79-19-6 CH_3N_3S F. Wt. 91.14	25 gm 100 gm	
479	THIOUREA, EXTRA PURE CAS 62-56-6 CH_4N_2S F. Wt. 76.12	500 gm 50 Kg	
4120	THIOUREA, AR GRADE CAS 62-56-6 $C_7H_7ClN_2S$ F. Wt. 186.66	500 gm	
4276	THORIN INDICATOR OCTAHYDRATE, AR GRADE CAS 3688-92-4 $C_{16}H_{11}AsN_2Na_2O_{10}S_2$ F. Wt. 576.3	5 gm 10 gm	
1131	L-THREONINE, 99%+ Crystalline CAS 72-19-5 $C_4H_9NO_3$ F. Wt. 119.12	5 gm 25 gm 500 gm	
4277	THROMBIN From Bovine plasma CAS 9002-04-4	1 gm	
961	THYMIDINE (2-Deoxyriboside Thymidine) CAS 50-89-5 $C_{10}H_{14}N_2O_5$ F. Wt. 242.22	1 gm 5 gm	
963	THYMINE (5-Methyluracil) CAS 65-71-4 $C_5H_4N_2O_2$ F. Wt. 126.11	5 gm 25 gm	
1668	THYMOL BLUE CAS 76-61-9 $C_{27}H_{30}O_5S$ F. Wt. 466.59	5 gm 25 gm	
4278	THYMOL BLUE SODIUM SALT CAS 62625-21-2 $C_{27}H_{29}NaO_5S$ F. Wt. 488.57	5 gm	
2343	THYMOL CRYSTAL, EXTRA PURE CAS 89-83-8 $C_{10}H_{14}O$ F. Wt. 150.22	100 gm 500 gm 50 Kg	
4123	THYMOLPHTHALEIN, AR GRADE CAS 125-20-2 $C_{28}H_{30}O_4$ F. Wt. 430.53	5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
964	THYMOLPHTHALEIN COMPLEXONE for E.D.T.A. Titration CAS 1913-93-5 $C_{38}H_{38}N_2Na_4O_{12}$ F. Wt. 808.71	1 gm 5 gm	
1667	THYMOLPHTHALEIN (pH Indicator) CAS 125-20-2 $C_{20}H_{30}O_4$ F. Wt. 430.53	5 gm 25 gm	
2752	THYMOL VIOLET (Indicator) CAS 7512-38-1 $C_{40}H_{39}N_2NaO_5S$ F. Wt. 746.8	5 gm 25 gm	
2988	TIN (METAL) GRANULATED CAS 7440-31-5 Sn F. Wt. 118.71	100 gm 500 gm	
760	TIN (METAL) POWDER CAS 7440-31-5 Sn F. Wt. 118.71	100 gm 500 gm	
2317	TIN (II) CHLORIDE DIHYDRATE, EXTRA PURE CAS 10025-69-1 $SnCl_2 \cdot 2H_2O$ F. Wt. 225.63	100 gm 250 gm	
2318	TIN (II) CHLORIDE DIHYDRATE, AR GRADE CAS 10025-69-1 $SnCl_2 \cdot 2H_2O$ F. Wt. 225.63	100 gm 250 gm	
4281	TIN (IV) CHLORIDE PENTAHYDRATE CAS 10026-06-9 $SnCl_4 \cdot 5H_2O$ F. Wt. 350.60	500 gm	
4282	TIN(IV) OXIDE, AR GRADE CAS 18282-10-5 SnO_2 F. Wt. 150.71	25 gm	
4283	TIRON, AR GRADE (4,5-Dihydroxy-1,3-benzenedisulfonic Acid Disodium Salt) CAS 149-45-1 $C_6H_4O_8S_2Na_2$ F. Wt. 314.2	10 gm 25 gm	
4284	TITANIUM (IV) BUTOXIDE CAS 5593-70-4 $C_{16}H_{36}O_4Ti$ F. Wt. 340.30	500 gm	
482	TITANIUM DIOXIDE, EXTRA PURE CAS 13463-67-7 TiO_2 F. Wt. 79.87	500 gm 5 Kg 25 Kg	
4285	TITANIUM METAL POWDER CAS 7440-32-6 Ti F. Wt. 47.87	100 gm	
483	TITANIUM DIOXIDE, AR GRADE CAS 13463-67-7 TiO_2 F. Wt. 79.87	100 gm	
1650	TITAN YELLOW CAS 1829-00-1 C.I. 19540 $C_{28}H_{19}N_5Na_2O_6S_4$ F. Wt. 695.71	10 gm 25 gm	
1666	O-TOLIDINE, AR GRADE (Reagent for Halogen & Gold) CAS 119-93-7 $C_{14}H_{16}N_2$ F. Wt. 212.29	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
2345	O-TOLIDINE DIHYDROCHLORIDE, AR GRADE CAS 612-82-8 $C_{14}H_{16}N_2 \cdot 2HCl$ F. Wt. 285.22	25 gm 100 gm	
4286	TOLUENE-3,4-DITHIOL, AR GRADE CAS 496-74-2 $C_7H_6S_2$ F. Wt. 156.27	5 gm	
4287	p-TOLUENESULPHONIC ACID MONOHYDRATE, AR GRADE CAS 6192-52-5 $C_7H_8SO_3 \cdot H_2O$ F. Wt. 190.22	100 gm 500 gm	
TR 009	O-TOLIDINE REAGENT (Free Chloride)	500 ml	
4288	p-TOLUIDINE, EXTRA PURE CAS 106-49-0 C_7H_9N F. Wt. 107.15	500 gm	
4290	TOMATO POWDER	250 gm	
1676	TOLUIDINE BLUE CAS 92-31-9 C.I. 52040 $(C_{15}H_{16}C_3N_3S)_2ZnCl_2$ F. Wt. 305.83	25 gm 100 gm	
TC 024	D-TREHALOSE DIHYDRATE CAS 6138-23-4 $C_{12}H_{22}O_{11} \cdot 2H_2O$ F. Wt. 378.33	5 gm 25 gm	
4703	D-TREHALOSE DIHYDRATE, ACS GRADE CAS 6138-23-4 $C_{12}H_{22}O_{11} \cdot 2H_2O$ F. Wt. 378.33	5 gm 25 gm	
4364	TRIACONTANOL (TRIA), 98% CAS 593-50-0 $C_{30}H_{62}O$ F. Wt. 438.81	500 mg 1 gm	
4124	1,2,4-TRIAZOLE (for Synthesis) CAS 288-88-0 $C_2H_3N_3$ F. Wt. 69.07	100 gm 500 gm	
851	TRIBUTYRIN (Glycerol Tributryrate) CAS 60-01-5 $C_{15}H_{26}O_6$ F. Wt. 302.36	100 ml 500 ml	
484	TRICHLORO ACETIC ACID, EXTRA PURE CAS 76-03-9 $C_2HCl_3O_2$ F. Wt. 163.39	100 gm 500 gm 25 Kg	
4292	2,4,6-TRICHLOROPHENOL CAS 88-06-2 $C_6H_3Cl_3O$ F. Wt. 197.45	500 gm	
3587	2,4,5-TRICHLOROPHENOXY ACETIC ACID CAS 93-76-5 $C_6H_2Cl_3O_3$ F. Wt. 255.49	25 gm 100 gm	
966	TRICINE (N-[Tris(Hydroxymethyl Methyl) Glycine]) CAS 5704-04-1 $C_6H_{13}NO_5$ F. Wt. 179.17	25 gm 100 gm	
4293	TRICLOSAN CAS 3380-34-5 $C_{12}H_7Cl_3O_2$ F. Wt. 289.54	100 gm	
2354	TRIETHANOLAMINE 85% CAS 102-71-6 $C_6H_{15}NO_3$ F. Wt. 149.19	500 ml 5 Ltr 50 Ltr	

CODE	PRODUCT NAME	PACK SIZE	
2356	TRIETHYLAMINE 99% CAS 121-44-8 $C_6H_{15}NO_3$ F. Wt. 101.19	500 ml 2.5 Ltr 25 Ltr	
4294	2,3,5-TRIIODOBENZOIC ACID (TIBA)	5 gm	
4125	TRIETHYLBENZYLAMMONIUM CHLORIDE, EXTRA PURE CAS 56-37-1 $C_{13}H_{22}ClN$ F. Wt. 227.77	500 gm	
2031	TRIMETHOPRIM CAS 738-70-5 $C_{14}H_{18}N_4O_3$ F. Wt. 290.30	5 gm	
4101	2,4,6-TRINITROPHENOL CAS 88-89-1 $C_6H_3N_3O_7$ F. Wt. 229.1	500 gm	
4129	TRIMETHYL SULPHOXONIUM IODIDE, Min 99% (for Synthesis) CAS 1774-47-6 C_3H_9IOS F. Wt. 220.07	100 gm 500 gm	
486	TRIS, EXTRA PURE (Tris Hydroxymethyl Aminomethane) CAS 77-86-1 $C_{14}H_{18}N_4O_3$ F. Wt. 121.14	100 gm 500 gm 25 Kg	
487	TRIS, 99.9% + PURITY (Biological Buffer), AR GRADE CAS 77-86-1 $C_4H_{11}NO_3$ F. Wt. 121.14	100 gm 500 gm 25 Kg	
4130	TRIS ULTRAPURE (for HPLC) CAS 77-86-1 $C_4H_{11}NO_3$ F. Wt. 121.14	500 gm	
672	TRIS-HYDROCHLORIDE, EXTRA PURE (Tris-Hydroxymethyl Aninomethane HCl) CAS 1185-53-1 $C_4H_{11}NO_3.HCl$ F. Wt. 157.60	100 gm 500 gm 25 Kg	
488	TRIS-HYDROCHLORIDE, AR GRADE CAS 1185-53-1 $C_4H_{11}NO_3.HCl$ F. Wt. 157.60	100 gm 500 gm 25 Kg	
4297	TRIS (HYDROXYMETHYL) AMINOMETHANE SULPHATE CAS 6992-38-7 $C_4H_{11}NO_3.0.5H_2SO$ F. Wt. 170.17	25 gm 100 gm	
4261	N-[TRIS (HYDROXYMETHYL) METHYL] -3-AMINOPROPANESULPHONIC ACID (TAPS BUFFER) CAS 29915-38-6 $C_7H_{17}NO_6S$ F. Wt. 243.28	25 gm 100 gm	

CODE	PRODUCT NAME	PACK SIZE	
669	TRITON X 100® CAS 9002-93-1 $C_{34}H_{62}O_{11}$ F. Wt. 646.86	100 ml 500 ml	
4298	TROPAEOLIN-O CAS 547-57-9 $C_{12}H_9N_2NaO_5S$ F. Wt. 316.27	25 gm	
1651	TRYPAN BLUE CAS 72-57-1 C.I. 23850 $C_{34}H_{28}N_6O_{14}S_4$ F. Wt. 872.88	25 gm 100 gm	
4300	TRYPSIN INHIBITOR * (From Soya Bean) CAS 9035-81-8	500 mg	
489	TRYPSIN 2000 U/G (0.2 anson U/G) * CAS 9002-07-7 $C_6H_{15}O_{12}P_3$ F. Wt. 372.10	100 gm 500 gm	
490	TRYPSIN 1: 250 * CAS 9002-07-7 $C_6H_{15}O_{12}P_3$ F. Wt. 372.10	25 gm 100 gm	
4301	TRYPSIN 1: 250 * (γ-Irradiated) CAS 9002-07-7 $C_6H_{15}O_{12}P_3$ F. Wt. 372.10	25 gm 100 gm 500 gm	
3514	TRYPTONE (Regular Grade) General Purpose for Bacteriology	500 gm	
1133	DL-TRYPTOPHAN, 99% + PURITY, CRYSTALLINE CAS 54-12-6 $C_{11}H_{12}N_2O_2$ F. Wt. 204.22	5 gm 25 gm	
1134	L-TRYPTOPHAN, 99% + PURITY, CRYSTALLINE CAS 73-22-3 $C_{11}H_{12}N_2O_2$ F. Wt. 204.23	5 gm 25 gm 500 gm	
4131	TUNGSTIC ACID 99%, AR GRADE CAS 7783-03-1 H_2O_7W F. Wt. 249.85	100 gm	
4302	DODECA-TUNGSTOPHOSPHORIC ACID HYDRATE, AR GRADE CAS 12501-23-4 $H_3PO_4.12WO_3.XH_2O$ F. Wt. 2880.05	100 gm 500 gm	
981	TUNGSTOPHOSPHORIC ACID, EXTRA PURE CAS 12501-23-4 $H_3PO_4.12WO_3.XH_2O$ F. Wt. 2880.05	100 gm 500 gm	

CODE	PRODUCT NAME	PACK SIZE	
4156	TWEEN 20, EXTRA PURE CAS 9005-64-5 $C_{58}H_{113}O_{26}$ F. Wt. 1226.5	500 ml 25 Ltr	
670	TWEEN 20 (Regd. Trademark ICI UK) CAS 9005-64-5 $C_{58}H_{113}O_{26}$ F. Wt. 1226.5	500 ml 25 Ltr	
4155	TWEEN 80, EXTRA PURE CAS 9005-65-6 $C_{32}H_{60}O_{10}$ F. Wt. 604.8	100 ml 500 ml	
673	TWEEN 80 (Regd. Trademark ICI UK) [®] CAS 9005-65-6 $C_{32}H_{60}O_{10}$ F. Wt. 604.8	100 ml 500 ml 25 Ltr	
4331	TWEEN 85 CAS 9005-70-3 CH_4 F. Wt. 16.04246	500 ml	
3234	DL-TYROSINE CAS 556-03-6 $C_9H_9NO_3$ F. Wt. 181.19	5 gm	
1135	L-TYROSINE, 99% + Crystalline CAS 60-18-4 $C_9H_9NO_3$ F. Wt. 181.19	25 gm 100 gm 500 gm	
4303	UNIVERSAL INDICATOR	100 ml 500 ml	
1663	UNIVERSAL INDICATOR POWDER * (pH 2 to 10)	25 gm 100 gm	
674	URACIL CAS 66-22-8 $C_4H_4N_2O_2$ F. Wt. 112.09	5 gm 25 gm	
4304	UREA (Meets USP/EP/JP/BP Testing specifications) CAS 57-13-6 N_2H_4CO F. Wt. 60.06	500 gm 1 Kg 5 Kg	
491	UREA, EXTRA PURE CAS 57-13-6 CH_4N_2O F. Wt. 60.06	500 gm 5 Kg 50 Kg	
492	UREA, AR GRADE CAS 57-13-6 CH_4N_2O F. Wt. 60.06	500 gm 5 Kg 50 Kg	
4305	UREA PHOSPHATE CAS 4861-19-12 $CN_2H_7PO_5$ F. Wt. 158.05	500 gm	
1933	UREASE (Jackbean Meal) CAS 9002-13-5	100 gm 500 gm 25 Kg	
2768	URIC ACID, AR GRADE CAS 69-93-2 $C_5H_4N_4O_3$ F. Wt. 168.0	25 gm 100 gm	
3451	URIDINE 99% (for Biochemistry) CAS 58-96-8 $C_9H_{12}N_2O_6$ F. Wt. 244.2	1 gm 5 gm 25 gm	
642	VACUUM GREASE SILICONE (High Vacuum Grease)	500 gm	
4306	D-VALINE CAS 640-68-6 $C_6H_{11}NO_2$ F. Wt. 117.15	1 gm 5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
4307	L-VALINE (From non-animal source) CAS 72-18-4 $C_6H_{11}NO_2$ F. Wt. 117.15	25 gm 100 gm 1 Kg	
1136	DL-VALINE CAS 516-06-3 $C_6H_{11}NO_2$ F. Wt. 117.15	25 gm 100 gm	
1137	L-VALINE, 99% + Purity, Crystalline CAS 72-18-4 $C_6H_{11}NO_2$ F. Wt. 117.15	25 gm 100 gm	
4308	VANADIUM (III) CHLORIDE CAS 7718-98-1 VCl_3 F. Wt. 157.3	10 gm	
2769	VANADIUM PENTOXIDE, EXTRA PURE CAS 1314-62-1 V_2O_5 F. Wt. 181.88	100 gm 500 gm	
4309	VANADYL SULPHATE HYDRATE CAS 27774-13-6 $VOSO_4 \cdot xH_2O$ F. Wt. 163.00	100 gm	
2770	VANILLIN (4-Hydroxy-3-Methoxy Benzaldehyde) CAS 121-33-5 $C_8H_8O_3$ F. Wt. 152.15	100 gm 500 gm	
1711	VITAMIN A ACETATE (Retinol acetate) * CAS 127-47-9 $C_{22}H_{32}O_2$ F. Wt. 328.5	25 gm 100 gm	
1702	VITAMIN B1 (Thiamine hydrochloride) * CAS 67-03-8 $C_{12}H_{17}ClN_4OS \cdot HCl_2$ F. Wt. 337.27	25 gm 100 gm 500 gm	
1703	VITAMIN B2 (Riboflavin), Vitamin G CAS 83-88-5 $C_{17}H_{20}N_4O_6$ F. Wt. 376.36	10 gm 25 gm 100 gm	
759	VITAMIN B3 NIACIN (Nicotinic Acid) CAS 59-67-6 $C_6H_5NO_2$ F. Wt. 123.11	25 gm 100 gm 1 Kg	
561	VITAMIN B4 (Adenine 99%) (6-Aminopurine) CAS 73-24-5 $C_5H_5N_6$ F. Wt. 135.13	5 gm 25 gm 1 Kg	
1701	VITAMIN B5 (Calcium-D-Pantothenate) * CAS 137-08-6 $C_{18}H_{32}CaN_2O_{10}$ F. Wt. 476.53	25 gm 100 gm 1 Kg	
1705	VITAMIN B6 (Pyridoxine Hydrochloride) * CAS 58-56-0 $C_8H_{11}NO_3 \cdot HCl$ F.Wt. 205.64	10 gm 25 gm 100 gm	
1709	VITAMIN B7 Coenzyme (D-Biotin) * CAS 58-85-5 $C_{10}H_{16}N_2O_3S$ F.Wt. 244.31	1 gm 10 gm 25 gm	
604	VITAMIN B8 (Inositol) CAS 61-19-8 $C_6H_{12}O_6$ F. Wt. 180.16	25 gm 100 gm 1 Kg	
335	VITAMIN B9 (Folic Acid) * CAS 59-30-3 $C_{19}H_{19}N_7O_6$ F. Wt. 441.40	5 gm 25 gm	

CODE	PRODUCT NAME	PACK SIZE	
233	VITAMIN B10 (PABA) CAS 150-13-0 $C_7H_7NO_2$ F. Wt. 137.14	100 gm 500 gm	
1704	VITAMIN B12 (Cyanocobalamin) * CAS 68-19-9 $C_{63}H_{88}CoN_{14}O_{14}P$ F. Wt. 1355.38	1 gm 10 gm	
1721	VITAMIN BT (L-Carnitine) * CAS 541-15-1 $C_7H_{15}NO_3$ F. Wt. 161.20	5 gm 25 gm	
1706	VITAMIN C (Ascorbic Acid) CAS 50-81-7 $C_6H_8O_6$ F. Wt. 176.12	100 gm 500 gm	
1715	VITAMIN C, AR GRADE CAS 50-81-7 $C_6H_8O_6$ F. Wt. 176.13	25 gm 100 gm 500 gm	
1707	VITAMIN D3 (Cholecalciferol) * CAS 67-97-0 $C_{27}H_{44}O$ F. Wt. 384.64	1 gm	
1708	VITAMIN E POWDER * (DL-Tocopherol Acetate) 50% Powder CAS 10191-41-0 $C_{29}H_{50}O_2$ F. Wt. 430.71	50 gm	
1709	VITAMIN H (D-Biotin) * CAS 58-85-5 $C_{10}H_{16}N_2O_3S$ F. Wt. 244.31	1 gm 10 gm 25 gm	
1716	VITAMIN K3 SODIUM BISULPHITE* (Menadione Sodium Bisulphite) Water soluble CAS 130-37-0 $C_{11}H_9NaO_5S$ F. Wt. 276.24	25 gm 100 gm	
1665	WRIGHT'S STAIN (Wright's Eosin Methylene Blue) CAS 68988-92-1	25 gm 100 gm	
3004	XYLENE CYANOL FF (ACID BLUE 147) CAS 2650-17-1 $C_{25}H_{27}N_3NaO_5S_2$ F. Wt. 538.61	5 gm 25 gm	
4310	XYLENOL ORANGE TETRASODIUM SALT, PRACTICAL GRADE CAS 3618-43-7 $C_{31}H_{28}N_2O_{13}SNa$ F. Wt. 760.58	5 gm 10 gm	
4311	XYLIDINE PONCEAU CAS 3761-53-3 $C_{18}H_{14}N_2Na_2O_5S_2$ F. Wt. 480.42	25 gm	
4312	XYLIDYL BLUE I SODIUM SALT, EXTRA PURE CAS 14936-97-1 $C_{26}H_{20}N_3NaO_5S$ F. Wt. 513.50	1 gm	
3093	X-GAL (5-BROMO-4-CHLORO-3 INDOLYL-b-D-GALACTOPYRANOSIDE) * CAS 7240-90-6 $C_{14}H_{15}BrClNO_6$ F. Wt. 408.63	100 mg 500 mg 1 gm	

CODE	PRODUCT NAME	PACK SIZE	
3095	X-GLUCOSIDE * CAS 15548-60-4 $C_{14}H_{15}BrClNO_6$ F. Wt. 408.63	100 mg	
4146	XANTHAN GUM, EXTRA PURE CAS 11138-66-2	100 gm 500 gm 25 Kg	
1890	XANTHAN GUM, AR GRADE CAS 11138-66-2	100 gm 500 gm 25 Kg	
3237	XANTHINE (2,4 Dihydroxypentane) CAS 69-89-6 $C_5H_8N_4O_2$ F. Wt. 152.11	5 gm 25 gm	
2772	XYLENOL ORANGE, AR GRADE CAS 1611-35-4 $C_{31}H_{32}N_2O_{13}S$ F. Wt. 672.67	5 gm 10 gm	
1312	D (+) XYLOSE CAS 58-86-6 $C_5H_{10}O_5$ F. Wt. 150.13	25 gm 100 gm 500 gm	
3515	YEAST EXTRACT PASTE for bacteriology	500 gm	
4204	YEAST EXTRACT POWDER, BACTO GRADE for general purpose for bacteriology	500 gm	
4132	YTTRIUM OXIDE 99.9% AR CAS 1314-36-9 Y_2O_3 F. Wt. 225.81	10 gm 25 gm	
516	ZEATIN (Plant Growth Regulator) * CAS 1637-39-4 $C_{10}H_{13}N_5O$ F. Wt. 219.25	100 mg 500 mg 1 gm	
4313	TRANS-ZEATIN RIBOSIDE CAS 6025-53-2 $C_{15}H_{21}N_5O_5$ F. Wt. 351.36	10 mg	
493	ZINC ACETATE (Dihydrate), EXTRA PURE CAS 5970-45-6 $C_4H_6O_4Zn.2H_2O$ F. Wt. 219.51	500 gm 50 Kg	
494	ZINC ACETATE (Dihydrate), AR GRADE CAS 5970-45-6 $C_4H_6O_4Zn.2H_2O$ F. Wt. 219.51	500 gm 50 Kg	
4314	ZINC BORATE, EXTRA PURE CAS 1332-07-6 $2ZnO.3B_2O_3.3.5H_2O$ F. Wt. 434.62	1 Kg	
677	ZINC CARBONATE BASIC CAS 5970-47-8 $3Zn(OH)_2.2ZnCO_3$ F. Wt. 549.01	500 gm	
2369	ZINC CHLORIDE (ANH.), EXTRA PURE CAS 7646-85-7 $ZnCl_2$ F. Wt. 136.30	500 gm 50 Kg	
4315	ZINC CHLORIDE (ANH) Meets USP/EP/JP/BP CAS 7646-85-7 $ZnCl_2$ F. Wt. 136.30	100 gm 500 gm 1 Kg	
2371	ZINC METAL (Granules) CAS 7440-66-6 Zn F. Wt. 65.37	500 gm 50 Kg	

Chemicals

Laboratory Chemicals (EP & AR Grades)

CODE	PRODUCT NAME	PACK SIZE	
700	ZINC METAL (Powder) (325 Mesh) CAS 7440-66-6 Zn F. Wt. 65.38	500 gm 50 Kg	
4133	ZINC (METAL) DUST 325 MESH, EXTRA PURE CAS 7440-66-6 Zn F. Wt. 65.38	500 gm	
4134	ZINC (METAL) DUST 325 MESH, AR GRADE CAS 7440-66-6 Zn F. Wt. 65.38	500 gm	
3507	ZINC CHLORIDE DRY, AR GRADE CAS 7646-85-7 ZnCl ₂ F. Wt. 136.28	500 gm	
4316	ZINC CITRATE TRIBASIC (Dihydrate) CAS 5990-32-9 C ₁₂ H ₁₀ O ₁₄ Zn ₃ ·2H ₂ O F. Wt. 610.4	500 gm	
4317	ZINC DIBENZYLTHIOCARBAMATE, AR GRADE CAS 14726-36-4 C ₃₀ H ₂₈ N ₂ S ₄ Zn F. Wt. 610.19	25 gm	
4318	ZINC GLUCONATE CAS 4468-02-4 C ₁₂ H ₂₂ O ₁₄ Zn F. Wt. 455.68	500 gm	
4136	ZINC NITRATE (Hexa), EXTRA PURE CAS 10196-18-6 Zn(NO ₃) ₂ ·6H ₂ O F. Wt. 297.49	500 gm	
4319	ZINC OXALATE (Hydra), EXTRA PURE CAS 547-68-2 ZnC ₂ O ₄ ·xH ₂ O F. Wt. 153.41	500 gm	

CODE	PRODUCT NAME	PACK SIZE	
495	ZINC OXIDE, EXTRA PURE CAS 1314-13-2 ZnO F. Wt. 81.38	500 gm 5 Kg 50 Kg	
496	ZINC OXIDE, AR GRADE CAS 1314-13-2 ZnO F. Wt. 81.38	500 gm	
2785	ZINC PICOLINATE CAS 17949-65-4 C ₁₂ H ₈ N ₂ O ₄ Zn ₂ F. Wt. 309.61	100 gm	
2445	ZINCON SODIUM, AR GRADE CAS 62625-22-3 C ₂₀ H ₁₅ N ₄ NaO ₆ S F. Wt. 480.43	1 gm 5 gm	
2777	ZINC STEARATE CAS 557-05-21 2(C ₁₈ H ₃₅ O ₂) ₂ ·Zn F. Wt. 632.33	500 gm	
497	ZINC SULPHATE (Hepta), EXTRA PURE CAS 7446-20-0 ZnSO ₄ ·7H ₂ O F. Wt. 287.55	500 gm 5 Kg 50 Kg	
3508	ZINC SULPHATE (Hepta), AR GRADE CAS 7446-20-0 ZnSO ₄ ·7H ₂ O F. Wt. 287.55	500 gm 5 Kg 50 Kg	
499	ZINC SULPHATE (Mono), PHARMA GRADE CAS 7446-19-7 ZnSO ₄ ·H ₂ O F. Wt. 179.48	500 gm 5 Kg 50 Kg	
4001	ZINC SULPHATE (Mono), AR GRADE CAS 7446-19-7 ZnSO ₄ ·H ₂ O F. Wt. 179.48	500 gm 5 Kg	

CODE	PRODUCT NAME	PACK SIZE	
569	ACETIC ACID, Glacial CAS 64-19-7 CH_3COOH F. Wt. 60.05	500 ml 2.5 Ltr 5 Ltr 25 Ltr	
568	ACETONE CAS 67-64-1 $(\text{CH}_3)_2\text{CO}$ F. Wt. 58.08	500 ml 2.5 Ltr 25 Ltr	
049	AMMONIA SOLUTION 25% CAS 1336-21-6 NH_3 F. Wt. 17.03	500 ml 2.5 Ltr 5 Ltr 25 Ltr	
133	BRIJ 35 (Polyoxyethylene Lauryl ether) CAS 9002-92-0 $\text{C}_{12}\text{H}_{25}\text{O}(\text{C}_2\text{H}_4\text{O})_n$	500 ml	
023	BENZALKONIUM CHLORIDE 50% CAS 63449-41-2	500 ml 5 Ltr 25 Ltr	
4169	BENZALKONIUM CHLORIDE 80% CAS 63449-41-2	25 Ltr 50 Ltr	
1861	ETHYL ACETATE CAS 141-78-6 $\text{CH}_3\text{COO}_2\text{H}_5$ F. Wt. 88.11	500 ml 2.5 Ltr 25 Ltr	
718	FORMALDEHYDE SOLUTION 37-41% CAS 50-00-0 CH_2O F. Wt. 30.03	500 ml 2.5 Ltr 5 Ltr 30 Ltr	
353	ISO PROPYL ALCOHOL CAS 67-63-0 $\text{C}_3\text{H}_8\text{O}$ F. Wt. 60.10	500 ml 2.5 Ltr 5 Ltr 25 Ltr	
2163	OLEIC ACID (Elainic Acid), EXTRA PURE (For Biochemistry) * CAS 112-80-1 $\text{C}_{18}\text{H}_{34}\text{O}_2$ F. Wt. 282.46	500 ml	
876	ORTHO-PHOSPHORIC ACID CAS 7664-38-2 H_3PO_4 F. Wt. 98.00	500 ml 2.5 Ltr 5 Ltr 25 Ltr	
1940	LACTIC ACID (Hydroxy Propionic Acid), EXTRA PURE CAS 50-21-5 $\text{C}_3\text{H}_6\text{O}_3$ F. Wt. 90.08	500 ml 2.5 Ltr 25 Ltr	
3200	POTASSIUM LACTATE 50% CAS 996-31-6 $\text{C}_3\text{H}_5\text{KO}_3$ F. Wt. 128.17	500 ml	

CODE	PRODUCT NAME	PACK SIZE	
4143	PROPIONIC ACID CAS 79-09-4 $\text{C}_3\text{H}_6\text{O}_2$ F. Wt. 74.08	500 ml 2.5 Ltr 25 Ltr	
2401	SODIUM HYPOCHLORITE 4% CAS 7681-52-9 NaClO F. Wt. 74.44	500 ml 1 Ltr 5 Ltr 25 Ltr	
2294	SODIUM LACTATE 60% CAS 72-17-3 $\text{C}_3\text{H}_5\text{NaO}_3$ F. Wt. 112.06	500 ml	
696	SOLUTION FOR LEISHMAN STAIN (S.F.L) CAS 67-56-1 CH_4O F. Wt. 32.04	500 ml 2.5 Ltr 25 Ltr	
665	TABOLENE GL (W/ Dispenser) for Glass & Computer Equipment Cleaning	500 ml 5 Ltr 50 Ltr	
666	TABOLENE GLW (DETERGENT TABOLENE) (W/ Neutral pH Phosphate Free Detergent) for Glassware Cleaning	500 ml 5 Ltr 50 Ltr	
667	TABOLENE GP General Purpose Detergent Lab Floor Cleaning, Like Teepol	500 ml 5 Ltr 50 Ltr	
2743	TERPINEOL (Mixture of Isomers), (PHARMA GRADE) CAS 8006-39-1 $\text{C}_{10}\text{H}_{18}\text{O}$ F. Wt. 154.25	500 ml 2.5 Ltr	
2354	TRIETHANOLAMINE 85% CAS 102-71-6 $\text{C}_6\text{H}_{15}\text{NO}_3$ F. Wt. 149.19	500 ml 5 Ltr 50 Ltr	
2356	TRIETHYLAMINE 99% CAS 121-44-8 $\text{C}_6\text{H}_{15}\text{NO}_3$ F. Wt. 101.19	500 ml 2.5 Ltr 25 Ltr	
670	TWEEN 20® (Regd. Trademark ICI UK) CAS 9005-64-5 $\text{C}_{58}\text{H}_{113}\text{O}_{26}$ F. Wt. 1226.5	500 ml 25 Ltr	
673	TWEEN 80® (Regd. Trademark ICI UK) CAS 9005-65-6 $\text{C}_{52}\text{H}_{99}\text{O}_{21}$ F. Wt. 604.8	100 ml 500 ml 25 Ltr	



Food & Water Chemicals

Tailored to meet the strict requirements set of the **Food and Beverage** industry, providing you with the assurance of quality, safety, and regulatory compliance.

Food Grades Chemicals

Enhancing food processing and production

- Preservatives
- Flavor Enhancers
- Antioxidants
- Emulsifiers
- Acidity Regulators

Water Industry Minerals

Elevate drinking water quality with our premium minerals, including:

- Calcium
- Potassium
- Magnesium
- Sodium



Certified



Certified



Certified



13485:2016



11133:2014



9001:2015

CODE	PRODUCT NAME	PACK SIZE	
5000	ACETIC ACID	25 Ltr	
5100	ACETIC ACID GLACIAL	25 Ltr	
5122	ACETIC ACID (9000)	25 Ltr	
5001	ACETONE	25 Ltr	
TBF 095	AGAR AGAR (TYPE-I)	25 Kg	
TBF 104	AGAR AGAR (TYPE-II)	25 Kg	
5002	AMMONIUM BICARBONATE, (Ammonium Hydrogen Carbonate)	50 Kg	
5003	AMMONIUM CHLORIDE	50 Kg	
5132	tri-AMMONIUM CITRATE	25 Kg	
5123	AMMONIUM DIHYDROGEN ORTHOPHOSPHATE	50 Kg	
5021	di-AMMONIUM HYDROGEN ORTHOPHOSPHATE (DAP)	50 Kg	
5004	AMMONIUM PERSULPHATE (For Distillery)	50 Kg	
5139	L-ARABINOSE	25 Kg	
TBF 137	CALCIUM LIQUID (VICTORY NA)	50 Kg	
5088	CALCIUM PEROXIDE	25 Kg	
5118	MONO CALCIUM PHOSPHATE (Monobasic) (MCP)	25 Kg	
5022	di-CALCIUM PHOSPHATE (Dihydrate) (DCP)	25 Kg	
5015	tri-CALCIUM PHOSPHATE (ANH.)	25 Kg	
5016	CALCIUM PHOSPHATE (TCP)	25 Kg	
5017	CALCIUM PROPIONATE	25 Kg	
TBF 128	CALCIUM PROPIONATE POWDER (VICTORY CAP)	25 Kg	
5109	CALCIUM SULPHATE (ANH.)	50 Kg	
5086	CALCIUM SULPHATE (Dihydrate)	50 Kg	
5019	CITRIC ACID (ANH.)	25 Kg	
5020	CITRIC ACID (Mono)	25 Kg	

CODE	PRODUCT NAME	PACK SIZE	
5113	L-ASCORBIC ACID	25 Kg	
5005	BENZOIC ACID	25 Kg	
5006	BUTYLATED HYDROXY ANISOLE, (B.H.A.)	25 Kg	
5007	BUTYLATED HYDROXY TOLUENE, (B.H.T.)	25 Kg	
5008	CAFFEINE (ANH.) (Food Grade)	25 Kg	
5009	CALCIUM ACETATE (Dried) (For Atorvastatin)	25 Kg	
5112	CALCIUM BUTYRATE	25 Kg	
5010	CALCIUM CARBONATE (Precipitated)	50 Kg	
5011	CALCIUM CHLORIDE (ANH.)	50 Kg	
5096	CALCIUM CHLORIDE (Dihydrate)	50 Kg	
5012	CALCIUM CITRATE (Tetra)	50 Kg	
TBF 155	CALCIUM DIACETATE	25 Kg	
5124	CALCIUM GLUCONATE (Mono)	50 Kg	
5013	CALCIUM HYDROXIDE (Purified)	50 Kg	
5014	CALCIUM LACTATE POWDER (Penta)	25 Kg	
5138	MANNOSE	25 Kg	
5127	METHYL PARABEN	10 Kg 25 Kg	
5035	METHYL PARABEN SODIUM	10 Kg 25 Kg	
5128	MICROSIL (Food Grade Sterilizer)	25 Ltr	
5036	NIACIN (Nicotinic Acid) (Vitamin B)	25 Kg	
5037	NIACINAMIDE (Nicotinamide) (Vitamin PP)	25 Kg	
TBF 224	NISIN	25 Kg	
5038	n-PROPYL GALLATE	25 Kg	
5039	ortho PHOSPHORIC ACID	25 Ltr	
5040	POLYVINYL PYRROLIDINE K-25 (PVP)	25 Kg	
5041	POLYVINYL PYRROLIDINE K-30 (PVP)	25 Kg	
5042	POTASSIUM ACETATE	25 Kg	
5101	POTASSIUM BENZOATE	25 Ltr	

Our all food products are produced under FSSC and FSSAI approved facilities

For more Laboratory Chemicals, please turn to page no. 280-326

CODE	PRODUCT NAME	PACK SIZE	
5098	COPPER SULPHATE (Penta)	25 Kg	
5125	DEGREASER	25 Kg	
TBF 016	DRIED YEAST INACTIVE (for Biscuits)	25 Kg	
5024	EDTA DISODIUM	25 Kg	
5025	FERRIC CHLORIDE (ANH.)	50 Kg	
5026	FERROUS SULPHATE (Hepta)	50 Kg	
5027	FOLIC ACID (Vitamin B)	25 Kg	
5137	GALACTOSE	25 Kg	
5028	GLYCEROL (Glycerine)	25 Ltr	
TBF 157	GELLAN GUM (Uses Milk, Juice Thicker & Smoother)	25 Kg	
TBF 158	GLYCINE (Amino Acetic Acid)	25 Kg	
5126	GLUCO DELTA LACTONE (For Cheese)	25 Kg	
5133	HYDROGEN PEROXIDE 35%	25 Kg	
5099	ISOPROPYL ALCOHOL	25 Ltr	
5030	LACTIC ACID (Hydroxy Propionic Acid)	25 Ltr	
5143	LACTIC ACID (Hydroxy Propionic Acid) (88%) (FOOD GRADE)	25 Ltr	
5031	MAGNESIUM CHLORIDE (Hexa)	50 Kg	
5032	MAGNESIUM SULPHATE (Hepta)	50 Kg	
5134	MANGANESE CHLORIDE	25 Kg	
5033	DL-MALIC ACID	25 Kg	
3527	MALT EXTRACT PASTE (Food Grade)	50 Kg	
3506	MALT EXTRACT POWDER (Food Grade)	25 Kg	
5034	MANGANESE (II) SULPHATE (Mono)	50 Kg	

CODE	PRODUCT NAME	PACK SIZE	
5043	POTASSIUM BICARBONATE	50 Kg	
5090	POTASSIUM BROMATE	50 Kg	
5044	POTASSIUM CARBONATE (ANH.)	50 Kg	
5045	POTASSIUM CHLORIDE	50 Kg	
5079	tri-POTASSIUM CITRATE (Mono)	25 Kg	
5046	POTASSIUM DIHYDROGEN ORTHOPHOSPHATE (Mono Potassium Phosphate) (ANH.)	50 Kg	
5141	POTASSIUM FLUORIDE	25 Kg	
5063	di-POTASSIUM HYDROGEN ORTHOPHOSPHATE (ANH.)	50 Kg	
5047	POTASSIUM HYDROXIDE FLAKES	50 Kg	
5048	POTASSIUM HYDROXIDE PELLETS	50 Kg	
5135	POTASSIUM HYDROXIDE POWDER	25 Kg	
5049	POTASSIUM IODATE	25 Kg	
5050	POTASSIUM IODIDE	25 Kg	
5051	POTASSIUM METABISULPHITE	50 Kg	
5052	POTASSIUM NITRATE	50 Kg	
5091	POTASSIUM PERSULPHATE	50 Kg	
5107	POTASSIUM PROPIONATE	25 Kg	
5092	POTASSIUM SODIUM TARTRATE	25 Kg	
TBF 197	POTASSIUM SORBATE (Granular)	25 Kg	
5054	POTASSIUM SULPHATE	50 Kg	
5055	PROPYLENE GLYCOL (1,2-Propanediol)	25 Ltr	
5103	PROPYLENE GLYCOL ALGINATE (PGA)	25 Kg	
5056	PROPYL PARABEN	25 Kg	
5057	PROPYL PARABEN SODIUM	25 Kg	

Our all food products are produced under FSSC and FSSAI approved facilities

For more Laboratory Chemicals, please turn to page no. 280-326

Qwikzee

The Smarter Way to Perfect Paneer!

A synergistic blend of Organic Acidulants

Benefits

- Highly soluble
- Retains fat & protein content
- Mild acidulants with a rapid drop in pH
- Helps spoilage control
- Provides additional shelf life
- Improves texture
- Maintains milk's natural sweetness
- Maintains softness
- Maintains sensory characteristics
- 5% Yield increase
- Long-term profitability



Chemicals

Food Grade Chemicals

CODE	PRODUCT NAME	PACK SIZE
5148	QWIKZEE (Paneer Coagulant)	1 Kg*20 25 Kg
5058	SODIUM ACETATE (ANH.)	50 Kg
5059	SODIUM ACETATE (Trihydrate)	50 Kg
5060	SODIUM ACID PYROPHOSPHATE (SAPP)	50 Kg
5114	SODIUM ALGINATE	50 Kg
5149	SODIUM ALGINATE HQ	25 Kg
5102	SODIUM ASCORBATE	25 Kg
5061	SODIUM BENZOATE	25 Kg
5062	SODIUM BICARBONATE	50 Kg
5064	SODIUM BUTYRATE	25 Kg
5065	SODIUM CARBONATE (ANH.)	50 Kg
5066	SODIUM CARBONATE (Mono)	50 Kg
TBF 105	SODIUM CASEINATE (Protein 90%)	25 Kg
5067	SODIUM CHLORIDE	50 Kg
5068	SODIUM CITRATE DIBASIC (di-Sodium Hydrogen Citrate)	50 Kg
5080	tri-SODIUM CITRATE (ANH.)	50 Kg
5115	tri-SODIUM CITRATE (Dihydrate)	50 Kg
TBF 006	SODIUM DIACETATE (Dry Vinegar)	25 Kg
5069	SODIUM DIHYDROGEN ORTHOPHOSPHATE (Mono Sodium Phosphate Dihydrate)	50 Kg
5142	SODIUM FLUORIDE	25 Kg
5023	di-SODIUM HYDROGEN ORTHOPHOSPHATE (Dihydrate)	50 Kg
5119	di-SODIUM HYDROGEN ORTHOPHOSPHATE (ANH.)	50 Kg
5085	SODIUM HYDROSULPHITE	50 Kg

CODE	PRODUCT NAME	PACK SIZE
5104	SODIUM HYPOCHLORITE 4%	25 Ltr
5070	SODIUM HYDROXIDE FLAKES	50 Kg
5071	SODIUM HYDROXIDE PELLETS	50 Kg
5117	SODIUM LAURYL SULPHATE	25 Kg
TBF 066	SOYA LECITHIN POWDER	25 Kg
5072	SODIUM METABISULPHITE	50 Kg
5073	SODIUM NITRATE	50 Kg
5074	SODIUM NITRITE	50 Kg
5121	SODIUM OLEATE	25 Kg
5075	SODIUM POTASSIUM TARTRATE (Tetra)	25 Kg
TBF 127	SODIUM PROPIONATE POWDER (VICTORY NAP)	25 Kg
5120	SODIUM SULPHATE (ANH.)	50 Kg
5076	SODIUM SULPHITE (ANH.)	50 Kg
5136	SODIUM TRIPOLYPHOSPHATE	25 Kg
TBF 2732	SORBITOL LIQUID	25 Kg
5110	SORBITOL LIQUID	25 Kg
5093	SORBITOL POWDER	25 Kg
5078	TITANIUM DIOXIDE	25 Kg
5082	TWEEN 20 (POLYSORBATE 20)	50 Kg
5083	TWEEN 80 (POLYSORBATE 80)	50 Kg
5108	UREA	50 Kg
5094	XANTHAN GUM	25 Kg
5129	YEAST EXTRACT	25 Kg
5130	ZINC CARBOANTE	50 Kg
5131	ZINC CHLORIDE	50 Kg
5084	ZINC SULPHATE (Hepta)	50 Kg
5105	ZINC SULPHATE (Mono)	50 Kg

Our all food products are produced under FSSC and FSSAI approved facilities

For more Laboratory Chemicals, please turn to page no. 280-326

Troubleshooting Guide

Issue	Possible Cause	Solutions
Dehydrated Culture Media		
Clumping or Agglomeration	<ul style="list-style-type: none"> Moisture exposure or improper storage 	<ul style="list-style-type: none"> Discard any clumps or agglomerated media Ensure proper storage in a cool and dry place Use fresh dehydrated media
Media Not Dissolving Properly	<ul style="list-style-type: none"> Insufficient mixing or improper water temperature 	<ul style="list-style-type: none"> Follow the instructions for media preparation Mix the media and water thoroughly using a sterile magnetic stirrer or vortex mixer Check the water temperature, as some media may require specific temperature conditions for dissolution
Media Contamination	<ul style="list-style-type: none"> Contaminated equipment or improper aseptic technique 	<ul style="list-style-type: none"> Sterilize all equipment, including containers, measuring tools, and water used for media preparation Practice proper aseptic technique, including working in a laminar flow hood or sterile environment Monitor the media for signs of contamination, such as discoloration, growth, or off-odours. Discard contaminated media
Inadequate Growth Results	<ul style="list-style-type: none"> Repeated Remelting pH Drift Storage at incorrect temperature 	<ul style="list-style-type: none"> Avoid repeated remelting Check the pH of distilled water prior to media preparation Store according to the instructions on product TD
Colour Issue in Prepared Media	<ul style="list-style-type: none"> Overheating 	<ul style="list-style-type: none"> Avoid overheating as it causes accumulation of carbohydrate causing the medium to become darker
Discoloration of Media	<ul style="list-style-type: none"> Repeated remelting 	<ul style="list-style-type: none"> Avoid repeated remelting, only reconstitute the amount needed
Improper Gelling	<ul style="list-style-type: none"> Low pH Excessive heating of the media 	<ul style="list-style-type: none"> Check the pH of the water before dissolution Do not over sterilize
Darkening	<ul style="list-style-type: none"> Overheating Incorrect weighing of the medium, excess amount of dehydrated powder Incorrect pH Repeated remelting 	<ul style="list-style-type: none"> Do not overheat Weigh the correct quantity Check pH Avoid remelting
Ready-to-Use Culture Media Plates		
Dry or Cracked Media Surface	<ul style="list-style-type: none"> Insufficient moisture during storage or drying out of the media 	<ul style="list-style-type: none"> Ensure proper storage conditions, including temperature and humidity Store plates in a sealed bag or container to maintain moisture If the plates are already dry or cracked, discard them and use fresh ones
Contamination on Media Surface	<ul style="list-style-type: none"> Improper storage Handling Aseptic technique 	<ul style="list-style-type: none"> Store plates in a sterile environment and avoid exposure to contaminants Follow aseptic techniques during handling and avoid touching the media surface Monitor plates for signs of contamination, and discard any contaminated plates
Inadequate Growth Results	<ul style="list-style-type: none"> Improper storage Handling Aseptic technique 	<ul style="list-style-type: none"> Store plates in a sterile environment and avoid exposure to contaminants Follow aseptic techniques during handling and avoid touching the media surface
Darkening of Medium on Storage	<ul style="list-style-type: none"> If media contains high concentration of sugars and exposed to continual heating or stored at high temperatures 	<ul style="list-style-type: none"> Do not overheat Avoid exposure to high temperatures for long durations

Issue	Possible Cause	Solutions
Ready-to-Use Culture Media Bottles		
Media pH Deviation	<ul style="list-style-type: none"> Contamination Exposure to air Improper storage 	<ul style="list-style-type: none"> Ensure the bottles are properly sealed and stored under appropriate conditions Check the expiration date and use media bottles within the specified period Perform pH testing of the media before use, and discard any bottles with deviated pH
Mold Growth or Contamination	<ul style="list-style-type: none"> Improper storage or contamination during handling 	<ul style="list-style-type: none"> Store media bottles in a cool and dry place, away from direct light Follow aseptic techniques during handling and avoid introducing contaminants Inspect the media visually before use, and discard any bottles showing signs of contamination
Inadequate Growth Results	<ul style="list-style-type: none"> Improper storage Handling Aseptic technique 	<ul style="list-style-type: none"> Store media in a sterile environment and avoid exposure to contaminants Follow aseptic techniques during handling and avoid touching the media surface
Darkening of Medium on Storage	<ul style="list-style-type: none"> If media contains high concentration of sugars and exposed to continual heating or stored at high temperatures 	<ul style="list-style-type: none"> Do not overheat Avoid exposure to high temperatures for long durations
Biological Media Bases		
Precipitation or Crystallization	<ul style="list-style-type: none"> Improper dissolution or reaction with other components 	<ul style="list-style-type: none"> Follow the instructions for media preparation Ensure thorough mixing and proper dissolution of the biological media base If precipitation or crystallization occurs, discard the media and prepare a fresh batch
Growth Inhibition or Poor Performance	<ul style="list-style-type: none"> Contaminated components, incorrect ratios, or expired materials 	<ul style="list-style-type: none"> Follow the instructions for media preparation Mix the media and water thoroughly using a sterile magnetic stirrer or vortex mixer Check the water temperature, as some media may require specific temperature conditions for dissolution
Incomplete Solubility	<ul style="list-style-type: none"> Insufficient mixing or improper water temperature 	<ul style="list-style-type: none"> Follow the instructions Mix the media base and water thoroughly using a sterile magnetic stirrer or vortex mixer Check the water temperature, as some media base may require specific temperature conditions for dissolution
Antibiotic Sensitivity Discs		
No Zones of Inhibition	<ul style="list-style-type: none"> Inoculum too heavy Incorrect antibiotic disc (e.g., wrong concentration, expired) Resistant organism Incorrect media (e.g. pH incorrect) Antibiotic degradation (e.g., improper storage, exposure to heat) Test performed incorrectly (e.g., discs not properly adhered, incubation conditions wrong) Contamination of the disc. 	<ul style="list-style-type: none"> Verify disc expiry date, antibiotic type, and concentration Perform additional tests to confirm resistance Correct agar pH Store discs as per instructions (usually refrigerated or frozen) Follow CLSI guidelines for disc placement, incubation time, and temperature Use sterile technique

Troubleshooting Guide

Issue	Possible Cause	Solutions
Small or Indistinct Zones of Inhibition	<ul style="list-style-type: none"> Inoculum too light Antibiotic disc potency reduced Reading zones prematurely Some slow growing organisms Some antibiotics diffuse slowly 	<ul style="list-style-type: none"> Verify disc expiry date and storage conditions Verify incubator temperature (usually 35°C) Read zones after the recommended incubation time (usually 16-18 hours) Extend the incubation time if the organism is known to be slow-growing Consult CLSI guidelines for specific antibiotic diffusion characteristics
Mold Growth or Contamination	<ul style="list-style-type: none"> Improper storage or contamination during handling 	<ul style="list-style-type: none"> Store media bottles in a cool and dry place, away from direct light Follow aseptic techniques during handling and avoid introducing contaminants Inspect the media visually before use, and discard any bottles showing signs of contamination
Large, Irregular Zones of Inhibition	<ul style="list-style-type: none"> Confluent growth (too many colonies) Moisture on the agar surface Discs placed too close together Spread of swarming organisms (e.g., <i>Proteus</i>) Some antibiotics may produce irregular zones naturally 	<ul style="list-style-type: none"> Ensure a lawn of confluent growth, not individual colonies Dry agar plates before inoculation Space discs according to CLSI guidelines Use appropriate swarming inhibition techniques Consult CLSI guidelines for specific antibiotic zone interpretations
Satellite Colonies within the Zone of Inhibition	<ul style="list-style-type: none"> Beta-lactamase production by the organism Mixed culture contamination Resistant mutants emerging 	<ul style="list-style-type: none"> Perform beta-lactamase testing Repeat the test with a pure culture Repeat the test and observe for consistent satellite colony growth
Inconsistent Zone Sizes Between Tests	<ul style="list-style-type: none"> Variations in inoculum density Variations in media preparation Variations in disc placement Variations in incubation conditions Variations in reading technique Disc degradation 	<ul style="list-style-type: none"> Use consistent media preparation procedures Use a disc dispenser to ensure consistent placement Monitor and maintain consistent incubation conditions Train personnel on proper zone reading techniques Verify disc expiration date and storage
Presence of Bubbles in the Agar	<ul style="list-style-type: none"> Improper pouring of the agar Aerobic bacteria producing gas 	<ul style="list-style-type: none"> Pour agar carefully to avoid bubbles Identify and address the gas-producing bacteria
Contamination of the Agar Plate	<ul style="list-style-type: none"> Non-sterile technique during preparation or inoculation Contaminated media Contaminated discs Airborne contamination 	<ul style="list-style-type: none"> Use aseptic technique throughout the process Prepare fresh media or use quality-controlled media Store discs as per instructions (usually refrigerated or frozen) Work in a clean environment or biosafety cabinet
Incorrect Reading of the Zone Size	<ul style="list-style-type: none"> Subjectivity of measurement Incorrect interpretation of zone endpoints Parallax error 	<ul style="list-style-type: none"> Use a calibrated ruler or automated zone reader Refer to CLSI guidelines for zone interpretation Read zones from directly above the plate
Staining Kits		
Weak or Faded Staining	<ul style="list-style-type: none"> Expired reagents Insufficient staining time Excessive washing Thick smear Fixation issues Incorrect pH Old cultures 	<ul style="list-style-type: none"> Check expiration dates Follow staining times Gentle washing Thinner smears Proper fixation Verify pH Use fresh cultures

Issue	Possible Cause	Solutions
Uneven Staining	<ul style="list-style-type: none"> • Uneven smear • Uneven reagent application • Air drying • Lipids in sample 	<ul style="list-style-type: none"> • Uniform smears • Even reagent application • Prevent drying • Degreasing step if necessary
Gram-Variable Staining	<ul style="list-style-type: none"> • Over-decolorization • Old cultures • Thick smears • Inadequate fixation • Natural variability 	<ul style="list-style-type: none"> • Careful decolorization • Fresh cultures • Thinner smears • Proper fixation • Consult reference material
Precipitate on Slide	<ul style="list-style-type: none"> • Contaminated reagents • Inadequate filtering • Dried reagents • Hard water 	<ul style="list-style-type: none"> • Filter reagents • Replace contaminated reagents • Thorough rinsing • Use distilled/deionized water
No Staining (Clear Smear)	<ul style="list-style-type: none"> • No bacterial smear • Excessive washing • Fixation failure • Missed steps 	<ul style="list-style-type: none"> • Ensure smear is present • Gentle washing • Repeat fixation/staining • Follow protocol
Background Staining	<ul style="list-style-type: none"> • Inadequate washing • Contaminated reagents • Excessive staining time • Dirty slides 	<ul style="list-style-type: none"> • Thorough washing • Replace contaminated reagents • Reduce staining times • Use clean slides
Water Microbiology Solutions		
No Color Change / No Reaction	<ul style="list-style-type: none"> • Expired or improperly stored reagents • Incorrect test procedure • Insufficient incubation time or incorrect temperature 	<ul style="list-style-type: none"> • Check the expiry date and ensure reagents are stored as per recommended conditions • Follow the test protocol carefully as per the kit instructions • Ensure proper incubation time and maintain recommended temperature conditions
Weak or Faint Color Development	<ul style="list-style-type: none"> • Low concentration of microorganisms • Insufficient sample volume • Improper mixing of reagents 	<ul style="list-style-type: none"> • Use the recommended sample volume for testing • Mix reagents thoroughly as per instructions • Allow full incubation time before interpreting results
Unexpected Positive / Negative Result	<ul style="list-style-type: none"> • Sample contamination during collection or handling • Use of non-sterile sampling containers • Improper incubation or storage conditions • Presence of disinfectants (e.g., chlorine) in the water sample 	<ul style="list-style-type: none"> • Use sterile containers and avoid touching inner surfaces • Follow proper aseptic handling during testing • Maintain recommended incubation temperature and time • Neutralize disinfectants using suitable reagents, if applicable
Turbidity or Sediment in Test Sample	<ul style="list-style-type: none"> • Highly turbid water sample • Particulate matter interfering with results 	<ul style="list-style-type: none"> • Allow the sample to settle before testing, if permitted by the method • Follow pre-treatment steps as specified in the kit instructions
Inconsistent Results Between Tests	<ul style="list-style-type: none"> • Variation in sample collection or handling • Inconsistent incubation conditions • Reagent degradation due to improper storage 	<ul style="list-style-type: none"> • Standardize sample collection and handling procedures • Maintain uniform incubation conditions for all tests • Store reagents in a cool, dry place away from direct sunlight

TM MEDIA



TITAN BIOTECH LTD.

Manufacturer & Exporter, Since 1992

Corporate Office: 903-909, 9th Floor, Big Jo's Tower, Netaji Subhash Place, Delhi-110034

www.tmmedia.in



marketing@titanbiotechltd.com

